

## 5. Modulation of autoISF aggressiveness. V 2.3

**Please note that with autoISF 3.0 you are in an early-dev. environment,** where the user interface is **not optimized for safety** of users who stray away from intended ways to use. Good safety features exist, but these are only as good as the development-oriented user understands and implements them. This is not a medical product, refer to disclaimer in [section 0](#)



Once the initial tuning according to [section 4](#), is done, you are ready to use autoISF for your fully automated meal management.

You will have three major *other* challenges to manage:

- recognize and manage (partial) occlusions, or other technical (CGM or BT related) obstacles (see [section 2](#) on pre-requisites of FCL, and related case studies)
- deal with times when insulin given by the loop must be restricted (e.g. a snack could be “misinterpreted” as a meal)
- deal with times when the loop should be set “milder” as a precaution (e.g. nights; or in an exercise context).

How big the remaining challenge is depends very much on your individual lifestyle. [Sections 5](#) and [6](#) discuss this in more detail.

**In order to run the loop fully automatically around the clock, the times outside the meal blocks must also be precisely analyzed, and solutions to problems (if any) must be sought.**

It is up to every user to decide where to draw the line.

- With a technically well functioning system, moderate meals, moderate or no exercise, moderate %TIR expectations and a bit of mindfulness it should be possible to go into Full Closed Loop 24/7, after working through, and observing, [sections 1-4](#).
- Especially if you are a bit shy of using the emulator for really detailed analysis, it is likely that you will not hit *one* real good system calibration ([section 4](#)) for your *entire range* of diets.

In that case you will occasionally run out of range, and your options to prevent, react, or improve are

- 35           ○ accepting a few % higher time outside range for that day (and, if feasible, in the
- 36           future avoiding what seemed to have caused it)
- 37           ○ taking a snack (whenever you tend to go low from the “tails” of insulin activity that
- 38           was required to fight a peak)
- 39           ○ doing a manual override (if you can think of one in time, to manage the problem
- 40           manually)
- 41           ○ temporarily resorting to the well-known hybrid closed loop.

42 Instead of accepting such instances, you could launch “improvement projects”

- 43       • that refine your initial tuning ([section 4.](#))
- 44       • that make you and your FCL loop fit to manage an increasing number of disturbances
- 45       either automatically, or via a user intervention ([sections 5-6](#)).

46 To tailor the loop’s response to disturbances *other-than* your major meals probably will require

47 specific **modulation of the aggressiveness** that you set according to [section 4](#) for your meal

48 management.

49 There are many avenues to achieve this. The main ones, that are also easy accessible via

50 Automations in AAPS, are:

- 51       • temporary shut-off SMBs (odd-numbered target)
- 52       • temporary reduce bgAccel\_ISF-weight
- 53       • temporary reduce iobTH
- 54       • temporary reduce set %profile
- 55       • temporary set higher TT (especially in connection with exercise mode)t

56

57 In setting up your FCL, you now have another difficult and time-consuming job at hand, to define

58 solutions for any of your „other“ situations (outside of meal management) that tend to drive glucose

59 outside of the desirable range.

- 60
- 61       • In [section 5.1](#) we explore avenues towards fully automated management that in daily life
- 62       will require no user intervention at all.
- 63       • In [section 5.2](#) and [5.3](#) we will look at solutions that involve an easy user interaction like a
- 64       data entry or button push.

65  
66  
67  
68  
69  
70  
71  
72  
73  
74  
75  
76  
77  
78  
79  
80  
81  
82  
83  
84  
85  
86  
87  
88  
89  
90  
91  
92  
93  
94  
95  
96  
97  
98  
99  
100  
101

## 5.1 Fully automatic modulation of FCL aggressiveness

The following subchapters describe set-ups you may want to use for allowing completely hands-off FCL in as many daily situations as possible.

### 5.1.1 autoISF generally switched off outside of meal-time windows

If, aside from having to bolus for meals, your hybrid closed loop was running pretty well *without* other interventions from your side, you could continue to run in that mode, and just focus your new autoISF FCL on management of meals.

In your initial transitioning phase this approach makes a lot of sense, and even by focusing autoISF on just a sub-set of them, like only dinners.

Also in the long run this avenue is taken by many FCL users for the night times, “hanging on” to their well performing hybrid closed loop with standard `oref(1)` SMB+UAM (A very good alternative is FCL with night time SMB shut-off, see next [section 5.1.2](#)).

For this, you define Automations

- that set meal time windows in which “Enable ISF adaptation by glucose behavior” (autoISF) is turned on in AAPS preferences/OpenAPS SMB
- or: that turn *all* autoISF’s ISF modulations (or just `bgAccel_ISF`) off in time windows in which surely no meal occurs. For instance, you can go for all nights back into your Hybrid Closed Loop, as you had before.

Other early DEV AAPS variants (see [section 13.3](#)) all work with meal-time windows. The window is either set by time of day in the settings, or it always must be „set“ by the user via giving a mandatory small pre-bolus before any meal starts. **Outside** of these time windows, these loops then runs with less aggressive SMBs like `oref(1)` SMB+UAM in AAPS Master. This mode is not really FCL, but an advance over traditional HCL that often achieves satisfying degrees of automation and performance.

The term **Meal Announcement** (MA) is often used to label this closed looping mode. Trigger to set a meal time window could also be a pre-bolus given by the user, a carb entry made, an EatingSoonTT set, or a meal announcement button pushed.

102 Note: Outside of the meal time windows you would be in hybrid closed loop. To the extent you  
103 rarely face disturbances (aside from meals), you could be looping in full automatic mode around  
104 the clock,

105 Your temp. autoISF shut-down is meant to prevent problems from the loop over-reacting to bumps  
106 in the glucose curve in times of day (night) when standard oref(1) performance is sufficient.

107

### 108 5.1.2 Odd-numbered *profile* targets used to block SMBs

109

110 An alternative route of preventing the FCL loop from over-reacting to bumps in the glucose curve  
111 would be to make use of the option to temporarily shut down SMBs

112

113 Ensure the even/odd logic in the settings is toggled on in Preferences> openAPS SMB>

114 autoISF settings> smb delivery settings>: "Enable alternative activation of SMB depending  
115 on profile target" ON.

116

117 In time blocks with an odd-numbered profile target you can prevent any SMBs being given by your  
118 loop. The (unchanged) aggressive settings then can only translate within the limits set by %TBR  
119 possible.

120

121 This will very much slow down any more insulin being given, and is an excellent solution for night  
122 times, especially if you occasionally experience compression lows.



123

124 Alternatively, you could use the new included options for Automation Conditions and  
125 temporarily tune your bgAccel\_ISF\_weight much lower ([section 5.](#) ).

126 Yet another alternative would be to go into hybrid closed loop for the night, with or without  
127 SMBs ([section 5.1.1](#)).

128

129 In case you occasionally do have nights that would benefit from a couple of SMBs (to treat temp.  
130 highs from a late fatty pizza, raclette and such): Define suitable Automations like the two „night“  
131 ones in this *example*:

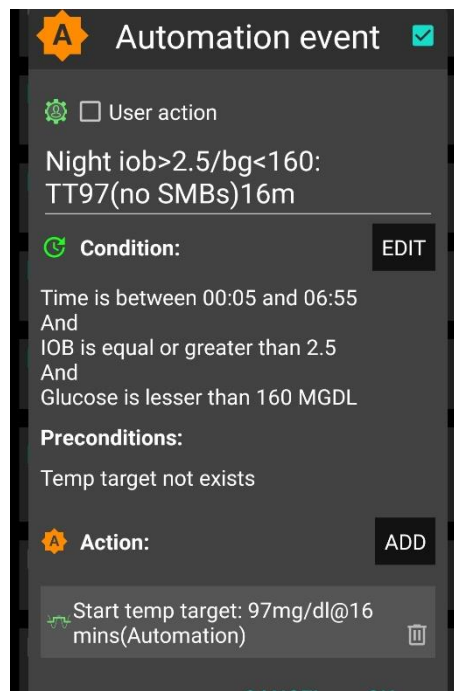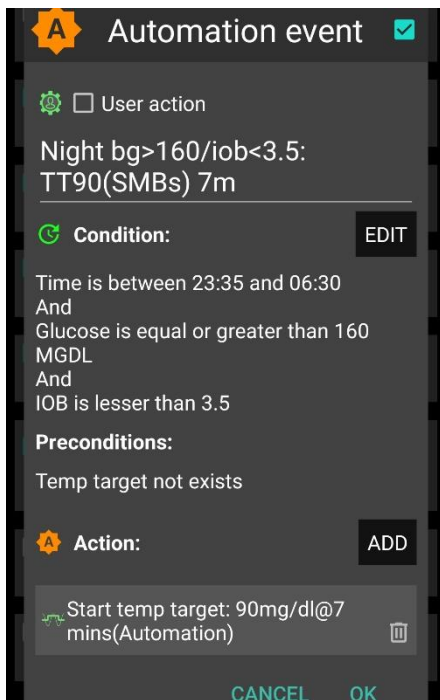

132  
133

134 Never underestimate the „trickyness“ of getting your Automations „right“.

135 With your thought-out Automations in place, night data need to be analyzed to see

- 136 • whether the bg and iob limits defined in the given example work sensibly four your data
- 137 pattern
- 138 • whether the TT duration is chosen appropriately
- 139 • how swapping the sequence in which the automations appear in the Automation list would
- 140 lead to different SMB impacts.

141

### 142 5.1.3 Odd-numbered *temp.* targets (TT) used to block SMBs

143

144 A widely used ACTION that strongly modifies how fast your FCL can add more iob is setting an  
145 **odd-numbered temp. glucose target** which makes the loop operate without giving any SMBs  
146 (%TBR modulation only).

147 Ensure the even/odd logic in the settings is toggled on in Preferences> openAPS SMB>

148 autoISF settings> smb delivery settings>: "Enable alternative activation of SMB depending  
149 on TempTarget" ON.

150

151 So, from patterns you find in YOUR data, at times where you want your loop act differently, you

152 need to carve out CONDITIONS that describe the respective situations (and either for how long it

153 typically lasts, or at which *other* CONDITIONS you want your loop get back to default FCL

154 operation).

155

156 An odd TT is often set for an *anti-hypo* snack or *sports* snack. In both instances, you do not want  
157 SMBs to quickly counter act.

158 In case of sweet “fun” snacks, this is entirely different -> [section, 5.2.1](#) or for regular snacks  
159 (e.g. at school break) see next [section 5.1.4](#)

160

#### 161 5.1.4 Automatic differentiation of FCL aggressiveness using Automations

162

163 **Personalized Automations** tailor the loop exactly to **YOUR** data so **fully automated handling** of  
164 situations with different aggressiveness of the loop can be made.

165

166 From, autoISF 3.0 onwards, also the following parameters are provided as CONDITION and/or as  
167 ACTION for defining YOUR Automations:

- 168 • Enable ISF adaptations by glucose behavior => Allows temp. ON/OFF for the key ISF  
169 modulation parts of autoISF (and, as a result, will usually decrease loop aggressiveness)
- 170 • Trigger/set iobTH percent => Keeps default aggressiveness, but only until a iob threshold  
171 (that your Automation modifies) is surpassed (which is when any further SMBs will be  
172 blocked blocked)
- 173 • Trigger/set bgAccel\_ISF\_weight => Modifies the default aggressiveness of just the  
174 acceleration component

175

176 To set up suitable Automations, you first must **analyze patterns** you find in **YOUR data**, at times  
177 (or geo-locationa, or bg and iob patterns that point to a problem ...) **where you want your loop**  
178 **act differently**, to carve out CONDITIONS that describe the respective situations (and either for  
179 how long it typically lasts, or at which *other* CONDITIONS you want your loop get back to default  
180 FCL operation).

181

182 A variant of this mode is to define several meal time windows in which autoISF aggressiveness  
183 (bgAccel\_ISF\_weight) and/or iobTH are set differently

- 184 • for different meal time slots of your day –

185 (*Breakfast at home, school lunches, school intermission snacks, dinners at home* could for  
186 example all deserve special settings regarding ISF\_weights and iobTH).

- 187 • or even for a geo-location etc –

188 (*School lunches, or mother-in-law visits, would be examples*).

189

190 Unless your meals differ vastly in size and in fast carb content all this may not be needed.

191 Still, personalized Automations might help ease your initial job of setting the various ISF\_weights,  
192 and a best-suitable iob\_threshold\_percent that would work “always”.

193

#### 194 5.1.5 Automatic adjustment of FCL aggressiveness via the Activity Monitor

195

196 With the autoISF variant of AAPS you can make use of your smartphone's **stepcounter** and use it  
197 to fully automatically adjust insulin sensitivity ratio to activity level in the past minutes to one hour  
198 time frame.

199

200 This is another little tuning opportunity, in which you study your body's response to light exercise  
201 (like walking) or to not moving at all (like desk, couch), and select appropriate settings which, in the  
202 future, will automatically adjust insulin delivery to suit activity state of the past minutes (up to 1  
203 hour).(AAPS Preferences/OpenAPS SMB/Activity modifies sensitivity/ -> set **two scaling factors**.)

204

205 This autoISF feature (new since V.3.0) is much quicker responding than Autosens or dynamicISF  
206 to adjust insulin sensitivity to your current „lifestyle state“.

207 For loopers who do not have huge variations in exercise levels in their everyday lives, this feature  
208 might fairly much close the gap towards being able to do a 24/7 hands-off FCL.

209

210 More see [sections 3.5](#) and [6.5](#)

211

#### 212 5.1.6 Pro/con completely hands-off Full Closed Loop

213

214 To stay 24/7 in a completely „hands-off“ FCL can be a realistic goal with autoISF 3.0 if besides  
215 meals also some special challenges, as discussed in this [section 5.1.](#), were analyzed and could be  
216 addressed.

217

218 Clearly it depends very much on your lifestyle, and how interested, willing, and capable you are to  
219 recognize, deal with, (and in the future avoid?) situations that get you outside of your desired %TIR  
220 on occasion.

221 So, this is also about what %TIR you are aiming at, and can accept, as it averages out for  
222 the week, for instance.

223

224 Everybody must weigh for her/himself

- 225 • how much **upfront effort** to put into getting it all 100% automatic
- 226 • **or** whether to take an **easier start, with a couple of situations left to take care of when**
- 227 **and as they arise in daily life**

228

229 Also, even if a principal capability for a fully automatic running FCL is given, this still  
230 means that

- 231 • the user should be knowledgeable about what exactly is going on, and
- 232 • has a principal capability to „nudge“, or to take over in a manual mode.

233

234 In the sections that immediately follow, we present the options to nudge or temporarily take over  
235 from the AAPS home screen that will be serving as your **FCL cockpit**:

- 236 • [Section 5.2](#) describes how you can use available “buttons” from your AAPS home screen,  
237 and how to complete it towards a suitable DIY FCL cockpit, for an even better FCL  
238 experience.

- 239 • [Section 5.3](#) describes how you might be able to manage “disturbances” even better (with  
240 more convenience) with an **improved** FCL cockpit **in the future**.

241

## 242 5.2 Modulating aggressiveness manually, from the DIY-FCL-Cockpit\*

243

244 \* Like in the airplane cockpit: Cruising in full auto mode should involve having an eye on the  
245 instruments, and on potential disturbances ahead in the environment.

246

247 In [section 4](#), we dealt with major meals. In [section 5.1](#) we looked into fully automatable manage-  
248 ment of other situations. Life in Full Closed Loop can become extremely easy then...

249

250 However: Other **disturbances** might come up, that:

- 251 • are not noticeable in-time, or foreseeable, by the loop (*e.g. your plan to start exercise in an*  
252 *hour or two*), but **that influence sensitivity dramatically** and therefore require temporary  
253 non-default settings in order to remain in-range, and/or
- 254 • **require** a different “starting point” regarding iob and bg, which translates into **a different**  
255 **iobTH** that should **temporarily** be set much lower (*in case of exercise*) or noticeably higher  
256 (*e.g. with very fast absorbing carbs in a sweet snack “sin”*) .

257

258 In [section 5.1](#) we looked into ways to automate also a modified loop response to *foreseeable* situa-  
259 tions (tied to a time of day, geo-location etc), or to those *the loop could recognize* (with enough  
260 time to react).

261

262 Other „disturbances“ might come up, and you must find an easy way to



- 263     • call up a pre-programmed routine for automatic management, with adjusted  
264         aggressiveness, or:
- 265     • manually tweak a setting or two, to temporarily adjust the aggressiveness
- 266     • There may also arise a desire to just exit the FCL mode, and be your own captain for  
267         mastering a special situation.

268 For peace of mind, to learn, and to stay informed (especially so in your initial tuning phase, or  
269 when your glucose curve goes in unexpected ways) we also must be able to

270     • find the key parameters that frame and drive the recent and upcoming loop decisions.

271

272 All this is facilitated within seconds right from the AAPS home screen, serving as a **FCL cockpit**  
273 after you built a couple of DIY cockpit features via Automations (as described below and in [case](#)  
274 [studies 5.2](#) and [6.2](#)):

275

276         Thoughts went also into [how to improve the cockpit in future releases](#), see [section 5.3](#)

277

#### 278 5.2.1 Status recognition

279

280 Before considering any manual interventions into the ongoing FCL, you should be aware what the  
281 current mode of action is, and hence which button eventually to fine-tune or lever to switch, in or-  
282 der to adjust to the disturbance you see coming up.

283

284 See [section 5.4](#)

285

#### 286 5.2.2 Manual interventions from the (DIY-) FCL cockpit

287

288 Trouble with all these is, not to forget to set back manually, too (=> better solutions in 5.3)

289

##### 290 5.2.2.1 Temporary tuning of FCL aggressiveness via temp. %profile or TT settings

291

292 The set **% profile** multiplies with both, the ISF resulting from autoISF, and also with the default  
293 iobTH you have set, so both are nicely modulated in a linear way with the % temporarily chosen

294

295         Just taking profile e.g. to 110% for an afternoon might be an easy way to explore whether  
296         you might benefit from 10% more “aggressiveness” in your core settings for lunches (like  
297         bgAccel\_ISF\_weight). Make sure, though, that the extra 10% are not cut away by set safety  
298         limits.

299

300 A lowered (relative to profile glucose target) temporary **bg target (TT)** signals lowered sensitivity  
301 (more insulin need), and an elevated TT (as often used with exercise) increases sensitivity and  
302 hence works in the direction of a lowered % profile to also reduce insulin given by the loop.

303

304 Moreover, the **exercise button** (top center on your AAPS home screen) can be activated (turns  
305 yellow, then). This will **further boost** how your set TT elevates the resulting ISF, and sharply  
306 lowers iobTH, as often desired for sports. See [section 6.1](#)).

307

#### 308 5.2.2.2 Making temporary changes in settings made in AAPS/preferences/Open APS SMB

309

310 Going into AAPS/preferences/Open APS SMB allows to:

- 311 - set milder or stronger ...\_ISF\_weights
- 312 - set different iob\_threshold\_percent (or iobMAX)
- 313 - elevate or lower the SMB\_delivers\_ratio
- 314 - limit or expand max. allowed SMB size
- 315 - change the the even <-> odd logic for SMB on/off

316

317 Doing temporary changes in AAPS/preferences should be the exception because

- 318 - they require multiple steps, including entering a password
- 319 - you will often forget to set everything back to default a couple of hours, or minutes, later

320

#### 321 5.2.2.3 Extra DIY cockpit buttons for pre-programmed “responses” via Automations with “User ac- 322 tion” trigger

323

324 Recognizing conditions for fully automatic handling by the loop may not be not possible, or come  
325 too late for the loop to act on. Examples would be

- 326 • *exercise*: Minimum an hour before starting “the loop should know” to be able to lower iob  
327 and elevate bg by the time exercise starts.
- 328 • *snacks*: High carb snacks, sweets, consuming ice cream or having a sweet drink comes  
329 with the problem of even steeper glucose rises, but overall a lesser insulin need, compared  
330 to major meals (for which we tuned our FCL according to [section 4](#)).

331

332 This not necessarily implies that snacks need different settings than a meal. After all, autoISF  
333 was designed to act to all available data, especially to where the developing glucose curve is

headed. So, depending on your effort to set parameters for a broad variety of meals (notably: how well you avoid to invariably bounce fast against your iobTH), you might be able to accommodate low carb, snack, and major meals with *one* set of settings.

In FCL autoISF, this is a bit more difficult than in HCL autoISF applications, because FCL involves revving up iob supply (largely via big bgAccel\_ISF-weights) often too much to be balanced by just a snack getting absorbed.

For that reason, or just for increased comfort and safety, you might want to differentiate, and make use of what follows for the *sweet snack* example.

### Tuning aggressiveness

Key is that a sweet snack likely benefits from even more aggressive initial FCL performance than the meals in your normal spectrum of diets require.

Therefore, you could set

- a higher **temp. profile%** and/or
- a temp.elevated **bgAccel\_ISF-weight** (see screenshot of my Automation).
- a **low temp. target** (76 for instance; this additionally helps maximize the first SMBs that will automatically be triggered at detection of acceleration)..

When first defining and testing this Automation, also check:

- that the safety limits as discussed in [section 2](#) will not block the intended elevated aggressiveness
- SMBs will not get outrageously big and iobTH sometimes exceeded by too much


Note that “the last SMB” is allowed to overshoot the effective iobTH by 30%

### Limiting iob

For “just a snack”, total insulin need will be lower than for a meal.

If you would just have your sweet drink, and your meal-oriented FCL would “attack”, iob likely would become too high, and a glucose rollercoaster would start, with you needing to consume more =>

If you just have a snack, or drink a glass of juice, you can lower the **iobTH\_percent** accordingly.

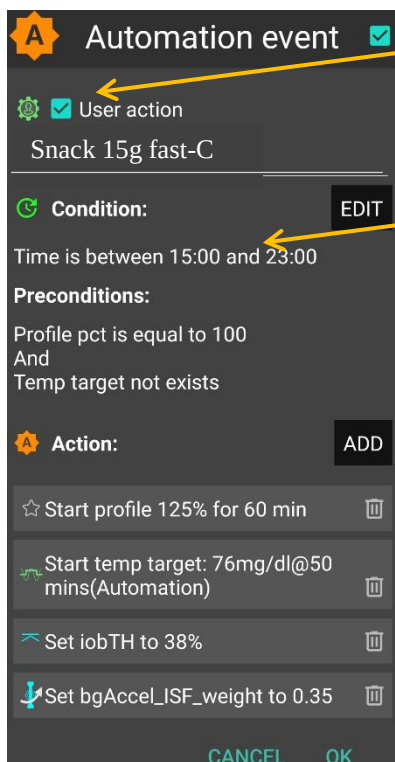

„User action“ is always ticked-on

This will, in the defined time space ..

..offer the “DIY cockpit” button..

..which I must press any time (~90...30 minutes) before my snack.

Note: Make sure that iobTH and bgAccel\_ISF\_weight revert to default afterwards



367

368

369 So, this can be a little extra “project” when setting up your FCL.

370 You need to research your snack habits (if any), and over time find out which settings in  
371 the snack-related Automations work well.

372

373 In everyday life you then just must press the related button in your cockpit (which is  
374 not time critical at all, except it should be clicked *latest* a couple of minutes after you took  
375 the drink or snack).

376

377 If you consume more, and also eat something with your sweet drink, this will more  
378 resemble a full meal... however, with unusual amounts of fast carbs.

379

380 Note: Pressing your snack button a *second time* would **not** help because the  
381 lowered iobTH does not allow iob going high enough. So you are better off just  
382 letting your normal FCL meal routine run, after your snack mode expired.

383

384 Other options when snacks keep extending would require a manual modulation  
385 regarding %profile and/or bgAccel\_ISF, but keeping the full default set  
386 iobTH\_percent, or even elevating it (refer to [section 5.2.3](#)). If that happens often,  
387 define for yourself an extra User action Automation for a bigger snack (= another  
388 grey DIY cockpit button).

## 389 Installing the DIY cockpit button

390 In the related Automation, just keep the “User action” box clicked at all times, and define in  
391 the Conditions when you want to see that button available for cockpit use (see screenshot  
392 above) => you will see that button offered.

393

394 Besides snacks, also any **other recurring special situations can be addressed via a**  
395 **DIY cockpit button, and receive different aggressiveness up to a suitable iobTH**  
396 **level.**

397

398 Over time you can have a big number of User action Automations, and keep them  
399 “shelved” rather invisibly (clicked in-active, top left box) in your long list of potential  
400 Automations. Even when active, they only show in your cockpit (bottom grey field of your  
401 AAPS home screen) in the time slot you assigned as potentially relevant.

402

403 In the future you might be able to set the stage for a snack and other “disturbances”  
404 also via an extended menu below the TT button on the AAPS home screen, see  
405 [section 5.3.3.1](#)

406

## 407 Discussion

408 In case you do have a snack habit and

- 409 • can not find settings, as in [section 4](#). defined for your meals, also suit your snacks
- 410 • can not pin a as in [section 5.1.4](#)

411 you minimum need a “snack announcement”

412

## 413 5.2.4 Temporarily exiting the FCL

414

415 The “last resort” alternative always is to **temporarily** leave the FCL mode, and handle any  
416 disturbance “the traditional way” in **hybrid closed loop**.

417

418 The suggested FCL cockpit user interface with an extra version of violet loop on the  
419 AAPS home screen ([section 5.3.1](#) ) would facilitate that, including automatic removal  
420 and re-appearance of the insulin button at the bottom of the APS home screen.

421

422 In case [this feature](#) is not yet available, you must:

423 Exit the FCL mode by going to AAPS/preferences/put in your password/OpenAPS SMB/scroll down  
424 to autoISF settings and switch “Enable ISF adaptation..” OFF.

425  
426 Unfortunately, there is no way yet for it to come automatically back on, after a selected time for  
427 instance. So do not forget to switch your autoISF fully back on, later.

428  
429 As this will often be forgotten, it may be worth doing a “User action” Automation, for a “temp.  
430 FCL OFF” grey button (see section 5.2.1).

431 Caution though, there is very limited experience with this brand new feature

432  
433 To recognize whether autoISF currently runs with ISF adaptation or not, you must consult the  
434 profile\_sens -> actual\_sens indicator below the Autosens%. However, this gets also modified  
435 by %profile switches or TT +/- exercise mode. So it is not as easy as it would be with the “violet  
436 loop” proposal mentioned already above.

437 Ultimately, you can of course study the SMB tab to find out what is going on.

438  
439  
440 **5.3 Modulating aggressiveness manually from the improved FCL-cockpit**

441  
442 autoISF 3.0 is an early dev variant of AAPS, and as user you are participating in an on-going  
443 development. Of note, autoISF 3.0 is first launched without many of the cockpit features that are  
444 presented below in this font color. (Only what is written in black is at this point of some relevance  
445 for using autoISF 3.0)

446  
447 For the time being, multi-step work-arounds may become necessary

- 448 • In many cases, going into AAPS Preferences and changing settings would be needed  
449 (...plus not forgetting to change these settings back, afterwards).

- 450 • Automations allow a DIY FCL cockpit, see [section 5.2](#) and [case studies 5.2](#) and [6.2](#)

451  
452 This is also an open invitation for you to contact us in case you could help program a  
453 module for one of the required user interface extras.

454 For future integration into AAPS Master, an eye should be kept also on the question which  
455 other modes (like FCL using Automations and others mentioned in [section 13](#).) might  
456 benefit from some of the extra features.

457  
458 Keep in mind, though, that the **goal should be to interfere with the loop as little as possible**.  
459 Under the described conditions it can run **fully automatically** without any user interaction ( = after

460 the initial tuning phase, and related settings made in AAPS /preferences/SMB/autoISF. See [section](#)  
461 [4.](#) and [5.1](#)).

462

463 However, just like in the airplane cockpit: Cruising in full auto mode should involve having an eye  
464 on the instruments, and on potential disturbances ahead in the environment.

465 *E.g.: storm ahead => instruct your plane to climb to another flight height.*

466 *Analogy: exercise ahead => setting an exercise TT, or => pressing a button that activates a*  
467 *sequence of instructions (some of them probably hinging on conditions, like actual iob) how*  
468 *to manage through that exercise situation).*

469

470 So, for the occasional „disturbance“ coming up, you should find an easy way to

- 471 • call up a pre-programmed routine for automatic management, with auto-adjusted  
472 aggressiveness, or:
- 473 • tweak a setting or two, to temporarily adjust the aggressiveness
- 474 • There may also arise a desire to just exit the FCL mode, and be your own captain for  
475 mastering a special situation.

476 All this is facilitated within seconds right from the AAPS home screen's **cockpit features to the**  
477 **extent they are already incorporated, or** to the extent you can build alike DIY cockpit features via  
478 Automations, as described in [section 4.1.3](#) and [case studies 5.2](#) and [6.2](#)):

479

- 480 • The button that is integrated into the **violet FCL icon** serves as emergency off button, to  
481 quickly stop FCL, or to at least to immediately stop any more SMBs (...just for a couple of  
482 minutes, or for the remaining meal time: pick from the options offered with just one  
483 keystroke).

484 Via the violet FCL icon on your AAPS home screen, you also can access a temp. switch-off  
485 button for SMBs (see section that next follows below).

486

- 487 • The **three top fields** (%profile, exercise, TT) provide access to temp. tuning of core  
488 parameters, and/or to some pre-programmed routines.

489 Taken together with some **new indicator fields** about your loop state, and the **grey DIY cockpit**  
490 **buttons** this makes the AAPS home screen your **cockpit** for Full Closed Looping.

491

492 Let us look on each of these cockpit elements in some detail:



493  
494 5.3.1 Violet FCL icon and underlying buttons


495  
496 Novices to FCL, or really anyone running into a very special situation, may appreciate that the new  
497 closed loop icon on the AAPS home screen in pink (for FCL) has buttons to quickly shut off getting  
498 more SMBs (1st row), or to enter other loop modes (second row).

499  
500 It functions very much as the other ones that you know from HCL already, and in fact you  
501 get offered some of the same options (for instance, to switch the (full) closed loop off for 15  
502 minutes for going to take a shower)

503  
504 Note that in FCL you leave all BG regulation, notably against meal spikes, to the loop. So, try not to  
505 disconnect in phases when your FCL must ramp up your job.

506 The required insulin would still be supplied *after* you reconnect. However, without the user  
507 pre-bolussing, the delay would be more of an issue in FCL than it had been in HCL.

508  
509 Just pressing on the FCL icon, a dialogue box comes up:

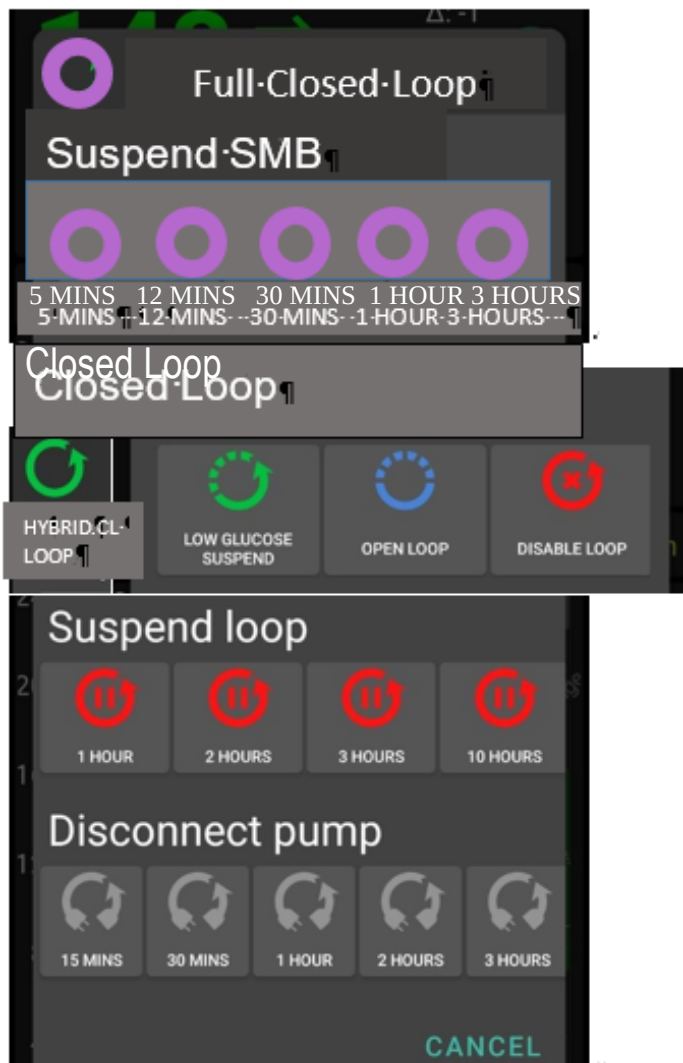

<-add-arrow, form-like-like-CL-Loop...¶

<-dotted+arrow, form-like-like-LGS...¶

**DEV:** when in Hybrid-Closed-Loop = top headline + green circle, there is no Suspend-SMB part but it starts right with „Closed-Loop“. In that screen. There, the 1st element must be option to pick ¶

.....pink circle, FULL-CL-LOOP ¶

where here the green HYBRID-CL-LOOP stands ¶

For loopers who did not set up FCL, a feedback must come up ~ „FCL not installed“ if they press on that button ¶

510



511 Pressing „**Suspend SMB**“ provides fast and easy „emergency braking“ regarding delivery of more  
512 SMBs:  
513 Select the one with the desired number of minutes: 5 or 12 for just blocking the potential next  
514 SMB(s), and up to 3 hours to manage the entire rest of this meal with %TBR from then on.  
515  
516 Whenever, and why-ever, your FCL is in „no SMBs allowed“ mode (e.g. automatically after  
517 surpassing an iobTH, or triggered by a set odd TT), the FCL icon will turn into a dotted one.  
518 Instead of remaining **duration to end time** it indicates in the middle „the condition“, „**iob**“ or „**TT**“  
519 Add an indication **if** suspend SMB comes from an Automation, e.g. add an „**(A)**“, **underneath** the  
520 #minutes, iob, or TT in the middle of the dotted violet field.  
521 So, as in other (already in HCL existing) cases, those icons show in the middle the minutes left that  
522 they will be running, or the condition which would have to go away for this temp. setting to stop.  
523 It always auto-reverts into the FCL state and FCL icon, when time (or other condition) has elapsed.  
524  
525 Pressing „**HYBRID CL. LOOP**“ or other buttons from the 2nd row provides fast and easy  
526 „emergency **exit**“ **into other modes**.  
527 This enables beginners an easy „temp. escape“ into their well-known HCL (green) at any  
528 point of time. bgAccel\_ISF\_weight is set to zero when going FCL->HCL. HCL can run with  
529 autoISF (for instance dura\_ISF) uninhibited otherwise. (check **implications** for HCL users of  
530 autoISF ?? ).  
531 Note: These options from row 2 have no time limit. Loop will **not** by itself go back to FCL. You see  
532 the different loop icon as a reminder to manually revert, when ready.  
533  
534  
535 5.3.2 Buttons „Insulin“, „Calculator“ etc at bottom of AAPS home screen  
536  
537 These buttons are **not useful any longer in FCL**, and automatically disappear whenever in FCL  
538 mode (also in Suspend SMB state), and re-appear when leaving FCL. This applies also when an  
539 Automation or technical system failure shut off FCL.  
540 Users who, maybe in the beginning phase, feel better having those buttons, can override  
541 the removal (of the insulin button, or any other) by going into /preferences/overview/buttons  
542 and forcing them on. They only remain on until the next re-entry into FCL mode, when auto-  
543 off happens again.  
544 The reason why we do this: It really is important to let the loop loop, and not interfere more  
545 than absolutely needed. Any bolus the user gives will sure distort the bg curve, on which  
546 autoISF, especially when aggressively tuned for FCL, builds a lot of its decisions!

547

548 5.3.3. Three top fields (%profile, exercise, TT)

549

550 Depending on the variedness of lifestyle, the desired %TIR, and the initial tuning effort put in, the  
551 user may want occasionally to „tweek“ the aggressiveness of her/his FCL.

552

553 The top 3 fields (grey in default mode, **yellow when temp. in mode with changed**  
554 **aggressiveness**) serve as quick and easy entry points to make temp. switches (as users will be  
555 used to for %profile switches, or for setting an EatingSoonTT in HCL, .. which they still can do in  
556 FCL ... but more:)

557

558 Expert FCL users might need this feature rarely, but probably at least to manage activity after  
559 meals: Each require opposite aggressiveness, and the switch has to come in a certain point in  
560 time that would be difficult to capture. (More see [section 6.4](#) )

561

562 5.3.3.1 TT dialogue field (( Currently not available in the pictured form and function ! )

563 *f extended design for FCL cockpit is already launched)*

564 The TT field (top right of AAPS home screen) is the primary daily interface, and a dialogue field

565 opens when pressing on it

Initially, the form can be entirely empty re. TT inputs and just show the default iobTH and bGAccel parameters.

Two lines appear for target and get labeled ES, AC or HY depending on what was pushed at the bottom buttons. (HY => red frame; evtl 1 empty => grey)

TT entries can be made or overridden.

iobTH calculates from Target AC and other settings shifts away from default set in /preferences (here 0.7) but could be overridden here 0!..200 % is allowed

bGAccel ISF\_weight can be modulated here, too. Note: it can change again if % profile is also changed. --0...200 % is allowed

CANCEL allows to start fresh (select one or two of the square buttons, ES, AC or HY)

OK needed to use the settings

Duration input is made in minutes. In the exceptional case that both, ES and AC targets are defined, the duration input is for AC and framed blue. (This is because the preceding AC mode is automatically determined in length by the loop observing when iobTH is exceeded)

566

567

568 This looks complicated but only because it allows 4 different modes of use. Each user will primarily

569 use her/his preferred one.

570 (1) Who is happy with the initially well tuned FCL and does not have huge variations in daily eating

571 and moving around, will **not use** the TT **at all**. FCL is possible without an intervention via the

572 TT button in your cockpit. Actually 4 of 8 modes (GGG ...YYY permutations, list see [section](#)  
573 [5.3.1](#)) are not making use of TT.

574 (2) Super easy is also, to just input **any odd-TT** (odd-numbered temporary target) that will shut out  
575 any SMBs for the set duration. *That can be a good idea when having a snack, for instance.*

576 Super quick access to stop SMBs is possible also via the loop icon ([section 5.2.1](#)).

577 Specifically, an **EatingSoon TT** can be activated here (limited relevance see [section 2.5](#)). It is  
578 time-un-critical, can be manually set, or come up via an Automation.

579 The cockpit enables you to set the iobTH differently (override) for the current meal.

580 Alternatively, iobTH can be temporarily changed in /preferences or using an Automation.  
581

582 Temp. iobTH will always revert to default when the TT expires. If another TT immediately  
583 follows, like in the example of the screen above, it will calculate, (then) show and use a new  
584 temp. iobTH.

585 (3) The third way is to **use the input mask** (if already ncluded in your software version see picture  
586 above) **to freely modulate the loop aggressiveness** for a declared number of minutes. Click  
587 the bottom big square(s): Either HYPO, or ACTIVITY, or EATING SOON, or ACTIVITY and  
588 EATING SOON (example in the pictured screen above). Make or override entries in the offered  
589 fields. Press OK.

590 (4) The fourth way is to exclusively use one of the 4+4+2 little buttons seen in the bottom part of  
591 the TT dialogue box (if already included in your software version). They provide a set of  
592 settings (as will immediately show in all input fields above) that the user has set up in  
593 Preferences/SMB/autoISF/FullLoop (refer to [section 6.3](#)), and can freely label there. For  
594 instance „hiC“ at high carb EatingSoon, „piz“ for Pizza/fatty meals, „grd“ for garden work,  
595 „mtb“ for mountain biking ...

596 **Capturing good settings for not-everyday situations in /preferences** (if already included)  
597 **allows calling them up within 1 second**, from your cockpit on the AAPS home screen (...and  
598 won't ruin the FCL experience at all , especially because in most cases it is not time-critical,  
599 how long before the intended exercise the buttons are pressed).

600 [Case study 6.2](#) demonstrates that nearly the same performance and comfort can be reached  
601 via the **DIY FCL cockpit** with the grey extra buttons appearing at the bottom of the AAPS home  
602 screen, based on Automations with User action (see also [section 4.2](#)).

603 The example picture given above, and also [case study 6.2](#), is the most complicated (but also most  
604 useful) case, **when exercise follows after a sizeable meal**. It is then that you need (a) aggressive

605 FCL initial performance at the meal, but, exactly when (!) a (for the intended sport already  
 606 temp.lowered) iobTH is exceeded, you need (b) to have SMBs automatically switched off and go  
 607 into the „milder“ mode, as defined for the exercise (with high instead of lowTT, that automatically  
 608 significantly reduces iobTH again, and insulin sensitivity(resistance) settings too).  
 609 Pressing exercise related buttons will automatically also light the **exercise button** on the main  
 610 screen yellow.


611  
 612 To summarize, the TT dialogue field offers easy but powerful ad-hoc modulation of loop  
 613 aggressiveness for FCL (if already included).

614  
 615  
 616 5.3.3.2 Exercise button (see more in section 6.)

617  
 618 The exercise button automatically lights yellow when exercise related TTs are activated in the TT  
 619 dialogue box. 4 of 8 modes are making use of the exercise button.

620 If pressing on the exercise button, a dialogue box appears (if extended design for FCL cockpit is  
 621 already launched) with info on exercise setting first (and opportunity to override), plus below the  
 622 activity monitor (experimental for auto-tracking of lighter movement during the day, and effects on  
 623 sensitivity that may have. See section 4.5).

624  
 625 So, first the exercise settings (as set under TT) are there to read. Example :

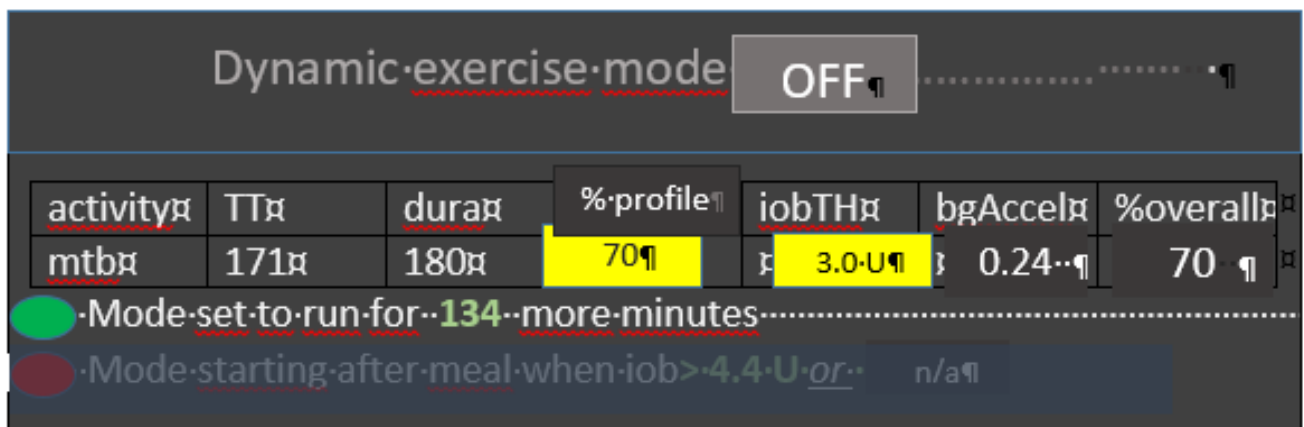

626  
 627  
 628 The exercise (here mtb) is selected in the dialogue box of the neighboring TT field, and there auto-  
 629 filled with settings made in the set-up and tuning stage by the user under preferences (see  
 630 above,...). They are reported also under the exercise button here, and TT, duration, and % sens  
 631 (which also shows active on the %-profile field on the left side of the exercise button) can be temp.  
 632 changed there. iobTH, bgAccel\_ISF and overall resulting sensitivity ratio is given in the other fields.  
 633 The **middle field** of the table, „% profile“ either picks up the % set under the %-profile button, or  
 634 an input can be made here, in the exercise button domain, which will:

- turn the neighboring %profile button on yellow and show that inputted % on it, too
- be multiplied with the result from the exercise mode settings per se, and change the % overall, accordingly.

So, if this middle field of above table (dialogue box of sports button) contains a figure other than 100, input field becomes yellow, and you are operating with a combination of traditional PLUS new exercise mode (with all three top buttons of your FCL cockpit yellow). This maximally will soften aggressiveness, for which you get an idea by the last calculated figure.

The mode is either running already (for another number of minutes, as probably also shown in the yellow TT field anyways). Or it is scheduled to run, after insulination for a started meal reaches iobTH (as in table). Or, no exercise is scheduled (both points red, no entries).

The lower part of the exercise dialogue box (not pictured above, but see in [section 6.5](#)) is dedicated to the Activity Monitor

#### 5.3.3.3 Profile button

The profile button can still be used to set a different profile, or profile%, for instance to adjust for days with sickness (as you are used to from hybrid closed looping). 4 of 8 modes are not making use of the profile button.

Any inputs made here will be used to modify profile\_ISF on which all further changes are made on (multiplied with).

The profile field remains grey if standard profile is applied.

It turns yellow, displaying a %figure relating to any altered loop overall aggressiveness:

- When no inputs (changes from 100% profile) are made here, but inputs in the TT field, e.g. for exercise, automatically lead to different insulin sensitivity ratio ~~that ratio is shown here~~
- when% is changed by input in the profile button itself, it will be multiplied with with profile\_ISF and be used in place of profile\_ISF *by the algorithm*.

However, for exercise (sports) you no longer must make an entry here, because reasonable %reductions should be automatically provided, driven by your set TT (and half-basal exercise target), see [section 6](#).

## 5.4 Recognizing your loop state in the AAPS home screen

### 5.4.1 Color scheme of top cockpit buttons tells kind of closed loop that is running

3 Buttons (%profile; exercise; TT) each in 2 states (yellow Y, or grey G) makes  $2 \times 2 \times 2 = 8$  possible combinations:

YYY = dynamic exercise mode

YGY = not-dynamic „traditional“ exercise mode (if <100%) or hypo mode (if >100%)

GYG = basic closed loop with Activity Monitor running

GGG = basic closed loop (FCL or HCL) without any altered sensitivities etc

YGG = basic closed loop but with a „long wave“ sensitivity shift (e.g. sickness)

GGY = temp. target like e.g. EatingSDoonTT is set; or Hypo mode

YYG = closed loop with „long wave“ sensitivity adjustment and Activity Monitor running

YYY = dynamic exercise mode in time with additional „long-waved“ sensitivity shift

### 5.4.2 Information printed on the top buttons

The yellow TT field shows the currently valid TT (and further duration):

(profile) stands for the abbreviation you labeled your selected running profile



In the special case of settings for meal preceding sports, the field will look slightly differently:



...and ...

... when iobTH is first time exceeded, this automatically switches to:



Likewise, if on the AAPS main screen just an **EatingSoonTT** is set (e.g. 72), this is entered with the desired duration. Afterwards, it automatically reverts to profile target and the display turns grey again there with e.g. 90 on it (and no time limit).

Without sports context, the middle field remains grey.





699  
700

701 Independently from setting a TT, the user can choose to set a **%profile in the left top field**, for an  
702 independent number of minutes, e.g. 70% in this screen example: Also, or additionally, this will  
703 influence the resulting ISF and sensitivity%

704



705  
706

707 The % might change and turn yellow also in context of making TT inputs in the related dialogue  
708 box (see chapter TT dialogue field, above). Still, the % (or the length of time the profile switch shall  
709 be active) can be independently overridden in the top left field, if so desired.

710

711 If an **Automation** sets a %profile, and/or a TT (e.g. *automatic detection of meal start at condition*  
712 *e.g. when delta >10*), this would automatically show in respective field(s) turning yellow and  
713 showing the temp. setting. To show the set parameter comes from an Automation, „ **(A)** „ is added  
714 in the end of button text.

715 Note that an Automation is **usually/ always (?)** only permitted to temp. change default profile  
716 settings, not other pre-existing temp. settings. This is for a good reason : Why should a  
717 sometimes in the past thought-out Automation supersede your - just for the occasion  
718 specified – temp.settings that you consciously activated for the day?

719 Advice: Try to stay away from Automations that also aim at temp. modifying  
720 aggressiveness. For the reason just given in above note, they often will not kick in anyways.  
721 Generally, it also is no good idea to double up sub-algorithms for tweaking loop behaviours.

722

723 Try to **keep things as simple and clear as possible.**

724

725 That said, a limited number of Automations can be of help in distinct scenarios (that differ in  
726 purpose and in applicable time of day).

727 A good one could be for night time, when your odd profile TT has SMBs shut off, but your  
728 experience after pizza nights tells you that, under certain condition patterns (bg, iob), an  
729 SMB or two should be „allowed in“. Another good example, if you go usually FCL without



730 any use of the TT button (which you could call a meal announcement of sorts), is to define  
731 an Automation that, after detecting a meal start, automatically sets a low TT to get  
732 maximally aggressive first SMBs.

733

#### 734 5.4.3 FCL related indicator fields in the AAPS home screen

735

736 In extra data fields of the AAPS main screen you can always see (not change) the key  
737 „aggressiveness“ parameters your loop currently operates with (see also home screen  
738 example below):

739 • how profile sensitivity (**ISF**) changes by the %profile input, by autoISF, and/or a set  
740 exerciseTT.

741 • next to current available iob number is an indication of your **valid iobTH** (the iob above  
742 which no more SMBs will be given)

743 • The AAPS home screen additionally shows, above the deltas, the current **acceleration**

744 Having a look at that can be valueable. For instance, when glucose is relatively low and still  
745 falling, a positive (and getting more positive) acceleration indicates that bg will swing back  
746 up, rather than crash low. This will give info about necessary snack size, and hence help  
747 avoid both, unnecessary calories, and going on a bg rioller coaster.

748  
749  
750 5.4.4 Overall home screen:  
751

Overall home screen:



752  
753  
754