# Tuning of autoISF settings for Full Closed Loop aided by the Emulator V.2.9 for 3.0.1

2 3

6

7

8

1

4 5

Please note that with autoISF you are in an early-dev. environment, where the user interface is **not optimized for safety** of users who stray away from intended ways to use. Good safety features exist, but these are only as good as the development-oriented user understands and implements them. This is not a medical product, refer to disclaimer in section 0



Available related case studies:

Based on older autoISF and older

Emulator versions, examples from

emulator use can be found in case

there), and case study 8.2

study 6.2, in case study 4.1 (last pages

Case study 10.1:

9

11

12

14

18

21

24

25

27

28

29

30

31

10 10.1 Installing the Emulator on your PC

10.1.1 Downloads for QPython

10.1.2 File structure on your PC

13 10.1.3 Start emulation button on desktop

10.1.4 .py scanning files

15 **10.2** Analyzing **loop decisions** in logfiles

16 10.2.1 **noChange**.vdf

17 10.2.2/3 Locate logfiles / prepare Emulator

10.2.4 Run emulation and inspect results

19 10.2.4.1 .txt (all SMB tab infos)

20 10.2.4.2 Tabular (.csv) presentation of all loop decisions

10.2.4.3 Manual extraction of key data in .xls or .odc

22 10.2.4. 4 .pdf chart

23 10.3 What-if analysis

10.3.1 Define (yourChanges).vdf

10.3.2 Run emulation

26 10.3.3 Inspect results

10.3.3.1 Logs (all SMB tab infos)

10.3.3.2 Tabular (.csv) presentation of all loop decisions

10.3.3.3 Semi-automated extraction of key data

10.3.3.4 .pdf chart

Rather than elaborating further, what to best do in data analysis, we should, over time, add case studies.

32 33

34 You can set up and tune the system for full closed loop as described in previous sections. Doing 35 this by frequently analyzing screenshots that must be taken in real-time of the AAPS SMB tab is

36 tedious, however.

37

38

More elegant and precise tuning can be done with a special evaluation software for the AAPS

39 logfiles, by using the **Emulator.**  In the Emulator, you can see in tabular and graphical form, which autoISF component, and other settings, contributed to SMB values that determined the glucose curve.


Note that iOS based variants of autoISF (Trio or iAPS) can not use the Emulator. Refer to section 11.3.

## 10.1 Installation of the emulator on your PC

# 10.1.1 Downloads

• First download **QPython 3L** onto your PC (from Google Playstore).

• Then download from: <a href="https://github.com/ga-zelle/APS-what-if/blob/A3.2.0.4\_ai3.0.1\_dev/software">https://github.com/ga-zelle/APS-what-if/blob/A3.2.0.4\_ai3.0.1\_dev/software</a>: each of these 5 py resp. config. files. To do this, you must press, for each one, here

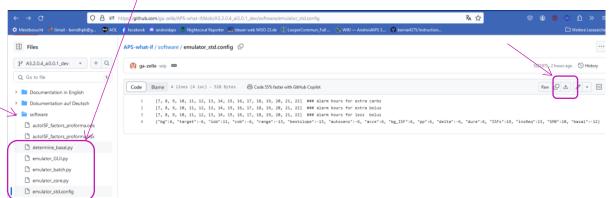

Always make sure you use the files from the branch with the same version number as your AAPS version (in the example above: These files will work with AAPS dev version 3.2.0.4 with autoISF version 3.0.1); (the dev connotation is temporary). Always keep your AAPS x autoISF and also the emulator related files up-to-date!


#### **10.1.2** Create your PC folder structure

Retrieve these 5 downloaded files on your PC (list of recent downloads), and shift them into a folder in

which, or neighboring to which, you also keep your logfiles.

(These you must copy-in about once a week from your phone )


My folder structure for Logfiles and Emulation on the PC has (three) monthly folders that I most likely will look at, plus one folder with data from previous months and years.

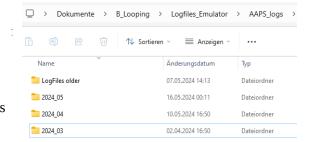

It is advisable to additionally store a pdf from **Nightscout Reporter** in the file for every month. From it, you can much easier find which days and times are of high interest to analyze with the Emulator.


Always factor in the time difference between your AAPS phone and the "Z" time used by AAPS.

Neighboring the AAPS\_logs is the Emulator file where the 5 downloaded files went:

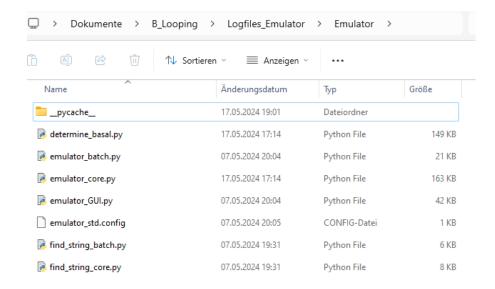

- **10.1.3** Create a "start emulation" button on your desktop
- One of these files is "emulator\_GUI.py"
- Create a link to it
  - Drag that link onto your **desktop**
- Name it something like Emulator start.

- 10.1.4 Fetch two more .py files
- Repeat steps 10.1.1. and step 10.1.2. for 2 more py files from <a href="https://github.com/ga-zelle/Scan-APS-">https://github.com/ga-zelle/Scan-APS-</a>

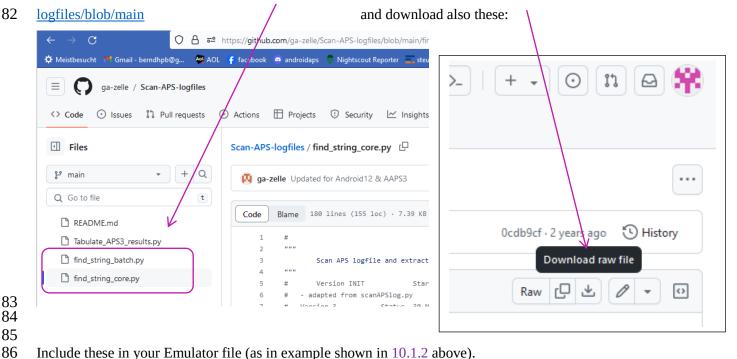

Include these in your Emulator file (as in example shown in 10.1.2 above).

## **10.2** Analyzing loop decisions in logfiles

89 Instead of making many screenshots every 5 minutes after meals, and analyzing them later, a much more

90 elegant and powerful way to analyze your loop decisions (and how you might want to influence them with

91 different settings, see section 10.3 for this), is to use the Emulator.

9293

88

10.2.1 Set up a "no change" .vdf file.

94 95

96

97

It is just empty in the lines that would define any change to be investigated.

Note: for "what-if" analysis, entries will be made (in a second .vdf later, see section 10.3)

To do this, just open **Notepad++** (from list of all programs on your PC). Store that in a file of your current

98 emulator project you are about to start (see my storage path in top line here)

99 The no change .vdf should look like something like this:



→ Dokumente → B\_Looping → Logfiles\_Emulator → Emulator\_Studien

autoISF\_factors\_ <- regarding this, see section 10.3.3.3 en

📋 🙆 🖻 🛈 🐧 Sortieren 🕆 🗏 Anzeigen 🗸 🚥

Studie\_1

Studie\_2

Studie\_3

noChange.vdf

Erase any entries after CR LF and also in lines 2 ff, if any

102103

100

101

Store this noChange.vdf in your Emul.-studies file (neighbors

 $104\,$   $\,$  your Emulator and Logfiles files). From that position, you

105 always make a copy and paste into each Study 1 ... n

106

107108

109

10.2.2 Locate relevant logfiles

110 Make sure you have the AAPS logfiles that you want to analyze in an immediate-neighboring file to

111 your Emulator-- and to your Emulator Studies- files (File structure as suggested and shown

112 above).

113 Copy (not: move!) your noChange.vdf (see above) also into your Study file (must be in all of them).

114

10.2.3 Prepare the Emulator

115116117

Now go to your PC desktop, and start the emulator by just pressing the button "Emulator


118 **start**" that you installed in step 10.1.3

119120

121

This opens a big dialogue box with 3 fields that you must fill in with the applicable path (without any

quotation marks "..") from your Windows Explorer file system, best done via (3x) Browse button:

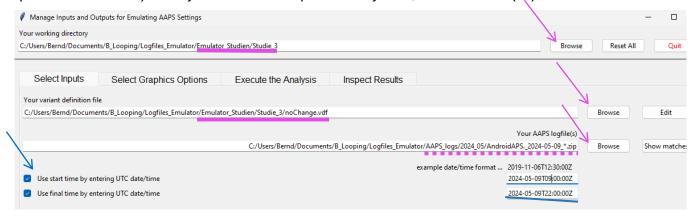

122123

a) The top box marks the path to your current emulator project ("Studie\_3" is where I want to store results)

125

129

130

131

132133

134

135

136

137

138139

140

141142

124

- b) The middle box marks the path to your current vdf (what kind of analysis; here:
- "...noChange.vdf" = read-only; see section 10.3 for what-if)

127 c) The third box marks the path to your AAPS logfiles you wish to look into. A good way to do this is:

- Browse in your Windows Explorer to any logfile from the desired day (2024-05-09 in above example)
- Replace the time with an asterix \* (this means you look at all-day data, in UTZ time).
   Check whether this will work by pressing Show matches.
   You should see all logfiles from that day in a pop-up info box.
- As I wanted to look at 11 am –midnight (for lunch and dinner related data), I:

A

- clicked the bottom left two boxes
- o copied the date 2024-05-09 over the default date in the bottom right two data fields

o after T (for time), I entered the desired time of analysis AFTER conversion into my local time (Central EU summer time minus 2 hours = UTZ; so to look at 11 to midnight of my AAPS screen, I must enter here 09.00:00Z, and below it 22:00:00Z).

Entries at the bottom are not mandatory, but when clicking these little boxes (bottom left) you can define a start and/or an end-point for analyzing, within the logfiles specified in the field above.

143

144

#### 146 10.2.4 Run emulation

147148

149

150

Now we are ready to go: Press "Run Emulation"

This produces sometimes an error message (e.g. if you have a syntax error, or incompatible software versions: => seek help, in the Github materials provided by ga-zelle, or in Discord/Full-Closed-

Looping/emulate-aaps here: https://discord.gg/n3tD5eXExC

151152153

After a short moment results should show up, which you can look into in a couple of ways. First you could have a quick look into the **.log** file to see whether the run had errors (see <u>section 3</u>.)

154155156


159

160

161

# 10.2.4.1 SMB tab contents in (date..) **noChange.txt** result file

- 157 This basically gives you "all the SMB tabs" without needing to make screenshots every 5 minutes.
- 158 Search options help find what lines are of interest to your analysis:
  - By using the search function you can jump, in that super long list, to all places that e.g. have "autoISF" in it or "script debug", or "SMB disabled" (if you want to analyze when that happened…)

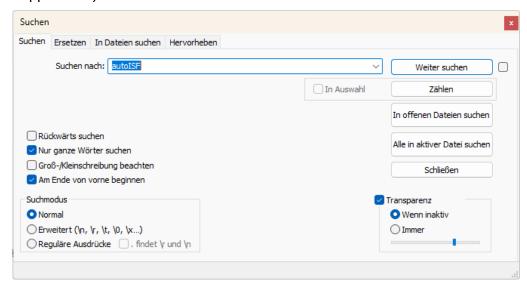

162 163 164


10.2.4.2 Table of results (...noChange.csv file)

- The .csv file in your project folder gives a tabular presentation of how parameters like bg, iob, iobTH, the various ISF contributors, bg target, insulinRequired etc. develop every 5 minutes, and what SMB size and %TBR resulted.
- 169 It is a vast table, so you may want to reduce it to something more "digestable", either after transfer
- to your standard calculation program (next <u>section 10.2.4.3</u>). You can also make settings to suppress
- information you are usually not interested in (or do not know how to interpret, anyways) under "Select
- 172 Graphics Options" when you open the emulator, before executing any analysis:
- 173 First, select your preferred way of outputting decimals (point or comma).

Then select whether you want "All" possible outputs in the graph, or "Most" = all except those you tick "off" in the boxes for each output parameter.

In case you would use "Some/just a few", you would have to tick those few you that do want to see, by ticking the corresponding boxes.

Recommendation is to look at (nearly) everything offered (as your default setting that you can leave untouched in most of your emulator runs):

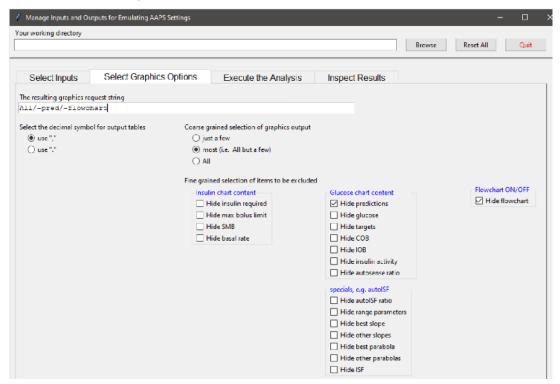

It might be easier, to not deal with customizing the csv file, and rather copy the data into your favorite calculation program:

10.2.4.3 Analysis of the **noChange.csv** table in Excel or LibreOffice calc.

Best copy the entire table into a new .xls or .ods sheet, where you can:

• add right next to the standard world time your corresponding "AAPS time"

For instance, adding +2/24 translates the UTZ column into central European summer time column next to it (where currently a row of Z stands). Likewise, subtract like -5/24 from UTZ for an US East Coast time scale.

Highlight all time fields (the entire columns), and switch from hh:mm:ss format to hh:mm. (While the seconds are important for the loop's calculations, for our comparison with Nightscout or other charts and data, it is much easier without the seconds attached)

• hide any column you find less important to look at for your intended analysis

That way, "boxes" (data fields) retain their original position in tables

196 197	Also, in case later you want to look into additional info, you can simply un-hide the relevant columns or time lines:
198	• hide lines (time segments) you find less important to look at for your intended analysis
199	
<ul><li>200</li><li>201</li><li>202</li><li>203</li></ul>	Usually you will color mark where relevant SMBs were given, which of the ISFs (and underlying weights) was strongly contributing (note that this can be good or not good). Also where iobTH was exceeded, whether an Automation kicked in e.g. setting a TT, when there were periods with zero insulinRequired.
204 205	In $\underline{\text{section } 10.3.4}$ we present an extra tool that does a standardized table reduction and color marking for you!
206	
207 208	You may be able to formulate a hypothesis or two, what settings (ISF_weights, iobTH%, SMB_range_extention, autoISFmax) should be changed for improvement (then go to <u>10.3</u> )
209 210 211 212	10.2.4.4 Graph noChange.pdf
<ul><li>213</li><li>214</li><li>215</li></ul>	After your emulation run, under Inspect Results, you can open the pdf file that is last in the results list offered.
<ul><li>216</li><li>217</li><li>218</li></ul>	This <b>noChange.pdf</b> is a chart that shows along the time axis (down), from right to left:  • Red: the bg curve
219 220	• Yellow: the bg target (note that I do no manual "EatingSoonTT" but for bg rises over +10 mg/dl I have an Automation that sets low TT for a couple of minutes)
221 222	• Light blue corridor: Left edge is set iobTH, and bandwidth +30% (would be +20% at elevated TT)
223	Dark blue line: iob (exceeding twice the iobTH, with temp. SMB shut-off
224 225	As bg did not convincingly come down enough, one could hypothesize that iobTH should be elevated. ((But, again, this would have to be confirmed also with other kinds of meals)).
226	Thin yellow line: Insulin activity
227	• Green dotted line: ISF as would result from AAPS w/Autosens
228	• Green scatter points: autoISF ISF no Chage (lighter points) or what-if (darker points)

229 •

Black line: Profile ISF

230

• Gray scatter points: ISF weakened (to the left of black line) or strengthened (to the right)


231

232 233

234

235236237238239240

• Orange line: cob=0 at all times (in FCL)

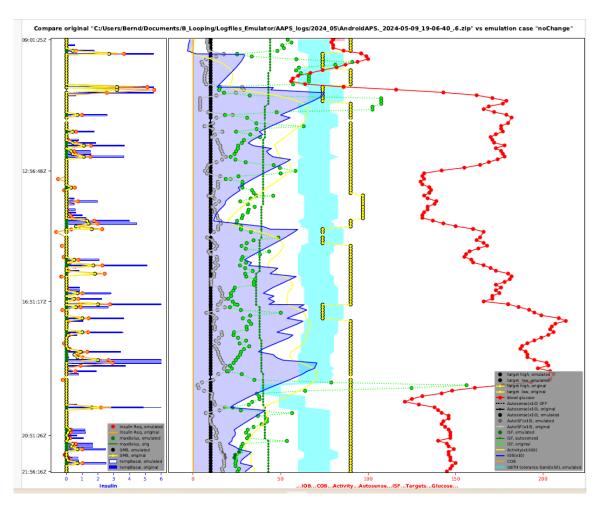

More see discussed together with (yourChanges).pdf in section 10.3.3.4

241 242 243	"What-if" analysis using the emulator
244	In the following you see an example how you can analyze a day of logfiles, and selecting the time span of
245	interest, for instance 11-24 h to look at how autoISF managed lunch and dinner.
246	
247	You have to go through the emulator twice. You already did it (section 10.2.) using the no change.vdf, and
248	now start over with the same emulator with your (desired changes).vdf. Repeat, if you have two or more such
249	vdf defined. (Just clear results before executing analysis each time. No worries, all results are automatically
250	captured for all runs, all in your selected study file).
251	
252	10.3.1 Define your investigated changes in a, or in several, (yourChanges).vdf
253	
254	Define for which one to max three parameters in your current profile settings you want to look into a
255	different setting. Recommendation is to use a factor, like for example $$ current setting * 0.9 $$ , or current
256	setting * 1.2, and use that in your naming for this vdf file, too.
257	
258	Within the same study, you can make several runs with several vdf files.
259	All results, like the csv results table, will appear then several times in your study file, only with different
260	name endings as in the underlying vdf.
261	
262	Example: I like to check in my actual data (they are in my noChange.vdf emulator run <b>), in which time</b>
263	points the following parameter changes would make a (how) big difference in the loop's decision:
264	• 20% higher bgAccel_ISF_weight to boost the first SMBs stronger: How would that tend to ramp
265	up early iob; and might that get too strong in other parts of the data? Or does it bounce into a
266	restriction (maxSMB size; autoISFmax; iobTH) that I might need to widen?
267	<ul> <li>Doubling my cautiously set bgBrake_ISF_weight shall give me insight into the workings of that</li> </ul>
268	parameter (and whether using a much smaller weight than for bgAccel_ISF_weight is really
269	what I should keep doing)
270	<ul> <li>As my bg came down from a persistent high quite slowly, I elevate the dura_ISF by 20%</li> </ul>
271	
272	Actually, it would make more sense to first find my "optimal", maybe indeed elevated,
273	bgAccel_ISF_weight. Then, do a noChange (!) run with that, plus a (yourChanges) run with the stronger
274	dura weight, investigated on that basis.
275	Reason: 1) As we always say, better do only one change at a time. 2) A better job with bg control via
276	bgAccel_ISF will reduce the peak height and provide a different (easier) scenario for dura_ISF to manage.
277	

Now, to **write** your **(yourChanges). vdf for the emulator** (this is same procedure as you did in section 10.2.1 for the noChange.vdf):

280281

284

285

286

287

288

289290

291

- just open Notepad++ (from list of all programs on your PC) to create a new vdf.
- Alternatively you can also take another pre-existing vdf file, and copy it into your current project giving it a new name (re-name it)
  - name your vdf (in our example: 1.2\_bgAccel\_2.0\_bgBrake\_1.2\_dura.vdf) ...
    - ... and store that in a file of your current emulator project you are about to start (see my storage path in top line here)

**Caution:** Make absolutely sure (best by looking it up in the SMB tab, down in the profile set section) to **spell each term exactly** as your loop uses it (probably w/ decimal points, not comma)

• ...when you make one line per parameter (separating entries with spacers->):

profile->(parameter) ->->profile['(parameter)']\*(factor)->->###(comment as you like)

292 293

The (yourChanges) .vdf should look like something like this:

294

```
1.2_bgAccel_2.0_bgBrake_1.2_dura.vdf AndroidAPS_2024-05-09_02_08-31_0.2rp.orig.bt 2 = 2024-05-09T09-00-00Z.1.2_bgAccel_2.0_bgBrake_1.2_dura.txt 2 = 2024-05-09T09-00-00Z.1.2_bgAccel_2.0_b
```

CR = LF= Erase any entries after CR LF and also any entries in lines below, if any

297298

295 296


10.3.2 Run the emulator with (yourChanges).vdf

299

300


301302

The "what-if" emulator run is done the same way as you did the noChanges run ( $\underbrace{section\ 10.2}$ ), however, now, the ( $\underbrace{yourChanges}$ ).vdf must be loaded into the  $2^{nd}$  input field, where formerly you had the noChange.vdf.:

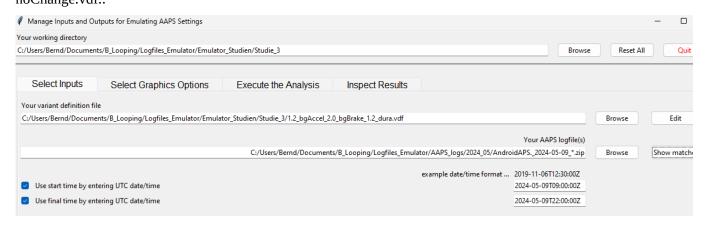


- 305 In the 3<sup>rd</sup> input field, give the path to your stored logfiles. A good way to do this is:
  - Browse in your Windows Explorer to any logfile from the desired day (2024-05-09 in above example)
  - Replace the time with an asterix \* (this means you look at all-day data, in UTZ time). Check
    whether this will work by pressing Show matches. You should see all logfiles from that day in a
    pop-up info box.
  - As I wanted to look at 11 am –midnight for lunch and dinner related data, I:
    - o clicked the bottom left two boxes
    - o copied the date 2024-05-09 over the default date in the bottom right two data fields
    - o after T (for time), I entered the desired time of analysis AFTER conversion into my local time (Central EU summer time minus 2 hours = UTZ; so to look at 11 to midnight of my AAPS screen, I must enter here 09.00:00Z, and below it 22:00:00Z).
  - After making these entries, press Execute the Analysis, (evtl also Clear old Data) and then press Run Emulation, you can look the results up under "Inspect Results". First you could have a quick look into the .log file to see whether the run had errors (see section 3.)

#### 10.3.3 Emulation results

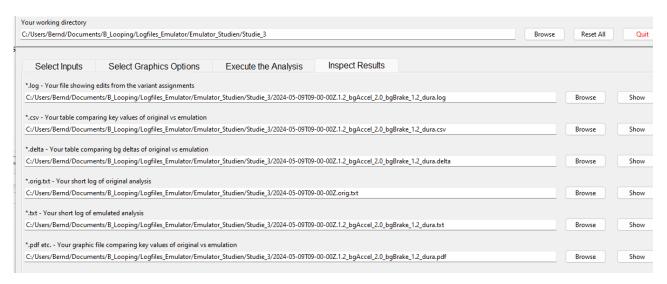

All results from your (yourChanges).vdf emulator go automatically where the noChange.vdf results are already stored, in our example into the "Studie 3" file, below:

Besides the 1.2\_bgAccel\_2.0\_bgBrake\_1.2\_dura.vdf case which I like to look into for the present high carb meal, I also prepared another vdf that investigates a factor 1.2 stronger pp\_ISF and a weaker, factor 0.8, bgAccel\_ISF (with the intention to test this, and a noChange, on a low carb meal later.

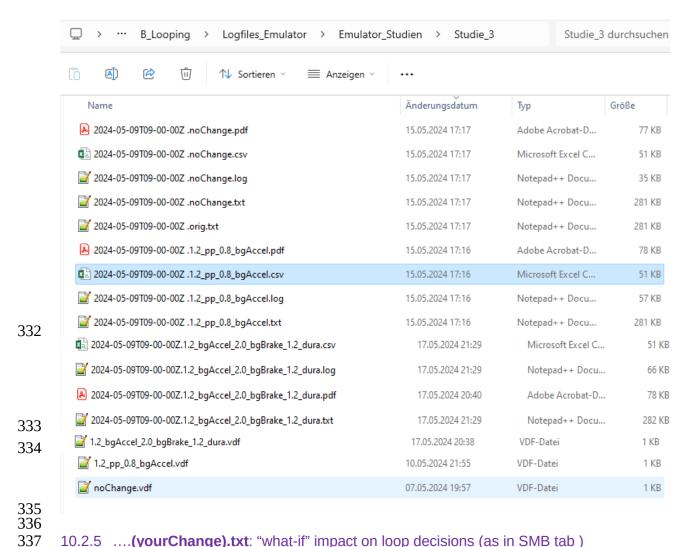

10.2.5 ....(yourChange).txt: "what-if" impact on loop decisions (as in SMB tab )

The **noChange.txt** has all the info your series of SMB tabs had that day.

How to search in this vast list is shown elsewhere (see section 10.2.4.3).

Likewise, the **(yourChanges).txt** gives *for each loop decision* in all detail how and why each single decision would have changed with the different parameter inputs you are checking out here

In the two (yourChanges) examples here, , it was a check on the difference

- a 20% stronger pp\_weight and 20% weaker bgAccel\_weight
- a 20% stronger weight for both, bgAccel\_ and dura\_ISF, and a doubling of bgBrake\_weight

347 would make.

338 339

340

341 342

343

344

345

346


348 349

350

351

Note that all these "what if" data can only give rough hints, notably about the first greater change that you would see with the investigated changed setting. So it works quite well for our main problem in FCL, investigating how to ramp up iob quickly after detection of acceleration.

352	Note that any relevant change would put your bg curve on a different trajectory, so that would influence <i>all</i>
353	following results. Therefore, what you get here is <b>not</b> a complete modelling how your bg would have
354	developed in the alternative scenario.
355	But you can investigate in which stages the parameter(s) you are looking at in your current "what-if" had big
356	influence, and in which direction the changes would go. (see also charts shown in section 10.3.3.4).
357	Analyzing how to safely come down from a high glucose plateau while limiting hypo danger towards the end
358	of digestion is also to some extent possible.
359	
360	A good other way to employ the what-if analysis is real time, on your smartphone, using speech synthesis
361	(see section 11): Then you get real-time info, as to exactly when a significantly different proposal would
362	emerge, and can decide (and watch!) real-time whether to follow the new idea and not was probably better.
363	
364	Observe that a setting change must work well for you
365	<ul> <li>not just in one point of time, and</li> </ul>
366	<ul> <li>not just for one kind of meal,</li> </ul>
367	but you must look at all time slots in the investigated meal, plus analyze with the same tool a totally different
368	meal within your usual spectrum, how things work out there
369	
370	10.3.3.2 Tabular results
371	
372	The <b>noChange.csv</b> table gives all relevant data. Besides development of bg and iob you see the calculated
373	insulinRequired in each loop decision, and how each of the autoISF categories contributed to the decision
374	(notably regarding SMB size).
375	
376	The <b>(your changes).csv</b> shows in detail how <b>every single</b> loop <b>decision</b> would be influenced by the different
377	settings you are investigating. To inspect that huge table, click on the Z behind the start UTC time entry
378	(see black box in the Z column of the following table
379	If you like to see the bg in each screen, too, go 3 or 4 columns farther to the right with your black
380	box.
381	Then, go to window/fix. Now you can scroll through the data and always see headline and time (or time and
382	bg level).
383	To further ease analysis, feel free to temporarily erase (hide) any columns that you (think you) do not
384	need for the intended analysis. More suggestions see in section 10.2.4.2
385	

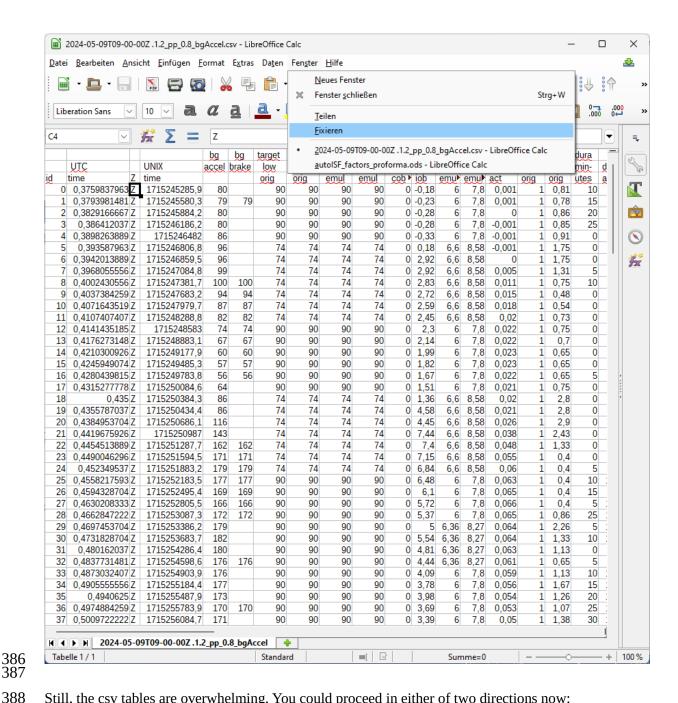

Still, the csv tables are overwhelming. You could proceed in either of two directions now:

389 390

391

392

393

394

395

396


397 398 A) Convert both (or all 3) csv files into one table in Excel or into Libre office calculator. Hide columns (and eventually also lines) that are of no particular interest for your analysis. Mark differences between noChange and (yourChanges) column data with color, add extra columns with additional calculations ...

This route is good to compare quantitative impacts of autoISF categories in critical time points.

B) For the core data relevant to assessing your autoISF settings, there is an extra tool for convenient analysis (see next section)

# 399 10.3.3.3 Automated extraction from tabular results

**autoISF\_factors\_proforma.ods** is provided as an **extra tool** that you download from here:

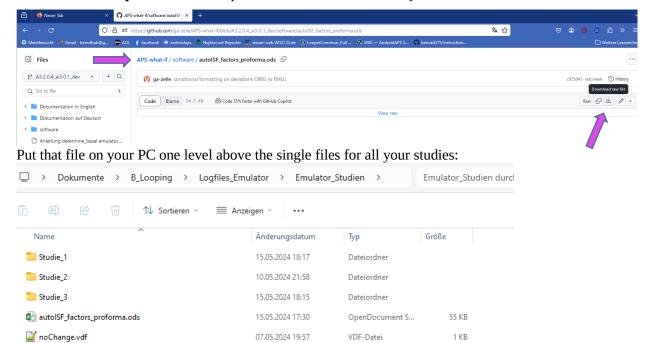

403 404


405

406

401 402

Now, if we want to use this tool on the two csv files of our Studie\_3 file, we must proceed as follows (for *each* of the two .csv files, *separately*):

- 1. Click on the .csv file and open in Libre office calculator.
- 407 2. Make sure the time column is set to US\_English:

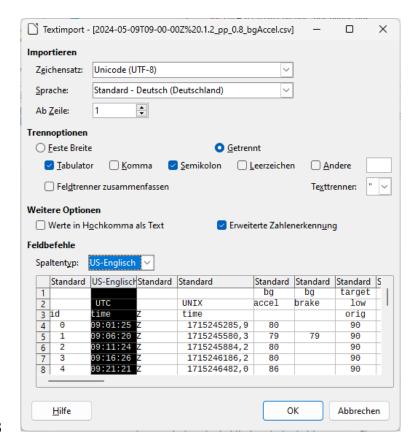

409 410


- 3. Now start, in Libre office calculator, the autoISF\_factors\_proforma.ods ...
- 411 This turns the first 30-some lines of your csv table (left side) into a form in which important effects are
- 412 highlighted in color, and formatting is improved:



413 414 415

Now, you want this for the entire table.


In the autoISF\_factors\_proforma table, highlight 20 or more lines (not including the first or last), and mouse right hand/insert above ...

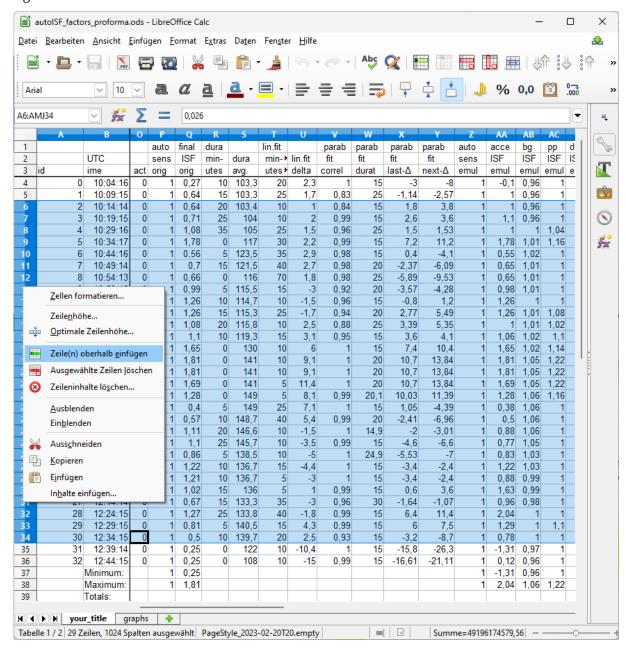

Do this as often as you need to create the number of lines that your emulated csv file comes with.

421 If you ended up with too many lines, erase the superfluous number (any four, in the example):

418

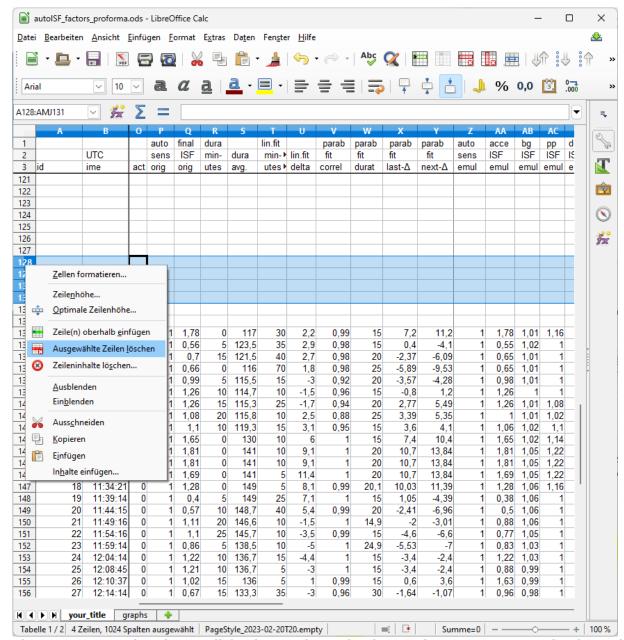

Then just copy it in, by selecting all data lines in the emulated csv, and pasting (paste special, values only)

into box A4 of your "elonged" autoISF\_factors\_performa.ods.

422

427 428

429

430

434

435


The bottom tab "your\_title" should be re-named by you, best with day of log you analyze, and your what-if parameters (so, the name of your csv file could be put in here)

Now you have a table with optimized lay-out that incorporates key data from both your no change AND of your investigated changes.csv files.

431 A super neat extra feature is already pre-programmed, which you can see if you click on the bottom **tab** 432 **"graphs".** 


433 The top graph is the bg curve (the actually seen bg).

Note that for the what-if no bg development over the time range is available. (The noChange one is also given there).

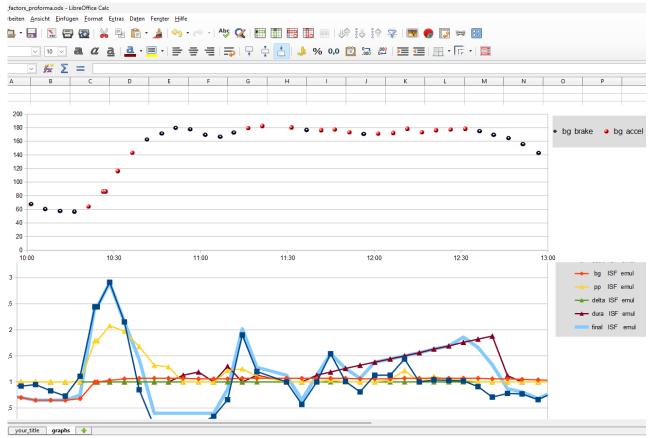

The bottom graph (do one for each, the noChange or the (yourChanges) case) shows the amplification factors coming from each autoISF category, and the overall resulting ISF amplification.

You probably have to widen the time scale (double click on the time axis, and type the desired time span (min and max UTC)(and spacing of data points, 00:30:00 or 00:15:00) into this box:

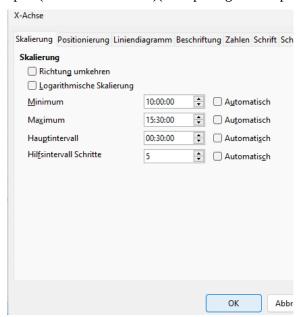

In the given example above, the 2.5 hours were not enough yet to analyze this 10:30 UTC (12:30 AAPS) lunch; we need to look until bg is near target (hopefully before dinner starts).

A similar graph is available on the (i-)phone if you use the autoISF dev variant of iAPS (and maybe of Trio, in the near future). (Later insert details here, or in <u>section 11</u>)

447	10.3.3.4 Chart coming with the Emulator
448	
449	In case you find the extra steps described in the preceding section "too much", also the emulator offers one
450 451	chart (the pdf offered at the bottom of the screen as shown below the "10.3.3 Emulaton results" headline).
452	First look at the initial bg rise in the noChange.pdf chart (emulation results from your noChange.vdf run),
453	and see how bgAccel_ISF and pp_ISF acted, or could have acted in improved ways.
454	Then look into in (yourChange).pdf to see potential effects (or what other change to try). (Actually, you
455	probably will have to go into a detrailed analysis of several lines and columns of the tables as discussed in
456	sections <u>10.3.3.2</u> and <u>10.3.3.3</u> ).
457	
458	Note that ideally we want FCL coverage of our entire "normal day" meal spectrum by one set of
459	settings. So, not much is gained if you put a lot of effort in optimizing FCL settings for one
460	meal.
461	
462	You will need iterations. Do such analysis for two or three very different meals that you wish the
463	algorithm to automatically handle. See $\underline{\text{section 4.2}}/\underline{4.3}$ on how meals with very different carb loads
464	might benefit or also suffer from too aggressive or to mild (category)_ISF_weights you could set.
465	
466	
467	The initial iob received might be limited by allowed SMB sizes, autoISFmax, or the (dynamic!) iobTH. You
468	will have to look into the data table to find out about this (a quick orientation - notably regarding the light
469	blue iobTH band, see next page - is also possible in the pdf result files you have in your project file (project
470	file example "Studie 3" in $2^{nd}$ chart under the $\underline{10.3.3}$ . headline).
471	
472	Only once you found OK weights for bgAccel- and pp_ISF_weights, does it make sense to go tune the
473	dura_ISF_weight. 12:00 – 12:45 UTC in above graph, the resulting effective ISF is dominated by dura_ISF.
474	Just judging from the picture, a stronger weight might be worth trying. However, we really need to see the
475	insulinRequired calculation and the further development because impatience about bringing bg values down
476	faster too often results in hypoglycemia later.
477	
478	The <b>noChange.pdf</b> is a chart that shows along the time axis (down), from right to left:
479	Red: the bg curve
480	• Yellow: the bg target (note that I do no manual "EatingSoonTT" but for bg rises over +10 mg/dl
481	I have an Automation that sets low TT for a couple of minutes)
482	• Light blue corridor: Left edge is set iobTH, and bandwidth +30% (would be +20% at elevated
483	TT)

• Dark blue line: iob (exceeding twice the iobTH, with temp. SMB shut-off


As bg did not convincingly come down enough, one could hypothesize that iobTH should be elevated. ((But, again, this would have to be confirmed also with other kinds of meals)).

• Thin yellow line: Insulin activity

- Green dotted line: ISF as would result from AAPS w/Autosens
- Green scatter points: autoISF ISF no Change (lighter points) or what-if (darker points)

Foreseeably, this is the strongest difference between our noChange (left) and 120% bgAcel\_ISF\_weight (right) in the picture below. (Note the red bg curve is *both times* the really seen bg, because the what-if case only looks at each single loop decision). The first ( -> ) time the dark green dot is far to the right, this *would* get the bg down, we *would start to see* a ( <- ) bg lowering effect, shifting the red curve to the left

- Black line: Profile ISF
- Gray scatter points: ISF weakened (to the left of black line) or strengthened (to the right)

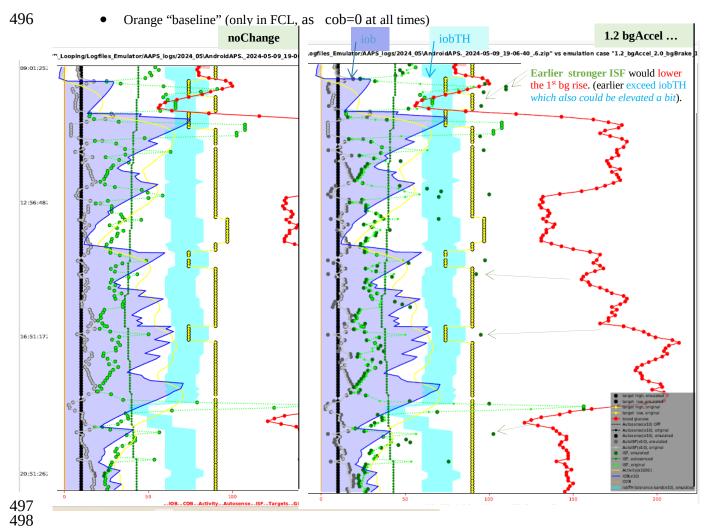

Regarding the other changed parameters: Stronger dura\_ISF would suggest more insulin towards the end of plateaus; this should have helped in the 1<sup>st</sup> plateau (red curve, top right quadrant of the picture). However,

same setting would have to work also on 2<sup>nd</sup> plateau; the chart cuts off there, so too early to see whether a hypo danger might result.

Effect from doubling the bgBrake\_ISF effect are hard to evaluate. Better probably to look in .csv tables, or run a separate emulation for that change only.

Always check for 2 or 3 kinds of your meals whether the "new" parameter settings really are on average better. (See negative example in <u>case study 8.2!</u>)

Part of both above shown charts (left side of each, with blue peaks) was cut out.....

(Unfinished / to be explained later) (...note: yourChanges = 1.2\_bgAccel\_2.0\_bgBrake\_1.2\_dura)

