

10. Tuning of autoISF settings for Full Closed Loop aided by the Emulator

V.2.8 for 3.0.1

Please note that with autoISF you are in an early-dev. environment, where the user interface is **not optimized for safety** of users who stray away from intended ways to use. Good safety features exist, but these are only as good as the development-oriented user understands and implements them. This is not a medical product, refer to disclaimer in [section 0](#)



10.1 Installing the Emulator on your PC

10.1.1 Downloads for QPython

10.1.2 File structure on your PC

10.1.3 Start emulation button on desktop

10.1.4 .py scanning files

10.2 Analyzing **loop decisions** in logfiles

10.2.1 **noChange.vdf**

10.2.2/3 Locate logfiles / prepare Emulator

10.2.4 Run emulation and inspect results

10.2.4.1 Logs (all SMB tab infos)

10.2.4.2 Tabular (.csv) presentation of all loop decisions

10.2.4.3 Manual extraction of key data in .xls or .odc

10.2.4.4 .pdf chart

10.3 What-if analysis

10.3.1 Define (**yourChanges**).vdf

10.3.2 Run emulation

10.3.3 Inspect results

10.3.3.1 Logs (all SMB tab infos)

10.3.3.2 Tabular (.csv) presentation of all loop decisions

10.3.3.3 Semi-automated extraction of key data

10.3.3.4 .pdf chart

Available related case studies:

Case study 10.1:

Based on older autoISF and older Emulator versions, examples from emulator use can be found in [case study 6.2](#), in [case study 4.1](#) (last pages there), and [case study 8.2](#)

Rather than elaborating further, what to best do in data analysis, we should, over time, add case studies.

You can set up and tune the system for full closed loop as described in previous sections. Doing this by frequently analyzing screenshots that must be taken in real-time of the AAPS **SMB tab** is tedious, however.

More elegant and precise tuning can be done with a special evaluation software for the AAPS logfiles, by using the **Emulator**.


In the Emulator, you can see in tabular and graphical form, which autoISF component, and other settings, contributed to SMB values that determined the glucose curve.

Note that iOS based variants of autoISF (Trio or iAPS) can not use the Emulator. Refer to [section 11.3.](#)

10.1 Installation of the emulator on your PC

10.1.1 Downloads

- First download **QPython 3L** onto your PC (from Google Playstore).
- Then download from: https://github.com/ga-zelle/APS-what-if/blob/A3.2.0.4_ai3.0.1_dev/software : each of these 5 py resp. config. files. To do this, you must press, for each one, here

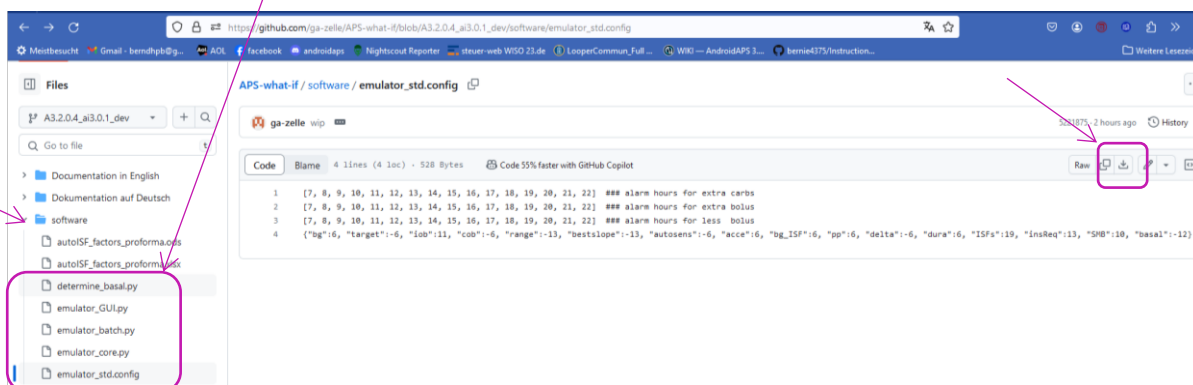

Always make sure you use the files from the branch with the same version number as your AAPS version (in the example above: These files will work with AAPS dev version 3.2.0.4 with autoISF version 3.0.1); (the dev connotation is temporary). Always keep your AAPS x autoISF and also the emulator related files up-to-date!

10.1.2 Create your PC folder structure

Retrieve these 5 downloaded files on your PC (list of recent downloads), and shift them into a folder in which, or neighboring to which, you also keep your logfiles.

(These you must copy-in about once a week from your phone) :

My folder structure for Logfiles and Emulation on the PC has (three) monthly folders that I most likely will look at, plus one folder with data from previous months and years.

Name	Änderungsdatum	Typ
LogFiles older	07.05.2024 14:13	Dateiordner
2024_05	16.05.2024 00:11	Dateiordner
2024_04	10.05.2024 16:50	Dateiordner
2024_03	02.04.2024 16:50	Dateiordner

It is advisable to additionally store a pdf from **Nightscout Reporter** in the file for every month. From it, you can much easier find which days and times are of high interest to analyze with the Emulator.

Always factor in the time difference between your AAPS phone and the “Z” time used by AAPS.

Neighboring the AAPS_logs is the Emulator file where the 5 downloaded files went:

Dokumente > B_Looping > Logfiles_Emulator > Emulator >			
Sortieren Anzeigen			
Name	Änderungsdatum	Typ	Größe
__pycache__	17.05.2024 19:01	Dateiordner	
determine_basal.py	17.05.2024 17:14	Python File	149 KB
emulator_batch.py	07.05.2024 20:04	Python File	21 KB
emulator_core.py	17.05.2024 17:14	Python File	163 KB
emulator_GUI.py	07.05.2024 20:04	Python File	42 KB
emulator_std.config	07.05.2024 20:05	CONFIG-Datei	1 KB
find_string_batch.py	07.05.2024 19:31	Python File	6 KB
find_string_core.py	07.05.2024 19:31	Python File	8 KB

10.1.3 Create a “start emulation” button on your desktop

One of these files is “**emulator_GUI.py**“

- Create a **link** to it
- Drag that link onto your **desktop**
- Name it something like **Emulator start**.

10.1.4 Fetch two more .py files

Repeat steps 10.1.1. and step 10.1.2. for 2 more py files from <https://github.com/ga-zelle/Scan-APS-logfiles/blob/main> and download also these:

Include these in your Emulator file (as in example shown in [10.1.2](#) above).

10.2 Analyzing loop decisions in logfiles

Instead of making many screenshots every 5 minutes after meals, and analyzing them later, a much more elegant and powerful way to analyze your loop decisions (and how you might want to influence them with different settings, see [section 10.3](#) for this), is to use the Emulator.


10.2.1 Set up a “no change” .vdf file.

It is just empty in the lines that would define any change to be investigated.

Note: for “what-if” analysis, entries will be made (in a second .vdf later, see [section 10.3](#))


To do this, just open **Notepad++** (from list of all programs on your PC). Store that in a file of your current emulator project you are about to start (see my storage path in top line here)

The no change .vdf should look like something like this:

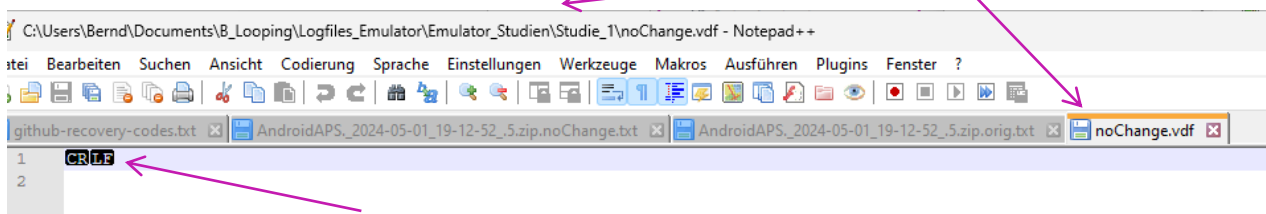

Erase any entries after CR LF and also in lines 2 ff, if any

Store this noChange.vdf in your Emul.-studies file (neighbors your Emulator and Logfiles files). From that position, you always make a copy and paste into each Studie_1 ..._n

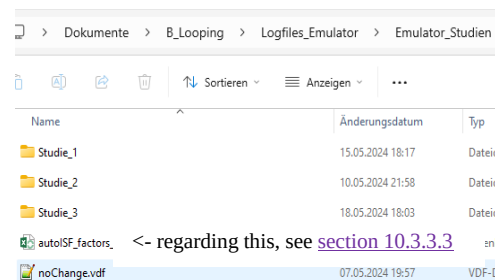

10.2.2 Locate relevant logfiles


Make sure you have the AAPS logfiles that you want to analyze in an immediate-neighboring file to your Emulator-- and to your Emulator_Studies- files (File structure as suggested and shown above).

Copy (not: move!) your noChange.vdf (see above) also into your Study file (must be in all of them).

10.2.3 Prepare the Emulator

Now go to your PC **desktop**, and start the emulator by **just pressing the button “Emulator start”** that you installed in step [10.1.3](#)

119
120 This opens a big dialogue box with 3 fields that you must fill in with the applicable path (*without* any
121 quotation marks “..”) from your Windows Explorer file system, best done via (3x) Browse button:

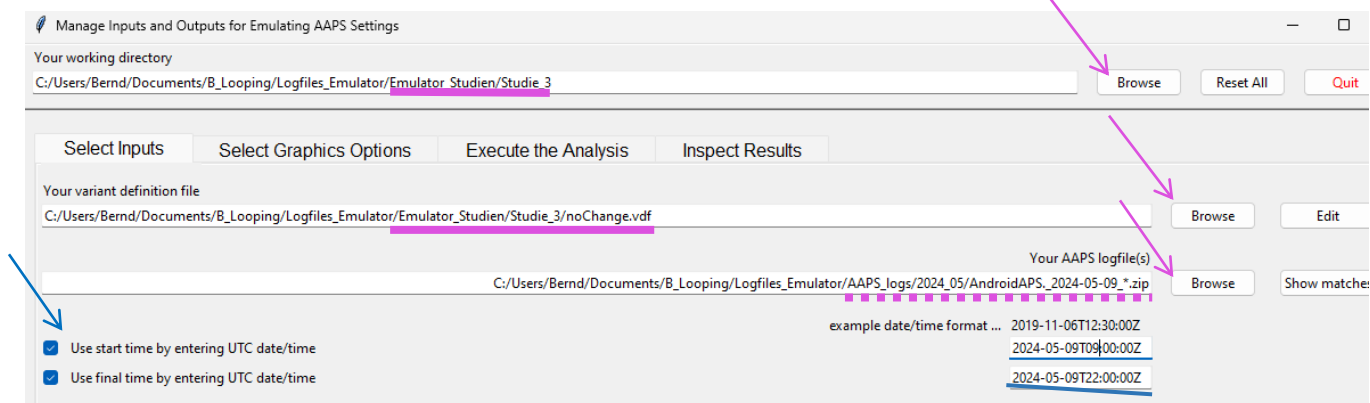

- 122
123 a) The top box marks the path to your current emulator project (“Studie_3” is where I want to
124 store results)
- 125 b) The middle box marks the path to your current vdf (what kind of analysis; here:
126 “...noChange.vdf” = *read-only*; see section 10.3 for *what-if*)
- 127 c) The third box marks the path to your AAPS logfiles you wish to look into. A good way to do
128 this is:
- 129 • Browse in your Windows Explorer to any logfile from the desired day (2024-05-09 in above ex-
130 ample)
 - 131 • Replace the time with an asterisk * (this means you look at **all-day** data, in UTZ time).
132 Check whether this will work by pressing Show matches .
133 You should see all logfiles from that day in a pop-up info box.
 - 134 • As I wanted to look at 11 am –midnight (for lunch and dinner related data), I :
135 ○ clicked the bottom left two boxes
136 ○ copied the date 2024-05-09 over the default date in the bottom right two data fields
137 ○ after T (for time), I entered the desired time of analysis AFTER conversion into my local
138 time (Central EU summer time minus 2 hours = UTZ; so to look at 11 to midnight of
139 my AAPS screen, I must enter here 09.00:00Z, and below it 22:00:00Z).

140 Entries at the bottom are not mandatory, but when clicking these little boxes (bottom left)
141 you can define a start and/or an end-point for analyzing, within the logfiles specified in the
142 field above.

143
144
145

10.2.4 Run emulation

Now we are ready to go: Press “Run Emulation”


This produces sometimes an error message (e.g. if you have a syntax error, or incompatible software versions: => seek help, in the Github materials provided by ga-zelle, or in Discord/Full-Closed-Looping/emulate-aaps here: <https://discord.gg/n3tD5eXExC>)

After a short moment results should show up, which you can look into in a couple of ways:

10.2.4.1 SMB tab contents in (date..) noChange.txt result file

This basically gives you “all the SMB tabs” without needing to make screenshots every 5 minutes. Search options help find what lines are of interest to your analysis:

By using the search function you can jump, in that super long list, to all places that e.g. have „autoISF“ in it or „script debug“, or „SMB disabled“ (if you want to analyze when that happened...)

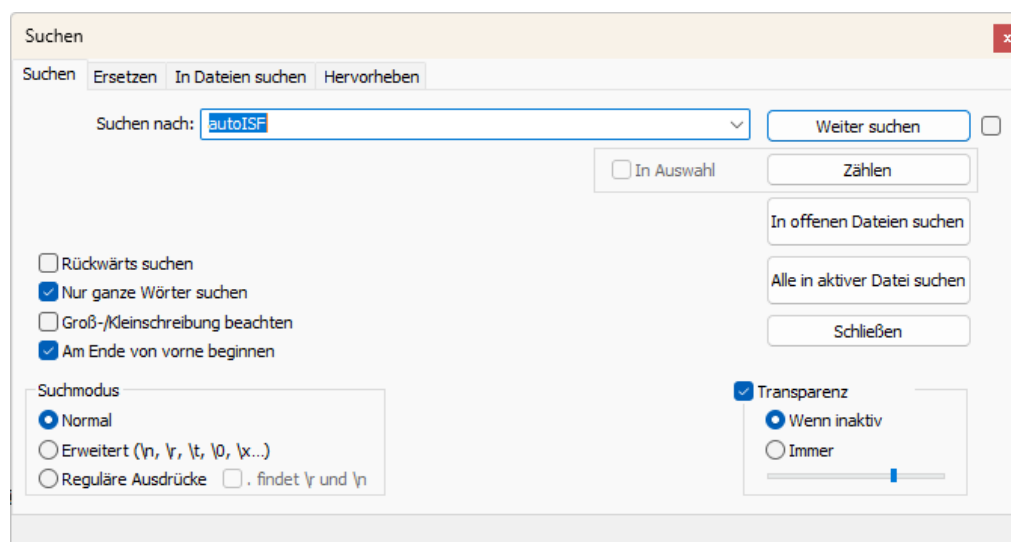

10.2.4.2 Table of results (...noChange.csv file)

The .csv file in your project folder gives a tabular presentation of how parameters like bg, iob, iobTH, the various ISF contributors, bg target, insulinRequired etc. develop every 5 minutes, and what SMB size and %TBR resulted.

It is a vast table, so you may want to reduce it to something more “digestable”:

10.2.4.2 Analysis of the (date..) noChange.csv table in Excel or LibreOffice calc.

Best copy the entire table into a new .xls or .ods sheet, where you can:

- add right next to the standard world time your corresponding “AAPS time”

174 For instance, adding +2/24 translates the UTZ column into central European summer time
175 column next to it (where currently a row of Z stands). Likewise, subtract like -5/24 from UTZ
176 for an US East Coast time scale.

177 Highlight all time fields (the entire columns), and switch from hh:mm:ss format to hh:mm.
178 (While the seconds are important for the loop's calculations, for our comparison with
179 Nightscout or other charts and data, it is much easier without the seconds attached)

180 • **hide** any column you find less important to look at for your intended analysis

181 That way, “boxes” (data fields) retain their original position in tables

182 Also, in case later you want to look into additional info, you can simply un-hide the relevant
183 columns ... or time lines:

184 • **hide** lines (time segments) you find less important to look at for your intended analysis

185

186 Usually you will color mark where relevant SMBs were given, which of the ISFs (and underlying
187 weights) was strongly contributing (note that this can be good or not good) . Also where iobTH was
188 exceeded, whether an Automation kicked in e.g. setting a TT, when there were periods with zero
189 insulinRequired.

190 In [section 10.3.4](#) we present an extra tool that does a standardized table reduction and color marking
191 for you!

192

193 You may be able to formulate a hypothesis or two, what settings (...ISF_weights, iobTH%,
194 SMB_range_extention, autoISFmax ...) should be changed for improvement (then go to [10.3](#))

195

196

197 [10.2.4.4..](#) Graph **noChange.pdf**

198


199 When you open the emulator, before executing any analysis, you can make settings under “Select Graphics
200 Options”.

201 Select your preferred way of outputting decimals (point or comma).

202 Then select whether you want “All” possible outputs in the graph, or “Most” = all except those you tick “off”
203 in the boxes for each output parameter.

204 If you would use “Some/just a few”, you would have to tick those few you that do want to see, by ticking the
205 corresponding boxes.

206 Recommendation is to look at (nearly) everything offered (as your default setting that you can leave
207 untouched in most of your emulator runs):

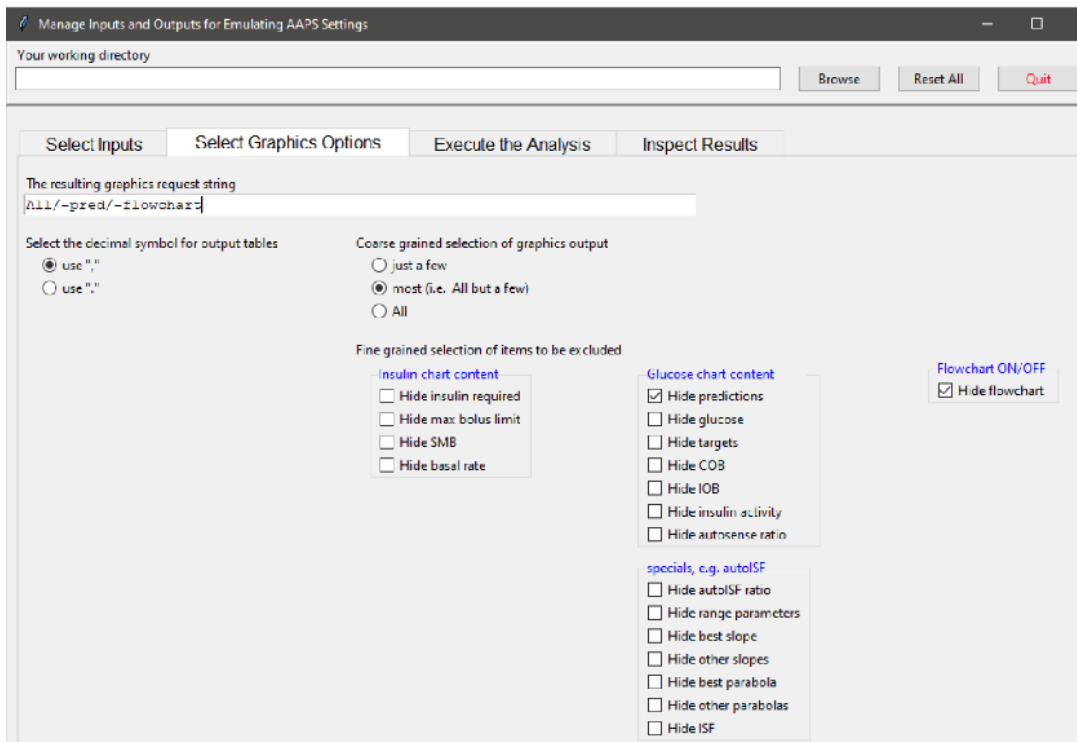

208
209
210 After your emulation run, under Inspect Results, you can open the pdf file that is last in the results list
211 offered.

212
213 This **noChange.pdf** is a chart that shows along the time axis (down), from right to left:


- 214
- Red: the bg curve
- 215
- Yellow: the bg target (note that I do no manual “EatingSoonTT” but for bg rises over +10 mg/dl
216 I have an Automation that sets low TT for a couple of minutes)
- 217
- Light blue corridor: Left edge is set iobTH, and bandwidth +30% (would be +20% at elevated
218 TT)

- 219
- Dark blue line: iob (exceeding twice the iobTH, with temp. SMB shut-off

220 As bg did not convincingly come down enough, one could hypothesize that iobTH should be
221 elevated. ((But, again, this would have to be confirmed also with other kinds of meals)).

- 222
- Thin yellow line: Insulin activity
- 223
- Green dotted line: ISF as would result from AAPS w/Autosens
- 224
- Green scatter points: autoISF ISF no Chage (lighter points) or what-if (darker points)
- 225
- Black line: Profile ISF

- 226
- 227
- Gray scatter points: ISF weakened (to the left of black line) or strengthened (to the right)
 - Orange line: cob=0 at all times (in FCL)

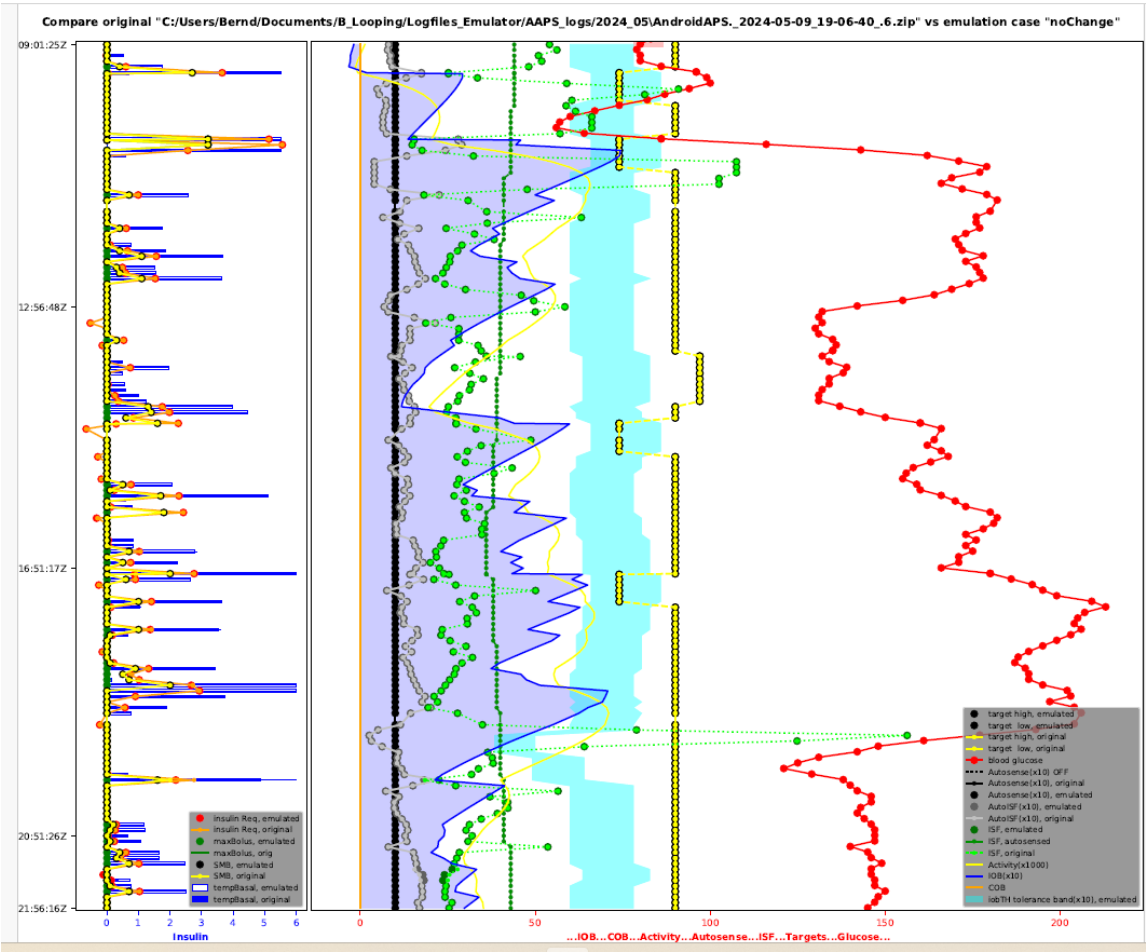

228

229

230

231

232

233

234

235

236

More see discussed together with (yourChanges).pdf in [section 10.3.3.4](#)

237 238 “What-if” analysis using the emulator 239

240 In the following you see an example how you can analyze a day of logfiles, and selecting the time span of
241 interest, for instance 11-24 h to look at how autoISF managed lunch and dinner.

242
243 You have to go through the emulator twice. You already did it ([section 10.2.](#)) using the no change.vdf, and
244 now start over with the same emulator with your (desired changes).vdf. Repeat, if you have two or more such
245 vdf defined. (Just clear results before executing analysis each time. No worries, all results are automatically
246 captured for all runs, all in your selected study file).

247 248 10.3.1 Define your investigated changes in a, or in several, (yourChanges).vdf 249

250 Define for which one to max three parameters in your current profile settings you want to look into a
251 different setting. Recommendation is to use a factor, like for example current setting * 0.9 , or current
252 setting * 1.2, and use that in your naming for this vdf file, too.

253
254 Within the same study, you can make several runs with several vdf files.

255 All results, like the csv results table, will appear then several times in your study file, only with different
256 name endings as in the underlying vdf.

257
258 Example: I like to check in my actual data (they are in my noChange.vdf emulator run), **in which time**
259 **points the following parameter changes would make a** (how) big **difference** in the loop’s decision:

- 260 • 20% higher bgAccel_ISF_weight to boost the first SMBs stronger: How would that tend to ramp
261 up early iob; and might that get too strong in other parts of the data? Or does it bounce into a
262 restriction (maxSMB size; autoISFmax; iobTH...) that I might need to widen?
- 263 • Doubling my cautiously set bgBrake_ISF_weight shall give me insight into the workings of that
264 parameter (and whether using a much smaller weight than for bgAccel_ISF_weight is really
265 what I should keep doing)
- 266 • As my bg came down from a persistent high quite slowly, I elevate the dura_ISF by 20%

267
268 Actually, it would make more sense to first find my “optimal”, maybe indeed elevated,
269 bgAccel_ISF_weight. *Then*, do a noChange (!) run **with that**, plus a (yourChanges) run with the stronger
270 dura weight, investigated on that basis.

271 Reason: 1) As we always say, better do only one change at a time. 2) A better job with bg control via
272 bgAccel_ISF will reduce the peak height and provide a different (easier) scenario for dura_ISF to manage.

273

274 Now, to **write** your **(yourChanges). vdf for the emulator** (this is same procedure as you did in section
275 10.2.1 for the noChange.vdf):

- 276
- 277 • just open Notepad++ (from list of all programs on your PC) to create a new vdf.

278 Alternatively you can also take another pre-existing vdf file, and copy it into your current project
279 giving it a new name (re-name it)

- 280 • name your vdf (in our example: 1.2_bgAccel_2.0_bgBrake_1.2_dura.vdf) ...


- 281 • ... and store that in a file of your current emulator project you are about to start (see my storage
282 path in top line here)

283 **Caution:** Make absolutely sure (best by looking it up in the SMB tab, down in the profile set
284 section) to **spell each term exactly** as your loop uses it (probably w/ decimal points, not comma)

- 285
- 286 • ...when you make one line per parameter (separating entries with spacers->):

287 profile->(parameter) ->->profile['(parameter)']* (factor)->->###(comment as you like)

288
289 The (yourChanges) .vdf should look like something like this:
290

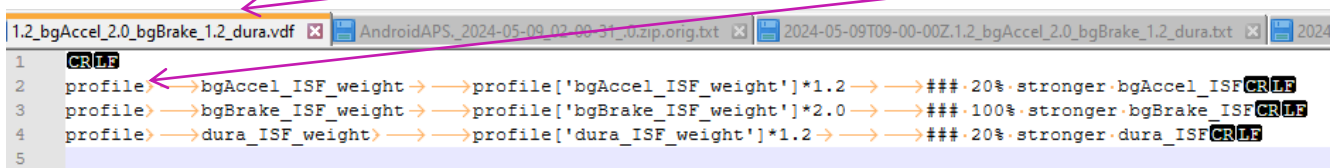


```
1 CR LF
2 profile ->bgAccel_ISF_weight ->->profile['bgAccel_ISF_weight']*1.2 ->->###.20%.stronger.bgAccel_ISF CR LF
3 profile ->bgBrake_ISF_weight ->->profile['bgBrake_ISF_weight']*2.0 ->->###.100%.stronger.bgBrake_ISF CR LF
4 profile ->dura_ISF_weight ->->profile['dura_ISF_weight']*1.2 ->->###.20%.stronger.dura_ISF CR LF
5
```

291
292 CR = LF= Erase any entries after CR LF and also any entries in lines below, if any

293 294 10.3.2 Run the emulator with (yourChanges).vdf

295
296 The “what-if” emulator run is done the same way as you did the noChanges run (section 10.2), however,
297 now , the **(yourChanges).vdf** must be loaded into the 2nd input field, where formerly you had the
298 noChange.vdf.:

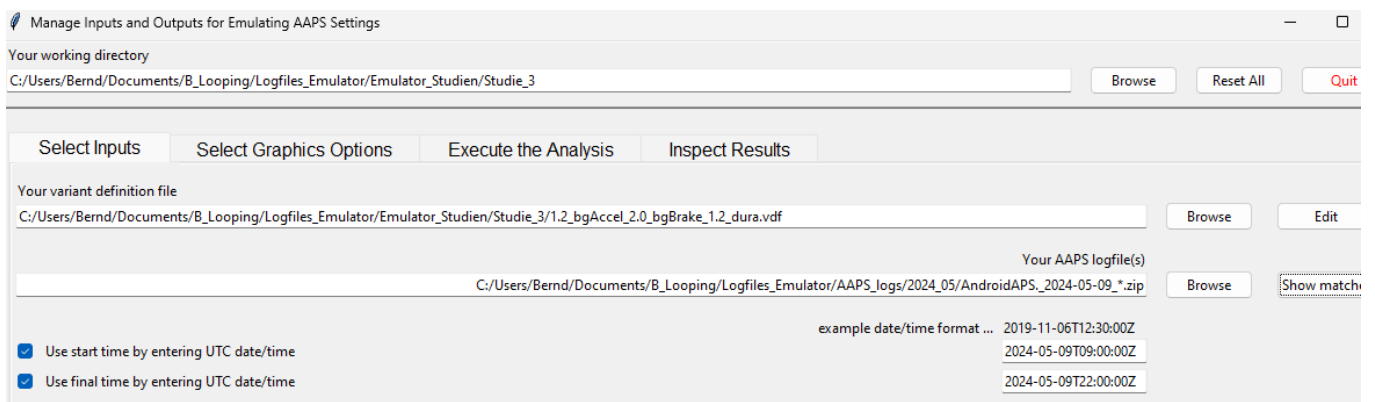

- 301 In the 3rd input field, give the path to your stored logfiles. A good way to do this is:
- 302 • Browse in your Windows Explorer to any logfile from the desired day (2024-05-09 in above ex-
- 303 ample)
- 304 • Replace the time with an asterix * (this means you look at all-day data, in UTZ time). Check
- 305 whether this will work by pressing Show matches . You should see all logfiles from that day in a
- 306 pop-up info box.
- 307 • As I wanted to look at 11 am –midnight for lunch and dinner related data, I :
- 308 ○ clicked the bottom left two boxes
- 309 ○ copied the date 2024-05-09 over the default date in the bottom right two data fields
- 310 ○ after T (for time), I entered the desired time of analysis AFTER conversion into my local
- 311 time (Central EU summer time minus 2 hours = UTZ; so to look at 11 to midnight of
- 312 my AAPS screen, I must enter here 09:00:00Z, and below it 22:00:00Z).

313 After making these entries, press Execute the Analysis, (evtl also Clear old Data) and then press Run

314 Emulation, I can look the results up under “Inspect Results”:

315

316 10.3.3 Emulation results

317

The screenshot shows a software window titled 'Your working directory' with a path 'C:/Users/Bernd/Documents/B_Looping/Logfiles_Emulator/Emulator_Studien/Studie_3'. Below this is a tabbed interface with four tabs: 'Select Inputs', 'Select Graphics Options', 'Execute the Analysis', and 'Inspect Results'. The 'Inspect Results' tab is active, showing a list of files with their paths and actions (Browse, Show).

File Type	File Path	Actions
*.log - Your file showing edits from the variant assignments	C:/Users/Bernd/Documents/B_Looping/Logfiles_Emulator/Emulator_Studien/Studie_3/2024-05-09T09-00-00Z.1.2_bgAccel_2.0_bgBrake_1.2_dura.log	Browse Show
*.csv - Your table comparing key values of original vs emulation	C:/Users/Bernd/Documents/B_Looping/Logfiles_Emulator/Emulator_Studien/Studie_3/2024-05-09T09-00-00Z.1.2_bgAccel_2.0_bgBrake_1.2_dura.csv	Browse Show
*.delta - Your table comparing bg deltas of original vs emulation	C:/Users/Bernd/Documents/B_Looping/Logfiles_Emulator/Emulator_Studien/Studie_3/2024-05-09T09-00-00Z.1.2_bgAccel_2.0_bgBrake_1.2_dura.delta	Browse Show
*.orig.txt - Your short log of original analysis	C:/Users/Bernd/Documents/B_Looping/Logfiles_Emulator/Emulator_Studien/Studie_3/2024-05-09T09-00-00Z.orig.txt	Browse Show
*.txt - Your short log of emulated analysis	C:/Users/Bernd/Documents/B_Looping/Logfiles_Emulator/Emulator_Studien/Studie_3/2024-05-09T09-00-00Z.1.2_bgAccel_2.0_bgBrake_1.2_dura.txt	Browse Show
*.pdf etc. - Your graphic file comparing key values of original vs emulation	C:/Users/Bernd/Documents/B_Looping/Logfiles_Emulator/Emulator_Studien/Studie_3/2024-05-09T09-00-00Z.1.2_bgAccel_2.0_bgBrake_1.2_dura.pdf	Browse Show

318

319 All results from your (yourChanges).vdf emulator go automatically where the noChange.vdf results are

320 already stored, in our example into the “Studie 3” file, below:

321

322 Besides the 1.2_bgAccel_2.0_bgBrake_1.2_dura.vdf case which I like to look into for the present high carb

323 meal, I also prepared another vdf that investigates a factor 1.2 stronger pp_ISF and a weaker, factor 0.8,

324 bgAccel_ISF (with the intention to test this, and a noChange, on a low carb meal later.

325

326

B_Looping > Logfiles_Emulator > Emulator_Studien > Studie_3				Studie_3 durchsuchen
Sortieren ▾ Anzeigen ▾ ...				
Name	Änderungsdatum	Typ	Größe	
2024-05-09T09-00-00Z .noChange.pdf	15.05.2024 17:17	Adobe Acrobat-D...	77 KB	
2024-05-09T09-00-00Z .noChange.csv	15.05.2024 17:17	Microsoft Excel C...	51 KB	
2024-05-09T09-00-00Z .noChange.log	15.05.2024 17:17	Notepad++ Docu...	35 KB	
2024-05-09T09-00-00Z .noChange.txt	15.05.2024 17:17	Notepad++ Docu...	281 KB	
2024-05-09T09-00-00Z .orig.txt	15.05.2024 17:17	Notepad++ Docu...	281 KB	
2024-05-09T09-00-00Z .1.2_pp_0.8_bgAccel.pdf	15.05.2024 17:16	Adobe Acrobat-D...	78 KB	
2024-05-09T09-00-00Z .1.2_pp_0.8_bgAccel.csv	15.05.2024 17:16	Microsoft Excel C...	51 KB	
2024-05-09T09-00-00Z .1.2_pp_0.8_bgAccel.log	15.05.2024 17:16	Notepad++ Docu...	57 KB	
2024-05-09T09-00-00Z .1.2_pp_0.8_bgAccel.txt	15.05.2024 17:16	Notepad++ Docu...	281 KB	327
2024-05-09T09-00-00Z.1.2_bgAccel_2.0_bgBrake_1.2_dura.csv	17.05.2024 21:29	Microsoft Excel C...	51 KB	
2024-05-09T09-00-00Z.1.2_bgAccel_2.0_bgBrake_1.2_dura.log	17.05.2024 21:29	Notepad++ Docu...	66 KB	
2024-05-09T09-00-00Z.1.2_bgAccel_2.0_bgBrake_1.2_dura.pdf	17.05.2024 20:40	Adobe Acrobat-D...	78 KB	328
2024-05-09T09-00-00Z.1.2_bgAccel_2.0_bgBrake_1.2_dura.txt	17.05.2024 21:29	Notepad++ Docu...	282 KB	
1.2_bgAccel_2.0_bgBrake_1.2_dura.vdf	17.05.2024 20:38	VDF-Datei	1 KB	329
1.2_pp_0.8_bgAccel.vdf	10.05.2024 21:55	VDF-Datei	1 KB	
noChange.vdf	07.05.2024 19:57	VDF-Datei	1 KB	

10.3.3.1 Logs check: log vs txt

The **noChange.log** has all the info your series of SMB tabs had that day.

How to search in this vast list is shown elsewhere (see [section 10.2.4.3](#)).

Likewise, the **(yourChanges).log** gives for each loop decision in all detail how and why that decision would have changed with the different parameter inputs you are checking out here

In the two examples here, , it was a check on the difference

- a 20% stronger pp_weight and 20% weaker bgAccel_weight
- a 20% stronger weight for both, bgAccel_ and dura_ISF, and a doubling of bgBrake_weight

would make.

Note that all these “what if” data can only give rough hints, notably about **the first** greater change that you would see with the investigated changed setting. So it works well for investigating how to ramp up iob quickly after detection of acceleration.

However, any relevant change would put your bg curve on a different trajectory, so that would influence all following results. Therefore, what you get here is **not** a complete modelling how your bg would have developed in the alternative scenario.

350 But you can investigate in which stages the parameter(s) you are looking at in your current “what-if” had big
351 influence, and in which direction the changes would go. (see also charts shown in [section 10.3.3.4](#)).
352 Analyzing how to safely come down from a high glucose plateau while limiting hypo danger towards the end
353 of digestion is also to some extent possible.

354

355 A good other way to employ the what-if analysis is real time, on your smartphone, using speech synthesis
356 (see [section 11](#)): Then you get real-time info, as to exactly when a significantly different proposal would
357 emerge, and can decide (and watch!) real-time whether to follow the new idea and not was probably better.

358

359 Observe that a setting change must work well for you

360 • not just in one point of time, and

361 • not just for one kind of meal,

362 but you must look at all time slots in the investigated meal, plus analyze with the same tool a totally different
363 meal within your usual spectrum, how things work out there

364

365 [10.3.3.2 Tabular results](#)

366

367 The **noChange.csv** table gives all relevant data. Besides development of bg and iob you see the calculated
368 insulinRequired in each loop decision, and how each of the autoISF categories contributed to the decision
369 (notably regarding SMB size).

370

371 The **(your changes).csv** shows in detail how **every single** loop **decision** would be influenced by the different
372 settings you are investigating. To inspect that huge table, click on the Z behind the start UTC time entry
373 (see black box in the Z column of the following table

374 If you like to see the bg in each screen, too, go 3 or 4 columns farther to the right with your black
375 box.

376 Then, go to window/fix. Now you can scroll through the data and always see headline and time (or time and
377 bg level).

378 To further ease analysis, feel free to temporarily erase (hide) any columns that you (think you) do not
379 need for the intended analysis. More suggestions see in [section 10.2.4.2](#)

380

id	UTC time	UNIX time	bg accel	bg brake	target low	orig	emul	emul	cob	job	emul	emul	act	orig	orig	dura	min	utes
0	0,3759837963	1715245285,9	80		90	90	90	90	0	-0,18	6	7,8	0,001	1	0,81	10		
1	0,3793981481	1715245580,3	79	79	90	90	90	90	0	-0,23	6	7,8	0,001	1	0,78	15		
2	0,3829166667	1715245884,2	80		90	90	90	90	0	-0,28	6	7,8	0	1	0,86	20		
3	0,386412037	1715246186,2	80		90	90	90	90	0	-0,28	6	7,8	-0,001	1	0,85	25		
4	0,3898263889	1715246482	86		90	90	90	90	0	-0,33	6	7,8	-0,001	1	0,91	0		
5	0,393587963	1715246806,8	96		74	74	74	74	0	0,18	6,6	8,58	-0,001	1	1,75	0		
6	0,3942013889	1715246859,5	96		74	74	74	74	0	2,92	6,6	8,58	0	1	1,75	0		
7	0,3968055556	1715247084,8	99		74	74	74	74	0	2,92	6,6	8,58	0,005	1	1,31	5		
8	0,4002430556	1715247381,7	100	100	74	74	74	74	0	2,83	6,6	8,58	0,011	1	0,75	10		
9	0,4037384259	1715247683,2	94	94	74	74	74	74	0	2,72	6,6	8,58	0,015	1	0,48	0		
10	0,4071643519	1715247979,7	87	87	74	74	74	74	0	2,59	6,6	8,58	0,018	1	0,54	0		
11	0,4107407407	1715248288,8	82	82	74	74	74	74	0	2,45	6,6	8,58	0,02	1	0,73	0		
12	0,4141435185	1715248583	74	74	90	90	90	90	0	2,3	6	7,8	0,022	1	0,75	0		
13	0,4176273148	1715248883,1	67	67	90	90	90	90	0	2,14	6	7,8	0,022	1	0,7	0		
14	0,4210300926	1715249177,9	60	60	90	90	90	90	0	1,99	6	7,8	0,023	1	0,65	0		
15	0,4245949074	1715249485,3	57	57	90	90	90	90	0	1,82	6	7,8	0,023	1	0,65	0		
16	0,4280439815	1715249783,8	56	56	90	90	90	90	0	1,67	6	7,8	0,022	1	0,65	5		
17	0,4315277778	1715250084,6	64		90	90	90	90	0	1,51	6	7,8	0,021	1	0,75	0		
18	0,435	1715250384,3	86		74	74	74	74	0	1,36	6,6	8,58	0,02	1	2,8	0		
19	0,4355787037	1715250434,4	86		74	74	74	74	0	4,58	6,6	8,58	0,021	1	2,8	0		
20	0,4384953704	1715250686,1	116		74	74	74	74	0	4,45	6,6	8,58	0,026	1	2,9	0		
21	0,4419675926	1715250987	143		74	74	74	74	0	7,44	6,6	8,58	0,038	1	2,43	0		
22	0,4454513889	1715251287,7	162	162	74	74	74	74	0	7,4	6,6	8,58	0,048	1	1,33	0		
23	0,4490046296	1715251594,5	171	171	74	74	74	74	0	7,15	6,6	8,58	0,055	1	0,4	0		
24	0,452349537	1715251883,2	179	179	74	74	74	74	0	6,84	6,6	8,58	0,06	1	0,4	5		
25	0,4558217593	1715252183,5	177	177	90	90	90	90	0	6,48	6	7,8	0,063	1	0,4	10		
26	0,4594328704	1715252495,4	169	169	90	90	90	90	0	6,1	6	7,8	0,065	1	0,4	15		
27	0,4630208333	1715252805,5	166	166	90	90	90	90	0	5,72	6	7,8	0,066	1	0,4	5		
28	0,4662847222	1715253087,3	172	172	90	90	90	90	0	5,37	6	7,8	0,065	1	0,86	25		
29	0,4697453704	1715253386,2	179		90	90	90	90	0	5	6,36	8,27	0,064	1	2,26	5		
30	0,4731828704	1715253683,7	182		90	90	90	90	0	5,54	6,36	8,27	0,064	1	1,33	10		
31	0,480162037	1715254286,4	180		90	90	90	90	0	4,81	6,36	8,27	0,063	1	1,13	0		
32	0,4837731481	1715254598,6	176	176	90	90	90	90	0	4,44	6,36	8,27	0,061	1	0,65	5		
33	0,4873032407	1715254903,9	176		90	90	90	90	0	4,09	6	7,8	0,059	1	1,13	10		
34	0,4905555556	1715255184,4	177		90	90	90	90	0	3,78	6	7,8	0,056	1	1,67	15		
35	0,4940625	1715255487,9	173		90	90	90	90	0	3,98	6	7,8	0,054	1	1,26	20		
36	0,4974884259	1715255783,9	170	170	90	90	90	90	0	3,69	6	7,8	0,053	1	1,07	25		
37	0,5009722222	1715256084,7	171		90	90	90	90	0	3,39	6	7,8	0,05	1	1,38	30		

381
382

383 Still, the csv tables are overwhelming. You could proceed in either of two directions now:

384

385 A) Convert both (or all 3) csv files into one table in Excel or into Libre office calculator. Hide columns
386 (and eventually also lines) that are of no particular interest for your analysis. Mark differences be-
387 tween noChange and (yourChanges) column data with color, add extra columns with additional cal-
388 culations ...

389 This route is good to compare quantitative impacts of autoISF categories in critical time points.


390 B) For the core data relevant to assessing your autoISF settings, there is an extra tool for convenient
391 analysis (see next section)

392


393

394 10.3.3.3 Automated extraction from tabular results

395 **autoISF_factors_proforma.ods** is provided as an **extra tool** that you download from here:

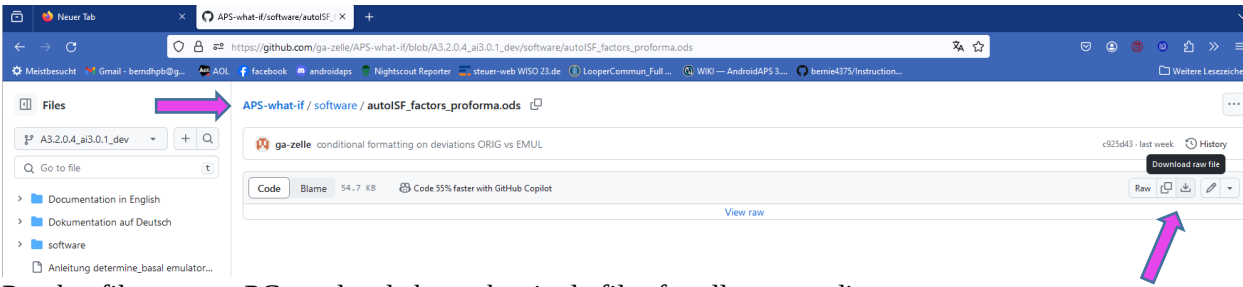


396 Put that file on your PC one level above the single files for all your studies:

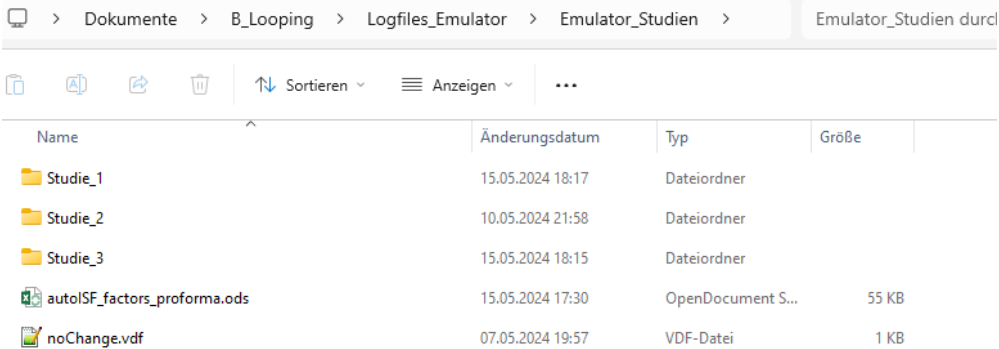

398 Now, if we want to use this tool on the two csv files of our Studie_3 file, we must proceed as follows (for
399 each of the two .csv files, *separately*):

- 401 1. Click on the .csv file and open in Libre office calculator.
- 402 2. Make sure the time column is set to US_English:

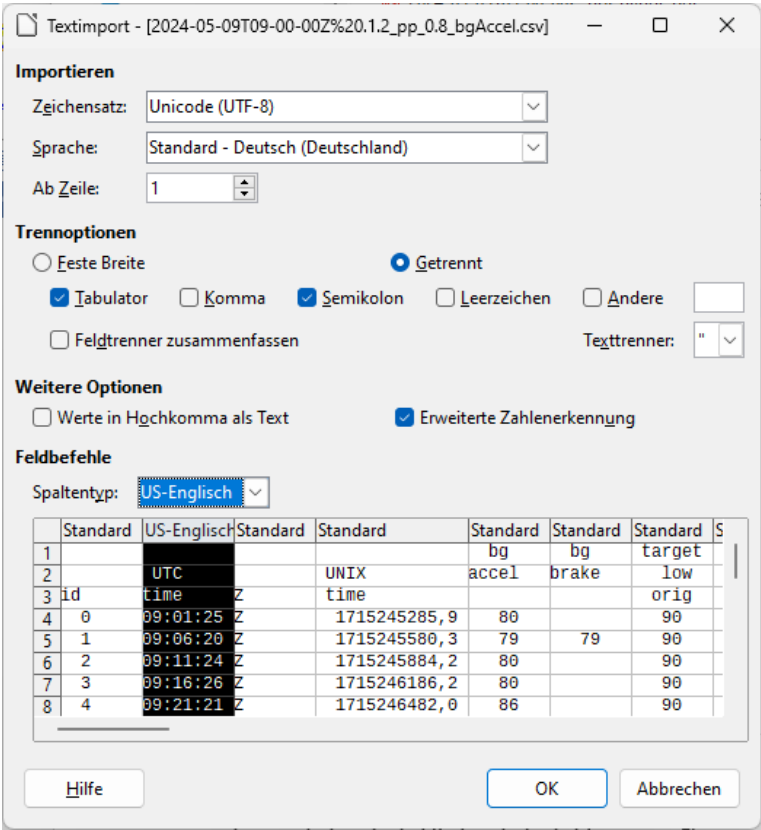

404
405 3. Now start, in Libre office calculator, the autoISF_factors_proforma.ods ...

406 This turns the first 30-some lines of your csv table (left side) into a form in which important effects are
407 highlighted in color, and formatting is improved:



408
409
410 Now, you want this for the entire table.

411 In the autoISF_factors_proforma table, highlight 20 or more lines (not including the first or last), and mouse
 412 right hand/insert above ...

The screenshot shows the LibreOffice Calc interface with the file 'autoISF_factors_proforma.ods' open. A range of cells is selected, and a context menu is displayed with the option 'Zeile(n) oberhalb einfügen' (Insert row(s) above) highlighted. The spreadsheet contains columns for various factors like 'auto', 'final', 'dura', 'min-utes', 'lin.fit', 'parab', etc. The status bar at the bottom indicates 'Tabelle 1 / 2 | 29 Zeilen, 1024 Spalten ausgewählt'.

413
 414
 415 Do this as often as you need to create the number of lines that your emulated csv file comes with.
 416 If you ended up with too many lines, erase the superfluous number (any four, in the example):

autolSF_factors_proforma.ods - LibreOffice Calc

Datei Bearbeiten Ansicht Einfügen Format Extras Daten Fenster Hilfe

Arial 10

A128:AMJ131

	A	B	O	P	Q	R	S	T	U	V	W	X	Y	Z	AA	AB	AC	AD
1				auto	final	dura		lin.fit		parab	parab	parab	parab	auto	acce	bg	pp	d
2		UTC		sens	ISF	min-	dura	min-	lin.fit	fit	fit	fit	fit	sens	ISF	ISF	ISF	ISF
3	id	ime	act	orig	orig	utes	avg.	utes	delta	correl	durat	last-Δ	next-Δ	emul	emul	emul	emul	e
121																		
122																		
123																		
124																		
125																		
126																		
127																		
128																		
129																		
130																		
131																		
132																		
133																		
134																		
135																		
136																		
137																		
138																		
139																		
140																		
141																		
142																		
143																		
144																		
145																		
146																		
147	18	11:34:21	0	1	1,28	0	149	5	8,1	0,99	20,1	10,03	11,39	1	1,28	1,06	1,16	
148	19	11:39:14	0	1	0,4	5	149	25	7,1	1	15	1,05	-4,39	1	0,38	1,06	1	
149	20	11:44:15	0	1	0,57	10	148,7	40	5,4	0,99	20	-2,41	-6,96	1	0,5	1,06	1	
150	21	11:49:16	0	1	1,11	20	146,6	10	-1,5	1	14,9	-2	-3,01	1	0,88	1,06	1	
151	22	11:54:16	0	1	1,1	25	145,7	10	-3,5	0,99	15	-4,6	-6,6	1	0,77	1,05	1	
152	23	11:59:14	0	1	0,86	5	138,5	10	-5	1	24,9	-5,53	-7	1	0,83	1,03	1	
153	24	12:04:14	0	1	1,22	10	136,7	15	-4,4	1	15	-3,4	-2,4	1	1,22	1,03	1	
154	25	12:08:45	0	1	1,21	10	136,7	5	-3	1	15	-3,4	-2,4	1	0,88	0,99	1	
155	26	12:10:37	0	1	1,02	15	136	5	1	0,99	15	0,6	3,6	1	1,63	0,99	1	
156	27	12:14:14	0	1	0,67	15	133,3	35	-3	0,96	30	-1,64	-1,07	1	0,96	0,98	1	

your_title graphs

Tabelle 1 / 2 4 Zeilen, 1024 Spalten ausgewählt PageStyle_2023-02-20T20.empty Summe=0 100 %

Then just copy it in, by selecting all data lines in the emulated csv, and pasting (paste special, values only) into box A4 of your “elonged” autoISF_factors_performa.ods.


The bottom tab “your_title” should be re-named by you, best with day of log you analyze, and your what-if parameters (so, the name of your csv file could be put in here)

Now you have a table with optimized lay-out that incorporates key data from both your no change AND of your investigated changes.csv files.

A super neat extra feature is already pre-programmed, which you can see if you click on the bottom tab “graphs”.

The top graph is the bg curve (the actually seen bg).

Note that for the what-if no bg development over the time range is available. (The noChange one is also given there).

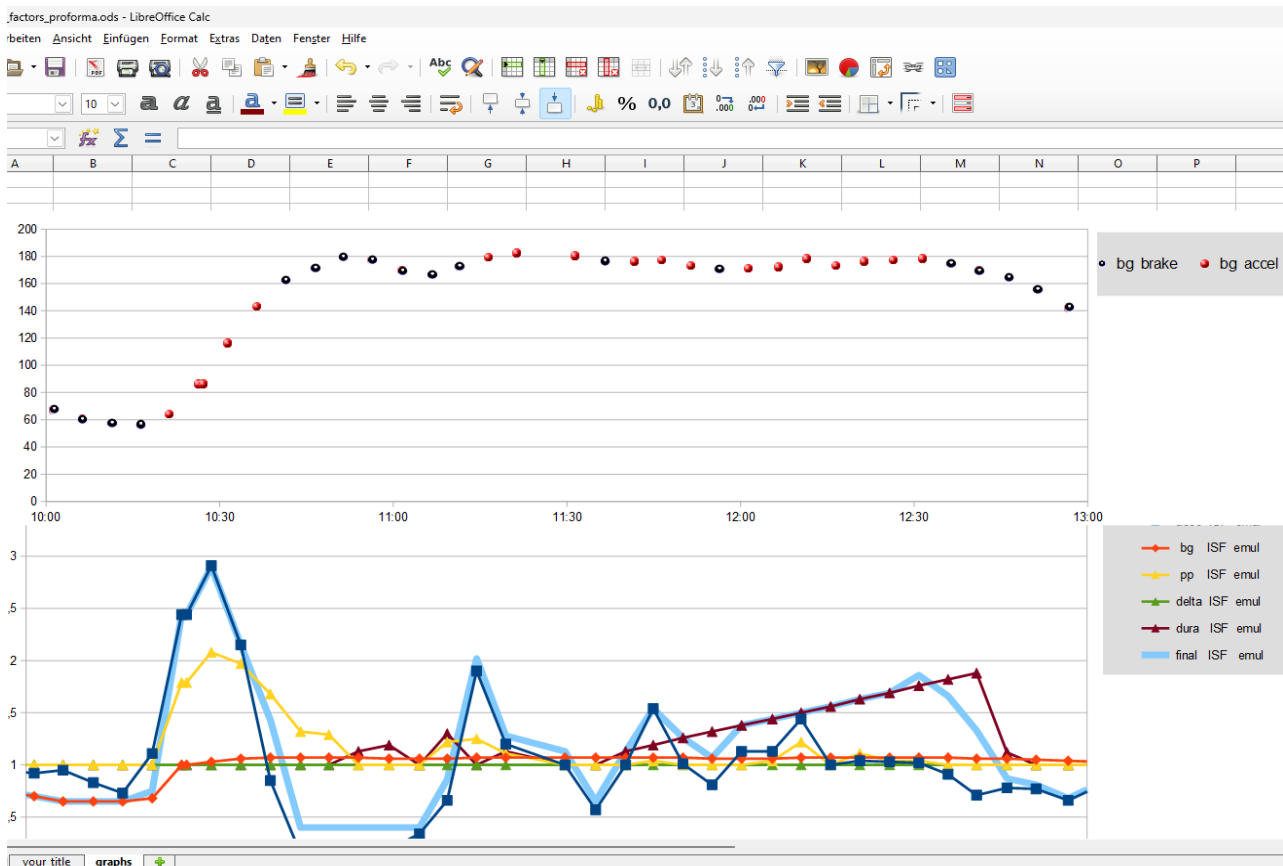

431
432 The bottom graph (do one for each, the noChange or the (yourChanges) case) shows the amplification factors
433 coming from each autoISF category, and the overall resulting ISF amplification.

434 You probably have to widen the time scale (double click on the time axis, and type the desired time
435 span (min and max UTC)(and spacing of data points, 00:30:00 or 00:15:00) into this box:

X-Achse

Skalierung Positionierung Liniendiagramm Beschriftung Zahlen Schrift Sch

Skalierung

☐ Richtung umkehren

☐ Logarithmische Skalierung

Minimum 10:00:00 ☐ Automatisch

Maximum 15:30:00 ☐ Automatisch

Hauptintervall 00:30:00 ☐ Automatisch

Hilfsintervall Schritte 5 ☐ Automatisch

OK Abbr

436
437 In the given example above, the 2.5 hours were not enough yet to analyze this 10:30 UTC (12:30 AAPS)
438 lunch; we need to look until bg is near target (hopefully before dinner starts).

439
440 A similar graph is available on the (i-)phone if you use the autoISF dev variant of iAPS (and maybe of Trio,
441 in the near future). (Later insert details here, or in [section 11](#))

10.3.3.4 Chart coming with the Emulator

443

444 In case you find the extra steps described in the preceding section “too much”, also the emulator offers one
445 chart (the pdf offered at the bottom of the screen as shown below the “[10.3.3](#) Emulaton results” headline).

446

447 First look at the initial bg rise in the noChange.pdf chart (emulation results from your noChange.vdf run),
448 and see how bgAccel_ISF and pp_ISF acted, or could have acted in improved ways.

449 Then look into in (yourChange).pdf to see potential effects (or what other change to try). (Actually, you
450 probably will have to go into a detrailed analysis of several lines and columns of the tables as discussed in
451 sections [10.3.3.2](#) and [10.3.3.3](#)).

452

453 Note that ideally we want FCL coverage of our entire “normal day” meal spectrum by one set of
454 settings. So, **not much is gained if you put a lot of effort in optimizing FCL settings for one**
455 **meal.**

456

457 You will need iterations. Do such analysis for **two or three very different meals** that you wish the
458 algorithm to automatically handle. See [section 4.2/4.3](#) on how meals with very different carb loads
459 might benefit or also suffer from too aggressive or to mild (category)_ISF_weights you could set.

460

461

462 The initial iob received might be limited by allowed SMB sizes, autoISFmax, or the (dynamic!) iobTH. You
463 will have to look into the data table to find out about this (a quick orientation - notably regarding the light
464 blue iobTH band, see next page - is also possible in the pdf result files you have in your project file (project
465 file example “Studie 3” in 2nd chart under the [10.3.3](#). headline).

466

467 Only once you found OK weights for bgAccel- and pp_ISF_weights, does it make sense to go tune the
468 dura_ISF_weight. 12:00 – 12:45 UTC in above graph, the resulting effective ISF is dominated by dura_ISF.
469 Just judging from the picture, a stronger weight might be worth trying. However, we really need to see the
470 insulinRequired calculation and the further development because impatience about bringing bg values down
471 faster too often results in hypoglycemia later.

472

473 The **noChange.pdf** is a chart that shows along the time axis (down), from right to left:

474

- Red: the bg curve
- Yellow: the bg target (note that I do no manual “EatingSoonTT” but for bg rises over +10 mg/dl I have an Automation that sets low TT for a couple of minutes)
- Light blue corridor: Left edge is set iobTH, and bandwidth +30% (would be +20% at elevated TT)

475

476

477

478

479 • Dark blue line: iob (exceeding twice the iobTH, with temp. SMB shut-off)

480 As bg did not convincingly come down enough, one could hypothesize that iobTH should be
481 elevated. ((But, again, this would have to be confirmed also with other kinds of meals)).

482 • Thin yellow line: Insulin activity

483 • Green dotted line: ISF as would result from AAPS w/Autosens


484 • **Green scatter points:** autoISF ISF no Change (lighter points) or **what-if (darker points)**

485 Foreseeably, this is the strongest difference between our noChange (left) and 120% bgAcel_ISF_weight
486 (right) in the picture below. (Note the red bg curve is *both times* the really seen bg, because the what-if
487 case only looks at each single loop decision). The first (->) time the dark green dot is far to the right, this
488 would get the bg down, we *would start to see* a (<-) bg lowering effect, shifting the red curve to the left

489 • Black line: Profile ISF

490 • Gray scatter points: ISF weakened (to the left of black line) or strengthened (to the right)

491 • Orange “baseline” (only in FCL, as cob=0 at all times)

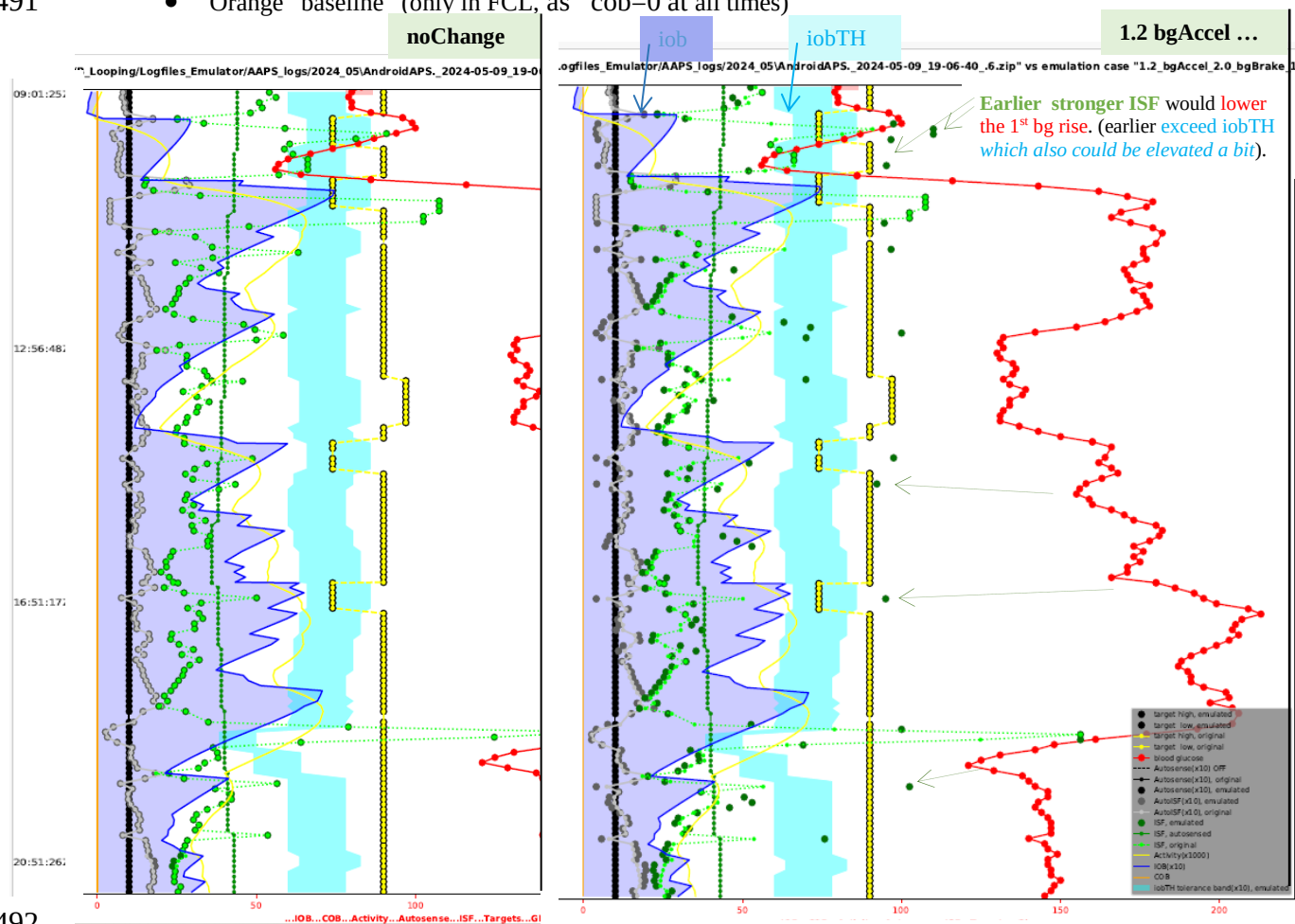

492
493 Regarding the other changed parameters: Stronger dura_ISF would suggest more insulin towards the end of
494 plateaus; this should have helped in the 1st plateau (red curve, top right quadrant of the picture). However,
495

496 same setting would have to work also on 2nd plateau; the chart cuts off there, so too early to see whether a
497 hypo danger might result.

498 Effect from doubling the bgBrake_ISF effect are hard to evaluate. Better probably to look in .csv tables, or
499 run a separate emulation for that change only.

500

501 Always check for 2 or 3 kinds of your meals whether the “new” parameter settings really are on average
502 better. (See negative example in [case study 8.2!](#))

503

504

505 Part of both above shown charts (left side of each, with blue peaks) was cut out.....

506

507 (Unfinished / to be explained later) (...note: yourChanges = 1.2_bgAccel_2.0_bgBrake_1.2_dura)

508

509

510

511

