

1. Pre-requisites for Full Closed Loop

V 3.2

Please note that with autoISF you are in an early-dev. environment, where the user interface is **not optimized for safety** of users who stray away from intended ways to use. Good safety features exist, but these are only as good as the development-oriented user understands and implements them. This is not a medical product, refer to disclaimer in [section 0](#)



- 1.1 Well tuned hybrid closed loop
- 1.2 Fast insulin
- 1.3 Reliable insulin delivery from pump and cannula
- 1.4 Excellent CGM
- 1.5 Meal-related limitations?
- 1.6 Lifestyle-related limitations?
- 1.7 Time required for setting-up

Available related case studies:

- Case study 1.1: Occlusion
- Case study 1.2: Comparing insulins for FCL
- Case study 1.3: Jumpy CGM
- Case study 1.4: Lost pump connection
- Case study 1.5: Overlapping 2 x G6
- Case study 1.6: Libre3 (1 minute) placeholder

1.1 Well-tuned hybrid closed loop (HCL)

It is advisable to first establish a well-tuned **hybrid closed loop** before considering the transition to FCL. Best if you achieved good HCL performance **without** using Autotune or dynamicISF (which can introduce, or cover up, problems that would get exposed in your transition to Full Closed Loop (FCL); more see at beginning of [section 4](#)).

There are two important reasons for starting on a solid basis (profile):

- The UAM full closed loop requires a highly personalized (individual) tuning of settings, so the loop will give insulin mimicking YOUR successful hybrid closed loop mode.
- The UAM full closed loop comes with new parameters to be set and tuned. It would be problematic to set and tune several new parameters before the basics were tuned "right". Errors could easily be balanced with counter-errors. This can work in single scenarios, but would create a highly unstable system. It would be very hard to re-calibrate better later, too.

1.2 Fast insulin (Lyumjev, Fiasp, Apidra?)

It should be clear without saying, but it is absolutely necessary to feed your loop with correct time-to-peak and DIA for your insulin, so it has any chance to know how the iob will get active in your body. See Insulin_DIA...pdf' in: <https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings>

If the user does not bolus for meals, clearly a very fast insulin is needed so, upon realization of a starting meal-related glucose rise, the loop has any chance to eventually keep glucose in range (by common definition, under 180 mg/dl (10 mmol/l))

40 A modelling study..

41 key findings are summarized in initial section of [case study 1.2](#); for more see:

42 [https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings/blob/HCL-.-settings-main-](https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings/blob/HCL-.-settings-main-repo-(pdf)/The%20Artificial%20Pancreas%20and%20Meal%20Control.pdf)
43 [repo-\(pdf\)/The%20Artificial%20Pancreas%20and%20Meal%20Control.pdf](https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings/blob/HCL-.-settings-main-repo-(pdf)/The%20Artificial%20Pancreas%20and%20Meal%20Control.pdf)

44 ...can show in quantitative terms, that **faster insulins**:

- 45 • will result in significantly **lower** glucose **peaks** than slower insulins
- 46 • **tolerate** a couple of minutes **delayed** first meal bolus while not incurring unacceptable
47 height of peaks
- 48 • minimize the effect on glucose peak from **different** carb load (**meal sizes**).

49 In conclusion, do not attempt FCL with other insulin than Lyumjev® or Fiasp®.

50 Potential exceptions:

- 51 ○ Being consistently on low carb, on a GLP-1 drug, or with gastroparesis -all of which ease the
52 job of a FCL significantly.
- 53 ○ According to [case study 1.2](#), Apidra® might work, too, but Humalog® would not work well
54 (with “normal diet”).

55

56 1.3 Reliable insulin delivery from the used pump/cannula/insulin system

57

58 **Good tolerance of Lyumjev (or Fiasp): Occlusions threaten the function of the full closed loop.**

59

60 It is very important to have an eye on the time a **cannula (or pod)** is in use (many find **48 hrs** to be
61 the **limit**), and whether hard-to-explain glucose rises happen at ever increasing „fake“ iob (even
62 before a 48 hr routine replacement). (See [case study 1.1](#): You easily lose 25% TIR that day)

63 It is absolutely contra-indicated to attempt FCL coming from leaking pods and associated erratic
64 sensitivity swings that may or may not have been somewhat controlled and tolerable by
65 dynamicISF or other measures when you were Hybrid Closed Looping,

66

67 Stable pump connection


68 In FCL you absolutely rely on your pump delivering,
69 without any further delay, the much needed insulin,
70 after any meal start.

71 Hence it is absolutely essential to avoid any
72 problems from a lost Bluetooth connection. In AAPS
73 preferences / Local alerts, switch alert on!

74 An Automation similar to the one pictured here →
75 could also help recognizing eventual problems.

76

See also [case study 1.4](#)

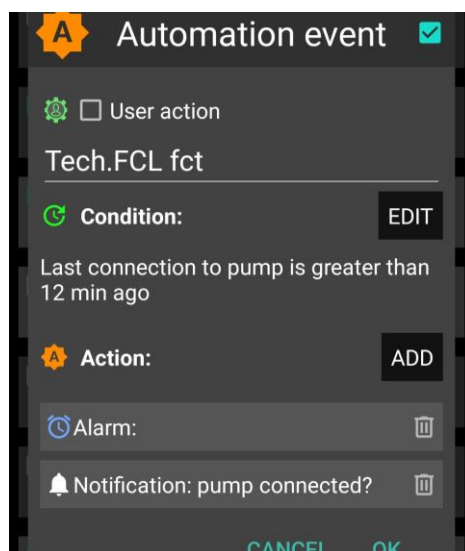

77 1.4 Excellent CGM

78

79 You do not give a meal size-related bolus any longer. That leaves all insulination jobs to the
80 algorithm! Around meals, a **stable Bluetooth** connectivity is absolutely essential, too, so CGM,
81 loop, and pump can do their job without losing more valuable time (see [case study 1.4](#)).

82

83 As glucose values are the very basis for your autoISF loop, please **inform yourself well about**
84 **your CGM:**

85

- 86 • How it principally performs (e.g. you absolutely must be “SMB always-ready” at cob=0)
- 87 • Whether you are using the best suited intermediate app that reports the “raw” value from
88 the CGM transmitter into AAPS
- 89 • Specifically, how and where any smoothing is done, and what this might imply for the ISF
90 boosting method you will be using See for instance here:
91 <https://androidaps.readthedocs.io/en/latest/Usage/Smoothing-Blood-Glucose-Data.html>
- 92 • Go through your data (in *all* day and also night *times*) to see whether your CGM produces
93 any artefacts (jumpy values; see [case study 1.3](#)) that the loop could **misinterpret** as sign of
94 a starting meal.
95 *For some* of these problems, e.g. “jumps” associated with nighttime compression lows,
96 there are options to mitigate (see [section 5.1.2./3.\)](#) . See also the User Action Automation
97 discussed about 2 pages below ([line 149ff](#)).
- 98 • In case your CGM requires calibrations: Note that calibrations often produce jumps. In that
99 case, be prepared to do, in the future, an extra handling step to protect from your FCL
100 reacting harshly.

101

102 autoISF has also a couple of in-built checks on the quality of the recent CGM values. Hence, a
103 CGM with more scatter will make the loop lose more time, and lead to higher peaks and
104 lower %TIR.

105

106 So, **if you are unhappy with a slow reaction of your loop it could be because the loop is**
107 **unhappy with your CGM.**

108 Consult the detail info given (at the time) in your SMB tab, or look it up later in the logfiles (using
109 the Emulator, [section 10](#), eventually).

110

111

112 1.4.1 Dexcom G6 and other 5-minute CGMs

113

114 The best proven way to stay out of trouble currently is to use Dexcom G5 or **G6**, and to ensure via
115 **overlapping** right and left arm sensor and transmitter utilization always good quality values, that
116 can be used by the Full Closed Loop ([case study 1.5](#)).

117

118 Other ways (using values from just *one* G6, or Dexcom ONE, G7, Libre2, or other new AAPS
119 integrated methods) are possible, but come with a lot of monitoring effort (best via watch), and
120 occasional time-outs for your FCL.

121

122 One safety feature in autoISF is a **blockage of SMB delivery whenever delta bg** (within the last
123 two 5 minute values) is **higher than 30% of that bg** (or higher +20%, at bg targets above 100
124 mg/dl).

125 *Example: From 74 mg/dl, a jump to 97 (+23 mg&dl = + 31% of 74) or more would not*
126 *receive SMB “response”. From 100 mg/dl to 131 mg/dl (+31) would neither.*

127

128 Check in your (HCL or FCL) data whether at meals or sweet drinks with rapid absorbing carbs you
129 could run into the problem that jumps are “too high” and much needed insulin will be blocked (only
130 come via very much smaller portions.

131 *For example, 400%TBR @ 0.6 U/h => 0.2 U in 5 minutes, instead of one ~3 U SMB. The*
132 *difference (of 2.8 U missed) translates @ ISF~ 40 mg/dl/U into up to + 112 mg/dl higher bg*
133 *peak! It will not become quite that bad, because the loop will catch up to the*
134 *insulinRequired with its next couple of decisions.*

135

136 Instead searching in old data, you can also just have an eye on instances where you think a first
137 SMB was due, *but blocked*. Confirm that (by looking in the SMB tab) and think about a solution that
138 would not require changing the 30% safety limit in the code.

139 *For instance, not drinking so much juice rapidly around meal start could be a “behavioral”*
140 *correction to get rid of the problem.*

141

142 This *blockage* (no SMBs) would likely last only 5 minutes (and go probably unnoticed). However,
143 not only would you lose 5 valueable minutes to get your iob substantially elevated, but also, all
144 following deltas are likely much smaller. As a consequence, if the >30% delta was in fact (largely)
145 due to carb absorption, you would, just when needed most, miss some of the boost sought from
146 bgAccel_ISF.

147

148 This example also underscores that the CGM in use cannot be allowed random scatter that leaves
149 no reasonable room for safe detection of (smaller and) bigger “truly carb related” deltas

If or when (like: first half day of a new sensor) you are not sure about sufficient CGM performance you might develop for yourself an Automation with User action ticked (along the lines as used for other purposes in [section 5.2.2.3](#)). It would “ask you” before giving a SMB whether you really want it delivered. That way you could a) have a look on your glucose curve b) .. and on the ai % (underneath Austosens%), which indicates the relative aggressiveness of ISF modulation from “what autoISF sees in your bg curve” c) think about what sense a SMB now makes with respect to your past meal, and the carbs to be still absorbed. Ultimately, you could also d) consult some of the detailed info given (every 5 minutes) in your SMB tab.

Such User action Automations need not be active at all times, but if you have it for your first half day of a new G6 sensor for instance, you could activate that Automation from your list of Automations; after the values have settled in, you can de-activate (“shelve”) it again.

For a brief period, and if you are tech savvy, another way to deal with uncertainty about CGM would be to employ the emulator method as presented in [section 11](#): Run a “too mildly” tuned FCL, and in parallel run a “what-if” with your more aggressive settings that you really would like to use once you are certain about your CGM.

However, I found it **easiest**, to **lay a solid groundwork by using** 1 Anubis, and **two overlapping G6**, to get rid of most problems (...that I keep seeing in my data, on the worse sensor of the two that often run for some days in parallel; see [case study 1.5](#)).

With a sensible **ioBTH** defined (see [section 2.4](#)), and your standard **alarms** for going towards a hypo not silenced, the worst consequence from any automatically “over-treated” glucose jump should be that you **need an unplanned snack** for the balance of “missing” carbs.

Like you should be used to from anti-hypo snacks also in your Hybrid Closed Looping past, you might also in this case here prevent SMBs for a short while by setting a (here: odd numbered) temporary glucose target (TT).

A disturbing and late sign of dealing with a too unreliable CGM (or too aggressive settings) could be a trend towards increasing TDD and body weight during your early FCL experience!

1.4.2 Libre 3 CGM with 1 minute values

Libre 3 showed promising performance too, and using it with **1-minute** values via Jugglucoo is integrated, starting with autoISF version 3.0.1.

188

189 autoISF automatically detects whether values come through Jugglucoo, and adjusts parabola fit
190 calculations to determine acceleration etc. - It is currently too early to tell, but the 1 minute values
191 applied to the parabolic bg curve fit *could* present an avenue to earlier/better acceleration detection
192 than the 5-minute CGMs enable.

193

194 When running on Libre3 / 1 minute, you get many more SMBs that each are, on average, much
195 smaller. This has **implications on** some related **settings** (notably, the smb_range_extention, see
196 [section 2.1](#)).

197

198 Tests done prior to introducing the 1 minute Libre3 option (from autoISF 3.0.1 onwards) showed
199 overall comparable but smoother insulin delivery.

200 *Call for a [case study 1.6](#) from a Libre3 user!.*

201

202 1.5 Meal-related limitations?

203

204 Setting up a full closed loop is relatively easy for people whose diet does not consist **mainly** of
205 components with rapid high effect on blood glucose (more see
206 <https://androidaps.readthedocs.io/en/latest/Usage/FullClosedLoop.html#meal-related-limitations>)

207

208 Meals do not have to be low on carb (provided you use a fast insulin for your FCL)
209 Fat or protein rich diets, or slow digestion/gastroparesis, make things easier rather than harder for
210 the full closed loop because late carbs nicely cover for inevitable “tails” of late action from SMBs
211 needed before or around peak time.

212

213 **Erratic consumption of smaller snacks with fast resorbing carbs can be a problem.**

214 In autoISF you can reduce this problem *to some extent* via one or two keystrokes from your
215 AAPS home screen. While certainly being a deviation from the FCL idea(I), this would be
216 one of the exceptional situations where you better do a quick “nudging” step from your “FCL
217 cockpit”. Details see in [section 5.2.1](#) and [case study 5.2](#)

218

219 **Really, there are no meal limitations.**

220 **The sketched more problematic options just force you to decide which extra “nudging”**
221 **efforts and/or worsened %TIR** (and/or “behavioral adjustments”, like less snacking) **you are**
222 **willing to accept**, occasionally.

223

224 225 1.6 Lifestyle-related limitations? 226

227 Providing a technically stable system 228

229 Full closed looping requires a 24/7 technically stable system, especially regarding

- 230 • reliable CGM signals
- 231 • Bluetooth stability with the pump (see [case study 1.4](#))
- 232 • keeping your phone in sufficient proximity at all times
- 233 • avoiding (or at least early recognition of) occlusion.

234 This requires a habit (or, unlikely, permanent attention to details) like keeping all components well
235 charged and in close proximity; making cannula (or pod) changes always early enough to lower the
236 risk of occlusion (see [case study 1.1](#)); having always potentially needed parts with you.

237 Depending on your system, your experience with it, but also on your acceptance and general
238 lifestyle, these aspects may or may not limit you.

239

240 Preparing for exercise 241

242 To prepare for exercise (sports, heavy work), the normal protocol with a pump or hybrid closed loop
243 is to take actions that reduce insulin on board prior to exercise

244 With your full closed loop, the algorithm is tuned to detect meals and to give you insulin to counter
245 glucose rises automatically. Setting a high temp. target and lower %profile right away (effective al-
246 ready around meal start) could be a problem.

247 Unusual activity levels therefore likely require **disciplined preparation** (especially **if you want to**
248 **keep the need to snack during sports low**)

249 In autoISF you can reduce this problem to some extent via two or three keystrokes on your
250 AAPS home screen. While certainly being a deviation from the FCL idea(I), this would be
251 one of the exceptional situations where you better “flick a lever” from your “FCL cockpit” to
252 keep iob low (example see [case study 6.2](#)).

253

254 Extra hurdles to establish FCL for kids 255

256 To establish and maintain a FCL for kids brings about some extra challenges if:

- 257 • Lyumjev is not available or well tolerated
- 258 • Hourly basal rate is very low, providing a poor basis for big SMBs

- Diet is rich in sweet components. With the typical low blood volume of a small body, strong tendency towards very high bg spikes!
- Going through marked changes of insulin sensitivity or of circadian pattern makes it difficult to keep the FCL appropriately tuned.

This problem is about the same in Hybrid Closed Looping. However, now you might expect miracles from the FCL. This is not going to happen. You still must be pro-active by setting suitable temporary % profiles. (These serve as a basis, also for your autoISF FCL).

- Discipline may be poor regarding keeping Bluetooth connectivity and infusion sites perfectly running
- Between kid and supervising parent it must be guaranteed, especially in the initial weeks, that an eye is kept on whether the FCL is working about as expected.

More see [section 7](#).

1.7 Time required for setting-up

Lastly, before enjoying a functioning full closed loop you need to have a period of a **some weeks** with some free time and „free head“ for set-up –. Can you get, in the time you are willing to invest, to a result that you consider good-enough is really the question. Depending on your „habits“, and which – if any - compromises (like doing cannula/pod changes more often, never starting meals when bg sits high and is not predicted to rush down soon ...) are you willing to make (and everyday able to stick to), for the ease of not having to deal with assessing meals and bolusing for them?

1.7.1 Recommended structured, step-by-step approach (see also beginning of [section 4](#))

Setting up your personal FCL using autoISF is a substantial project, for which you should **follow the sequence as described** in this e-book.

But there is **no need to implement it fully in one step**.

- There is nothing wrong to go in your well running Hybrid Closed Loop mostly, while switching to FCL only for dinners, for instance, or only for weekend lunches, as a start.
- Once you found feasible settings, you can expand to other meal times...
- and lastly towards figuring out your best strategies for challenges aside from meal management, as we shall discuss in [sections 5](#). and [6](#).

296 There are alternatives to using autoISF for FCL, as well. See [sections 7](#) and [13](#). for more info.

297

298 Notably FCL using AAPS Master and Automations (see in [section 13.1](#)).could be a much
299 easier and more error-tolerant way of stepping into FCL. In a clinical study with 16
300 participants about 80% TIR was achieved without much tuning effort

301

302 To close the circle to where we had started; A very time consuming pre-requisite might actually be
303 to **first sort out your Hybrid Closed Loop** ([section 1.1](#)), so your profile parameters are set „right“,
304 and your “old” data really can serve as a blueprint for what, now, you would like *your loop* to do in
305 FCL mode

306 If you feel you better do your homework first there, I highly recommend some of the material in the
307 neighboring HCL repo: <https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings>

308

309 1.7.2 “Trial and error” fast track alternative

310

311 Note that if you had used dynamic parameters or special Automations („loops inside the loop“) this
312 might have balanced some principal errors, but leaves you now *without a good starting point*, as
313 you must get rid of these over-patches (see also warnings at start of [section 4](#))..

314

315 Nevertheless, you will find FCL success stories also from loopers who continue(d) in that spirit, and
316 just jump(ed) into using more powerful tools, in kind of a trial and error mode, frequently adding the
317 latest add-on, or self-constructed patch (often in form of an Automation), to counter-balance
318 encountered problems.

319

320 Resulting solutions may be good-enough.

321 But they tend to be unstable and not well-understood. That is a poor basis for managing arising
322 problems, and for temporarily adjusting to special situations.

323

324 Nevertheless, it is an alternative avenue for the impatient, less analytically, and more adventurous
325 inclined.

326 Note though, that it is hard to consult (help) someone who, over time, constructed his/her own
327 complicated maze of constructs.

328 Still, I welcome everyone's creativity. “Not macht erfinderisch” (necessity is the mother of
329 invention) we say in German: Something amazing of broader interest could emerge, and
330 help push innovation for us all.

331

332 However, be prepared to eventually needing **a complete fresh re-start, if your trial-and-error**
333 **got you lost**.. Depending on your knowledge level and experience, this easy can happen on a
334 “fast track” route.

335

336 1.7.3 Safety first

337

338 Regardless which route you choose, PLEASE always observe the **safety** settings/instructions
339 coming with the DIY dev- variant of FCL software you select.

340

341 All FCL methods come with boosted SMBs. So a key safety measure every user going towards
342 any FCL should have in place is to set an **iob threshold** (iobTH; size a bit below what you used as
343 a bolus for bigger meals in HCL) above which no more SMBs can be given by your FCL.

344

- 345 • iobTH is an integrated feature of *autoISF* (see [section 2.4](#)).
- 346 • *Other AAPS-based FCL methods* may require to set up an Automation for a temporary iob
347 threshold that blocks SMBs from being delivered, see e.g. here for AAPS FCL
348 w/Automations; :<https://androidaps.readthedocs.io/de/latest/Usage/FullClosedLoop.html#iob-threshold> .
- 349 • *In case there are other methods than autoISF* for FCL also on the *iAPS* or *Open iAPS* platforms, you
350 may have to rely on an adjusted iob_max border, or watch the iob development, and intervene with a
351 SMB shut-off, or by opening the loop, when deemed necessary.

352

353
354 Also, make use of the easy-to-use feature of **SMB shut-off at odd** profile or temporary **target**. This
355 can **at any time easily** be done manually, via the top right “TT” field in your AAPS screen (set, and
356 time, an odd-numbered target; [section 5.1.3](#)), and can be of enormous help **to temporarily**
357 **safeguard you from aggressive loop actions** (i.e. further growth of iob, no matter how close you
358 already may be to the iobTH)-

359 This is the same concept that you already know from your HCL times, when you wanted to “tame your loop”
360 so it does not “fight” your anti-Hypo Snack with a SMB (An elevated TT> 100 mg/dl was then used, to shut
361 off SMBs in HCL for a while).

362

363 Lastly, make sure you “train”, for your set-up weeks, how/when to **transition between FCL and**
364 **your prior HCL** mode (refer to [section 5.1.1](#) and [5.2.3](#)).

365 A new User Interface has been suggested to ease this transition via a modified loop button in the AAPS main
366 screen (developers: see [section 5.3.1](#)).