

5. Modulation of autoISF aggressiveness. V 2.2

Please note that with autoISF 3.0 you are in an early-dev. environment, where the user interface is **not optimized for safety** of users who stray away from intended ways to use. Good safety features exist, but these are only as good as the development-oriented user understands and implements them. This is not a medical product, refer to disclaimer in [section 0](#)



Once the initial tuning according to [section 4](#) is done, you are ready to use autoISF for your automated meal management.

You will have three major other challenges to manage:

- recognize and manage (partial) occlusions, or other technical (CGM or BT related) obstacles (see [section 2](#) on pre-requisites of FCL, and related case studies)
- deal with times when the loop should be set "milder" as a precaution (e.g. nights; or in an exercise context)
- deal with times when insulin given by the loop must be restricted (e.g. a snack could be "misinterpreted" as a meal)

How big the remaining challenge is depends very much on your individual lifestyle. [Sections 5 and 6](#) discuss this in more detail.

In order to run the loop fully automatically around the clock, the times outside the meal blocks must also be precisely analyzed, and solutions to problems (if any) must be sought.

It is up to every user to decide where to draw the line.

- With a technically well functioning system, moderate meals, moderate or no exercise, moderate %TIR expectations and a bit of mindfulness it should be possible to go into Full Closed Loop 24/7, after working through, and observing, [sections 1-4](#).
- Especially if you are a bit shy of using the emulator for really detailed analysis, it is likely that you will not hit a real good system calibration ([section 4](#)) for an entire range of diets.
- In that case you will occasionally run out of range, and your options to prevent, react, or improve are
 - accepting a few % higher time outside range for that day (and avoiding what seemed to have caused it in the future)

- taking a snack (whenever you tend to go low from the “tails” of insulin activity that was required to fight a peak)
- doing a manual override (if you can think of one in time, to manage the problem manually)
- temporarily resorting to the well-known hybrid closed loop.
- Instead of accepting such instances, you could launch “improvement projects”
 - that refine your initial tuning ([section 4.](#))
 - that make you and your FCL loop fit to manage an increasing number of disturbances either automatically, or via a user intervention ([sections 5-6](#)).

To tailor the loop’s response to disturbances may require specific **modulation of the aggressiveness** as set for your meal management.

There are many avenues to achieve this. The main ones, that are also easy accessible via Automations in AAPS, are:

- temp. shut-off SMBs (odd-numbered target)
- temp. reduce bgAccel_ISF-weigh
- temp. reduce iobTH
- temp. reduce set %profile
- temp. set higher TT (especially in connection with exercise mode)t

5.1 Permanent fully automatic modulation of FCL aggressiveness

In setting up your FCL, you now have another difficult and time-consuming job at hand, to define fully automated solutions (that require **no user intervention** at all) for any of your „other“ situations, outside of meal management.

In [section 5.2](#) and [5.3](#) we will look at options *to ease that job by “allowing” a 1-button push or data entry* intervention, like for a snack or exercise announcement.

5.1.1 autoISF generally switched off outside of meal-time windows

If, aside from meal management, you were rather happy in hybrid closed loop, you could continue to run in that mode, and just focus your new autoISF FCL on management of meals (on all meals,

70 or only on a sub-set of them, like only dinners – which might make sense especially in your initial
71 transitioning phase).

72

73 For this, you define Automations

74

- 75 • that set **meal time windows** in which autoISF gets fully turned on
- 76 • or: that turn *all* autoISF's ISF modulations (*or just bgAccel_ISF*) off in time windows in
77 which surely no meal occurs. For instance, you can go for all nights back into your Hybrid
78 Closed Loop, as you had before.

79

80 Other early DEV AAPS variants (see [section 13.3](#)) all work with meal-time windows. The
81 window is either set by time of day in the settings, or it always must be „set“ by the user via
82 giving a mandatory small pre-bolus before any meal starts. **Outside** of these time windows,
83 these loops then runs with less aggressive SMBs like oref(1) SMB+UAM in AAPS Master.
84 This mode is not really FCL, but an advance over traditional HCL that often achieves satisfying
85 degrees of automation and performance.

86 The term **Meal Announcement** (MA) is often used to label this closed looping mode. Trigger to
87 set a meal time window could also be a pre-bolus given by the user, a carb entry made, an
88 EatingSoonTT set or a meal announcement button pushed.

89

90 Note: Outside of the meal time windows you would be in hybrid closed loop. To the extent you
91 rarely face disturbances (aside from meals), you could be looping in full automatic mode around
92 the clock. Your temp. autoISF shut-down is only meant to prevent problems from the loop over-
93 reacting to bumps in the glucose curve in times of day (night) when standard oref(1) performance
94 is sufficient.

95

96 5.1.2 Odd-numbered profile targets used to block SMBs

97

98 Ensure the even/odd logic in the settings is toggled on in Preferences> openAPS SMB>
99 autoISF settings> smb delivery settings>: "Enable alternative activation of SMB depending
100 on profile target" ON.

101

102 In time blocks with an odd-numbered profile target you can prevent any SMBs being given by your
103 loop. The (unchanged) aggressive settings can only translate within the limits set by %TBR
104 possible.

105



106 This will very much slow down any more insulin being given, and is an excellent solution for night
107 times, especially if you occasionally experience compression lows.

108
109 Notably it is an excellent solution to prevent getting too much insulin because of jumpy CGM
110 values, like after a compression low. This is therefore a good solution for night time.

111
112 Alternatively, you could use the new included options for Automation Conditions and
113 temporarily tune your bgAccel_ISF_weight much lower.

114 Yet another alternative would be to go into hybrid closed loop for the night, with or without
115 SMBs ([section 5.1.1](#)).

116
117 In case you occasionally do have nights that would benefit from a couple of SMBs (to treat temp.
118 highs from a late fatty pizza, raclette and such): Define suitable Automations like the two „night“
119 ones here:

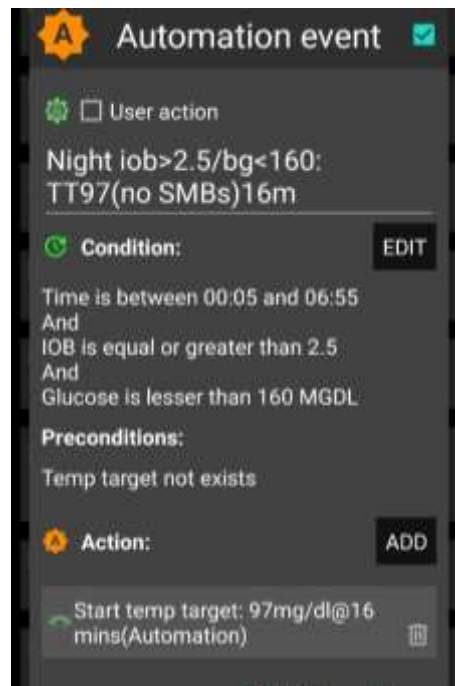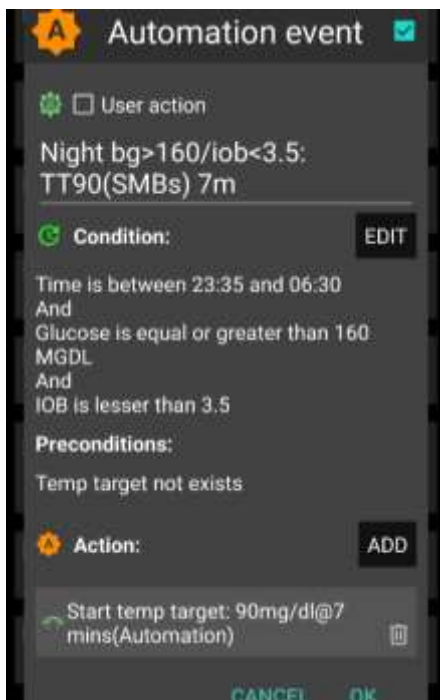

120
121
122 **Never underestimate the „trickyness“ of getting your Automations „right“.**
123 With your thought-out Automations in place, night data need to be analyzed to see
124 • whether the bg and iob limits defined in the given example work sensibly
125 • whether the TT duration is chosen appropriately
126 • how swapping the sequence in which the automations appear in the Automation list would
127 lead to different SMB impacts.

128
129
130

131 5.1.3 Odd-numbered temp. targets used to block SMBs

132

133 A widely used ACTION that strongly modifies how fast your FCL can add more iob is setting an
134 **odd-numbered temp. glucose target** which makes the loop operate without giving any SMBs
135 (%TBR modulation only).

136 Ensure the even/odd logic in the settings is toggled on in Preferences> openAPS SMB>
137 autoISF settings> smb delivery settings>: "Enable alternative activation of SMB depending
138 on TempTarget" ON.

139

140 So, from patterns you find in YOUR data, at times where you want your loop act differently, you
141 need to carve out CONDITIONS that describe the respective situations (and either for how long it
142 typically lasts, or at which *other* CONDITIONS you want your loop get back to default FCL
143 operation).

144

145 From, autoISF 3.0 onwards, also the following parameters are provided as CONDITION and/or as
146 ACTION for defining YOUR Automations:

- 147 • Enable ISF weights / Disable ISF weights => Allows temp. ON/OFF for the key ISF
148 modulation parts of autoISF
- 149 • Trigger/set iobTH percent => Keeps default aggressiveness, but only until a modified iob
150 threshold is surpassed
- 151 • Trigger/set bgAccel_ISF_weight => Modifies the default aggressiveness of just the
152 acceleration component

153 An odd TT is often set for an *anti-hypo* snack or *sports* snack. In both instances, you do not want
154 SMBs to quickly counter act.

155

156 In case of *sweet “fun”* snacks, this is entirely different -> next [section, 5.1.4](#)

157

158 5.1.4 Automatic differentiation of FCL aggressiveness using Automations

159

160 **Personalized Automations** tailor the loop exactly to **your** data so **fully automated handling** of
161 situations with different aggressiveness of the loop can be made.

162

163 So, you first must analyze patterns you find in YOUR data, at times where you want your loop act
164 differently, to carve out CONDITIONS that describe the respective situations (and either for how
165 long it typically lasts, or at which *other* CONDITIONS you want your loop get back to default FCL
166 operation).

167

168 A variant of this mode is to define **several meal time windows** in which autoISF aggressiveness
169 (bgAccel_ISF_weight) and/or iobTH are **set differently** for different meal time slots of your day (or
170 even for a geo-location etc).

171 Unless your meals differ vastly in size and in fast carb content this may not be needed, but it can
172 make initial tuning (setting the various ISF_weights, iobTH and size limits) a lot easier

173

174 5.1.5 Automatic adjustment of FCL aggressiveness via the Activity Monitor

175

176 If you choose to make use of your smartphone's **stepcounter**, you can (automatically) adjust
177 insulin sensitivity ratio to activity level in the past minutes to one hour time frame.

178 This is another little tuning opportunity, in which you study your body's response to light exercise
179 (like walking) or to not moving at all (like desk, couch), and select appropriate settings which, in the
180 future, will automatically adjust insulin delivery to suit activity state of the past minutes (up to 1
181 hour).(AAPS Preferences/OpenAPS SMB/Activity modifies sensitivity/ -> set **two scaling factors**.)

182

183 This autoISF feature (new since V.3.0) is much quicker responding than Autosens or dynamicISF
184 to adjust insulin sensitivity to your current „lifestyle state“.

185

186 More see [sections 3.5](#) and [6.5](#)

187

188 5.1.6 Pro/con completely hands-off Full Closed Loop

189

190 Remaining 24/7 in a completely „hands-off“ FCL can be a realistic goal with autoISF 3.0 if besides
191 meals also some special challenges as discussed in [section 5.1.](#) were analyzed and addressed.

192

193 Clearly it depends very much on your lifestyle, and how interested, willing, and capable you are to
194 recognize, deal with, (and in the future avoid?) situations that get you outside of your desired %TIR
195 on occasion.

196 So, this is also about what %TIR you are aiming at, and can accept, as it averages out for
197 the week, for instance.

198

199 Everybody must weigh for her/himself how much **upfront effort** to put into getting it all 100%
200 automatic, **or** whether to take an **easier start, with a couple of situations left to take care of**
201 **when and as they arise in daily life**

202

203 Even if a principal capability for a fully automatic running FCL is given, this still
204 means that

205 • the user should be knowledgeable about what exactly is going on, and

206 • has a principal capability to „nudge“, or to take over in a manual mode.

207
208 In the sections that immediately follow, we present the options to nudge or temporarily take over
209 from the AAPS home screen serving as your **FCL cockpit**:

210 • [Section 5.2](#) describes how you can build your own DIY cockpit

211 • [Section 5.3](#) describes how a FCL cockpit might look like in the future

212 213 5.2 Modulating aggressiveness manually, from the DIY-FCL-Cockpit* 214

215 * Like in the airplane cockpit: Cruising in full auto mode should involve having an eye on the
216 instruments, and on potential disturbances ahead in the environment.

217
218 Life in Full Closed Loop is easy if your lifestyle largely consists of real meals, and not much other
219 disturbances like from snacking in between. In [section 4](#). we dealt with major meals.

220
221 In [section 5.1](#) we looked into ways to automate also a modified loop response to foreseeable situa-
222 tions, or to those the loop could recognize (with enough time to react).

223
224 Other „disturbances“ might come up, and you must find an easy way to

225 • call up a pre-programmed routine for automatic management, with auto-adjusted
226 aggressiveness, or:

227 • tweak a setting or two, to temporarily adjust the aggressiveness

228 • There may also arise a desire to just exit the FCL mode, and be your own captain for
229 mastering a special situation.

230 For peace of mind, to learn, and to stay informed (especially so in your initial tuning phase, or
231 when your glucose curve goes in unexpected ways) we also must be able to

232 • find the key parameters that frame and drive the recent and upcoming loop decisions.

233
234 All this is facilitated within seconds right from the AAPS home screen, serving as a **FCL cockpit**
235 after you built a couple of DIY cockpit features via Automations (as described below and in [case](#)
236 [studies 5.2](#) and [6.2](#)):

237 Thoughts went also into [how to improve the cockpit in future releases](#), see [section 5.3](#)

238

239

240

241 5.2.1 Triggered Automations: Grey extra DIY cockpit buttons for pre-programmed “responses”

242

243 Recognizing conditions for fully automatic handling by the loop may not be not possible, or come
244 too late for the loop to act on. Examples would be exercise, where minimum an hour before start-
245 ing “the loop should know”.

246

247 Another example: snacks: High carb snacks, sweets, consuming ice cream or having a sweet drink
248 comes with the problem of even steeper glucose rises, but overall a lesser insulin need, compared
249 to major meals (for which we tuned our FCL according to [section 4](#)).

250

251 This not necessarily implies that snacks need different settings than a meal. After all, autoISF was
252 designed to act to all available data, especially to where the developing glucose curve is headed.
253 So, depending on your effort to set parameters for a broad variety of meals (notably: how well you
254 avoid to invariably bounce fast against your iobTH), you might be able to accommodate low carb,
255 snack, and major meals with *one* set of settings.

256 In FCL autoISF, this is a bit more difficult than in HCL autoISF applications, because FCL
257 involves revving up iob supply (largely via big bgAccel_ISF-weights) often too much to be
258 balanced by just a snack getting absorbed.

259 For that reason, or just for increased comfort and safety, you might want to differentiate, and make
260 use of what follows.

261

262 Tuning aggressiveness

263 Key is that a sweet snack likely benefits from even more aggressive initial FCL
264 performance than the meals in your normal spectrum of diets require.

265 Therefore, you could set

- 266
- a higher **temp. profile%** and/or
 - a temp.elevated **bgAccel_ISF-weight** (see screenshot of my Automation).
 - a **low temp. target** (76 for instance) additionally helps maximize the first SMBs that
269 will automatically be triggered at detection of acceleration.

270

271 When first defining and testing this Automation, also check:

- 272
- that the safety limits as discussed in [section 2](#) will not block the intended elevated
273 aggressiveness

- 274
- SMBs will not get outrageously big and iobTH sometimes exceeded by too much


275 Note that “the last SMB” is allowed to overshoot the valid iobTH by 30%

276 Limiting iob

277 For “just a snack”, total insulin need will be lower than for a meal.

278 If you would just have your sweet drink, and your meal-oriented FCL would “attack”,
279 iob likely would become too high, and a glucose rollercoaster would start, with you
280 needing to consume more =>

281 If you just have a snack, or drink a glass of juice, you can lower the **iobTH_percent**
282 accordingly.

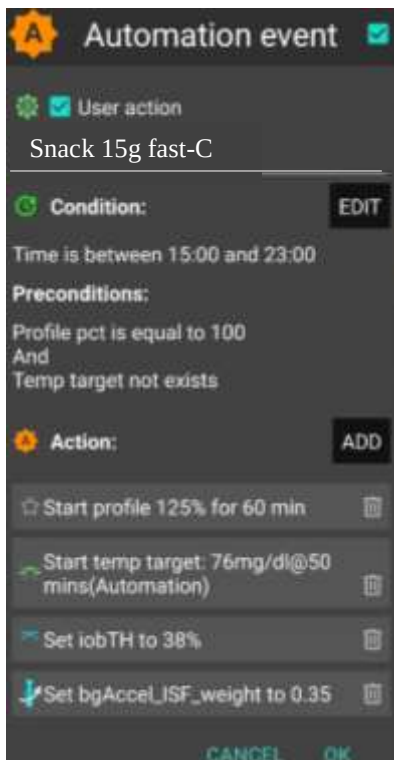

„User action“ is always ticked-on

This will, in the defined time space ..

..offer the “DIY cockpit” button..

..which I must press any time (~90...30 minutes) before my snack.

Note: Make sure that iobTH and bgAccel_ISF_weight revert to default afterwards



283

284

285 So, this is a little **extra “project” when setting up your FCL**. You need to research your
286 snack habits (if any), and over time find out which settings in the snack-related Automation
287 work well.

288

289 **In everyday life you then just must press the related button in your cockpit** (which is
290 not time critical at all, except it should be clicked *latest* a couple of minutes after you took
291 the drink or snack).

292

293 If you consume more, and also eat something with your sweet drink, this will more
294 resemble a full meal... however, with unusual amounts of fast carbs.

295

Note: Pressing your snack button *a second time* would **not** help because the lowered iobTH does not allow iob going high enough. So you are better off just letting your normal FCL meal routine run, after your snack mode expired.

Other options when snacks keep extending would require a manual modulation regarding %profile and/or bgAccel_ISF, but keeping the full default set iobTH_percent, or even elevating it (refer to [section 5.2.3](#)). If that happens often, define for yourself an extra User action Automation for a bigger snack (grey DIY cockpit button).

Installing the DIY cockpit button

In the related Automation, just keep the “User action” box clicked at all times, and define in the Conditions when you want to see that button available for cockpit use (see screenshot above) => you will see that button offered.

Besides snacks, also any other recurring special situations can be addressed via a cockpit button, and receive different aggressiveness up to a suitable iobTH level.

In the future you might be able to set the stage for a snack and other “disturbances” also via an extended menu below the TT button on the AAPS home screen, see [section 5.3.3.1](#)

Discussion

If I had *regular* snacking habits in certain parts of day, I might take an *alternative* route and **modify my FCL settings in those time slots** to run automatically upon acceleration detection.

Yet another **alternative** would be to temporarily leave the FCL mode and handle the sweet snack or drink “the traditional way” in **hybrid closed loop**.

The suggested FCL cockpit user interface with an extra version of violet loop on the AAPS home screen facilitates that, including automatic removal and re-appearance of the insulin button at the bottom of the APS home screen.

331 As mentioned in section “Limiting iob” above, it is essential though to either avoid snacks,
332 or select one of the discussed easy ways to deal with them in everyday life.

333

334 5.2.2 Status recognition

335

336 Before considering any manual interventions into the ongoing FCL, you should be aware what the
337 current mode of action is, and hence which button eventually to fine-tune or lever to switch, in or-
338 der to adjust to the disturbance you see coming up.

339

340 See [section 5.4](#)

341

342 5.2.3 Manual interventions from the (DIY-) FCL cockpit

343

344 Trouble with all these is, not to forget to set back manually, too (=> better solutions in 5.3)

345

346 5.2.3.1 Tuning aggressiveness via temp. %profile or TT settings

347

348 The set **% profile** multiplies with both, the ISF resulting from autoISF, and also with the default
349 iobTH you have set, so both are nicely modulated in a linear way with the % temporarily chosen
350

351 A lowered (relative to profile glucose target) temporary **bg target (TT)** signals lowered sensitivity
352 (more insulin need), and an elevated TT (as often used with exercise) increases sensitivity and
353 hence works in the direction of a lowered % profile to also reduce insulin given by the loop.

354

355 Moreover, the **exercise button** (top center on your AAPS home screen) can be activated (turns
356 yellow, then). This will **further boost** how your set TT elevates the resulting ISF, and sharply
357 lowers iobTH, as often desired for sports. See [section 6.1](#)).

358

359

360 5.5.3.2 Making temp. changes in settings made in /preferences

361

362 Going into AAPS/preferences/Open APS SMB allows to:

363 - Limit/expand SMB sizes

364 - set different iobTH

365 - set milder or stronger ..._ISF_weights

366 - switch OFF the even <-> odd SMB on/off

367 5.2.4 Temporarily exiting the FCL

368 One **alternative** always is to temporarily leave the FCL mode, and handle any disturbance
369 “the traditional way” in **hybrid closed loop**.

370

371 The suggested FCL cockpit user interface with an extra version of violet loop on the
372 AAPS home screen ([section 5.3.1](#)) would facilitate that, including automatic removal
373 and re-appearance of the insulin button at the bottom of the APS home screen.

374 In case this feature is not yet available, you must:

375

376 Exit the FCL mode by going to AAPS/preferences/put in your password/OpenAPS SMB/scroll down
377 to autoISF settings and switch “Enable ISF adaptation..” OFF.

378

379 Unfortunately, there is no way yet for it to come automatically back on, after a selected time for
380 instance. So do not forget to switch your autoISF fully back on, later.

381

382 As this will often be forgotten, it may be worth doing a “User action” Automation, for a “temp.
383 FCL OFF” grey button (see section 5.2.1).

384 Caution though, there is very limited experience with this brand new feature, and it might be
385 accidentally activated . – To make things worse:

386

387 To recognize whether autoISF currently runs with ISF adaptation or not, you must consult the
388 profile_sens -> actual_sens indicator below the Autosens%. However, this gets also modified
389 by %profile switches or TT +/- exercise mode. So it is not as easy as it would be with the “violet
390 loop” proposal mentioned already above.

391 Ultimately, you can of course study the SMB tab to find out what is going on.

392

393

394 5.3 Modulating aggressiveness manually from the improved FCL-cockpit

395

396

397 autoISF 3.0 is an early dev variant of AAPS, and as user you are participating in an on-going
398 development. Of note, autoISF 3.0 is first launched without many of the cockpit features that are
399 presented below in this font color. (Only what is written in black is at this point of some relevance
400 for using autoISF 3.0)

401

402 For the time being, multi-step work-arounds may become necessary

403 • In many cases, going into AAPS Preferences and changing settings would be needed
404 (...plus not forgetting to change these settings back, afterwards).

405 • Automations allow a DIY FCL cockpit, see [section 5.2](#) and [case studies 5.2](#) and [6.2](#)

406
407 This is also an open invitation for you to contact us in case you could help program a
408 module for one of the required user interface extras.

409 For future integration into AAPS Master, an eye should be kept also on the question which
410 other modes (like FCL using Automations and others mentioned in [section 13](#).) might benefit
411 from some of the extra features.

412
413 Keep in mind, though, that the **goal should be to interfere with the loop as little as possible**.

414 Under the described conditions it can run **fully automatically** without any user interaction (= after
415 the initial tuning phase, and related settings made in AAPS /preferences/SMB/autoISF. See [section](#)
416 [4](#). and [5.1](#)).

417
418 However, just like in the airplane cockpit: Cruising in full auto mode should involve having an eye
419 on the instruments, and on potential disturbances ahead in the environment.

420 *E.g.: storm ahead => instruct your plane to climb to another flight height.*

421 *Analogy: exercise ahead => setting an exercise TT, or => pressing a button that activates a*
422 *sequence of instructions (some of them probably hinging on conditions, like actual iob) how*
423 *to manage through that exercise situation).*

424

425 So, for the occasional „disturbance“ coming up, you should find an easy way to

426 • call up a pre-programmed routine for automatic management, with auto-adjusted
427 aggressiveness, or:

428 • tweak a setting or two, to temporarily adjust the aggressiveness

429 • There may also arise a desire to just exit the FCL mode, and be your own captain for
430 mastering a special situation.

431 All this is facilitated within seconds right from the AAPS home screen's **cockpit features to the**
432 **extent they are already incorporated, or** to the extent you can build alike DIY cockpit features via
433 Automations, as described in [section 4.1.3](#) and [case studies 5.2](#) and [6.2](#)):

434

435 • The button that is integrated into the **violet FCL icon** serves as emergency off button, to
436 quickly stop FCL, or to at least to immediately stop any more SMBs (...just for a couple of
437 minutes, or for the remaining meal time: pick from the options offered with just one
438 keystroke).

439 Via the violet FCL icon on your AAPS home screen, you also can access a temp. switch-off
440 button for SMBs (see section that next follows below).

441

- 442 • The **three top fields** (%profile, exercise, TT) provide access to temp. tuning of core
443 parameters, and/or to some pre-programmed routines.

444 Taken together with some **new indicator fields** about your loop state, and the **grey DIY cockpit**
445 **buttons** this makes the AAPS home screen your **cockpit** for Full Closed Looping.

446

447 Let us look on each of these cockpit elements in some detail:

448

449 5.3.1 Violet FCL icon and underlying buttons

450

451 Novices to FCL, or really anyone running into a very special situation, may appreciate that the new
452 closed loop icon on the AAPS home screen in pink (for FCL) has buttons to quickly shut off getting
453 more SMBs (1st row), or to enter other loop modes (second row).

454

455 It functions very much as the other ones that you know from HCL already, and in fact you
456 get offered some of the same options (for instance, to switch the (full) closed loop off for 15
457 minutes for going to take a shower)


458

459 Note that in FCL you leave all BG regulation, notably against meal spikes, to the loop. So, try not to
460 disconnect in phases when your FCL must ramp up your iob.

461 The required insulin would still be supplied *after* you reconnect. However, without the user
462 pre-bolussing, the delay would be more of an issue in FCL than it had been in HCL.

463

464 Just pressing on the FCL icon, a dialogue box comes up:

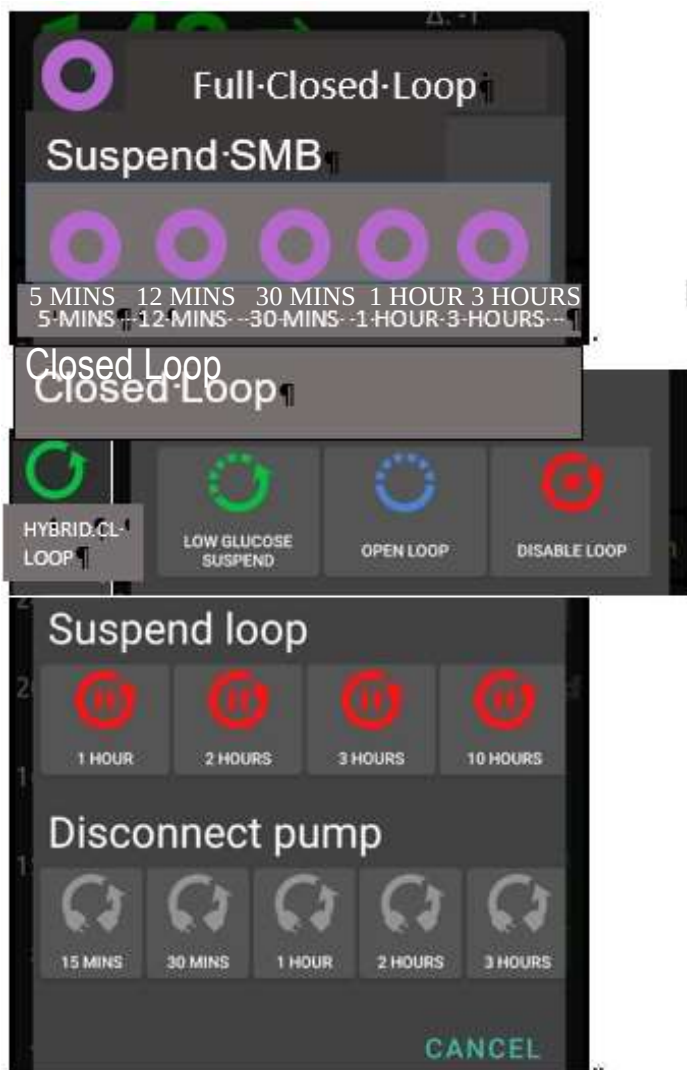

<-add-arrow, form-like-like-CL-Loop...¶

<-dotted++arrow, form-like-like-LGS-¶

DEV: when in Hybrid-Closed-Loop = top headline++green-circle), there is no Suspend SMB part but it starts right with „Closed Loop“. In that screen. There, the 1st element must be option to pick.¶

.....pink-circle, FULL-CL-LOOP-¶

where here the green HYBRID-CL-LOOP stands.¶

For loopers who did not set up FCL, a feedback must come up ~ „FCL not installed“ if they press on that button.¶

465

466 Pressing „**Suspend SMB**“ provides fast and easy „emergency braking“ regarding delivery of more
467 SMBs:

468 Select the one with the desired number of minutes: 5 or 12 for just blocking the potential next
469 SMB(s), and up to 3 hours to manage the entire rest of this meal with %TBR from then on.

470

471 Whenever, and why-ever, your FCL is in „no SMBs allowed“ mode (e.g. automatically after
472 surpassing an iobTH, or triggered by a set odd TT), the FCL icon will turn into a dotted one.

473 Instead of remaining **duration to end time** it indicates in the middle „the condition“, „**iob**“ or „**TT**“

474 Add an indication if suspend SMB comes from an Automation, e.g. add an „**(A)**“, **underneath** the
475 #minutes, iob, or TT in the middle of the dotted violet field.

476 So, as in other (already in HCL existing) cases, those icons show in the middle the minutes left that
477 they will be running, or the condition which would have to go away for this temp. setting to stop.

478 It always auto-reverts into the FCL state and FCL icon, when time (or other condition) has elapsed.

479

480 Pressing „**HYBRID CL. LOOP**“ or other buttons from the 2nd row provides fast and easy

481 „emergency exit“ into other modes.

482 This enables beginners an easy „temp. escape“ into their well-known HCL (green) at any
483 point of time. bgAccel_ISF_weight is set to zero when going FCL->HCL. HCL can run with
484 autoISF (for instance dura_ISF) uninhibited otherwise. (check implications for HCL users of
485 autoISF ??).

486 Note: These options from row 2 have no time limit. Loop will **not** by itself go back to FCL. You see
487 the different loop icon as a reminder to manually revert, when ready.

488

489

490 5.3.2 Buttons „Insulin“, „Calculator“ etc at bottom of AAPS home screen

491

492 These buttons are **not useful any longer in FCL**, and automatically disappear whenever in FCL
493 mode (also in Suspend SMB state), and re-appear when leaving FCL. This applies also when an
494 Automation or technical system failure shut off FCL.

495 Users who, maybe in the beginning phase, feel better having those buttons, can override
496 the removal (of the insulin button, or any other) by going into /preferences/overview/buttons
497 and forcing them on. They only remain on until the next re-entry into FCL mode, when auto-
498 off happens again.

499 The reason why we do this: It really is important to let the loop loop, and not interfere more
500 than absolutely needed. Any bolus the user gives will sure distort the bg curve, on which
501 autoISF, especially when aggressively tuned for FCL, builds a lot of its decisions!

502

503 5.3.3. Three top fields (%profile, exercise, TT)

504

505 Depending on the variedness of lifestyle, the desired %TIR, and the initial tuning effort put in, the
506 user may want occasionally to **„tweek“ the aggressiveness of her/his FCL**.

507

508 The top 3 fields (grey in default mode, **yellow when temp. in mode with changed**
509 **aggressiveness**) serve as quick and easy entry points to make temp. switches (as users will be
510 used to for %profile switches, or for setting an EatingSoonTT in HCL, .. which they still can do in
511 FCL ... but more:)

512

513 Expert FCL users might need this feature rarely, but probably at least to manage activity after
514 meals: Each require opposite aggressiveness, and the switch has to come in a certain point in
515 time that would be difficult to capture. (More see [section 6.4](#))

516
 517 5.3.3.1 TT dialogue field ((Currently not available in the pictured form and function !)
 518 *f extended design for FCL cockpit is already launched)*
 519 The TT field (top right of AAPS home screen) is the primary daily interface, and a dialogue field
 520 opens when pressing on it

Initially, the form can be entirely empty re. TT inputs and just show the default iobTH and bGAccel parameters.

Two lines appear for target and get labeled ES, AC or HY depending on what was pushed at the bottom buttons. (HY => red frame; evtl. empty => grey)

TT entries can be made or overridden.

iobTH calculates from Target AC and other settings shifts away from default set in /preferences (here 0.7) but could be overridden here. 0!..200 % is allowed

bGAccel ISF_weight can be modulated here, too. Note: it can change again if % profile is also changed. --0...200 % is allowed

CANCEL allows to start fresh (select one or two of the square buttons, ES, AC or HY)

OK needed to use the settings

Duration input is made in minutes. In the exceptional case that both, ES and AC targets are defined, the duration input is for AC and framed blue. (This is because the preceding AC mode is automatically determined in length by the loop observing when iobTH is exceeded)

521
 522
 523 This looks complicated but only because it allows 4 different modes of use. Each user will primarily
 524 use her/his preferred one.
 525 (1) Who is happy with the initially well tuned FCL and does not have huge variations in daily eating
 526 and moving around, will **not use** the TT **at all**. FCL is possible without an intervention via the

527 TT button in your cockpit. Actually 4 of 8 modes (GGG ...YYY permutations, list see [section](#)
528 [5.3.1](#)) are not making use of TT.

529 (2) Super easy is also, to just input **any odd-TT** (odd-numbered temporary target) that will shut out
530 any SMBs for the set duration. *That can be a good idea when having a snack, for instance.*

531 *Super quick access to stop SMBs is possible also via the loop icon ([section 5.2.1](#)).*

532 Specifically, an **EatingSoon TT** can be activated here (limited relevance see [section 2.5](#)). It is
533 time-un-critical, can be manually set, or come up via an Automation.

534 *The cockpit enables you to set the iobTH differently (override) for the current meal.*

535 *Alternatively, iobTH can be temporarily changed in /preferences or using an Automation.*
536

537 *Temp. iobTH will always revert to default when the TT expires. If another TT immediately*
538 *follows, like in the example of the screen above, it will calculate, (then) show and use a new*
539 *temp. iobTH.*

540 (3) The third way is to **use the input mask** (*if already ncluded in your software version see picture*
541 *above*) **to freely modulate the loop aggressiveness** for a declared number of minutes. Click
542 the bottom big square(s): Either HYPO, or ACTIVITY, or EATING SOON, or ACTIVITY and
543 EATING SOON (*example in the pictured screen above*). Make or override entries in the offered
544 fields. Press OK.

545 (4) The fourth way is to exclusively use one of the 4+4+2 little buttons seen in the bottom part of
546 the TT dialogue box (*if already included in your software version*). They provide a set of
547 settings (as will immediately show in all input fields above) that the user has set up in
548 Preferences/SMB/autoISF/FullLoop (refer to [section 6.3](#)), and can freely label there. *For*
549 *instance „hiC“ at high carb EatingSoon, „piz“ for Pizza/fatty meals, „grd“ for garden work,*
550 *„mtb“ for mountain biking ...*

551 **Capturing good settings for not-everyday situations in /preferences** (*if already included*)
552 **allows calling them up within 1 second**, from your cockpit on the AAPS home screen (...and
553 won't ruin the FCL experience at all , especially because in most cases it is not time-critical,
554 how long before the intended exercise the buttons are pressed).

555 [Case study 6.2](#) demonstrates that nearly the same performance and comfort can be reached
556 via the **DIY FCL cockpit** with the grey extra buttons appearing at the bottom of the AAPS home
557 screen, based on Automations with User action (see also [section 4.2](#)).

558 *The example picture given above, and also [case study 6.2](#), is the most complicated (but also most*
559 *useful) case, **when exercise follows after a sizeable meal**. It is then that you need (a) aggressive*

560 FCL initial performance at the meal, but, exactly when (!) a (for the intended sport already
 561 temp.lowered) iobTH is exceeded, you need (b) to have SMBs automatically switched off and go
 562 into the „milder“ mode, as defined for the exercise (with high instead of lowTT, that automatically
 563 significantly reduces iobTH again, and insulin sensitivity(resistance) settings too).
 564 Pressing exercise related buttons will automatically also light the **exercise button** on the main
 565 screen yellow.

566
 567 To summarize, the TT dialogue field offers easy but powerful ad-hoc modulation of loop
 568 aggressiveness for FCL (if already included).

569
 570
 571 5.3.3.2 Exercise button (see more in [section 6.](#))

572
 573 The exercise button automatically lights yellow when exercise related TTs are activated in the TT
 574 dialogue box. 4 of 8 modes are making use of the exercise button.

575 If pressing on the exercise button, a dialogue box appears (*if extended design for FCL cockpit is*
 576 *already launched*) with info on exercise setting first (and opportunity to override), plus below the
 577 activity monitor (experimental for auto-tracking of lighter movement during the day, and effects on
 578 sensitivity that may have. See [section 4.5](#)).

579
 580 So, first the exercise settings (as set under TT) are there to read. Example :

| activity | TT | dura | %-profile | iobTH | bgAccel | %overall |
|----------|-----|------|-----------|-------|---------|----------|
| mtb | 171 | 180 | 70 | 3.0-U | 0.24 | 70 |

Mode set to run for 134 more minutes

Mode starting after meal when iob > 4.4 U or n/a

581
 582
 583 The exercise (here mtb) is selected in the dialogue box of the neighboring TT field, and there auto-
 584 filled with settings made in the set-up and tuning stage by the user under preferences (see
 585 above,...). They are reported also under the exercise button here, and TT, duration, and % sens
 586 (which also shows active on the %-profile field on the left side of the exercise button) can be temp.
 587 changed there. iobTH, bgAccel_ISF and overall resulting sensitivity ratio is given in the other fields.
 588 The **middle field** of the table, „% profile“ either picks up the % set under the %-profile button, or
 589 an input can be made here, in the exercise button domain, which will:

- turn the neighboring %profile button on yellow and show that inputted % on it, too
- be multiplied with the result from the exercise mode settings per se, and change the % overall, accordingly.

So, if this middle field of above table (dialogue box of sports button) contains a figure other than 100, input field becomes yellow, and you are operating with a combination of traditional PLUS new exercise mode (with all three top buttons of your FCL cockpit yellow). This maximally will soften aggressiveness, for which you get an idea by the last calculated figure.

The mode is either running already (for another number of minutes, as probably also shown in the yellow TT field anyways). Or it is scheduled to run, after insulination for a started meal reaches iobTH (as in table). Or, no exercise is scheduled (both points red, no entries).

The lower part of the exercise dialogue box (not pictured above, but see in [section 6.5](#)) is dedicated to the Activity Monitor

5.3.3.3 Profile button

The profile button can still be used to set a different profile, or profile%, for instance to adjust for days with sickness (as you are used to from hybrid closed looping). 4 of 8 modes are not making use of the profile button.

Any inputs made here will be used to modify profile_ISF on which all further changes are made on (multiplied with).

The profile field remains grey if standard profile is applied.

It turns yellow, displaying a %figure relating to any altered loop overall aggressiveness:

- When no inputs (changes from 100% profile) are made here, but inputs in the TT field, e.g. for exercise, automatically lead to different insulin sensitivity ratio ~~that ratio is shown here~~
- when% is changed by input in the profile button itself, it will be multiplied with with profile_ISF and be used in place of profile_ISF *by the algorithm*.

However, for exercise (sports) you no longer must make an entry here, because reasonable %reductions should be automatically provided, driven by your set TT (and half-basal exercise target), see [section 6](#).

5.4 Recognizing your loop state in the AAPS home screen

627

5.4.1 Color scheme of top cockpit buttons tells kind of closed loop that is running

629 3 Buttons (%profile; exercise; TT) each in 2 states (yellow Y, or grey G) makes $2 \times 2 \times 2 = 8$ possible
630 combinations:

631 GYY = dynamic exercise mode

632 YGY = not-dynamic „traditional“ exercise mode (if <100%) or hypo mode (if >100%)

633 GYG = basic closed loop with Activity Monitor running

634 GGG = basic closed loop (FCL or HCL) without any altered sensitivities etc

635 YGG = basic closed loop but with a „long wave“ sensitivity shift (e.g. sickness)

636 GGY = temp. target like e.g. EatingSDoonTT is set; or Hypo mode

637 YYG = closed loop with „long wave“ sensitivity adjustment and Activity Monitor running

638 YYY = dynamic exercise mode in time with additional „long-waved“ sensitivity shift

639

5.4.2 Information printed on the top buttons


641

642 The yellow TT field shows the currently valid TT (and further duration):

643 (profile) stands for the abbreviation you labeled your selected running profile




645 In the special case of settings for meal preceding sports, the field will look slightly differently:

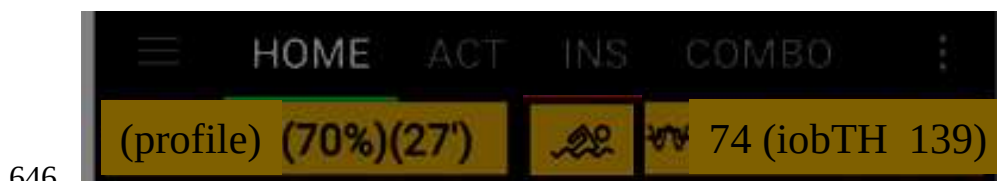

...and ...

647 ... when iobTH is first time exceeded, this automatically switches to:

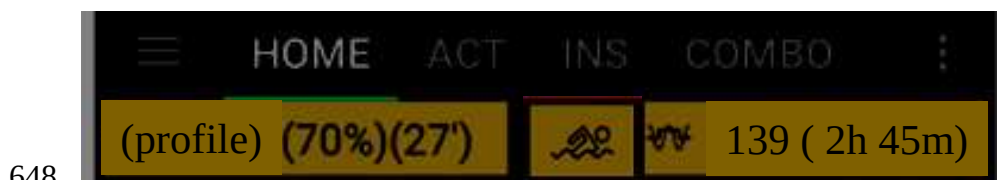

649

650 Likewise, if on the AAPS main screen just an **EatingSoonTT** is set (e.g.72), this is entered with the
651 desired duration. Afterwards, it automatically reverts to profile target and the display turns grey

652 again there with e.g. 90 on it (and no time limit).

653 Without sports context, the middle field remains grey.



654
655

656 Independently from setting a TT, the user can choose to set a **%profile in the left top field**, for an
657 independent number of minutes, e.g. 70% in this screen example: Also, or additionally, this will
658 influence the resulting ISF and sensitivity%

659



660
661

662 The % might change and turn yellow also in context of making TT inputs in the related dialogue
663 box (see chapter TT dialogue field, above). Still, the % (or the length of time the profile switch shall
664 be active) can be independently overridden in the top left field, if so desired.

665

666 If an **Automation** sets a %profile, and/or a TT (e.g. *automatic detection of meal start at condition*
667 *e.g. when delta >10*), this would automatically show in respective field(s) turning yellow and
668 showing the temp. setting. To show the set parameter comes from an Automation, „ **(A)** „ is added
669 in the end of button text.

670 Note that an Automation is **usually/ always (?)** only permitted to temp. change default profile
671 settings, not other pre-existing temp. settings. This is for a good reason : Why should a
672 sometimes in the past thought-out Automation supersede your - just for the occasion
673 specified – temp.settings that you consciously activated for the day?

674 Advice: Try to stay away from Automations that also aim at temp. modifying
675 aggressiveness. For the reason just given in above note, they often will not kick in anyways.
676 Generally, it also is no good idea to double up sub-algorithms for tweaking loop behaviours.

677

678 Try to **keep things as simple and clear as possible.**

679

680 That said, a limited number of Automations can be of help in distinct scenarios (that differ in
681 purpose and in applicable time of day).

682 A good one could be for night time, when your odd profile TT has SMBs shut off, but your
683 experience after pizza nights tells you that, under certain condition patterns (bg, iob), an
684 SMB or two should be „allowed in“. Another good example, if you go usually FCL without

685 any use of the TT button (which you could call a meal announcement of sorts), is to define
686 an Automation that, after detecting a meal start, automatically sets a low TT to get
687 maximally aggressive first SMBs.

688

689 5.4.3 FCL related indicator fields in the AAPS home screen

690

691 In extra data fields of the AAPS main screen you can always see (not change) the key
692 „aggressiveness“ parameters your loop currently operates with (see also home screen
693 example below):

694 • how profile sensitivity (**ISF**) changes by the %profile input, by autoISF, and/or a set
695 exerciseTT.


696 • next to current available iob number is an indication of your **valid iobTH** (the iob above
697 which no more SMBs will be given)

698 • The AAPS home screen additionally shows, above the deltas, the current **acceleration**

699 Having a look at that can be valueable. For instance, when glucose is relatively low and still
700 falling, a positive (and getting more positive) acceleration indicates that bg will swing back
701 up, rather than crash low. This will give info about necessary snack size, and hence help
702 avoid both, unnecessary calories, and going on a bg rioller coaster.

703
704
705 5.4.4 Overall home screen:
706

Overall home screen:

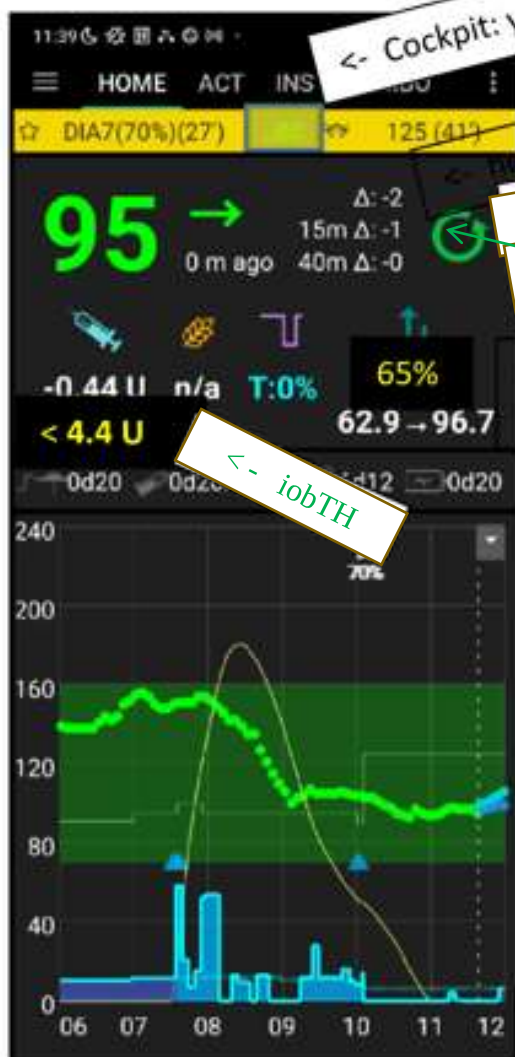

<- Cockpit: yellow fields=>temp. modulated sens.

are additionally: acceleration-factor

<- Violet <-> green circle for FCL <-> HCL dotted if SMB off
in the middle: minutes counting down if temp.set;
below: (A) if coming from an Automation that is running

<- iobTH

<- % reduced insulin supply... because of resulting...
<- higher (weaker) ISF

<- buttons „bolus“ „carbs“ etc. eliminated
(auto- re-appearing when violet -> green loop)

707
708
709