

## 5. Modulation of autoISF aggressiveness. V 2.5-

**Please note that with autoISF 3.0 you are in an early-dev. environment,** where the user interface is **not optimized for safety** of users who stray away from intended ways to use. Good safety features exist, but these are only as good as the development-oriented user understands and implements them. This is not a medical product, refer to disclaimer in [section 0](#)



Once the initial tuning according to [section 4](#) is done, you are ready to use autoISF for your fully automated meal management.

You will have three major *other* challenges to manage:

- recognize and manage (partial) occlusions, or other technical (CGM or BT related) obstacles (see [section 2](#) on pre-requisites of FCL, and related case studies)
- deal with times when insulin given by the loop must be restricted (e.g. a snack could be “misinterpreted” as a meal)
- deal with times when the loop should be set “milder” as a precaution (e.g. nights; or in an exercise context).

How big the remaining challenge is depends very much on your individual lifestyle.

[Sections 5](#) and [6](#) discuss this in more detail.

**In order to run the loop fully automatically around the clock, the times *outside* the meal blocks must also be precisely analyzed, and solutions to problems (if any) must be sought.**

It is up to every user to decide where to draw the line.

- With a technically well functioning system, moderate meals, moderate or no exercise, moderate %TIR expectations and a bit of mindfulness it should be possible to go into Full Closed Loop 24/7, after working through, and observing, [sections 1-4](#).
- Especially if you are a bit shy of using the emulator for really detailed analysis, it is likely that you will not hit *one* real good system calibration ([section 4](#)) for your *entire range* of diets.  
In that case you will occasionally run out of range, and your options to prevent, react, or improve are

- 35           ○ accepting a few % higher time outside range for that day (and, if feasible, in the
- 36           future avoiding what seemed to have caused it)
- 37           ○ taking a snack (whenever you tend to go low from the “tails” of insulin activity that
- 38           was required to fight a peak)
- 39           ○ doing a manual override (if you can think of one in time, to manage the problem
- 40           manually)
- 41           ○ temporarily resorting to the well-known hybrid closed loop.

42 Instead of accepting such instances, you could launch “improvement projects”

- 43       • that refine your initial tuning ([section 4.](#))
- 44       • that make you and your FCL loop fit to manage an increasing number of disturbances
- 45       either automatically, or via a user intervention ([sections 5-6](#)).

46 To tailor the loop’s response to disturbances *other-than* your major meals probably will require

47 specific **modulation of the aggressiveness** that you set according to [section 4](#) for your meal

48 management.

49 There are many avenues to achieve this. The main ones, that are also easy accessible via

50 Automations in AAPS, are:

- 51       • temporary shut-off SMBs (odd-numbered target)
- 52       • temporary reduce bgAccel\_ISF-weight
- 53       • temporary reduce iobTH
- 54       • temporary reduce set %profile
- 55       • temporary set higher TT (especially in connection with exercise mode)t

56

57 In setting up your FCL, you now have another difficult and time-consuming job at hand, to define

58 solutions for any of your „other“ situations (outside of meal management) that tend to drive glucose

59 outside of the desirable range.

- 60
- 61       • In [section 5.1](#) we explore avenues towards *fully automated* management that in daily life
- 62       will require no user intervention at all.
- 63       • In [section 5.2](#) and [5.3](#) we will look at solutions that involve an easy *user interaction like a*
- 64       *data entry or button push.*

65  
66  
67  
68  
69  
70  
71  
72  
73  
74  
75  
76  
77  
78  
79  
80  
81  
82  
83  
84  
85  
86  
87  
88  
89  
90  
91  
92  
93  
94  
95  
96  
97  
98  
99  
100

## 5.1 Fully automatic modulation of FCL aggressiveness

The following subchapters describe set-ups you may want to use for allowing **completely hands-off FCL in as many daily situations as possible.**

### 5.1.1 autoISF generally switched off outside of meal-time windows

If, aside from having to bolus for meals, your hybrid closed loop was running pretty well *without* other interventions from your side, you could continue to run in that mode, and just focus your new autoISF FCL on management of meals.

In your initial transitioning phase this approach makes a lot of sense, and even by focusing autoISF on just a sub-set of them, like only dinners.

Also in the long run this avenue is taken by many FCL users for the night times, “hanging on” to their well performing hybrid closed loop with standardoref(1) SMB+UAM

For this, you define Automations

- that set meal time windows in which “Enable ISF adaptation by glucose behavior” (autoISF) is turned on in AAPS preferences/OpenAPS SMB
- or: that turn *all* autoISF’s ISF modulations (*or just bgAccel\_ISF*) off in time windows in which surely no meal occurs. For instance, you can go for all nights back into your Hybrid Closed Loop, as you had before.

Other early DEV AAPS variants (see [section 13.3](#)) all work with meal-time windows. The window is either set by time of day in the settings, or it always must be „set“ by the user via giving a mandatory small pre-bolus before any meal starts. **Outside** of these time windows, these loops then runs with less aggressive SMBs likeoref(1) SMB+UAM in AAPS Master. This mode is not really FCL, but an advance over traditional HCL that often achieves satisfying degrees of automation and performance.

The term **Meal Announcement** (MA) is often used to label this closed looping mode. Trigger to set a meal time window could also be a pre-bolus given by the user, a carb entry made, an EatingSoonTT set, or a meal announcement button pushed.

101 Note: Outside of the meal time windows you would be in hybrid closed loop. To the extent you  
102 rarely face disturbances (aside from meals), you could be looping in full automatic mode around  
103 the clock,

104

105 Your temp. "autoISF shut-down" (exiting autoISF FCL = shutting off "Enable ISF adaptation by  
106 glucose behaviour") is meant to prevent problems from the loop *over-reacting* to bumps in the  
107 glucose curve in times of day (night) when standard oref(1) performance is sufficient.

108

109 A very good alternative to fully resorting to Hybrid Closed Loop is "taming" the FCL via a night time  
110 SMB shut-off (see next [section 5.1.2](#)).

111

### 112 5.1.2 Odd-numbered *profile* targets used to block SMBs

113

114 An alternative route of preventing the FCL loop from over-reacting to bumps in the glucose curve  
115 would be to make use of the option to temporarily shut down SMBs

116

117       Ensure the even/odd logic in the settings is toggled on in Preferences> openAPS SMB>  
118       autoISF settings> smb delivery settings>: "Enable alternative activation of SMB depending  
119       on profile target" ON.

120

121 In time blocks with an odd-numbered profile target you can prevent any SMBs being given by your  
122 loop. The (unchanged) aggressive settings then can only translate within the limits set by %TBR  
123 possible.

124

125 This will very much slow down any more insulin being given, and is an excellent solution for night  
126 times, especially if you occasionally experience compression lows.

127

128       Alternatively, you could use the new included options for Automation Conditions and  
129       temporarily tune your bgAccel\_ISF\_weight much lower ([section 5.1.4](#)).

130

131       Yet another alternative was already presented ([section 5.1.1](#)) = to go into hybrid closed  
132       loop for the night.


133       That is possible to do with SMBs available (without them getting boosted via autoISF), and,  
134       for a long time, was the author's favoured solution for the nights.

135

136       But, my current favourite builds on the method of this section (5.1.2, odd profile target  
137       provides SMB shut off), but then allowing some, automatically triggered when needed:

138

139 In case you occasionally do have nights that would benefit from a couple of SMBs (to treat temp.  
 140 highs from a late fatty pizza, raclette and such): Define suitable Automations like the two „night“  
 141 ones in this *example*:

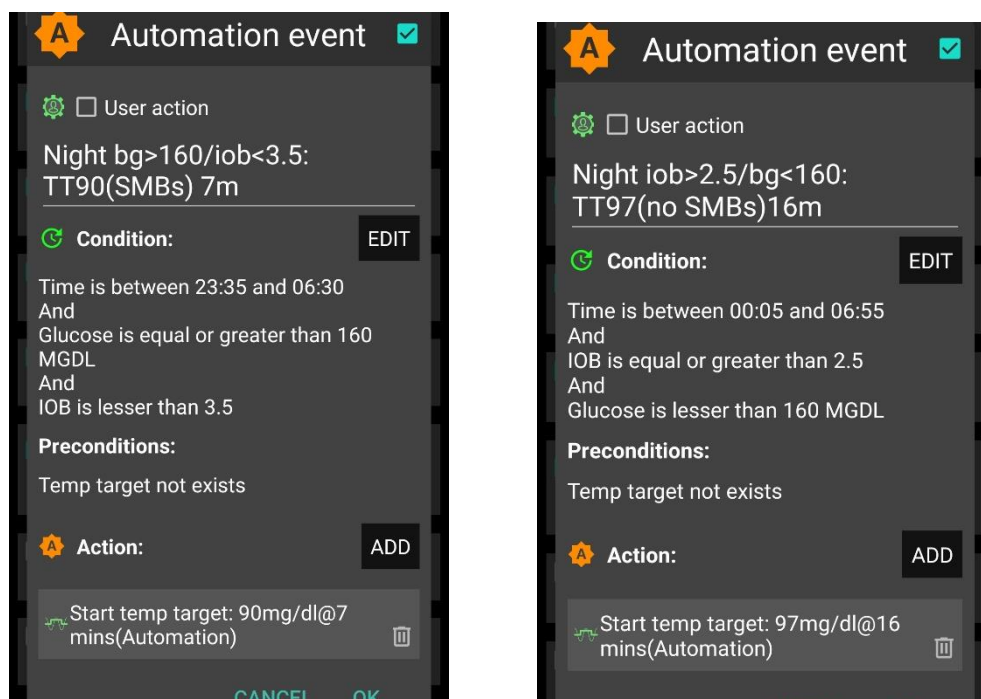

142  
 143  
 144 Never underestimate the „trickyness“ of getting your Automations „right“.  
 145 With your thought-out Automations in place, night data need to be analyzed to see  
 146 • whether the bg and iob limits defined in the given example work sensibly four your data  
 147 pattern  
 148 • whether the TT duration is chosen appropriately  
 149 • how swapping the sequence in which the automations appear in the Automation list would  
 150 lead to different SMB impacts.

### 152 5.1.3 Odd-numbered *temp.* targets (TT) set via Automation used to block SMBs

153  
 154 A widely used Action that strongly modifies how fast your FCL can add more iob is setting an **odd-**  
 155 numbered **temp. glucose target** which makes the loop operate without giving any SMBs (%TBR  
 156 modulation only).

157 Ensure the even/odd logic in the settings is toggled on in Preferences> openAPS SMB>  
 158 autoISF settings> smb delivery settings>: "Enable alternative activation of SMB depending  
 159 on TempTarget" ON.

160  
 161 So, from patterns you find in YOUR data, at times where you want your loop act differently, you  
 162 need to carve out Conditions that describe the respective situations (and either *for how long* it  
 163 typically lasts, or at which *other* Conditions you want your loop get back to default FCL operation).

164 An odd TT is often set for an *anti-hypo* snack or *sports* snack. In both instances, you do not want  
165 SMBs to quickly counter act.

166 In case of *sweet “fun”* snacks, this is entirely different -> [section, 5.2.1](#) or for regular snacks  
167 (*e.g. at school break*) see next [section 5.1.4](#)

168

#### 169 5.1.4 Automatic differentiation of FCL aggressiveness using Automations

170

171 **Personalized Automations** tailor the loop exactly to **YOUR** data so **fully automated** handling of  
172 situations with **different aggressiveness** of the loop can be made.

173

174 From, autoISF 3.0 onwards, also the following parameters are provided as Condition and/or as  
175 Action for defining YOUR Automations:

- 176 • Enable ISF adaptations by glucose behavior => Allows temp. ON/OFF for the key ISF  
177 modulation parts of autoISF (and, as a result, will usually decrease loop aggressiveness)
- 178 • Trigger/set iobTH percent => Keeps default aggressiveness, but only until a iob threshold  
179 (that your Automation modifies) is surpassed (which is when any further SMBs will be  
180 blocked blocked)
- 181 • Trigger/set bgAccel\_ISF\_weight => Modifies the default aggressiveness of just the  
182 acceleration component

183

184 To set up suitable Automations, you first must **analyze patterns** you find in **YOUR data**, at times  
185 (or geo-locationa, or bg and iob patterns that point to a problem ...) **where you want your loop**  
186 **act differently**, to carve out Conditions that describe the respective situations (and either for how  
187 long it typically lasts, or at which *other* Conditions you want your loop get back to default FCL  
188 operation).

189

190 A variant of this mode is to define several windows in which autoISF aggressiveness  
191 (bgAccel\_ISF\_weight) and/or iobTH are automatically set differently

- 192 • for different meal time slots of your day –

193 (*Breakfast at home, school lunches, school intermission snacks, dinners at home* could for  
194 example all deserve special settings regarding ISF\_weights and iobTH).

- 195 • or even for a geo-location etc –

196 (*School lunches, or mother-in-law visits, would be examples*).

197


198 Unless your meals differ vastly in size and in fast carb content all this may not be needed.

199 Still, personalized Automations might help ease your initial job of setting the various ISF\_weights,  
200 and a best-suitable iob\_threshold\_percent that would work “always”.

201

202 **Caution:** If (as in autoISF 3.0) setting a different iobTH or bgAccel\_ISF\_weight can not be done  
203 temporarily (i.e. with a duration attached) you **must** define a suitable **additional Automation**, that  
204 must be active in tandem, that **restores the default**  
205 **set iobTH or bgAccel-ISF\_weight again**. Else, once  
206 your Automation set in, it will *forever* shift this  
207 important parameter setting!

208 If for instance you have several Automations that, in  
209 combination with a set elevated TT also set a lower  
210 iobTH: Don't be fooled, the duration only applies to  
211 the TT. You need an extra Automation for all of them.  
212 I picked out the highest of the altered iobTH values  
213 that these Automations can set (45\_percent), and  
214 then I can automatically restore my default desired  
215 60% via this one Automation (see screenshot - - > )

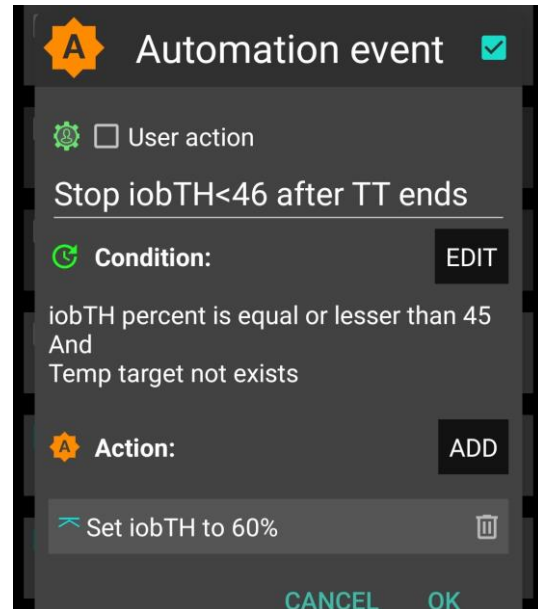

216

#### 217 5.1.5 Automatic adjustment of FCL aggressiveness via the Activity Monitor

218

219 With the autoISF variant of AAPS you can make use of your smartphone's **stepcounter** and use it  
220 to fully automatically adjust insulin sensitivity ratio to **activity level in the past minutes to one**  
221 **hour** time frame.

222

223 This feature comes with yet another little tuning opportunity, in which you study your body's  
224 response to light exercise (like walking) or to not moving at all (like desk, couch), and select  
225 appropriate settings which, in the future, will automatically adjust insulin delivery to suit activity  
226 state of the past minutes (up to 1 hour). (AAPS Preferences/OpenAPS SMB/Activity modifies  
227 sensitivity/ -> set **two scaling factors**.)

228

229 This autoISF feature (new since V.3.0) is much quicker responding than Autosens or dynamicISF  
230 to adjust insulin sensitivity to your current „lifestyle state“.

231

232 For loopers who do not have huge variations in exercise levels in their everyday lives, this feature  
233 might fairly much close the gap towards being able to do a 24/7 hands-off FCL.

234

235 [Sections 3.5](#) and [6.5](#) describe the Activity monitor in more detail.

236



## 237 5.1.6 Pro/con completely hands-off Full Closed Loop

238

239 To stay 24/7 in a completely „hands-off“ FCL can be a realistic goal with autoISF 3.0 if besides  
240 meals also some special challenges, as discussed in this [section 5.1](#), were analyzed and could be  
241 addressed.

242

243 Clearly it depends very much on your lifestyle, and how interested, willing, and capable you are to  
244 recognize, deal with, (and in the future avoid?) situations that get you outside of your desired %TIR  
245 on occasion.

246 So, this is also about what %TIR you are aiming at, and can accept, as it averages out for  
247 the week, for instance.

248

249 Everybody must weigh for her/himself

- 250 • how much **upfront effort** to put into the setting up process for getting it all 100% automatic
- 251 • **or whether to take an easier start, with a couple of situations left to take care of when**
- 252 **and as they arise in daily life**

253

254 Also, even if a principal capability for a fully automatic running FCL is given, this still  
255 means that

- 256 • the user should be knowledgeable about what exactly is going on, and
- 257 • have a principal capability to „nudge“, or to take over in a manual mode.

258

259 In the sections that immediately follow, we present the options to nudge or temporarily take over  
260 from the AAPS home screen which will be serving as your **FCL cockpit**:

- 261 • [Section 5.2](#) describes how you can use available “buttons” from your AAPS home screen,  
262 and how to complete it towards a suitable DIY FCL cockpit, for an even better FCL  
263 experience.
- 264 • [Section 5.3](#) describes how you might be able to manage “disturbances” even better (with  
265 more convenience) with an **improved FCL cockpit in the future. (No need to read any of**  
266 **the green lines, unless you are interested in contributing to define/design/program further**  
267 **improvements)**

268

## 269 5.2 Modulating aggressiveness manually, from the DIY-FCL-Cockpit\*

270



\* Like in the airplane cockpit: Cruising in full auto mode should involve having an eye on the instruments, and on potential disturbances ahead in the environment.

In [section 4](#), we dealt with major meals. In [section 5.1](#) we looked into fully automatable management of other situations. Life in Full Closed Loop can become extremely easy then...

However: Other **disturbances** might come up, that:

- are not noticeable in-time, or foreseeable, by the loop (*e.g. your plan to start exercise in an hour or two*), but **that influence sensitivity dramatically** and therefore require temporary non-default settings in order to remain in-range, and/or
- **require** a different “starting point” regarding iob and bg, which translates into **a different iobTH** that should **temporarily** be set much lower (*in case of exercise*) or noticeably higher (*e.g. with very fast absorbing carbs in a sweet snack “sin”*).

In [section 5.1](#) we looked into ways to automate also a modified loop response to *foreseeable* situations (tied to a time of day, geo-location etc), or to those *the loop could recognize* (with enough time to react).

Other „disturbances“ might come up, and you must find an easy way to

- call up a pre-programmed routine for automatic management, with adjusted aggressiveness, or:
- manually tweak a setting or two, to temporarily adjust the aggressiveness
- There may also arise a desire to just exit the FCL mode, and be your own captain for mastering a special situation.

For peace of mind, to learn, and to stay informed (especially so in your initial tuning phase, or when your glucose curve goes in unexpected ways) we also must be able to

- find the key parameters that frame and drive the recent and upcoming loop decisions.

All this is facilitated within seconds right from the AAPS home screen, serving as a **FCL cockpit** after you built a couple of DIY cockpit features via Automations (as described below and in [case studies 5.2](#) and [6.2](#)):

Thoughts went also into [how to improve the cockpit in future releases](#), see [section 5.3](#)

## 306 5.2.1 Status recognition

307

308 Before considering any manual interventions into the ongoing FCL, you should be aware what the  
309 current mode of action is, and hence which button eventually to fine-tune or lever to switch, in or-  
310 der to adjust to the disturbance you see coming up.

311

312 See [section 5.4](#)

313

## 314 5.2.2 Manual interventions from the (DIY-) FCL cockpit

315

316 Trouble with all these is, not to forget to set back manually, too (=> better solutions in 5.3)

317

### 318 5.2.2.1 Temporary tuning of FCL aggressiveness via temp. %profile or TT settings

319

320 The set **% profile** multiplies with both, the ISF resulting from autoISF, and also with the default  
321 iobTH you have set, so both are nicely modulated in a linear way with the % temporarily chosen

322

323 Just taking profile e.g. to 110% for an afternoon might be an easy way to explore whether  
324 you might benefit from 10% more “aggressiveness” in your core settings for lunches (like  
325 bgAccel\_ISF\_weight). Make sure, though, that the extra 10% are not cut away by set safety  
326 limits.

327

328 A lowered (relative to profile glucose target) temporary **bg target (TT)** signals lowered sensitivity  
329 (more insulin need), and an elevated TT (as often used with exercise) increases sensitivity and  
330 hence works in the direction of a lowered % profile to also reduce insulin given by the loop.

331

332 Moreover, the **exercise button** (top center on your AAPS home screen) can be activated (turns  
333 yellow, then). This will **further boost** how your set TT elevates the resulting ISF, and sharply  
334 lowers iobTH, as often desired for sports. See [section 6.1](#)).

335

### 336 5.2.2.2 Making temporary changes in settings made in AAPS/preferences/Open APS SMB

337

338 Going into AAPS/preferences/Open APS SMB allows to:

- 339 - set milder or stronger ...\_ISF\_weights
- 340 - set different iob\_threshold\_percent (or iobMAX)
- 341 - elevate or lower the SMB\_delivers\_ratio

- limit or expand max. allowed SMB size
- change the the even <-> odd logic for SMB on/off

344

345 Doing temporary changes in AAPS/preferences should be the exception because

- they require multiple steps, including entering a password
- you will often forget to set everything back to default a couple of hours, or minutes, later

348

### 349 5.2.2.3 Triggered Automations: Grey extra DIY cockpit buttons for pre-programmed “responses”

350

351 Recognizing conditions for fully automatic handling by the loop may not be not possible, or come  
352 too late for the loop to act on. Examples would be

- *exercise*: Minimum an hour before starting “the loop should know” to be able to lower iob and elevate bg by the time exercise starts.
- *snacks*: High carb snacks, sweets, consuming ice cream or having a sweet drink comes with the problem of even steeper glucose rises, but overall a lesser insulin need, compared to major meals (for which we tuned our FCL according to [section 4](#)).

358

359 This not necessarily implies that snacks need different settings than a meal. After all, autoISF  
360 was designed to act to all available data, especially to where the developing glucose curve is  
361 headed. So, depending on your effort to set parameters for a broad variety of meals (notably:  
362 how well you avoid to invariably bounce fast against your iobTH), you might be able to accom-  
363 modate low carb, snack, and major meals with *one* set of settings.

364

365 In FCL autoISF, this is a bit more difficult than in HCL autoISF applications, because FCL  
366 involves revving up iob supply (largely via big bgAccel\_ISF-weights) often too much to be  
367 balanced by just a snack getting absorbed.

368

369 For that reason, or just for increased comfort and safety, you might want to differentiate, and make  
370 use of what follows for the *sweet snack* example.

371

### 372 Tuning aggressiveness

373 Key is that a sweet snack likely benefits from even more aggressive initial FCL  
374 performance than the meals in your normal spectrum of diets require.

375 Therefore, you could set

- a higher **temp. profile%** and/or
- a temp.elevated **bgAccel\_ISF-weight** (see screenshot of my Automation).

- a **low temp. target** (76 for instance; this additionally helps maximize the first SMBs that will automatically be triggered at detection of acceleration)..

When first defining and testing this Automation, also check:

- that the safety limits as discussed in [section 2](#) will not block the intended elevated aggressiveness
- SMBs will not get outrageously big and iobTH sometimes exceeded by too much


Note that “the last SMB” is allowed to overshoot the effective iobTH by 30%

### Limiting iob

For “just a snack”, total insulin need will be lower than for a meal.

If you would just have your sweet drink, and your meal-oriented FCL would “attack”, iob likely would become too high, and a glucose rollercoaster would start, with you needing to consume more =>

If you just have a snack, or drink a glass of juice, you can lower the **iobTH\_percent** accordingly.

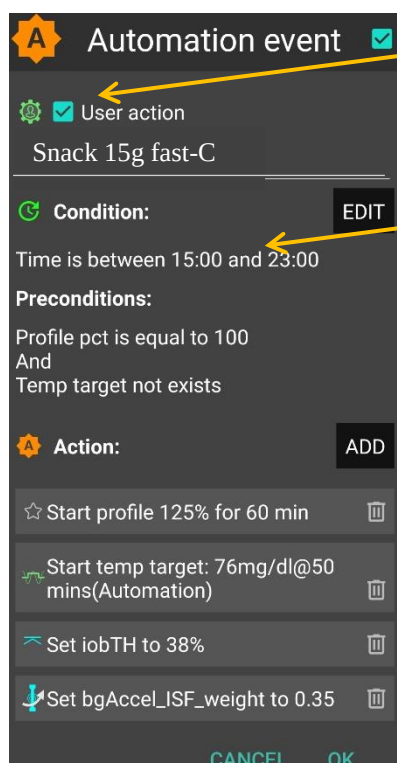

„User action“ is always ticked-on

This will, in the defined time space ..

..offer the “DIY cockpit” button..

..which I must press any time (~90...30 minutes) before my snack.

Note: Make sure that iobTH and bgAccel\_ISF\_weight revert to default afterwards



399

400 **In everyday life** you **then just** must **press the related button in your cockpit** (which is  
401 not time critical at all, except it should be clicked *latest* a couple of minutes after you took  
402 the drink or snack).

403

404 If you consume more, and also eat something with your sweet drink, this will more  
405 resemble a full meal... however, with unusual amounts of fast carbs.

406

407 Note: Pressing your snack button *a second time* would **not** help because the  
408 lowered iobTH does not allow iob going high enough. So you are better off just  
409 letting your normal FCL meal routine run, after your snack mode expired.

410

411 Other options when snacks keep extending would require a manual modulation  
412 regarding %profile and/or bgAccel\_ISF, but keeping the full default set  
413 iobTH\_percent, or even elevating it (refer to [section 5.2.3](#)). If that happens often,  
414 define for yourself an extra User action Automation for a bigger snack (= another  
415 grey DIY cockpit button).

416

417 **Caution:** If (as in autoISF 3.0) setting a different iobTH or bgAccel\_ISF\_weight can not be done  
418 temporarily (i.e. with a duration attached) you **must** define a suitable **additional Automation**, that  
419 must be active in tandem, that **restores the default**  
420 **set iobTH or bgAccel-ISF\_weight again**. Else, once  
421 your Automation set in, it will *forever* shift this  
422 important parameter setting!


423 If for instance you have several Automations that, in  
424 combination with a set elevated TT also set a lower  
425 iobTH: Don't be fooled, the duration only applies to  
426 the TT. You need an extra Automation for all of them.  
427 I picked out the highest of the altered iobTH values  
428 that these Automations can set (45\_percent), and  
429 then I can automatically restore my default desired  
430 60% via this one Automation (see screenshot - - > )

431

432

433

434

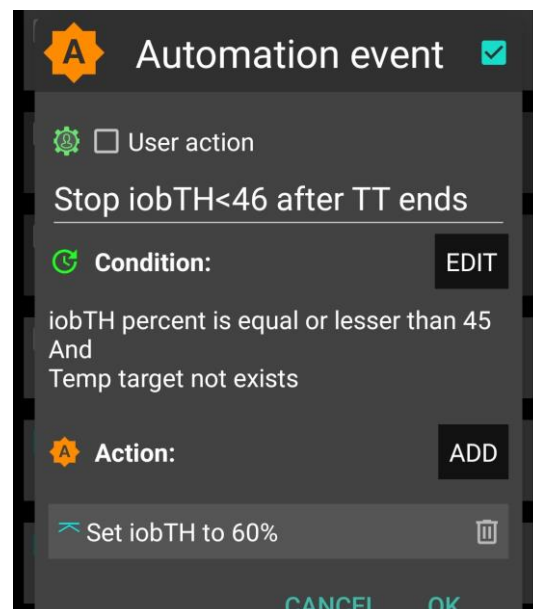

## 435 Installing the DIY cockpit button

436 In the related Automation, just keep the “User action” box clicked at all times, and define in  
437 the Conditions when you want to see that button available for cockpit use (see screenshot  
438 above) => you will see that button offered.

439

440 Besides snacks, also any **other recurring special situations can be addressed via a**  
441 **DIY cockpit button, and receive different aggressiveness up to a suitable iobTH**  
442 **level.**

443

444 Over time you can have a big number of User action Automations, and keep them  
445 “shelved” rather invisibly (clicked in-active, top left box) in your long list of potential  
446 Automations. Even when active, they only show in your cockpit (bottom grey field of your  
447 AAPS home screen) in the time slot you assigned as potentially relevant.

448

449 In the future you might be able to set the stage for a snack and other “disturbances”  
450 also via an extended menu behind the TT button on the AAPS home screen, see  
451 [section 5.3.3.1](#)

452

## 453 Discussion

454 In case you do have a snack habit and

- 455 • can not find settings, as in [section 4](#). defined for your meals, also suit your snacks
- 456 • can not pin a as in [section 5.1.4](#)

457 you minimum need a “snack announcement”

458

## 459 5.2.3 Temporarily exiting the FCL

460

461 The “last resort” alternative always is to **temporarily** leave the FCL mode, and handle any  
462 disturbance “the traditional way” in **hybrid closed loop**. For this, we switch the automatic  
463 aggressive adaptations of ISF to the bg curve off that are only needed in FCL ....

464 (if in hybrid closed loop you like e.g. the dura\_ISF adapation still, you alternatively

465 could elect to just set bgAccel\_ISF\_weight temp. to zero, instead)

466 ... and now around meal starts giving a bolus will be necessary.

467

468 The suggested improved FCL cockpit user interface with an extra version of violet loop on  
469 the AAPS home screen ([section 5.3.1](#) ) would facilitate this transition FCL < - > HCL,  
470 including automatic removal and re-appearance of the insulin button at the bottom of the  
471 APS home screen.

472

473 In case [this feature](#) is not yet available, you must:

474 Exit the FCL mode by going to AAPS/preferences/put in your password/OpenAPS SMB/scroll down  
475 to autoISF settings and switch “Enable ISF adaptation..” OFF

476 (or, alternatively, set bgAccel\_ISF\_weight to zero).

477

478 **Caution:** Unfortunately, there is no way yet for your full closed loop with ISF adaptations to come  
479 automatically back on, after a selected time for instance. So do not forget to switch your autoISF  
480 fully back on, later.

481

482 As this will often be forgotten, it may be worth doing a “User action” Automation, for a “temp.  
483 FCL OFF” grey button (see [section 5.2.2.3](#)).

484 Caution though, there is very limited experience with this brand new feature. Make sure your  
485 Automation definition really applies a duration (or other condition) that will automatically  
486 terminate all non-default settings it made. As we have seen e.g. in [section 5.1.4](#), this is not  
487 always the case.

488

489 To recognize whether autoISF currently runs with ISF adaptation or not, you must consult the  
490 profile\_sens -> actual\_sens indicator below the Autosens%. However, this gets also modified  
491 by %profile switches or TT +/- exercise mode. So it is not as easy as it would be with the “violet  
492 loop” proposal mentioned already above.

493 Ultimately, you can of course study the SMB tab to find out what is going on.

494

495

## 496 5.3 Modulating aggressiveness manually from the improved FCL-cockpit

497

498

499 autoISF 3.0 is an early dev variant of AAPS, and as user you are participating in an on-going  
500 development. Of note, autoISF 3.0 is first launched [without many](#) of the cockpit features that are  
501 presented below [in this font color](#).

502



503 Only what is written in black is at this point of some relevance for using autoISF 3.0.  
504 **No need to read** any of the green lines, unless you are interested in contributing to  
505 define/design/program further improvements.  
506 This is also an open invitation for you to contact us in case you could help program a  
507 module for one of the suggested user interface extras.  
508 For future integration into AAPS Master, an eye should be kept also on the question which  
509 *other* modes (like FCL using Automations and others mentioned in [section 13](#); and maybe  
510 also HCL) might benefit from some of the extra features.

511 For the time being, multi-step work-arounds may become necessary

- 512 • In many cases, going into AAPS Preferences and changing settings would be needed  
513 (...plus not forgetting to change these settings back, afterwards).
- 514 • Automations allow a DIY FCL cockpit, see [section 5.2](#) and [case studies 5.2](#) and [6.2](#)

515  
516 Keep in mind, though, that the **goal should be to interfere with the loop as little as possible**.  
517 Under the described conditions it can run **fully automatically** without any user interaction (= after  
518 the initial tuning phase, and related settings made in AAPS /preferences/SMB/autoISF. See [section](#)  
519 [4](#). and [5.1](#)).

520  
521 However, just like in the airplane cockpit: Cruising in full auto mode should involve having an eye  
522 on the instruments, and on potential disturbances ahead in the environment.

523 *E.g.: storm ahead => instruct your plane to climb to another flight height.*  
524 *Analogy: exercise ahead => setting an exercise TT, or => pressing a button that activates a*  
525 *sequence of instructions (some of them probably hinging on conditions, like actual iob) how*  
526 *to manage through that exercise situation).*

527  
528 So, for the occasional „disturbance“ coming up, you should find an easy way to

- 529 • call up a pre-programmed routine for automatic management, with auto-adjusted  
530 aggressiveness, or:
- 531 • tweak a setting or two, to temporarily adjust the aggressiveness
- 532 • There may also arise a desire to just exit the FCL mode, and be your own captain for  
533 mastering a special situation.

534 All this is facilitated within seconds right from the AAPS home screen's **cockpit features** to the  
535 **extent they are already incorporated**, or to the extent you can build alike DIY cockpit features via  
536 Automations, as described in [section 4.1.3](#) and [case studies 5.2](#) and [6.2](#)):  
537

538 • The button that is integrated into the **violet FCL icon** serves as emergency off button, to  
539 quickly stop FCL, or to at least to immediately stop any more SMBs (...just for a couple of  
540 minutes, or for the remaining meal time: pick from the options offered with just one  
541 keystroke).

542 Via the violet FCL icon on your AAPS home screen, you also can access a temp. switch-off  
543 button for SMBs (see section that next follows below).

544  
545 • The **three top fields** (%profile, exercise, TT) provide access to temp. tuning of core  
546 parameters, and/or to some pre-programmed routines.

547 Taken together with some **new indicator fields** about your loop state ([section 5.4.3 – 5.4.4](#)), and  
548 the **grey DIY cockpit buttons** ([section 5.2.2.3](#)) this makes the AAPS home screen your **cockpit**  
549 for Full Closed Looping.

550

551 Let us look on each of these cockpit elements in some detail:

552  
553 5.3.1 Violet FCL icon and underlying buttons

554

555 Novices to FCL, or really anyone running into a very special situation, may appreciate that the new  
556 closed loop icon on the AAPS home screen in pink (for FCL) has buttons to quickly shut off getting  
557 more SMBs (1st row), or to enter other loop modes (second row).

558

559 It functions very much as the other ones that you know from HCL already, and in fact you  
560 get offered some of the same options (for instance, to switch the (full) closed loop off for 15  
561 minutes for going to take a shower)


562

563 Note that in FCL you leave all BG regulation, notably against meal spikes, to the loop. So, try not to  
564 disconnect in phases when your FCL must ramp up your job.

565 The required insulin would still be supplied *after* you reconnect. However, without the user  
566 pre-bolussing, the delay would be more of an issue in FCL than it had been in HCL.

567

568 Just pressing on the FCL icon, a dialogue box comes up:

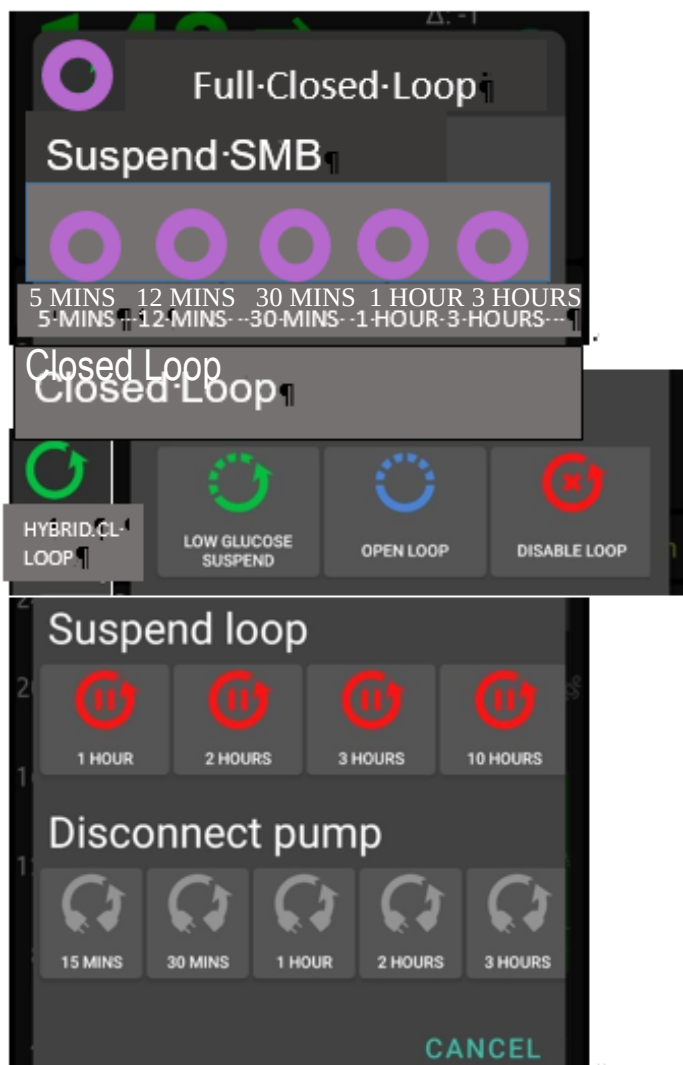

<-add-arrow,·form-like·like·CL·Loop··¶

<-dotted+·arrow,·form-like·like·LGS·¶



**DEV:**·when·in·Hybrid·Closed·Loop·=·top·  
headline+·green·circle),·there·is·no·Suspend·  
SMB·part·but·it·starts·right·with·„Closed·  
Loop“.·In·that·screen·There·the·1st·element·  
must·be·option·to·pick·¶

·····pink·circle·FULL·CL·LOOP·¶

where·here·the·green·HYBRID·CL·LOOP·  
stands·¶

For·loopers·who·did·not·set·up·FCL,·a·  
feedback·must·come·up·~·„FCL·not·installed“.·  
if·they·press·on·that·button·¶

569

570 Pressing „**Suspend SMB**“ provides fast and easy „emergency braking“ regarding delivery of more  
571 SMBs:  
572 Select the one with the desired number of minutes: 5 or 12 for just blocking the potential next  
573 SMB(s), and up to 3 hours to manage the entire rest of this meal with %TBR from then on.  
574  
575 Whenever, and why-ever, your FCL is in „no SMBs allowed“ mode (e.g. automatically after  
576 surpassing an iobTH, or triggered by a set odd TT), the FCL icon will turn into a dotted one.  
577 Instead of remaining **duration to end time** it indicates in the middle „the condition“, „**iob**“ or „**TT**“  
578 Add an indication **if** suspend SMB comes from an Automation, e.g. add an „**(A)**“, **underneath** the  
579 #minutes, iob, or TT in the middle of the dotted violet field.  
580 So, as in other (already in HCL existing) cases, those icons show in the middle the minutes left that  
581 they will be running, or the condition which would have to go away for this temp. setting to stop.  
582 It always auto-reverts into the FCL state and FCL icon, when time (or other condition) has elapsed.  
583  
584 Pressing „**HYBRID CL. LOOP**“ or other buttons from the 2nd row provides fast and easy  
585 „emergency **exit**“ **into other modes**.  
586 This enables beginners an easy „temp. escape“ into their well-known HCL (green) at any  
587 point of time. bgAccel\_ISF\_weight is set to zero when going FCL->HCL. HCL can run with  
588 autoISF (for instance dura\_ISF) uninhibited otherwise. (check **implications** for HCL users of  
589 autoISF ?? ).  
590 Note: These options from row 2 have no time limit. Loop will **not** by itself go back to FCL. You see  
591 the different loop icon as a reminder to manually revert, when ready.  
592  
593  
594 5.3.2 Buttons „Insulin“, „Calculator“ etc at bottom of AAPS home screen  
595  
596 These buttons are **not useful any longer in FCL**, and automatically disappear whenever in FCL  
597 mode (also in Suspend SMB state), and re-appear when leaving FCL. This applies also when an  
598 Automation or technical system failure shut off FCL.  
599 Users who, maybe in the beginning phase, feel better having those buttons, can override  
600 the removal (of the insulin button, or any other) by going into /preferences/overview/buttons  
601 and forcing them on. They only remain on until the next re-entry into FCL mode, when auto-  
602 off happens again.  
603 The reason why we do this: It really is important to let the loop loop, and not interfere more  
604 than absolutely needed. Any bolus the user gives will sure distort the bg curve, on which  
605 autoISF, especially when aggressively tuned for FCL, builds a lot of its decisions!

606

607 5.3.3. Three top fields (%profile, exercise, TT)

608

609 Depending on the variedness of lifestyle, the desired %TIR, and the initial tuning effort put in, the  
610 user may want occasionally to „tweek“ the aggressiveness of her/his FCL.

611

612 The top 3 fields (grey in default mode, **yellow when temp. in mode with changed**  
613 **aggressiveness**) serve as quick and easy entry points to make temp. switches (as users will be  
614 used to for %profile switches, or for setting an EatingSoonTT in HCL, .. which they still can do in  
615 FCL ... but more:)

616

617 Expert FCL users might need this feature rarely, but probably at least to manage activity after  
618 meals: Each require opposite aggressiveness, and the switch has to come in a certain point in  
619 time that would be difficult to capture. (More see [section 6.4](#) )

620  
 621 5.3.3.1 TT dialogue field ( Currently not available in the pictured form and function ! )  
 622  
 623 The TT field (top right of AAPS home screen) is a primary daily interface, and a dialogue field  
 624 opens when pressing on it

Initially, the form can be entirely empty re. TT inputs and just show the default iobTH and bGAccel parameters.

Two lines appear for target and get labeled ES, AC or HY depending on what was pushed at the bottom buttons. (HY => red frame; evtl. 1. empty => grey)

TT entries can be made or overridden.

iobTH calculates from Target AC and other settings shifts away from default set in /preferences (here 0.7) but could be overridden here 0!..200% is allowed

bGAccel ISF\_weight can be modulated here, too. Note: it can change again if % profile is also changed. -- 0...200% is allowed

CANCEL allows to start fresh (select one or two of the square buttons, ES, AC or HY)

OK needed to use the settings

Duration input is made in minutes. In the exceptional case that both, ES and AC targets are defined, the duration input is pre-AC and framed blue. (This is because the preceding AC mode is automatically determined in length by the loop observing when iobTH is exceeded)

625  
 626  
 627 This looks complicated but only because it allows 4 different modes of use. Each user will primarily  
 628 use her/his preferred one.  
 629 (1) Who is happy with the initially well tuned FCL and does not have huge variations in daily eating  
 630 and moving around, will **not use** the TT **at all**. FCL is possible without an intervention via the

631 TT button in your cockpit. Actually 4 of 8 modes (GGG ...YYY permutations, list see [section](#)  
632 [5.4.1](#)) are not making use of TT.

633 (2) Super easy is also, to just input **any odd-TT** (odd-numbered temporary target) that will shut out  
634 any SMBs for the set duration. *That can be a good idea when having a small snack, for*  
635 *instance.*

636 Super quick access to stop SMBs is possible also via the loop icon ([section 5.3.1](#)).

637 Specifically, an **EatingSoon TT** can be activated here (limited relevance see [section 2.5](#)). It is  
638 time-un-critical, can be manually set, or come up via an Automation.

639 The cockpit enables you to set the iobTH differently (override) for the current meal.

640 Alternatively, iobTH can be temporarily changed in /preferences or using an Automation.  
641

642 Temp. iobTH will always revert to default when the TT expires. If another TT immediately  
643 follows, like in the example of the screen above, it will calculate, (then) show and use a new  
644 temp. iobTH.

645 (3) The third way is to **use the input mask** (if already ncluded in your software version see picture  
646 above) **to freely modulate the loop aggressiveness** for a declared number of minutes. Click  
647 the bottom big square(s): Either HYPO, or ACTIVITY, or EATING SOON, or ACTIVITY and  
648 EATING SOON (example in the pictured screen above). Make or override entries in the offered  
649 fields. Press OK.

650 (4) The fourth way is to exclusively use one of the 4+4+2 little buttons seen in the bottom part of  
651 the TT dialogue box (if already included in your software version). They provide a set of  
652 settings (as will immediately show in all input fields above) that the user has set up in  
653 Preferences/SMB/autoISF/FullLoop (refer to [section 6.3](#)), and can freely label there. *For*  
654 *instance „hiC“ at high carb EatingSoon, „piz“ for Pizza/fatty meals, „grd“ for garden work,*  
655 *„mtb“ for mountain biking ...*

656 **Capturing good settings for not-everyday situations in /preferences** (if already included)  
657 **allows calling them up within 1 second**, from your cockpit on the AAPS home screen (...and  
658 won't ruin the FCL experience at all , especially because in most cases it is not time-critical,  
659 how long before the intended exercise the buttons are pressed).

660 [Case study 6.2](#) demonstrates that nearly the same performance and comfort can be reached  
661 via the **DIY FCL cockpit** with the grey extra buttons appearing at the bottom of the AAPS home  
662 screen, based on Automations with User action (see also [section 5.2.2.3](#)).



663 The example picture given above, and also [case study 6.2](#), is the most complicated (but also most  
 664 useful) case, **when exercise follows after a sizeable meal**. It is then that you need (a) aggressive  
 665 FCL initial performance at the meal, but, *exactly when* (!) a (for the intended sport already  
 666 temp.lowered) *iobTH* is exceeded, you need (b) to have SMBs automatically switched off and go  
 667 into the „milder“ mode, as defined for the exercise (with *high* instead of the immediately prior  
 668 *low*TT, that automatically significantly reduces iobTH again, and insulin sensitivity(resistance)  
 669 settings too).

670  
 671 Pressing exercise related buttons will automatically also light the **exercise button** on the main  
 672 screen yellow.

673  
 674 To summarize, the TT dialogue field offers easy but powerful ad-hoc [modulation of loop](#)  
 675 [aggressiveness](#) for FCL (if already included).

676

677

678 5.3.3.2 Exercise button (see more in [section 6](#).)

679

680 The exercise button automatically lights yellow when exercise related TTs are activated [in the TT](#)  
 681 [dialogue box](#).


682 4 of 8 principal FCL modes ([section 5.4.1](#)) are making use of the exercise button.

683

684 If pressing on the exercise button, a dialogue box appears (*if extended design for FCL cockpit is*  
 685 *already launched*) with info on exercise setting first (and opportunity to override), plus below the  
 686 activity monitor (experimental for auto-tracking of lighter movement during the day, and effects on  
 687 sensitivity that may have. See [section 4.5](#)).

688

689 So, first the exercise settings (as set under TT) are there to read. Example :

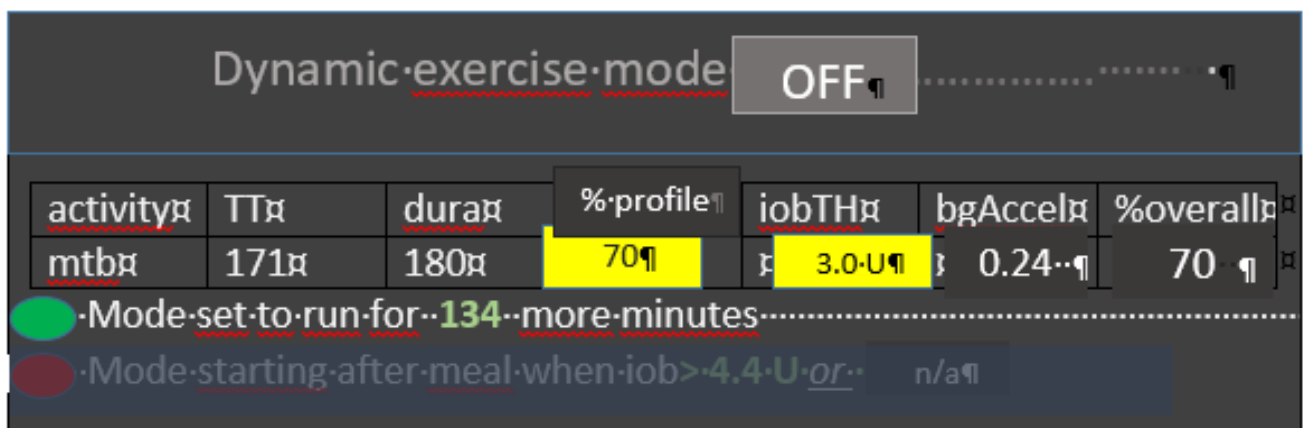

690

691

692 The exercise (here mtb) is selected in the dialogue box of the neighboring TT field, and there auto-  
693 filled with settings made in the set-up and tuning stage by the user under preferences (see  
694 above,...). They are reported also under the exercise button here, and TT, duration, and % sens  
695 (which also shows active on the %profile field on the left side of the exercise button) can be temp.  
696 changed there. iobTH, bgAccel\_ISF and overall resulting sensitivity ratio is given in the other fields.  
697 The **middle field** of the table, „% **profile**“ either picks up the % set under the %profile button, or  
698 an input can be made here, in the exercise button domain, which will:

- 699 • turn the neighboring %profile button on yellow and show that inputted % on it, too
- 700 • be multiplied with the result from the exercise mode settings per se, and change the %  
701 overall, accordingly.

702 So, if this middle field of above table (dialogue box of sports button) contains a figure other than  
703 100, input field becomes yellow, and you are operating with a combination of traditional PLUS new  
704 exercise mode (with all three top buttons of your FCL cockpit yellow). This maximally will soften  
705 aggressiveness, for which you get an idea by the last calculated figure.

706

707 The mode is either running already (for another number of minutes, as probably also shown in the  
708 yellow TT field anyways). Or it is scheduled to run, after insulination for a started meal reaches  
709 iobTH (as in table). Or, no exercise is scheduled (both points red, no entries).

710

711 The lower part of the exercise dialogue box (not pictured above, but see in [section 6.5](#)) is  
712 dedicated to the Activity Monitor

713

#### 714 5.3.3.3 Profile button

715 The profile button can still be used to set a different profile, or profile%, for instance to adjust for  
716 days with sickness (as you are used to from hybrid closed looping). 4 of 8 modes are not making  
717 use of the profile button.

718

719 Any inputs made here will be used to modify profile\_ISF on which all further changes are made on  
720 (multiplied with).

721

722 The profile field remains grey if standard profile is applied.

723 It turns yellow, displaying a %figure relating to any altered loop overall aggressiveness:

- 724 • When no inputs (changes from 100% profile) are made here, but inputs in the TT field,  
725 e.g. for exercise, automatically lead to different insulin sensitivity ratio ~~that ratio is shown~~  
726 ~~here~~
- 727 • when% is changed by input in the profile button itself, it will be multiplied with with  
728 profile\_ISF and be used in place of profile\_ISF *by the algorithm*.

729 However, for exercise (sports) you no longer must make an entry here, because  
730 reasonable %reductions should be automatically provided, driven by your set TT (and half-basal  
731 exercise target), see [section 6](#).

732

## 733 5.4 Recognizing your loop state in the AAPS home screen

734

### 735 5.4.1 Color scheme of top cockpit buttons tells kind of closed loop that is running

736

737 3 Buttons (%profile; exercise; TT) each in 2 states (yellow Y, or grey G) make  $2 \times 3 =$   
738 **eight principal FCL modes** possible:

739

740 GYY = dynamic exercise mode

741 YGY = not-dynamic „traditional“ exercise mode (if <100%) or hypo mode (if >100%)

742 GYG = basic closed loop with Activity Monitor running

743 GGG = basic closed loop (FCL or HCL) without any altered sensitivities etc

744 YGG = basic closed loop but with a „long wave“ sensitivity shift (e.g. sickness)

745 GGY = temp. target like e.g. EatingSDoonTT is set; or Hypo mode

746 YYG = closed loop with „long wave“ sensitivity adjustment and Activity Monitor running

747 YYY = dynamic exercise mode in time with additional „long-waved“ sensitivity shift

748

### 749 5.4.2 Information printed on the top buttons

750

751 The yellow TT field shows **the currently valid TT** (and further duration):

752 (profile) stands for the abbreviation you labeled your selected running profile



754 In the special case of settings for meal preceding sports, the field will look slightly differently:



...and ...

756 ... when iobTH is first time exceeded, this automatically switches to:



758

759 Likewise, if on the AAPS main screen just an **EatingSoonTT** is set (e.g.72), this is entered with the  
760 desired duration. Afterwards, it automatically reverts to profile target and the display turns grey  
761 again there with e.g. 90 on it (and no time limit).  
762 Without sports context, the middle field remains grey.



763

764

765 Independently from setting a TT, the user can choose to set a **%profile in the left top field**, for an  
766 independent number of minutes, e.g. 70% in this screen example: Also, or additionally, this will  
767 influence the resulting ISF and sensitivity%

768



769

770

771 The % might change and turn yellow also in context of making TT inputs in the related dialogue  
772 box (see chapter TT dialogue field, above). Still, the % (or the length of time the profile switch shall  
773 be active) can be independently overridden in the top left field, if so desired.

774

775 If an **Automation** sets a %profile, and/or a TT (e.g. *automatic detection of meal start at condition*  
776 *e.g. when delta >10*) , this would automatically show in respective field(s) turning yellow and  
777 showing the temp. setting. To show the set parameter comes from an Automation, „ **(A)** „ is added  
778 in the end of button text.

779 Note that an Automation is **usually/ always (?)** only permitted to temp. change default profile  
780 settings, not other pre-existing temp. settings. This is for a good reason : Why should a  
781 sometimes in the past thought-out Automation supersede your - just for the occasion  
782 specified – temp.settings that you consciously activated for the day?

783 Advice: Try to stay away from Automations that also aim at temp. modifying  
784 aggressiveness. For the reason just given in above note, they often will not kick in anyways.  
785 Generally, it also is no good idea to double up sub-algorithms for tweaking loop behaviours.

786

787 Also, as pointed to above twice already:

788

789 Try to **keep things as simple and clear as possible.**

790

791 That said, a limited number of Automations can be of help in distinct scenarios (that differ in  
792 purpose and in applicable time of day).

793 A good one could be for night time, when your odd profile TT has SMBs shut off, but your  
794 experience after pizza nights tells you that, under certain condition patterns (bg, iob), an  
795 SMB or two should be „allowed in“ (see example given in [section 5.1.2](#) ).

796 Another good example, if you go usually FCL without any use of the TT button (which would  
797 be a meal announcement of sorts), is to define an Automation that, after detecting a meal  
798 start, automatically sets a low TT to get maximally aggressive first SMBs (as is the author's  
799 preferred way, mentioned already in [section 2.5](#) ).

800

#### 801 5.4.3 FCL related indicator fields in the AAPS home screen

802

803 In extra data fields of the AAPS main screen you can always see (not change) the key  
804 „aggressiveness“ parameters your loop currently operates with (see also home screen  
805 example below):

806 • how profile sensitivity (**ISF**) changes by the %profile input, by autoISF, and/or a set  
807 exerciseTT.

808 • next to current available iob number is an indication of your **valid iobTH** (the iob above  
809 which no more SMBs will be given)

810 • The AAPS home screen additionally shows, above the deltas, the current **acceleration**

811 Having a look at that can be valueable. For instance, when glucose is relatively low and still  
812 falling, a positive (and getting more positive) acceleration indicates that bg will swing back  
813 up, rather than crash low. This will give info about necessary snack size, and hence help  
814 avoid both, unnecessary calories, and going on a bg roller coaster.

815  
816  
817 5.4.4 Overall home screen:  
818

Overall home screen:



819  
820  
821