

1. Pre-requisites for Full Closed Loop

V 2.2

Please note that with autoISF 3.0 you are in an early-dev. environment, where the user interface is **not optimized for safety** of users who stray away from intended ways to use. Good safety features exist, but these are only as good as the development-oriented user understands and implements them. This is not a medical product, refer to disclaimer in [section 0](#)



1.1 Well-tuned hybrid closed loop

It is advisable to first establish a well-tuned hybrid closed loop before considering the transition to FCL. There are two important reasons for this:

- The UAM full closed loop requires a highly personalized (individual) tuning of settings, so the loop will give insulin mimicking YOUR successful hybrid closed loop mode.
- The UAM full closed loop comes with new parameters to be set and tuned. It would be problematic to set and tune several new parameters before the basics were tuned “right”. Errors could easily be balanced with counter-errors. This can work in single scenarios, but would create a highly unstable system, hard to re-calibrate better later.

1.2 Fast insulin (Lyumjev, Fiasp, Apidra?)

If the user does not bolus for meals, clearly a very fast insulin is needed so, upon realization of a starting meal-related glucose rise, the loop has any chance to eventually keep glucose in range (by common definition, under 180 mg/dl (10 mmol/l))

A modelling study (details see

<https://androidaps.readthedocs.io/en/latest/Usage/FullClosedLoop.html#fast-insulin-lyumjev-fiasp>)

can show in quantitative terms that **faster insulins**

- will result in significantly **lower** glucose **peaks** than slower insulins
- **tolerate** a couple of minutes **delayed** first meal bolus while not incurring unacceptable height of peaks
- minimize the effect on glucose peak from **different** carb load (**meal sizes**).

In conclusion, do not attempt FCL with other insulin than Lyumjev® or Fiasp®, unless, maybe, if you are on a very moderate to low carb diet. (According to [case study 1.2](#), Apidra® might work, too, but Humalog® would not work well).

1.3 Reliable insulin delivery from the used pump/cannula/insulin system

Good Tolerance of Lyumjev (or Fiasp): Occlusions threaten the function of the full closed loop.

It is very important to have an eye on the time a **cannula (or pod)** is in use (many find **48 hrs** to be the **limit**), and whether hard-to-explain glucose rises happen at ever increasing „fake“ iob (even before a 48 hr routine replacement). (See [case study 1.1](#): You easily lose 25% TIR that day)

It is absolutely contra-indicated to attempt FCL coming from leaking pods and associated erratic sensitivity swings that may or may not have been somewhat controlled and tolerable by dynamicISF or other measures when you were Hybrid Closed Looping,


Stable pump connection

In FCL you absolutely rely on your pump delivering, without any further delay, the much needed insulin, after any meal start.

Hence it is absolutely essential not to “miss” any problems from a lost Bluetooth connection.

An Automation similar to the one pictured here → should help recognizing eventual problems.

See also [case study 1.4](#)

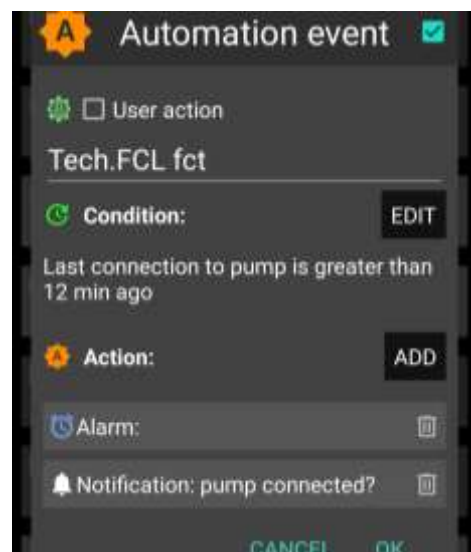

1.4 Excellent CGM

You do not give a meal size-related bolus any longer. That leaves all insulination jobs to the algorithm! As glucose values are the very basis for this, please **inform yourself well about** how **your CGM** 1) principally performs 2) whether and how this may depend on data flow and intermediate apps you use 3) specifically, how and where any smoothing is done, and what this might imply for the ISF boosting method you will be using See for instance here:

<https://androidaps.readthedocs.io/en/latest/Usage/Smoothing-Blood-Glucose-Data-in-xDrip.html>

In AAPS Preferences/OpenAPS SMB/autoISF/Glucose source data for parabola fit, you must select between 4 options relating to your CGM (1 or 5 minute values, raw or smoothened).

Around meals, a stable Bluetooth connectivity is absolutely essential, too, so CGM, loop, and pump can do their job without losing more valuable time (see [case study 1.4](#)).

Then, but even more importantly in *all other* day and night *times*, the CGM should not produce any artefacts (jumpy values; see [case study 1.3](#)) that the loop could **misinterpret** as sign of a starting meal. Note that also calibrations could produce jumps.

76 But, in any case, a CGM with more scatter will make the loop lose more time, and lead to higher
77 peaks and lower %TIR.

78

79 The best way currently is to use Dexcom G5 or **G6**, and to ensure via **overlapping** right and left
80 arm sensor and transmitter utilization, that always good quality values can be used by the loop.
81 Other ways are possible, but come with a lot of monitoring effort (via watch) and occasional time-
82 outs for the loop.

83

84 One safety feature in autoISF is a **blockage of SMB delivery whenever delta bg** (within the last
85 two 5 minute values) is **higher than 30% of that bg**. *So from 74 mg/dl, a jump to 97 (+23) or more
86 would not receive SMB "response", or from 100 mg/dl to 131 mg/dl (+31) would neither.*

87

88 Check in your (HCL or FCL) data whether at meals or sweet drinks with rapid absorbing
89 carbs you could run into the problem that jumps are "too high" and much needed insulin will
90 be blocked (only come via very much smaller portions: *e.g. 400%TBR @ 0.6 U/h => 0.2 U in
91 5 minutes, instead of one ~3 U SMB. The difference of 2.8 U missed translates @ ISF~ 40
92 mg/dl/U into up to + 112 mg/dl higher bg peak! It will not become quite that bad, because
93 the loop will catch up to the insulinRequired with it's next couple of decisions*).

94

95 Instead searching in old data, you can also just have an eye on instances where you think a
96 first SMB was due, but blocked. Confirm that by looking in the SMB tab, and think about a
97 solution that would not require changing the 30% safety limit in the code. *For instance, not
98 drinking so much juice rapidly around meal start could be a likely "behavioral" correction to
99 get rid of the problem.*

100

101 This blockage (no SMBs) would likely last only 5 minutes (and go probably unnoticed - **in a
102 suggested improved User Interface it would be indicated for about 5 minutes via a dotted violet full
103 closed loop logo/button on the AAPS home screen**). However, not only would you lose 5 valueable
104 minutes to get your iob substantially elevated; all following deltas are likely much smaller, and, as a
105 consequence, you will miss some of the boost sought from bgAccel_ISF if the >30% delta was in
106 fact (largely) due to carb absorption.

107

108 This example also underscores that the CGM in use cannot be allowed random scatter that leaves
109 no reasonable room for safe detection of (smaller and) bigger "truly carb related" deltas

110

111 Besides Dexcom, also FS **Libre 3** is useable. Observe info in the general section about autoISF
112 (<https://github.com/ga-zelle/autoISF>) regarding how it (or maybe by then a couple of others, too)
113 should be implemented.

Hints for users with 1-minute CGM (Libre 3) *(status of 05 Dec. 2023)*

In AAPS 3.2.0.2 the glucose history only uses data at 5 minute intervals and ignores the ones in between. That has awkward consequences for the graph display and for the exponential smoothing method: Each minute a new reading is received the 5-minute pattern moves forward by 1 minute and a completely fresh and new subset of glucose values and times is used.

Alternatives of determining the glucose acceleration based on the 1 minute data are still under evaluation. Once a promising method can be found and validated then an interim release will be provided.

As a libre user you have 3 alternatives:

- AAPS gets values from xDrip+ where you smooth and reduce the 1-minute data to 5-minute data before sending to AAPS
- you should stick to the 5 minute mode in your Libre setup if possible
- AAPS gets values from Juggluco every minute but AAPS uses only the 5-minute subset for looping and for parabola fitting

1.5 Meal-related limitations

Setting up a full closed loop is relatively easy for people whose diet does not consist **mainly** of components with rapid high effect on blood glucose (more see <https://androidaps.readthedocs.io/en/latest/Usage/FullClosedLoop.html#meal-related-limitations>)

Meals do not have to be low on carb (provided you use a fast insulin for your FCL)

Fat or protein rich diets, or slow digestion/gastroparesis, make things easier rather than harder for the full closed loop because late carbs nicely cover for inevitable “tails” of late action from SMBs needed around peak time.

Erratic consumption of snacks with fast resorbing carbs can be a problem.

In autoISF you can reduce this problem to some extent via one or two keystrokes from your AAPS home screen. While certainly being a deviation from the FCL idea(I), this would be one of the exceptional situations where you better do a quick “nudging” step from your “FCL cockpit”. Details see in [section 5.2.1](#) and [5.3.3.1 \(4\)](#)

152 1.6 Lifestyle-related limitations

153
154 Technically stable system
155
156 Full closed looping requires a 24/7 technically stable system, especially regarding

- 157 • reliable CGM signals
- 158 • Bluetooth stability with the pump (see [case study 1.4](#))
- 159 • keeping your phone in sufficient proximity at all times
- 160 • avoiding (or at least early recognition of) occlusion.

161 This requires a habit (or, unlikely, permanent attention to details) like keeping all components well
162 charged and in close proximity; making cannula (or pod) changes always early enough to lower the
163 risk of occlusion (see [case study 1.1](#)); having always potentially needed parts with you.

164 **Depending on your system, your experience with it, but also on your acceptance and general**
165 **lifestyle, these aspects may or may not limit you.**

166 167 Preparing for exercise

168
169 To prepare for exercise (sports, heavy work), the normal protocol with a pump or hybrid closed loop
170 is to take actions that reduce insulin on board prior to exercise

171 With your full closed loop, the algorithm is tuned to detect meals and to give you insulin to counter
172 glucose rises automatically. Setting a high temp. target and lower %profile right away (effective al-
173 ready around meal start) could be a problem.

174 Unusual activity levels therefore likely require **disciplined preparation** (especially **if you want to**
175 **keep the need to snack during sports low**)

176 In autoISF you can reduce this problem to some extent via two or three keystrokes on your
177 AAPS home screen. While certainly being a deviation from the FCL idea(l), this would be
178 one of the exceptional situations where you better “flick a lever” from your “FCL cockpit” to
179 have temporarily adjusted settings for the planned exercise. Details see in [section 6.2-6.3](#)

180

181 Extra hurdles to establish FCL for kids

182
183 To establish and maintain a FCL for kids brings about some extra challenges if:

- 184 • Lyumjev is not available or well tolerated
- 185 • Hourly basal rate is very low, providing a poor basis for big SMBs
- 186 • Diet is rich in sweet components. With the typical low blood volume of a small body, strong
187 tendency towards very high bg spikes!
- 188 • Going through marked changes of insulin sensitivity or of circadian pattern makes it difficult
189 to keep the FCL appropriately tuned.

190 This problem is about the same in Hybrid Closed Looping. However, now you might
191 expect miracles from the FCL. This is not going to happen. You still should try to set
192 appropriate temp. changed profiles, that serve also as a basis for your autoISF FCL.

- 193 • Discipline is poor regarding keeping Bluetooth connectivity and infusion sites perfectly run-
194 ning
- 195 • Between kid and supervising parent it must be guaranteed, especially in the initial weeks,
196 that an eye is kept on whether the FCL is working about as to be expected.

197 More see [section 7](#).

198 199 200 1.7 Time required for setting-up

201
202 Lastly, before enjoying a functioning full closed loop you need to have a period of a some weeks
203 with some free time and „free head“ for set-up –. Can you get, in the time you are willing to invest,
204 to a result that you consider good-enough is really the question. Depending on your „habits“, and
205 which – if any - compromises (like doing cannula/pod changes more often, never starting meals
206 when bg sits high ...) are you willing to make (and everyday able to stick to), for the ease of not
207 having to deal with assessing meals and bolussing for them?

208
209 While setting up your personal FCL using autoISF is a substantial project, there is no need
210 to implement it fully in one step. There is nothing wrong to go in your well running Hybrid
211 Closed Loop mostly, while switching to FCL only for dinners, for instance, or only for
212 weekend lunches, as a start. Once you found feasible settings, you can expand to other
213 meal times, and lastly towards figuring out your best strategies for challenges outside of
214 meal windows, as we shall discuss in [sections 5. and 6](#).

215
216
217

218 There are alternatives to using autoISF for FCL, as well. See [section 13](#). for more info.
219
220 Notably [13.1](#) FCL using AAPS Master and Automations could be a much easier and more
221 error-tolerant way of stepping into FCL. In a clinical study with 16 participants about 80%
222 TIR was achieved without much tuning effort (source: see in [section 13.1](#)).
223
224 To close the circle to where we had started ([section 1.1](#)): A very time consuming pre-requisite might
225 actually be to first sort out your Hybrid Closed Loop, so your profile parameters are set „right“, and
226 your “old” data really can serve as a blueprint for what, now, you would like *your loop* to do in FCL
227 mode
228 Note that if you had used dynamic parameters or special Automations („loops inside the loop“) this
229 might have balanced some principal errors, but leaves you now without a good starting point as
230 you must get rid of these over-patches (see also warnings at start of [section 4](#))..
231
232 You will see also success stories of loopers who just jump into using more powerful tools, in
233 kind of a trial and error mode, and frequently add the latest add-on, or self-constructed
234 patch (often in form of an Automation) to counter-balance problems.
235 So, yes, you can also continue in that spirit. Resulting solutions may be good-enough. But
236 they tend to be unstable and not well-understood. That is a poor basis for managing arising
237 problems (-> fine tuning), and for adjusting to special situations (-> which setting to
238 temporarily change). But it certainly is an alternative avenue for the impatient, less
239 analytically, and more adventurous inclined.
240 In any case, PLEASE always observe the safety settings/instructions coming with the DIY
241 dev- variant of software you select.
242
243 One key safety measure every AAPS user going towards FCL should have in place is to set an **iob**
244 **threshold** (iobTH; size a bit below what you used as a bolus for bigger meals in HCL) above which
245 no more SMBs can be given by your FCL.
246 This is an integrated feature of autoISF, from 3.0 version onwards (see [section 2.4](#)).