

5. Modulation of autoISF aggressiveness.

V 3.9

Please note that with autoISF you are in an early-dev. environment, where the user interface is **not optimized for safety** of users who stray away from intended ways to use. Good safety features exist, but these are only as good as the development-oriented user understands and implements them. This is not a medical product. Refer to disclaimer in [section 0/readme.md](#)



5.1 Automatic modulation of loop aggressiveness

- 5.1.1 “autoISF off” outside of meal times
- 5.1.2 SMB off @ odd profile target
- 5.1.3 SMB off @ odd temp. target
- 5.1.4 diff. of FCL aggressiveness via Automations
- 5.1.5 diff. of FCL aggressiveness via Activity Monitor
- 5.1.6 Pro/con completely hands-off FCL

5.2 Manual modulation of FCL aggressiveness (DIY cockpit)

- 5.2.1 Status recognition
- 5.2.2 Manual interventions from DIY cockpit
 - 5.2.2.1 Temp. %profile or TT settings
 - 5.2.2.2 Temp. settings in /preferences
 - 5.2.2.3 Grey DIY cockpit buttons for FCL responses
- 5.2.3 Temporary exit from FCL

5.3 Manual modulation via improved cockpit

- 5.3.1 Violet FCL icon and underlying buttons
- 5.3.2 Bottom buttons “insulin” etc.
- 5.3.3 Top three fields
 - 5.3.3.1 TT dialogue field
 - 5.3.3.2 Exercise button / dialogue field
 - 5.3.3.3 Profile dialogue field

5.4 Recognizing loop state from the AAPS home screen

- 5.4.1 Modulated loop aggressiveness using the 3 top buttons
- 5.4.2 Color scheme of the top 3 fields
- 5.4.3 Info on the top 3 fields (profile, exercise, TT)
- 5.4.4 FCL related indicator fields
- 5.4.5 Overall AAPS home screen
- 5.4.6 Info given every 5 minutes in the SMB tab
- 5.4.7 SMB tab info when operating 1-minute/Libre3

[Available \(related\) case studies:](#)

Case study 5.2: Sweet snack.

[Skip what is in green writing:](#)

= Drafted fragments or not implemented ideas.
Please contribute, or wait for update with the missing info

Once the initial tuning according to [section 4.](#) is done, you are ready to use autoISF for your fully automated meal management.

You will have up to four major *other* challenges to manage:

- recognize and manage (partial) occlusions, or other technical (CGM or BT related) obstacles (see [section 1](#) on pre-requisites of FCL, and related case studies).
- deal with times when insulin given by the loop must be restricted (e.g. a snack should not be “treated” as a meal)
- deal with times when the loop should be set “milder” as a precaution (e.g. nights; or in an exercise context)

48 • deal with times when the loop should be set "more aggressive" (e.g. tempting hotel
49 breakfast buffet)

50 How big the remaining challenge really is, depends very much on your individual lifestyle.

51 • This [section 5](#) discusses this in more detail

52 • [Section 6](#) will extend this discussion regarding how to deal with exercise.

53

54 **In order to run the FCL around the clock** (preferably fully automatically, which can be possible, see
55 [case study 4.3](#)), **also the times *outside* the meal blocks must be precisely analyzed, and**
56 **solutions to problems (if any) must be sought.**

57

58 It is up to every user to decide where to draw the line:

59

60 • With a technically well-functioning system, moderate meals, moderate or no exercise,
61 moderate %TIR expectations and a bit of mindfulness it should be possible to go into **Full**
62 **Closed Loop 24/7**, after working through, and observing, [sections 1-4](#). See [case study 4.3](#).

63 • Especially if you are a bit shy of using the Emulator ([section 10](#) and [11](#)) for really detailed
64 analysis, it can well be that you will not hit *one* real good system calibration ([section 4](#)) for
65 your *entire range* of diets.

66 In that case you will occasionally run out of range (bg =70...180 mg/dl), and your options to
67 prevent, react, or improve are:

68 ○ accepting a few % higher time outside of range for that day (and, if feasible, in the
69 future avoiding what seemed to have caused it)

70 ○ taking a snack (whenever you tend to go low from the "tails" of insulin activity that
71 was required to fight a peak)

72 ○ doing a manual "tweak" (if you can think of one in time), to manage the problem
73 manually. For example, briefly going into an odd TT (=temp. blocking more SMBs)
74 can be a very easy-to-handle remedy sometimes.

75 ○ temporarily resorting to "your old" hybrid closed loop.

76

77 Instead of accepting such instances, you could launch "improvement projects"

78

79 • that refine your initial tuning ([section 4](#). and [sections 8](#) and [9](#))

Note, though, that it could be near-impossible to fine-tune *if your basics never were “right”* and you got lost in a maze of errors and counter-errors. Then only a fresh start might convincingly help.

- that make you and your FCL loop fit **to manage an increasing number of disturbances either automatically, or via an “informed”, maybe pre-programmed, user intervention** (notably, an exercise “announcement”) ([sections 5 and 6](#)).

To tailor the loop’s response to disturbances *other-than* your usual major meals probably will require specific **modulation of the aggressiveness** (which you have set according to [section 4](#) for your *usual meal* spectrum).

There are many **avenues to achieve this**.

The main ones, that are also easy accessible via Automations in AAPS, are:

- temporary shut-off SMBs (odd-numbered target)
- temporary change bgAccel_ISF-weight
- temporary change iobTH_percent
- temporary change the set %profile
- temporary set different bg target (especially in connection with exercise mode)t

After set up of your core FCL for fully automatic meal management according to [section 4](#), you now can progress to define solutions for any of your „other“ situations (outside of meal management) that tend to drive glucose outside of the desirable range.

- In [section 5.1](#) we explore avenues towards **fully automated** management that in daily life will require no user intervention at all.
- In [section 5.2](#) and [5.3](#) we will look at solutions that involve an easy **user interaction** like a *data entry or button push*.

116 5.1 Fully automatic modulation of FCL aggressiveness

117

118 The following subchapters describe set-ups you may want to use for allowing **completely hands-**
119 **off FCL** in as many daily situations as possible (and potentially all the time, as in [case study 4.3](#))

120

121 5.1.1 All autoISF ISF adaptations switched off outside of meal-time windows

122

123 If, aside from having to bolus for meals, your hybrid closed loop was running pretty well *without*
124 other interventions from your side, you could continue to run in that mode, and just focus your new
125 autoISF FCL on management of meals.

126

127 In your initial transitioning phase, this approach makes a lot of sense, and even by focusing
128 autoISF on just a sub-set of meals, like only dinners.

129

130 Also in the long run this avenue is taken by many FCL users for the night times, “hanging
131 on” to their well performing hybrid closed loop with standard oref(1) SMB+UAM

132

133 For this, you define Automations

134

135 • that set meal time windows in which “Enable ISF adaptation by glucose behavior” (autoISF) is
136 turned on in AAPS preferences/OpenAPS SMB

137 • or: that turn *all* autoISF's ISF modulations (*or just bgAccel_ISF*) off in time windows in which
138 surely no meal occurs.

139 For instance, you can go *for all nights* back into your Hybrid Closed Loop, as you had before.

140

141 Your temp. “autoISF shut-down” (exiting autoISF FCL = shutting off “Enable ISF adaptation by
142 glucose behavior” in AAPS Preferences / SMB) is meant to prevent problems from the loop over-
143 reacting to bumps in the glucose curve in times of day (night), when standard oref(1) performance
144 is sufficient.

145

146 A very good alternative to defining “meal-time windows”, or to resorting to night-time Hybrid Closed
147 Loop, is letting the autoISF FCL run 24/7, and “taming” the FCL via a night time SMB shut-off (see
148 next [section 5.1.2](#)).

149

150 *Off topic: Other “early dev” AAPS variants* (see [section 13.3](#)), do require working with meal-time
151 windows. The window is either set by time of day in the settings, or it always must be „set“ by the user.
152 Trigger to set a meal time window could be a pre-bolus given by the user, a carb entry made, an

EatingSoonTT set, or a meal-announcement-button pushed (none of these things are required in autoISF FCL).
Outside of these time windows, these loops then run with less aggressive SMBs, just like oref(1) SMB+UAM in AAPS Master, with or without (modified versions of) dynamicISF. So, outside of the meal time windows, you would be in “your normal” hybrid closed loop.
The term **Meal Announcement** (MA) is often used to label this closed looping mode.
It is not really FCL, but an advance over traditional HCL.

5.1.2 Odd-numbered profile targets, to block SMBs

An alternative route of preventing the FCL loop from over-reacting to bumps in the glucose curve would be to make use of the option to temporarily shut down SMBs

Ensure the even/odd logic in the settings is toggled on in Preferences> openAPS SMB> autoISF settings> smb delivery settings>: "Enable alternative activation of SMB depending on bg target": ON.

In time blocks with an odd-numbered profile target you can prevent any SMBs being given by your loop. The (unchanged) aggressive settings then can only translate within the limits set by %TBR possible.

This will very much slow down any more insulin being given, and is an excellent solution for night times, especially if you occasionally experience compression lows.

Alternatively, you could use the new included options for Automation Conditions and temporarily tune your bgAccel_ISF_weight much lower (see [section 5.1.4](#)).

The same situation can be achieved if you generally operate with a mild bgAccel_ISF, and make your autoISF only really aggressive for meal-time slots (if you have similar enough times every day, or also can “employ” geo-fencing in your Automation (or middleware, in iAPS) conditions).

In these cases you would not need to have night profiles that disable SMBs: - Which is the better way would depend on a lot of personal factors relating to how high-carb the diet is, regularity of meals, snacking habit, CGM quality and incidence of compression lows, and probably more. - I would try both routes, or, as this is fairly complex to tune, just one, and stick with what is working good enough.

Yet another alternative was already presented ([section 5.1.1](#)) = to go into hybrid closed loop for the night.



That is possible, with SMBs available (without them getting boosted via autoISF), and, for a long time, was the author’s favored solution for the nights.

This solution is evidently similar to the prior discussed one, of having a mildly tuned autoISF 24/7, boosted to high aggressiveness only in meal-time slots.

Adjunct Automations to allow a few SMBs, in nights with odd profile target SMB shut-out

My current **favorite** builds on the **method** of this section (5.1.2, odd profile target provides SMB shut out), but then allowing some, automatically triggered SMBs, when needed:

In case you occasionally do have nights that would benefit from a couple of SMBs (to treat temp. highs from a late fatty pizza, raclette and such): Define suitable Automations like the two „night“ ones in this *example*:

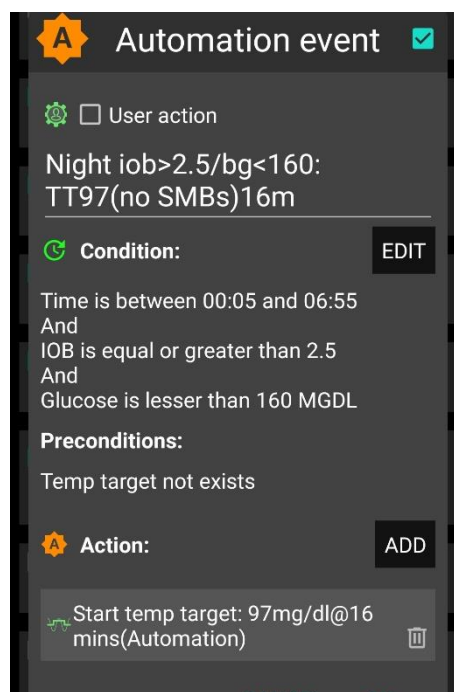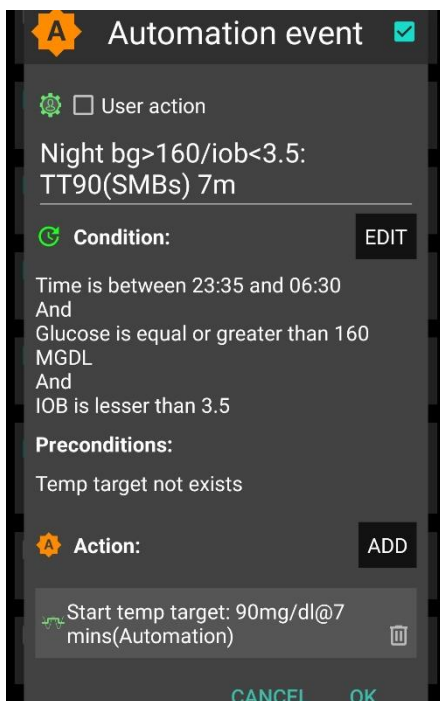

Caution: Never underestimate the „trickiness“ of getting your Automations „right“.

With your thought-out Automations in place, night data need to be analyzed to see

- whether the bg and iob limits, as defined in the given example, work sensibly four your data pattern
- whether the TT duration is chosen appropriately
- how swapping the sequence in which the automations appear in the Automation list would lead to different SMB impacts.

217 5.1.3 Odd-numbered temp. targets (TT) set via Automation, to block SMBs

218

219 A widely used Action that strongly modifies how fast your FCL can add more iob is setting an **odd-**
220 numbered **temp. glucose target** which makes the loop operate without giving any SMBs (%TBR
221 modulation only).

222 Ensure the even/odd logic in the settings is toggled on in Preferences> openAPS SMB>

223 autoISF settings> smb delivery settings>: "Enable alternative activation of SMB depending
224 on bg target": ON.

225

226 Then, from patterns you find in *your* data, at times where you want your loop act differently, you
227 need to carve out Conditions that describe the respective situations (and either *for how long* it
228 typically lasts, or at which *other* Conditions you want your loop get back to default FCL operation).

229

230 An odd TT is often set for an *anti-hypo* snack or *sports* snack. In both instances, you do not want
231 SMBs to quickly counter-act.

232

233 In case of sweet "*fun*" snacks, this is entirely different -> [section, 5.2.1](#) or for regular snacks
234 (e.g. *at school break*) see next [section 5.1.4](#)

235

236 5.1.4 Automatic differentiation of FCL aggressiveness using Automations (or middleware)

237

238 **Personalized Automations** tailor the loop exactly to **your** data so **fully automated** handling of
239 situations with **different aggressiveness** of the loop can be made.

240

241 Automations are an integrated and very easy-to-use feature in AAPS.

242 (The i-Phone platforms Trio or iAPS lack this feature. However, so-called **middleware** has been
243 developed as add-in to your code, see: <https://github.com/macconnellk/RoboSurfer/tree/main>)

244 .

245 From, autoISF 3.0 onwards, also the following parameters are provided as Condition and/or as
246 Action for defining YOUR Automations:

- 247 • Enable ISF adaptations by glucose behavior => Allows temp. ON/OFF for the key ISF
248 modulation parts of autoISF (and, as a result, will usually decrease loop aggressiveness)
- 249 • Trigger/set iobTH percent => Keeps default aggressiveness, but only until a iob threshold
250 (that your Automation modifies) is surpassed (which is when any further SMBs will be
251 blocked blocked)
- 252 • Trigger/set bgAccel_ISF_weight => Modifies the aggressiveness of just the acceleration
253 component

254 To set up suitable Automations, you first must **analyze patterns** you find **in your data**, at times (or
255 geo-location, or bg and iob patterns that point to a problem ...) **where you want your loop act**
256 **differently**, to carve out Conditions that describe the respective situations (and either for how long
257 it typically lasts, or at which *other* Conditions you want your loop get back to default FCL
258 operation).

259

260 A variant of this mode is to define several windows in which autoISF aggressiveness
261 (bgAccel_ISF_weight) and/or iobTH are automatically set differently

262 • for **different meal time slots** of your day –

263 *(Breakfast at home, school lunches, school intermission snacks, dinners at home could for*
264 *example all deserve special settings regarding ISF_weights and iobTH).*

265 Note: Circadian differences in insulin sensitivity between meal times are included via your
266 ISF profile and should not be a reason for different _weights needed between meals!

267 • or even for a geo-location etc –

268 *(School lunches, or mother-in-law visits, would be examples).*

269 An example for this was given in section 3 already:

Here is an example set of automations to alternate between two values of iobTH:

I use two different values of *iob_threshold_percent* during a normal day. It is 40% for lunch time and 60% for dinner time. I have these two rules to switch by time of day and only if the current value equals the value from the earlier shift. Any other value is treated as a manual override for special occasions until I manually set it to its regular value. The time windows for switching are long enough to catch an opportunity to be processed and do not need to be actioned half a day each.


270

271 Unless your meals differ vastly in size and in fast carb content **all this may not be needed**.

272

273 Still, personalized Automations might help ease your initial job of setting the various ISF_weights,
274 and a best-suitable iob_threshold_percent that would work “always”.


275 Note: Only the main two parameters (bgAccel_ISF_weight for “initial aggressiveness”, and
276 iobTH_percent for “where SMBs stop”) are available in Automations. So, finding your
277 parameter sets *for each of* the time slots, will not be trivial. => **Spending more effort to set**
278 **the .._weights so they accommodate just one, broader spectrum** ([section 4.](#)) **should be**
279 **the first, and standard, approach.**

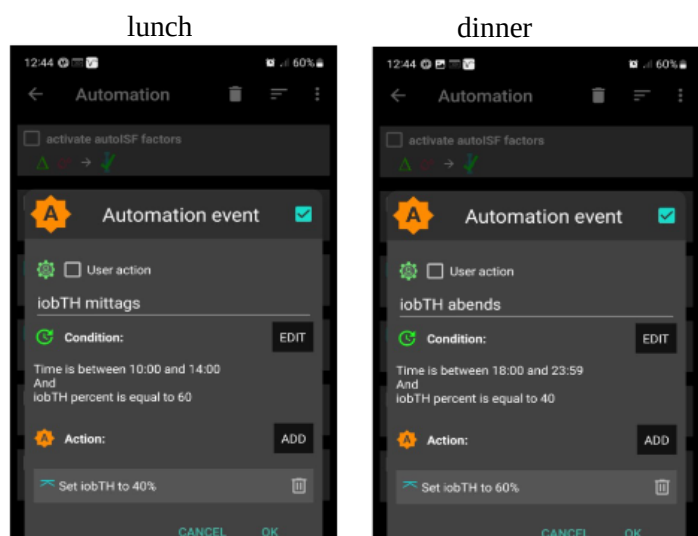

An intermediate (maybe only temporary) approach could be to use a profile switch (for low carb meal, or eating half, setting %profile to 60% for instance, and only for the brief, less than an hour, initial meal period). See [section 5.2.2](#) and [5.4.5](#)

Caution: Setting a different iobTH% or bgAccel_ISF_weight can probably not be done with a duration attached. Then you **must** define a suitable **additional Automation** that must be active in tandem, to **restore the values you had set in /Preferences for your iobTH% or bgAccel-ISF_weight**. Else, once your Automation set in, it will *forever* shift these important parameter settings!

If for instance you have several Automations that, in combination with a set elevated TT also set a lower iobTH: Don't be fooled, the duration only applies to the TT. You need an extra Automation for all of them. I picked out the highest of the altered iobTH values that these Automations can set (45_percent), and then I can automatically restore my default desired 60% via this one Automation (see screenshot - - >)

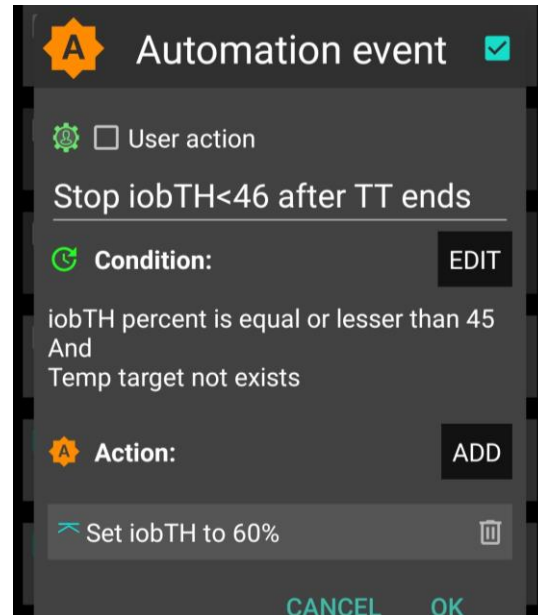

5.1.5 Different FCL aggressiveness set by the Activity Monitor

With the autoISF variant of AAPS you can make use of your smartphone's **stepcounter** and use it to fully automatically adjust insulin sensitivity ratio to **activity level in the past minutes to one hour** time frame.

This feature comes with yet another little tuning opportunity, in which you study your body's response to light exercise (like walking) or to not moving at all (like desk, couch), and select appropriate settings (in AAPS Preferences/OpenAPS SMB/Activity modifies sensitivity/ -> set **two scaling factors**).

This **will automatically adjust insulin delivery** (basal, ISF, and iobTH; see 1st screen of AAPS SMB tab (example in [section 5.4.5](#))), to suit activity state of the past minutes (up to 1 hour).

This autoISF feature (new since V.3.0) is much quicker responding than Autosens or dynamicISF to adjust insulin sensitivity to your current „lifestyle state“.

For loopers who do not have huge variations in exercise levels in their everyday lives, this feature might be a superior replacement for using Autosens (and also for dynamicISF, which, however, is

anyways contra-indicated in autoISF), and fairly much close the gap towards being able to do a 24/7 hands-off FCL.

autoISF Quick guide p.9 ([Section 3.5](#)) and [section 6.6](#) describe the Activity monitor in more detail.

While the Activity Monitor takes automatic care of light deviations from an average activity level, exercise enthusiasts, or heavy workers, should make use of the “heavier tools” (which automatically supercede (shut off) Activity Monitor; discussed in [section 6](#).)

5.1.6 Pro/con completely hands-off Full Closed Loop

To stay **24/7** in a **completely „hands-off“ FCL** can be a realistic goal with autoISF, if besides meals also some special challenges, as discussed in this [section 5.1](#), were analyzed and could be addressed. Example see in [Case study 5.3](#).

Clearly it depends very much on your lifestyle, and how interested, willing, and capable you are to recognize, deal with, (and in the future avoid?) situations that get you outside of your desired %TIR on occasion.

So, this is also about what %TIR you are aiming at, and can accept, as it averages out for the week, for instance.

Everybody must weigh for her/himself

- how much **upfront effort** to put into the setting up process, for getting it all 100% automatic
- **or whether to take an easier start, with a couple of situations left to take care of when and as they arise in daily life**

Even if a principal capability for a fully automatic “hands-off” running FCL is given, this still means that

- the user should be knowledgeable about what exactly is going on, and
- have a principal capability to „nudge“, or even to completely take over, in a manual mode.

In the sections that immediately follow, we present the options to nudge or temporarily take over from the AAPS home screen which will be serving as your **FCL cockpit**:

- [Section 5.2](#) describes how you can use available “buttons” from your AAPS home screen, and how to complete it towards a suitable DIY FCL cockpit, for an even better FCL experience.

- [Section 5.3](#) describes how you might be able to manage “disturbances” even better (with more convenience) with an improved FCL cockpit in the future. (No need to read any of the green lines, unless you are interested in contributing to define/design/program further improvements)

-

5.2 Modulating aggressiveness manually, from the DIY-FCL-Cockpit*

* Like in the airplane cockpit: Cruising in full auto mode should involve having an eye on the instruments, and on potential disturbances ahead in the environment.

In [section 4](#), we dealt with major meals.

In [section 5.1](#) we looked into ways to automate also a modified loop response to *foreseeable* situations (tied to a time of day, geo-location etc), or to those *the loop could recognize* (with enough time to react).

However: *Other disturbances* might come up, that:

- are not noticeable in-time, or foreseeable, by the loop (*e.g. your plan to start exercise in an hour or two*), but **that influence sensitivity dramatically**, and therefore require temporary modified settings in order to remain in-range, and/or
- **require** a different “starting point” regarding iob and bg, which translates into **a different iobTH** that should **temporarily** be set much lower (*in case of exercise*) or noticeably higher (*e.g. with very fast absorbing carbs in a sweet snack “sin”*) .

In these scenarios, you must find an easy way to

- call up a pre-programmed routine for automatic management, with adjusted aggressiveness, or:
- manually tweak a setting or two, to temporarily adjust the aggressiveness
- There may also arise a desire to just exit the FCL mode, and “be your own captain” for mastering a special situation.

Lastly, for peace of mind, to learn, and to stay informed (especially so in your initial tuning phase, or when your glucose curve goes in unexpected ways) we also must be able to...

- ...find the key parameters that frame and drive the recent and upcoming loop decisions.

All this is facilitated within seconds right from the AAPS home screen, serving as a **FCL cockpit** after you built a couple of DIY cockpit features via Automations (as described below and in [case studies 5.2](#) and [6.2](#)):

Thoughts went also into how to improve the cockpit in future releases, see [section 5.3](#)

5.2.1 Status recognition

Before considering any manual interventions into the ongoing FCL, you should be aware what the current mode of action is (refer to [section 5.4](#)), and hence how you might be able to “nudge” your loop in order to adjust to the disturbance that you see coming up.

5.2.2 Manual interventions from the (DIY-) FCL cockpit

Trouble with most of these is, not to forget to set back manually, too (=> better solutions in 5.3)

5.2.2.1 Temporary tuning of FCL aggressiveness via temp. %profile or TT settings

The set **% profile** multiplies with both, the ISF resulting from autoISF, and also with the default iobTH you have set, so both are nicely modulated in a linear way with the % temporarily chosen

Just taking profile e.g. to 110% for an afternoon might be an easy way to explore whether you might benefit from 10% more “aggressiveness” in your core settings for lunches (like bgAccel_ISF_weight).

Make sure, though, that the extra 10% are not cut away by set safety limits.

A lowered (relative to profile glucose target) temporary **bg target (TT)** signals lowered sensitivity (more insulin need), and

An elevated TT (as often used with exercise) increases sensitivity and hence works in the direction of a lowered % profile to also reduce insulin given by the loop.

Moreover, the **exercise button** (top center on your AAPS home screen) can be activated (turns yellow, then). This will **further boost** how your set TT elevates the resulting ISF, and sharply lowers iobTH, as often desired for sports. (See [section 6.1](#)).

425 5.2.2.2 Making temporary changes in settings made in AAPS/preferences/Open APS SMB

426

427 Going into AAPS/preferences/Open APS SMB allows to:

- 428 - set milder or stronger ..._ISF_weights
- 429 - set different iob_threshold_percent (or iobMAX)
- 430 - elevate or lower the SMB_delivery_ratio
- 431 - limit or expand max. allowed SMB size
- 432 - change the the even <-> odd logic for SMB on/off

433 Doing temporary changes *in AAPS/preferences* should be the exception because

- 434 - they require multiple steps, including entering a password
- 435 - you will often forget to set everything *back to original* settings, a couple of hours (or already
- 436 minutes) later.

437

438 5.2.2.3 Triggered Automations: Grey extra DIY cockpit buttons for pre-programmed “responses”

439

440 Recognizing conditions for fully automatic handling by the loop may not be not possible, or come
441 too late for the loop to act on. Examples would be

442

- 443 • *exercise*: Minimum an hour *before* starting exercise, “the loop should know” to be able to
444 lower iob and elevate bg by the time exercise starts.
- 445
- 446 • *snacks*: High carb snacks, sweets, consuming ice cream or having a sweet drink comes
447 with the problem of even steeper glucose rises, but overall a lesser insulin need, compared
448 to major meals (for which we tuned our FCL according to [section 4](#)).

449

450 This not necessarily implies that snacks *need* different settings than a meal. After all, autoISF
451 was designed to react to all available data, especially to where the developing glucose curve is
452 headed. So, depending on your effort to set parameters for a broad variety of meals (notably:
453 how well you avoid to invariably bounce fast against your iobTH), you might be able to accom-
454 modate low carb, snack, and major meals with *one* set of settings.

455

456 In FCL autoISF, this is a bit more difficult than in HCL autoISF applications, because FCL
457 involves revving up iob supply (largely via big bgAccel_ISF-weights) sometimes too much,
458 to be balanced by just a snack getting absorbed.

In case a snack did trigger a “full meal response”: (1) You probably must continue snacking to prevent a hypo from your initial FCL over-reaction. (2) For future days, analyze your data (and snacking habit) to define how to prevent this from happening often.

For increased comfort and safety, you might have to differentiate, and make use of what follows for the *sweet snack* example, [case study 5.2](#).

Note that in the iPhone versions of autoISF (**Trio** and **iAPS**) there are no Automations . Instead you need so-called **Middleware**, like for instance suggested for %sensitivity (profile ISF) adaptation by one user here: <https://discord.com/channels/953929437894803478/1025731124615458848/1238099464531611668>

Tuning aggressiveness

A sweet snack likely benefits from even more aggressive initial FCL performance than set for the meals in your normal spectrum of diets.

Therefore, you could set

- a higher **temp. profile%** and/or
- a temp.elevated **bgAccel_ISF-weight** (see screenshot of my Automation).
- a **low temp. target** (76 mg/dl for instance; this additionally helps maximize the first SMBs that will automatically be triggered at detection of acceleration)..

When first defining and testing this Automation, also check:

- that the safety limits as discussed in [section 2](#) will not block the intended elevated aggressiveness
- SMBs will not get outrageously big, and iobTH sometimes exceeded by too much.


Note that “the last SMB” is allowed to overshoot the effective iobTH by up to 30%, where it will be cut (or by up to 20% at even target> 100 mg/dl).

Limiting iob

For “just a snack”, total insulin need will be lower than for a meal.

If you would just have your sweet drink, and your meal-oriented FCL would “attack”, iob likely would become too high, and a glucose rollercoaster would start, with you needing to consume more =>

495 If you just have a snack, or drink a small glass of juice, you can lower the **iobTH_percent**
496 accordingly.

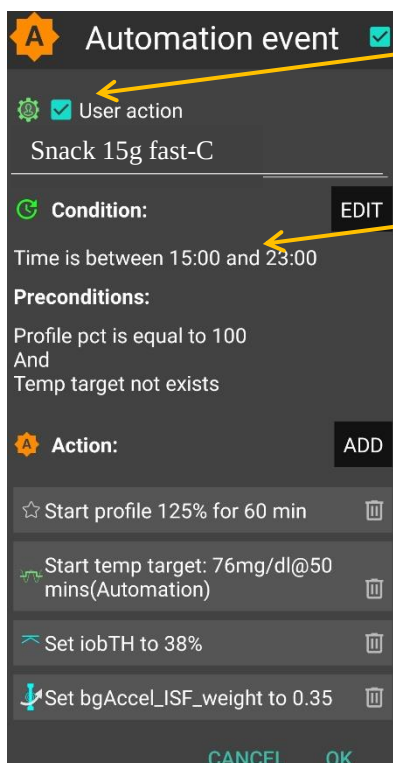


„User action“ is always ticked-on

This will, in the defined time space ..

..offer the “DIY cockpit” button..

..which I must press any time (~90...30 minutes) before my snack.

Note: Make sure that iobTH and bgAccel_ISF_weight revert to default afterwards

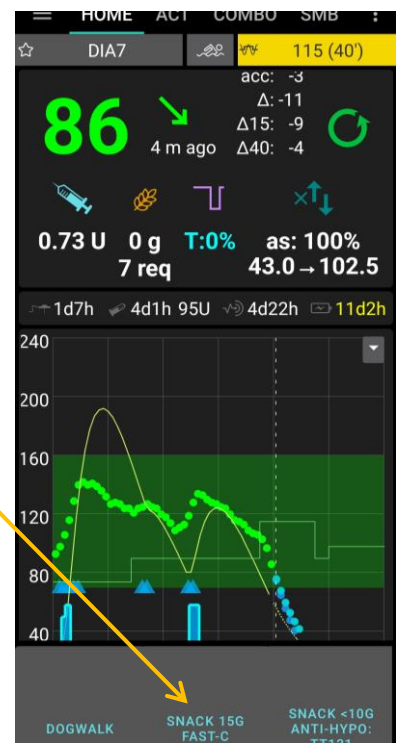

497
498

499 So, this can be a little **extra “project”** when setting up your FCL.

500 You need to research your snack habits (if any), and over time find out which settings in
501 the snack-related Automations work well.

502

503 **In everyday life** you **then just** must **press the related button in your cockpit** (which is
504 not time critical at all, except it should be clicked *latest* a couple of minutes after you took
505 the drink or snack).

506

507 If you consume more, and also eat something with your sweet drink, this will more
508 resemble a full meal... however, with unusual amounts of fast carbs.

509

510 **Caution:** Pressing your snack button a *second time* would **not** help because the lowered
511 iobTH does not allow iob going high enough. So you are better off just letting your *normal*
512 FCL meal routine run, after your snack mode expired.

513


514 Question to developers: Do we already have, or can we get, the option to automatically block out
515 using the same User Action Automation a second time within, say, 2 hours?

516

517 Other options (when you just can't stop snacking) would require a manual modulation
518 regarding %profile and/or bgAccel_ISF, but keeping the full default set iobTH_percent, or even
519 elevating it (refer to [section 5.2.3](#)). If that happens often, define for yourself an extra User action
520 Automation for a bigger snack (= another grey DIY cockpit button).

521
522 **Caution:** Setting a different iobTH or bgAccel_ISF_weight can not be done with a duration
523 attached. Hence you **must** define a suitable **additional Automation, that** must be active
524 in tandem, and **restores the iobTH or bgAccel-**
525 **ISF_weight** in AAPS/Preferences. Else, once your
526 Automation set in, it will *forever* shift these important
527 parameter settings!

528
529 If for instance you have several Automations that, in
530 combination with a set elevated TT also set a lower iobTH:
531 Don't be fooled, the duration only applies to the TT. You
532 need an extra Automation for all of them.
533 I picked out the highest of the altered iobTH values that
534 these Automations can set (45_percent), and then I can
535 automatically restore my default desired 60% via this one
Automation (see screenshot - - >)

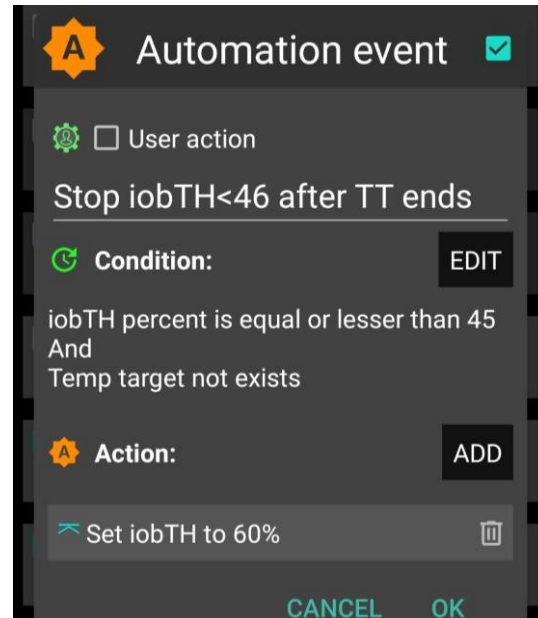

536
537 [Installing the DIY cockpit button](#)

538
539 In the related Automation, just keep the "User action" box clicked at all times, and define in the
540 Conditions when you want to see that button available for cockpit use (see screenshot above) =>
541 you will see that button offered.

542
543 Besides snacks, also any **other recurring special situations can be addressed via a DIY**
544 **cockpit button, and receive different aggressiveness up to a suitable iobTH level.**

545
546 Over time you can have a big number of User action Automations, and keep them "shelved" rather
547 invisibly (clicked in-active, top left box) in your long list of potential Automations. Even when active,
548 they only show in your cockpit (bottom grey field of your AAPS home screen) *in the time slot you*
549 *assigned* as potentially relevant.

550
551 [In the future you might be able to set the stage for a snack and other "disturbances" also via](#)
552 [an extended menu behind the TT button on the AAPS home screen, see \[section 5.3.3.1\]\(#\)](#)

553
554

555 Discussion

556

557 In case you do have a snack habit and ...

- 558 • can not find settings, as in [section 4](#). defined for your meals, also suit your snacks
- 559 • can not pin a time slot or other Condition to it for programming an Automation response
- 560 as in [section 5.1.4](#)

561 .. then you minimum need a “snack announcement” for which the extra button in your DIY cockpit
562 provides a time-uncritical 1-button-push solution.

563

564 This could be a good solution for kids in kindergarten, too. Make sure caregivers
565 understand to use it *only once* for *one* snack. Continued snacking would require iob as for a
566 meals. This is what the FCL loop takes care of automatically; **using the snack button**
567 **several times in a row would limit iobTH at a too-low level!**

568 In a software update, we might try to automatically block usage of that type of
569 Automation for 2 hours, after it was once used.

570

571 5.2.3 Temporarily exiting the FCL

572

573 The “last resort” alternative always is to **temporarily** leave the FCL mode, and handle any
574 disturbance “the traditional way” in **hybrid closed loop**.

575

576 For this, we switch the automatic aggressive adaptations of ISF to the bg curve OFF that are only
577 needed in FCL

578 (if in hybrid closed loop you like e.g. the dura_ISF adaption still, you alternatively could elect
579 to just set bgAccel_ISF_weight temp. to zero, instead)

580

581 Do not forget that, before meal starts, giving a bolus will then be necessary again.

582

583 The suggested improved FCL cockpit user interface with an extra version of violet loop on the
584 AAPS home screen ([section 5.3.1](#)) would facilitate this transition FCL < - > HCL, including
585 automatic removal and re-appearance of the insulin button at the bottom of the APS home screen.

586

587 In case this feature is not yet available, you must:

588 Exit the FCL mode by going to AAPS/preferences/put in your password/OpenAPS SMB/scroll down
589 to autoISF settings and switch “Enable ISF adaptation..” OFF

590 (or, alternatively, set bgAccel_ISF_weight to zero).

591

592 **Caution:** Unfortunately, there is no way yet for your full closed loop with ISF adaptations to come
593 automatically back on, after a selected time for instance. So **do not forget to switch your**
594 **autoISF fully back on, later.**

595

596 As this will often be forgotten, it may be worth doing a “User action” Automation, for a “temp.
597 FCL OFF” grey button (see [section 5.2.2.3](#)).

598 Caution though, there is very limited experience with this brand new feature. Make sure your
599 Automation definition really applies a duration (or other condition) that will automatically
600 terminate all settings changes it made. As we have seen e.g. in [section 5.1.4](#), this is not always
601 the case.

602

603 To recognize whether autoISF currently runs with **ISF adaptation** or not, please consult the
604 **“ai: %” indicator** below the Autosens% **on the AAPS home screen.**

605

606 From autoISF 3.0.1 onwards, there is also a very easy way to see effective ISF and effective iobTH
607 in the 1st screen of the **SMB tab**. At the same time, there you see the adaptation of sensitivity to:

- 608 • a set %profile change (or effect from Autosens, in case you have that activated)
- 609 • a set temporary target
- 610 • the Activity Monitor
- 611 • +/- exercise mode

612 So, occasionally (especially in your early set-up phase, after starting of a meal) it is a great idea to study the
613 SMB tab to find out what is going on. See example given in [section 5.4.5](#)

614

615

5.3 Modulating aggressiveness manually from the improved FCL-cockpit

Skip this section 5.3 (next 7-8 pages) unless you are deeper interested in discussing further user interface upgrades. Actually, some suggestions made are probably an “over-design”. After trying a lot of options for refinements out, the author returned pretty much to a “keeping-it-simple” route.

My main suggestion is to get that violet loop button (sections 5.3.1-5.3.2), something I think many would use - very handy certainly in the setting-up stage, too, for easy switching between the “old” HCL, and new territory in FCL.

autoISF is an early dev variant of AAPS, and as user you are participating in an on-going development. Of note, autoISF 3.0.x was launched without many of the cockpit features that are suggested below in green font color.

Only what is written in black is at this point of some relevance for using autoISF.

No need to read any of the green lines, unless you are interested in contributing to define/design/program further improvements.

This is also an open invitation for you to contact us in case you could help program a module for one of the suggested user interface extras.

For future integration into AAPS Master, an eye should be kept also on the question which *other* modes (like FCL using Automations and others mentioned in [section 13](#); and maybe also HCL) might benefit from some of the extra features.

For the time being, multi-step work-arounds may become necessary

- In many cases, going into AAPS Preferences and changing settings would be needed (...plus not forgetting to change these settings back, afterwards).
- Automations allow a DIY FCL cockpit, see [section 5.2](#) and [case studies 5.2](#) and [6.2](#)

Keep in mind, though, that the **goal should be to interfere with the loop as little as possible**.

Under certain conditions, it can run **fully automatically** without any user interaction, as described in the preceding [section 4](#). + [section 5.1](#).

Just like in the airplane **cockpit**, also cruising in full auto mode should involve having an eye on the instruments, and on potential disturbances ahead in the environment.

E.g.: storm ahead => instruct your plane to climb to another flight height.

Analogy: exercise ahead => setting an exercise TT, or => pressing a button that activates a sequence of instructions (some of them probably hinging on conditions, like actual iob), how to manage through that exercise situation).

656 So, **for the occasional „disturbance“ coming up**, you should find an easy way to
657 • call up a pre-programmed routine for automatic management, with auto-adjusted
658 aggressiveness, or:
659 • tweak a setting or two, to temporarily adjust the aggressiveness
660 • There may also arise a desire to just exit the FCL mode, and “be your own captain” for
661 mastering a special situation.

662 All this is facilitated within seconds right from the AAPS home screen’s **cockpit features to the**
663 **extent they are already incorporated**, or to the extent you can build **alike** DIY cockpit features via
664 Automations, as described in [section 4.1.3](#) and [case studies 5.2](#) and [6.2](#)):

665
666 • The button that is integrated into the **violet FCL icon** serves as emergency off button, to
667 quickly stop FCL, or to at least to immediately stop any more SMBs (...just for a couple of
668 minutes, or for the remaining meal time: pick from the options offered with just one
669 keystroke).

670 Via the violet FCL icon on your AAPS home screen, you also can access a temp. switch-off
671 button for SMBs (see section that next follows below).

672
673 • The **three top fields** (%profile, exercise, TT) provide access to temp. tuning of core
674 parameters, and/or to some pre-programmed routines.

675 Taken together with some **new indicator fields** about your loop state ([section 5.4.3](#) and [5.4.4](#)),
676 and the **grey DIY cockpit buttons** ([section 5.2.2.3](#)) this makes the AAPS home screen your
677 **cockpit** for Full Closed Looping.

678
679 Let us look on each of these **suggested** cockpit elements in some detail:

680
681 **5.3.1 Violet FCL icon and underlying buttons**

682
683 Novices to FCL, or really anyone running into a very special situation, may appreciate that the new
684 closed loop icon on the AAPS home screen in pink (for FCL) has buttons to quickly shut off getting
685 more SMBs (1st row), or to enter other loop modes (second row).

686
687 It functions very much as the other ones that you know from HCL already, and in fact you
688 get offered some of the same options (for instance, to switch the (full) closed loop off for 15
689 minutes for going to take a shower)


690

691 Note that in FCL you leave all BG regulation, notably against meal spikes, to the loop. So, try not to
692 disconnect in phases when your FCL must ramp up your job.

693 The required insulin would still be supplied *after* you reconnect. However, without the user
694 pre-bolussing, the delay would be more of an issue in FCL than it had been in HCL.

695

696 Just pressing on the FCL icon, a dialogue box comes up:

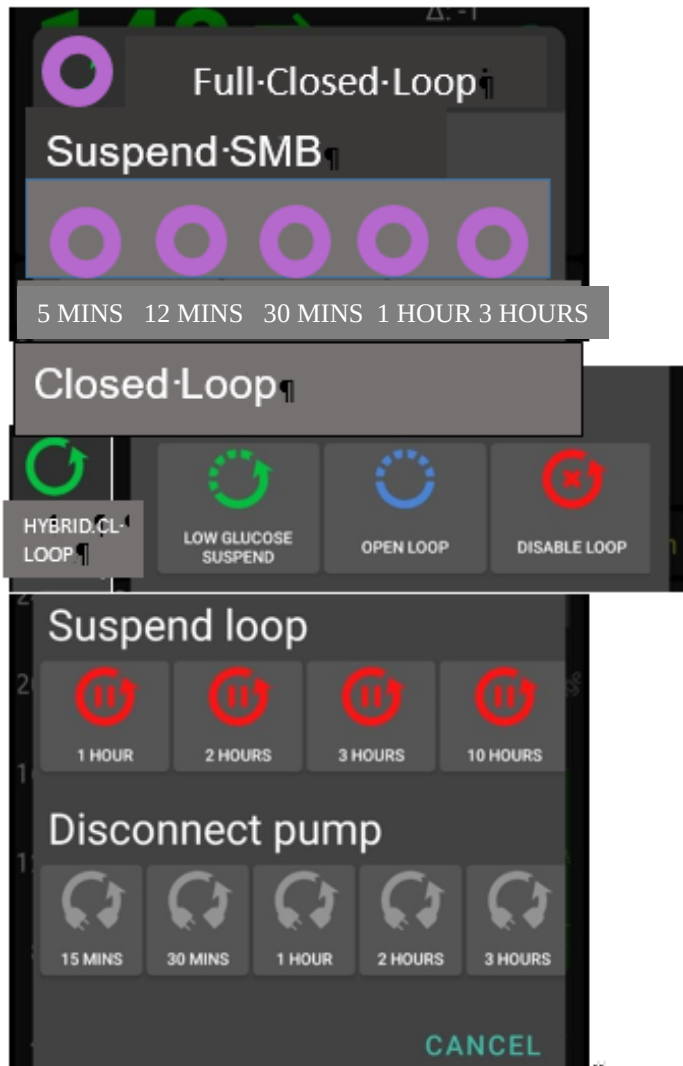

<-add-arrow, form-like-like-CL-Loop...¶

<-dotted+arrow, form-like-like-LGS...¶



DEV: when-in-Hybrid-Closed-Loop= top headline+green-circle), there is no Suspend SMB part but it starts right with „Closed Loop“. In that screen, There, the 1st element must be option to pick ¶

.....pink-circle, FULL-CL-LOOP ¶

where here the green-HYBRID-CL-LOOP stands ¶

For loopers who did not set up FCL, a feedback must come up ~ „FCL not installed“ if they press on that button ¶

Skip:

= This User Interface is a not implemented idea only.
Please contribute, or wait for update with the missing info

697
698 Pressing „Suspend SMB“ provides fast and easy „emergency braking“ regarding delivery of more
699 SMBs:

700 Select the one with the desired number of minutes: 5 or 12 for just blocking the potential next
701 SMB(s), and up to 3 hours to manage the entire rest of this meal with %TBR from then on.

702

703 Whenever, and why-ever, your FCL is in „no SMBs allowed“ mode (e.g. automatically after
704 surpassing an iobTH, or triggered by a set odd TT), the FCL icon will turn into a dotted one.

705 Instead of remaining **duration to end time** it indicates in the middle „the condition“, „**iob**“ or „**TT**“

706 Add an indication if suspend SMB comes from an Automation, e.g. add an „**(A)**“, **underneath** the
707 #minutes, iob, or TT in the middle of the dotted violet field.

708 So, as in other (already in HCL existing) cases, those icons show in the middle the minutes left that
709 they will be running, or the condition which would have to go away for this temp. setting to stop.
710 It always auto-reverts into the FCL state and FCL icon, when time (or other condition) has elapsed.

711
712 Pressing „**HYBRID CL. LOOP**“ or other buttons from the 2nd row provides fast and easy
713 „emergency exit“ into other modes.

714 This enables beginners an easy „temp. escape“ into their well-known HCL (green) at any
715 point of time. bgAccel_ISF_weight is set to zero when going FCL->HCL. HCL can run with
716 autoISF (for instance dura_ISF) uninhibited otherwise. (check implications for HCL users of
717 autoISF ??).

718 Note: These options from row 2 have no time limit. Loop will not by itself go back to FCL. You see
719 the different loop icon as a reminder to manually revert, when ready.

720

721

722 5.3.2 Buttons „Insulin“, „Calculator“ etc at bottom of AAPS home screen

723

724 These buttons are **not useful any longer in FCL**, and automatically disappear whenever in FCL
725 mode (also in Suspend SMB state), and re-appear when leaving FCL. This applies also when an
726 Automation or technical system failure shut off FCL.

727 Users who, maybe in the beginning phase, feel better having those buttons, can override
728 the removal (of the insulin button, or any other) by going into /preferences/overview/buttons
729 and forcing them on. They only remain on until the next re-entry into FCL mode, when auto-
730 off happens again.

731 The reason why we do this: It really is important to let the loop loop, and not interfere more
732 than absolutely needed. Any bolus the user gives will sure distort the bg curve, on which
733 autoISF, especially when aggressively tuned for FCL, builds a lot of its decisions!

734

735

736 5.3.3. Three top fields (%profile, exercise, TT)

737

738 Depending on the variedness of lifestyle, the desired %TIR, and the initial tuning effort put in, the
739 user may want occasionally to „tweak“ the aggressiveness of her/his FCL.

740

741 The top 3 fields (grey in default mode, **yellow when temp. in mode with changed**
742 **aggressiveness**) serve as quick and easy entry points to make temp. switches (as users will be
743 used to for %profile switches, or for setting an EatingSoonTT in HCL, .. which they still can do in
744 FCL ... but more:)

745

746 Expert FCL users might need this feature rarely, but probably at least to manage activity after
747 meals: Each require opposite aggressiveness, and the switch has to come in a certain point in
748 time that would be difficult to capture. (More see [section 6.4](#))

752 The TT field (top right of AAPS home screen) is a primary daily interface, and a dialogue field
753 opens when pressing on it

755

758 (1) Who is happy with the initially well tuned FCL and does not have huge variations in daily eating
759 and moving around, will **not use** the TT **at all**. FCL is possible without an intervention via the

760 TT button in your cockpit. Actually 4 of 8 modes (GGG ...YYY permutations, list see [section](#)
761 [5.4.1](#)) are not making use of TT.

762 (2) Super easy is also, to just input **any odd-TT** (odd-numbered temporary target) that will shut out
763 any SMBs for the set duration. *That can be a good idea when having a small snack, for*
764 *instance.*

765 Super quick access to stop SMBs is possible also via the loop icon ([section 5.3.1](#)).

766 Specifically, an **EatingSoon TT** can be activated here (limited relevance see [section 2.5](#)). It is
767 time-un-critical, can be manually set, or come up via an Automation.

768 The cockpit enables you to set the iobTH differently (override) for the current meal.

769 Alternatively, iobTH can be temporarily changed in /preferences or using an Automation.
770

771 Temp. iobTH will always revert to default when the TT expires. If another TT immediately
772 follows, like in the example of the screen above, it will calculate, (then) show and use a new
773 temp. iobTH.

774 (3) The third way is to **use the input mask** (if already ncluded in your software version see picture
775 above) **to freely modulate the loop aggressiveness** for a declared number of minutes. Click
776 the bottom big square(s): Either HYPO, or ACTIVITY, or EATING SOON, or ACTIVITY and
777 EATING SOON (example in the pictured screen above). Make or override entries in the offered
778 fields. Press OK.

779 (4) The fourth way is to exclusively use one of the 4+4+2 little buttons seen in the bottom part of
780 the TT dialogue box (if already included in your software version). They provide a set of
781 settings (as will immediately show in all input fields above) that the user has set up in
782 Preferences/SMB/autoISF/FullLoop (refer to [section 6.3](#)), and can freely label there. *For*
783 *instance „hiC“ at high carb EatingSoon, „piz“ for Pizza/fatty meals, „grd“ for garden work,*
784 *„mtb“ for mountain biking ...*

785 **Capturing good settings for not-everyday situations in /preferences** (if already included)
786 **allows calling them up within 1 second**, from your cockpit on the AAPS home screen (...and
787 won't ruin the FCL experience at all , especially because in most cases it is not time-critical,
788 how long before the intended exercise the buttons are pressed).

789

790 [Case study 6.2](#) demonstrates that nearly the same performance and comfort can be reached
791 via the **DIY FCL cockpit** with the grey extra buttons appearing at the bottom of the AAPS home
792 screen, based on Automations with User action (see also [section 5.2.2.3](#)).

793 The example picture given above, and also [case study 6.2](#), is the most complicated (but also most
 794 useful) case, **when exercise follows after a sizeable meal**. It is then that you need (a) aggressive
 795 FCL initial performance at the meal, but, *exactly when (!)* a (for the intended sport already
 796 temp.lowered) *iobTH* is exceeded, you need (b) to have SMBs automatically switched off and go
 797 into the „milder“ mode, as defined for the exercise (with *high* instead of the immediately prior
 798 *low*TT, that automatically significantly reduces iobTH again, and insulin sensitivity(resistance)
 799 settings too).

800
 801 Pressing exercise related buttons will automatically also light the **exercise button** on the main
 802 screen yellow.

803
 804 To summarize, the TT dialogue field offers easy but powerful ad-hoc [modulation of loop](#)
 805 [aggressiveness](#) for FCL (if already included).


806
 807
 808 5.3.3.2 Exercise button (see more in [section 6](#).)

809
 810 The exercise button automatically lights yellow when exercise related TTs are activated [in the TT](#)
 811 [dialogue box](#).

812 4 of 8 principal FCL modes ([section 5.4.1](#)) are making use of the exercise button.

813
 814 If pressing on the exercise button, a dialogue box appears (*if extended design for FCL cockpit is*
 815 *already launched*) with info on exercise setting first (and opportunity to override), plus below the
 816 activity monitor (experimental for auto-tracking of lighter movement during the day, and effects on
 817 sensitivity that may have. See [section 4.5](#)).

818
 819 So, first the exercise settings (as set under TT) are there to read. Example :

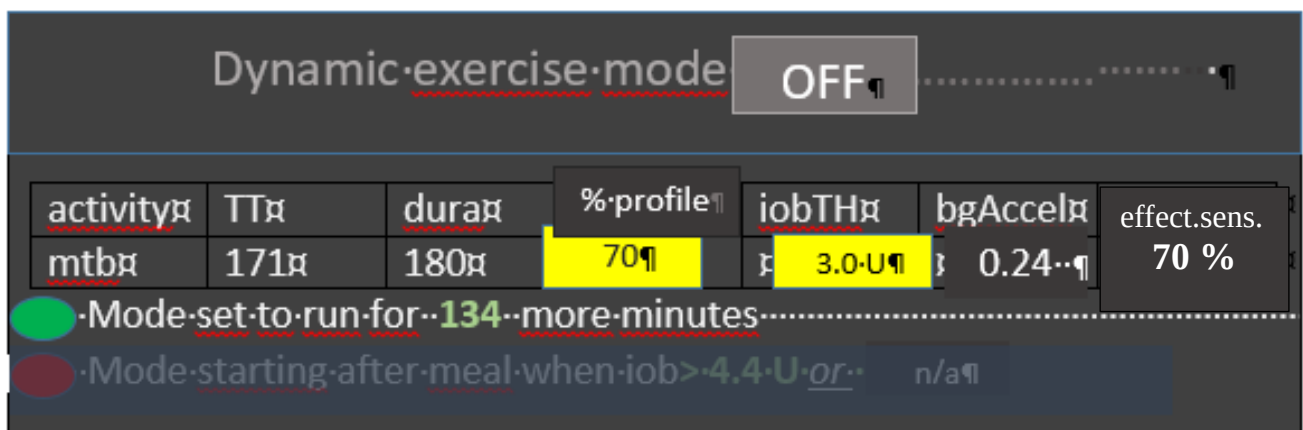

820
 821

822 The exercise (here mtb) is selected in the dialogue box of the neighboring TT field, and there auto-
823 filled with settings made in the set-up and tuning stage by the user under preferences. They are
824 reported also under the exercise button here, and TT, duration, and % sens (the temp. profile
825 sensitivity that also shows on the %profile field on the left side of the exercise button) can be temp.
826 changed there.

827 iobTH, bgAccel_ISF and overall resulting effective sensitivity ratio (effect.sens. %) is given in the
828 other fields.

829 The **middle field** of the table, „% **profile**“ either picks up the % set under the %profile button, or
830 an input can be made here, in the exercise button domain, which will:

- 831 • turn the neighboring %profile button on yellow and show that inputted % on it, too
- 832 • be multiplied with the result from the exercise mode settings per se, and change the %
833 overall, accordingly.

834 So, if this middle field of above table (dialogue box of sports button) contains a figure other than
835 100, input field becomes yellow, and you are operating with a combination of traditional PLUS new
836 exercise mode (with all three top buttons of your FCL cockpit yellow). This maximally will soften
837 aggressiveness, for which you get an idea by the last calculated figure.

838

839 The mode is either running already (for another number of minutes, as probably also shown in the
840 yellow TT field anyways). Or it is scheduled to run, after insulination for a started meal reaches
841 iobTH (as in table). Or, no exercise is scheduled (both points red, no entries.

842

843 The lower part of the exercise dialogue box (not pictured above, but see in [section 6.5](#)) is
844 dedicated to the Activity Monitor

845

846 5.3.3.3 Profile button

847 The profile button can still be used to set a different profile, or profile%, for instance to adjust for
848 days with sickness (as you are used to from hybrid closed looping). 4 of 8 modes are not making
849 use of the profile button.

850

851 Any inputs made here will be used to modify profile_ISF on which all further changes are made on
852 (multiplied with).

853

854 The profile field remains grey if standard profile is applied.

855 It turns yellow, displaying a %figure relating to any altered loop overall aggressiveness:

- 856 • When no inputs (changes from 100% profile) are made here, but inputs in the TT field,
857 e.g. for exercise, automatically lead to different effective sensitivity ratio

858 • when% is changed by input in the profile button itself, it will be multiplied with with
859 profile_ISF and be used in place of profile_ISF by the algorithm.

860 However, for exercise (sports) you no longer must make an entry here, because
861 reasonable %reductions should be automatically provided, driven by your set TT (and half-basal
862 exercise target), see [section 6](#).

863

864 5.4 Recognizing your loop state in the AAPS home screen

865

866 5.4.1 Modulated loop aggressiveness via the 3 top buttons in the AAPS home screen

867

868 3 Buttons (%profile; exercise; TT) allow temporary modulation of aggressiveness (sensitivity) as
869 may be needed e.g.

870 • going into meals (EatingSoonTT => providing lower bg “starting point” and more pos. iob),

871 • or doing exercise (Exercise mode, limiting basal, iob etc for hypo prevention).

872 ,

873 All buttons grey (GGG) indicates: The loop is running using un-modulated profile parameters

874 • Note that features like autoISF might still build on profile (regardless whether original, or a
875 temp. modulated profile) for further modifications (See calculation cascade for used
876 sensitivity (sens) in the SMB tab. Example see FCL e-book [section](#)).

877 Any **yellow** button indicates you are running with a **modified profile sensitivity**.

878 *Actually, the middle button could be yellow without being active, see table below*

879

	Less aggressive loop = more sensitive user needs temp. less insulin	Standard loop using the set profile	More aggressive loop = less sensitive user needs temp. more insulin
% profile	under 100 % (y)	= 100% (G)	above 100% (Y)
Exercise button	ON (y) Note that it can	OFF (G) <i>If ON (y), it is <u>not</u> active but ready to automatically become active whenever any TT above target is set (e.g. by any Automation)</i>	OFF (G)
TT	above profile target (y)	= profile bg target (G)	below profile target (Y)
Example	Exercise mode with weakened %profile, yyy		EatingSoonTT, GGY
..used e.g. in :	Case study 6.2 (yyy)		Case study 6.2 (G_Y via Automation) Case study 5.2 (YGY)

- 880 Case study 6.2 uses extremes from both ends of this sensitivity modulation spectrum:
- 881 • yyy for reducing basal and iob, and getting weaker ISF, during exercise
- 882 • GGY for temporary "EatingSoonTT" boost when, in FCL, a high carb meal had
- 883 • The even stronger YGY case is contained in Case study 5.2, where for a high carb snack an extra boost,
- 884 using >100% profile sensitivity, is temp. added

885

886 5.4.2 Color scheme of top cockpit buttons tells kind of closed loop that is running

887

888 3 Buttons (%profile; exercise; TT), each in 2 states (yellow Y, or grey G)

889

890 GGG is the loop running with set profile. Any yellow setting modulates sensitivity

891 or attempts to modulate, whenever another required condition may become true (e.g. also via an

892 Automation):.

893

894 Overall, 2 exp 3 = **eight principal FCL modes** are possible. They need to be **differentiated**

895 **further**, based on whether the modulation goes into the **more aggressive (Y, higher % and/or**

896 lower TT), or the *less aggressive (y, lower % and/or higher TT and evtl exercise)* direction. Note

897 the Exercise button works only one way, to make the loop less aggressive.

898 **DEV idea: Should we make the Y to violet or orange (more aggressive) and keep the y**

899 **yellow (less aggressive than profile)? Then beginners would easier recognize whether they**

900 **go more or less aggressive ...**

Color combinations	Lower aggressiveness	Higher aggressiveness
GGG		profile	
Gyy GyY	dynamic exercise mode (if TT > profile target)		EatingSoonTT (as TT < profile target: middle y inactive, GyY = GGY 2)
yGy-yGY- YGy-YGY	„traditional“ exercise mode (if <100% sens AND TT > profile target)	mixed cases 1)	boosted EatingSoonTT (if >100% sens AND TT < profile target)
GyG		GyG = GGG 2)	
yGG YGG	Loop adjusted for period of lower insulin need (e.g. 80% for night after exercise)		Loop adjusted for period of increased insulin need (e.g. 120% during sickness)
GGy GGY	Loop running with temp elevated target e.g. for enhanced safety against going low (HypoTT)		EatingSoonTT (if TT < profile target)
yyG YyG	same as yGG	y inactive 2)	same as YGG

yyy - Yyy- yyY- YyY	dynamic exercise mode further modulated by <100% profile	mixed cases 3) yyY = yGY, 2)	Boosted EatingSoonTT (if >100% sens AND TT < profile target)(YyY = YGY, 2)
---------------------	--	---------------------------------------	---

901

902 1) **yGY** would be the combination of setting a milder % profile, yet orienting the loop to a lower
 903 target. Usually we do not use this combination, but it would be OK for an instance where we
 904 want to start the next meal at a low target, but we want our loop to be careful and not rush
 905 getting there.

906 **YGY** would be the combination of setting a more aggressive % profile, yet orienting the loop to
 907 a higher target. Usually we do not use this combination, but it would be OK if we believe we are
 908 in a state of enhanced insulin need, and want to get to a set elevated TT (which could be what
 909 we want to go max down to to have room for exercise, or just to be sure not to go into a
 910 hypoglycemia if we misjudged the symptoms about enhanced insulin need.

911

912 2) A yellow middle button y works really only in the combination ..yy. Any combination ending yG
 913 or yY will work the same as if it were ending GG or GY. Although the yellow middle button
 914 does not make a difference in these cases, it might be “justified” by the loop showing its
 915 readiness to employ the exercise mode, whenever a TT > profile target comes around, e.g. via
 916 an Automation.

917 Same as GGG, only *if* (e.g. an Automation) sets a *TT>profile target*, automatically a softer
 918 response due to additional exercise mode kicking in *then*.

919

920 3) **yyY** would essentially working like **yGY** (discussed under 1), because the middle button,
 921 whether y or G, makes no difference with an elevated TT.

922 Again, this the combination of setting a milder % profile, yet orienting the loop to a lower target.
 923 Usually we do not use this combination, but it would be OK for an instance where we want to
 924 start the next meal at a low target, but we want our loop to be careful and not rush getting
 925 there.

926 **Yyy** would be the combination of setting a more aggressive % profile, yet orienting the loop to a
 927 higher target. As opposed to the **YGY** case (discussed under 1), in **Yyy** the Exercise mode is
 928 active and strongly reduces aggressiveness. It is hard to imagine, why this combination this
 929 combination of settings could make sense; maybe if a “long waved” state of reduced insulin
 930 sensitivity meets a case where normally Gyy would be sought (see there).

931

932

933 5.4.3 Information printed on the top buttons

934

935 The yellow TT field shows the currently valid TT (and further duration):

936 (profile) stands for the abbreviation you labeled your selected running profile

937



938 In the special case of settings for meal preceding sports, the field will look slightly differently:

939



...and ...

940 ... when iobTH is first time exceeded, this automatically switches to:

941



942

943 Likewise, if on the AAPS main screen just an **EatingSoonTT** is set (e.g.72), this is entered with the

944 desired duration. Afterwards, it automatically reverts to profile target and the display turns grey

945 again (can also turn green first, for a short transition period), with the profile target (e.g. 90) on it.

946 (No time limit, then, for the profile value as set in preferences).

947 Without sports context, the middle field remains grey.

948



949

950 Independently from setting a TT, the user can choose to set a **%profile in the left top field**, for an

951 independent number of minutes, e.g. 70% in this screen example: Also, or additionally, this will

952 influence the resulting ISF and sensitivity%

953

954



955

956 The % might change and turn yellow also in context of making TT inputs in the related dialogue
957 box (see chapter TT dialogue field, above). Still, the % (or the length of time the profile switch shall
958 be active) can be independently overridden in the top left field, if so desired.

959

960 If an **Automation** sets a %profile, and/or a TT (e.g. *automatic detection of meal start at condition*
961 *e.g. when delta >10*), this would automatically show in respective field(s) turning yellow and
962 showing the temp. setting. To show the set parameter comes from an Automation, „ (A) „ is added
963 in the end of button text.

964 Note that an **Automation might not be permitted** to change settings by “killing” another
965 **still running Automation** (always consider that, when putting the duration into your
966 Automations!). For instance, you cannot switch from 130% profile to 110%. Either the 130%
967 times out, or you **need an extra “in-between” Automation that terminates** the 130%
968 under described conditions (example see around line 100 in [Case study 6.2](#)). – This
969 “design” is for a good reason: The assumption here is, that your 1st Automation (the 130%
970 in the example) is designed well and runs for a reason. It should either “get finished” when
971 the job might be done (and kick in again, if not), or, in exceptional cases, it should be
972 consciously terminated by another well thought through 2nd Automation (describing the
973 conditions in which you would find that other Automation more important than “finishing up”
974 the one that was already running). That “in-between” Automation makes the loop return to
975 base profile, which is a signal *to all Automations*, to now check whether any conditions
976 exist, to activate a 3rd Automation (as in example of [Case study 6.2](#)).

977

978 Advice: Try to stay away from Automations that also aim at temp. modifying aggressiveness (e.g.
979 temp. setting different bgAccel_ISF_weight). For the reason just given in above note, they often will
980 not kick in anyways. Generally, it also is no good idea to double up sub-algorithms for tweaking
981 loop behaviors (“loop inside a loop”).

982

983 Summary

984

985 A lot of avenues were shown that could help you or your loop manœuvre through a variety of
986 “disturbances”. You should not have to try out many of them, and (like the author), you should find
987 a way to narrow it down to what really helps in **your** everyday T1D management.

988

989 Try to **keep things as simple and clear as possible**.

990

991 That said, a limited number of Automations can be of help in distinct scenarios (that differ in
992 purpose and in applicable time of day).

993

A good one could be for night time, when your odd profile TT has SMBs shut off, but your experience after pizza nights tells you that, under certain condition patterns (bg, iob), an SMB or two should be „allowed in“ (see example given in [section 5.1.2](#); used also in [case study 4.3](#)).

Another good example, if you go usually FCL without any use of the TT button (which would be a meal announcement of sorts), is to define an Automation that, after detecting a meal start, automatically sets a low TT to get maximally aggressive first SMBs (as is the author's preferred way, mentioned already in [section 2.5](#), used also in [case study 4.3](#)).

1005 5.4.4 FCL related indicator fields in the AAPS home screen

1006

1007 In extra data fields of the AAPS main screen you can always see (not change) the key
1008 „aggressiveness“ parameters your loop currently operates with (see also home screen
1009 example below):

- 1010 • how profile sensitivity (**ISF**) adjusts by the %profile input, by autoISF, and/or a set
1011 exerciseTT, resulting in an effective sensitivity (ISF that is used to determine
1012 insulinRequired. Details for every loop decision see result/debug section of the SMB tab).
- 1013 • next to current available iob number is an indication of your **valid iobTH** (the iob above
1014 which no more SMBs will be given)
- 1015 • The AAPS home screen additionally shows, above the deltas, the current **acceleration**
1016 Having a look at that can be valueable. For instance, when glucose is relatively low and still
1017 falling, a positive (and getting more positive) acceleration indicates that bg will swing back
1018 up, rather than crash low. This will give info about necessary snack size, and hence help
1019 avoid both, unnecessary calories, and going on a bg roller coaster.

1022

(auto- re-appearing when violet \rightarrow green loop)

1025

1026 5.4.6 Info given every 5 minutes in the SMB tab

1027 When clicking on the SMB tab, you see how your standard and temporary settings, as well as the latest bg and iob status, influenced the last decision of your FCL.

Example 1: A 80% temp. profile modulates 60% iobTH to 80% of 60% = 48%:



The profile ISF of 41 mg/dl/U got modified by the set 80% temp. profile to $41/0.8 = 51.3$ mg/dl/U, called “ISF unchanged” (before “start autoISF”).

autoISF applies the 4 sub-categories (acce, bg, pp and dura_ISF), and depending on the bg curve form suggests various ISF changes.

The final resulting factor “sens” (see flowcharts in [section 3](#)) is 1.11 (in our case, driven by bgAccel_ISF). This changes the 51.3 “unchanged” ISF to $51.3/1.11 = 42.4$ mg/dl/U

(If **tab** not shown -> Config.builder/APS: tick box on the right)

```

COMBOV2 SMB AUTO NSI
Last run : 3/30/24 22:18

Result

Script debug : d:
  Activity monitor disabled
  inactivity detection: sleeping
  hours
  Autosens ratio: 1;
  Basal unchanged: 0.48;
  ISF unchanged: 51.3
  CR: 10

-----
start autoISF 3.0.1
-----
User setting iobTH=60%
modulated to 48% or 4.8U
due to profile % and/or
exercise/activity mode
SMB enabled; current target
90 is even number
Loop allows maximum power
acce_ISF adaptation is 1.11
bg_ISF adaptation is 1
pp_ISF adaptation is 1.09
dura_ISF adaptation is 1.03
because ISF 51.3 did not do it
for 15 m
final ISF factor is 1.11

-----
end autoISF
-----
currenttemp: 0.34
lastTempAge: 0 m
tempModulus: 28 m

```

Further down in the SMB tab, you can see how this ISF is applied to define the SMB size to be given, and whether any limitations – notably by autoISFmax, max possible SMB size, or maxIOB – cut the amount.

Message	Condition	What does it affect?
Loop allows maximum power	even target < 100	increase in bg limited to 30%, otherwise no SMB; actual SMB delivery ratio is max of fixed smb_delivery_ratio and linearly growing ratio
Loop allows medium power	even target >= 100	increase in bg limited to 20%, the AAPS default, otherwise no SMB; actual SMB delivery ratio is either fixed smb_delivery_ratio or linearly growing ratio
Loop allows minimal power	odd target	no SMB, only TBR available for action
Loop power level temporarily capped	IOB > effective iobTH	Last SMB capped to stay below iob threshold + 30% overrun; IOB getting above user defined iobTH, potentially modulated by exercise mode, activity monitor and profile percent
Loop allows AAPS power level	no even/odd target option active	SMB enabled/disabled according to standard AAPS rules and settings; no iobTH threshold is active

Note that in the SMB tab you can– in “real time” - capture and analyze *only one* decision.

Refer to [section 11](#) for an option that enables extended analysis of the on-going ISF modulations from autoISF. (To do this on your loop phone requires QPython and a logfile emulator).

5.4.7 SMB tab info when operating in 1-minute mode with Libre3

Users: Please supply text and screenshot