

1. Pre-requisites for Full Closed Loop

V 2.52

Please note that with autoISF 3.0 you are in an early-dev. environment, where the user interface is **not optimized for safety** of users who stray away from intended ways to use. Good safety features exist, but these are only as good as the development-oriented user understands and implements them. This is not a medical product, refer to disclaimer in [section 0](#)



- 1.1 Well tuned hybrid closed loop
- 1.2 Fast insulin
- 1.3 Reliable insulin delivery from pump and cannula
- 1.4 Excellent CGM
- 1.5 Meal-related limitations
- 1.6 Lifestyle-related limitations
- 1.7 Time required for setting-up

Available related case studies:

- Case study 1.1: Occlusion
- Case study 1.2: Comparing insulins for FCL
- Case study 1.3: Jumpy CGM
- Case study 1.4: Lost pump connection
- Case study 1.5: Overlapping 2 x G6

1.1 Well-tuned hybrid closed loop

It is advisable to first establish a well-tuned hybrid closed loop before considering the transition to FCL. There are two important reasons for this:

- The UAM full closed loop requires a highly personalized (individual) tuning of settings, so the loop will give insulin mimicking YOUR successful hybrid closed loop mode.
- The UAM full closed loop comes with new parameters to be set and tuned. It would be problematic to set and tune several new parameters before the basics were tuned “right”. Errors could easily be balanced with counter-errors. This can work in single scenarios, but would create a highly unstable system, hard to re-calibrate better later.

1.2 Fast insulin (Lyumjev, Fiasp, Apidra?)

If the user does not bolus for meals, clearly a very fast insulin is needed so, upon realization of a starting meal-related glucose rise, the loop has any chance to eventually keep glucose in range (by common definition, under 180 mg/dl (10 mmol/l))

A modelling study (key findings are summarized in initial section of [case study 1.2](https://androidaps.readthedocs.io/en/latest/Usage/FullClosedLoop.html#fast-insulin-lyumjev-fiasp); for more see: <https://androidaps.readthedocs.io/en/latest/Usage/FullClosedLoop.html#fast-insulin-lyumjev-fiasp>) can show in quantitative terms that **faster insulins**

- will result in significantly **lower** glucose **peaks** than slower insulins
- **tolerate** a couple of minutes **delayed** first meal bolus while not incurring unacceptable height of peaks
- minimize the effect on glucose peak from **different** carb load (**meal sizes**).

39 In conclusion, do not attempt FCL with other insulin than Lyumjev® or Fiasp®, unless, maybe, if
40 you are on a very moderate to low carb diet.

41 According to [case study 1.2](#), Apidra® might work, too, but Humalog® would not work well.

42

43 1.3 Reliable insulin delivery from the used pump/cannula/insulin system

44

45 **Good Tolerance of Lyumjev (or Fiasp): Occlusions threaten the function of the full closed loop.**

46

47 It is very important to have an eye on the time a **cannula (or pod)** is in use (many find **48 hrs** to be
48 the **limit**), and whether hard-to-explain glucose rises happen at ever increasing „fake“ iob (even
49 before a 48 hr routine replacement). (See [case study 1.1](#): You easily lose 25% TIR that day)

50 It is absolutely contra-indicated to attempt FCL coming from leaking pods and associated erratic
51 sensitivity swings that may or may not have been somewhat controlled and tolerable by
52 dynamicISF or other measures when you were Hybrid Closed Looping,

53

54 Stable pump connection

55 In FCL you absolutely rely on your pump delivering,
56 without any further delay, the much needed insulin,
57 after any meal start.

58 Hence it is absolutely essential to avoid any
59 problems from a lost Bluetooth connection. In AAPS
60 preferences / Local alerts, switch alert on!

61 An Automation similar to the one pictured here →
62 could also help recognizing eventual problems.

63

64 See also [case study 1.4](#)

65

66 1.4 Excellent CGM

67


68 You do not give a meal size-related bolus any longer. That leaves all insulination jobs to the
69 algorithm! As glucose values are the very basis for this, please **inform yourself well about** how
70 **your CGM** 1) principally performs 2) whether and how this may depend on data flow and
71 intermediate apps you use 3) specifically, how and where any smoothing is done, and what this
72 might imply for the ISF boosting method you will be using See for instance here:

73 <https://androidaps.readthedocs.io/en/latest/Usage/Smoothing-Blood-Glucose-Data-in-xDrip.html>

74

75 Around meals, a stable Bluetooth connectivity is absolutely essential, too, so CGM, loop, and
76 pump can do their job without losing more valuable time (see [case study 1.4](#)).

77

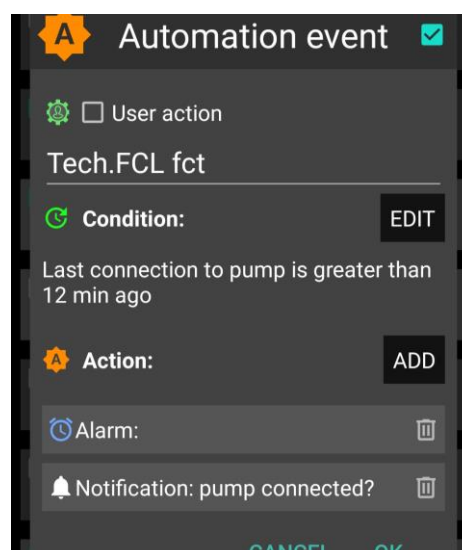

78 Then, but even more importantly in *all other* day and night *times*, the CGM should not produce any
79 artefacts (jumpy values; see [case study 1.3](#)) that the loop could **misinterpret** as sign of a starting
80 meal. Note that also calibrations could produce jumps.

81

82 autoISF has also a couple of in-built checks on the quality of the recent CGM values. Hence, a
83 CGM with more scatter will make the loop lose more time, and lead to higher peaks and
84 lower %TIR.

85 So, if you are unhappy with a slow reaction of your loop it could be because the loop is unhappy
86 with your CGM. Consult the detail info given (at the time) in your SMB tab, or look it up later in the
87 logfiles (using the emulator eventually).

88

89 The best proven way to stay out of trouble currently is to use Dexcom G5 or **G6**, and to ensure via
90 **overlapping** right and left arm sensor and transmitter utilization for always good quality values that
91 can be used by the Full Closed Loop ([case study 1.5](#)).

92 Other ways (making use also of day-1 sensor values, G7, Libre2/3...) are possible, but come with
93 a lot of monitoring effort (best via watch) and occasional time-outs for the FCL.

94

95 One safety feature in autoISF is a **blockage of SMB delivery whenever delta bg** (within the last
96 two 5 minute values) is **higher than 30% of that bg**. *So from 74 mg/dl, a jump to 97 (+23) or more*
97 *would not receive SMB "response", or from 100 mg/dl to 131 mg/dl (+31) would neither.*

98

99 Check in your (HCL or FCL) data whether at meals or sweet drinks with rapid absorbing
100 carbs you could run into the problem that jumps are "too high" and much needed insulin will
101 be blocked (only come via very much smaller portions: *e.g. 400%TBR @ 0.6 U/h => 0.2 U in*
102 *5 minutes, instead of one ~3 U SMB. The difference of 2.8 U missed translates @ ISF~ 40*
103 *mg/dl/U into up to + 112 mg/dl higher bg peak! It will not become quite that bad, because*
104 *the loop will catch up to the insulinRequired with it's next couple of decisions).*

105

106 Instead searching in old data, you can also just have an eye on instances where you think a
107 first SMB was due, but blocked. Confirm that by looking in the SMB tab, and think about a
108 solution that would not require changing the 30% safety limit in the code. *For instance, not*
109 *drinking so much juice rapidly around meal start could be a likely "behavioral" correction to*
110 *get rid of the problem.*

111

112 This blockage (no SMBs) would likely last only 5 minutes (and go probably unnoticed - **in a**
113 **suggested improved User Interface it would be indicated for about 5 minutes via a dotted violet full**
114 **closed loop logo/button on the AAPS home screen**). However, not only would you lose 5 valueable
115 minutes to get your job substantially elevated; all following deltas are likely much smaller, and, as a

116 consequence, you will miss some of the boost sought from bgAccel_ISF if the >30% delta was in
117 fact (largely) due to carb absorption.

118

119 This example also underscores that the CGM in use cannot be allowed random scatter that leaves
120 no reasonable room for safe detection of (smaller and) bigger “truly carb related” deltas

121

122 If or when (like: first half day of a new sensor) you are not sure about sufficient CGM
123 performance you might develop for yourself an Automation with User action ticked (along
124 the lines as used for other purposes in [section 5.2.2.3](#)). It would “ask you” before giving a
125 SMB whether you really want it delivered. That way you can a) have a look on your glucose
126 curve b) on the delta and acc data underneath the TT field of your AAPS home screen c)
127 think about what sense a SMB now makes with respect to your last meal, and the carbs to
128 be still absorbed. Ultimately, you could also d) consult some of the detailed info given
129 (every 5 minutes) in your SMB tab.

130 Such User action Automations need not be active at all times, but if you have it for
131 your first half day of a new G6 sensor for instance, you could activate that
132 Automation from your list of Automations; after the values have settled in, you can
133 deactivate (“shelve”) it again.

134 For a brief period, and if you are tech savvy, another way to deal with uncertainty about
135 CGM would be to employ the emulator method as presented in [section 11](#): Run a “too
136 mildly” tuned FCL, and in parallel run a “what-if” with your more aggressive settings that you
137 really would like to use once you are certain about your CGM.

138 However, I found it easiest to lay a solid groundwork by using 1 Anubis, and 2 overlapping
139 G6 to get rid of most problems that I saw (and keep seeing, on the worse sensor of the two
140 running for some days often in parallel) in my data.

141

142 With a sensible iobTH defined, and your standard alarms for going towards a hypo not silenced,
143 the worst consequence from any automatically “over-treated” glucose jump should be that you
144 need an unplanned snack for the balance of “missing” carbs.

145

146 Also FS Libre 3 is useable in the autoISF FCL context. Observe info in the general section about
147 autoISF (<https://github.com/ga-zelle/autoISF>) regarding established ways to use with autoISF.

148 As of Dec.2023, there is still development work going around what to best make of the 1
149 minute values (it might help recognize acceleration at meal start on average 2 or 3 minutes
150 earlier).

151

152

153

154 1.5 Meal-related limitations

155

156 Setting up a full closed loop is relatively easy for people whose diet does not consist **mainly** of
157 components with rapid high effect on blood glucose (more see

158 <https://androidaps.readthedocs.io/en/latest/Usage/FullClosedLoop.html#meal-related-limitations>)

159

160 Meals do not have to be low on carb (provided you use a fast insulin for your FCL)

161 Fat or protein rich diets, or slow digestion/gastroparesis, make things easier rather than harder for
162 the full closed loop because late carbs nicely cover for inevitable “tails” of late action from SMBs
163 needed around peak time.

164

165 Erratic consumption of snacks with fast resorbing carbs can be a problem.

166 In autoISF you can reduce this problem to some extent via one or two keystrokes from your
167 AAPS home screen. While certainly being a deviation from the FCL idea(l), this would be
168 one of the exceptional situations where you better do a quick “nudging” step from your “FCL
169 cockpit”. Details see in [section 5.2.1](#) and [5.3.3.1 \(4\)](#) and [case study 5.2](#)

170

171

172 1.6 Lifestyle-related limitations

173

174

175 Technically stable system

176

177 Full closed looping requires a 24/7 technically stable system, especially regarding

- 178
- reliable CGM signals
 - 179 • Bluetooth stability with the pump (see [case study 1.4](#))
 - 180 • keeping your phone in sufficient proximity at all times
 - 181 • avoiding (or at least early recognition of) occlusion.

182 This requires a habit (or, unlikely, permanent attention to details) like keeping all components well
183 charged and in close proximity; making cannula (or pod) changes always early enough to lower the
184 risk of occlusion (see [case study 1.1](#)); having always potentially needed parts with you.

185 **Depending on your system, your experience with it, but also on your acceptance and general**
186 **lifestyle, these aspects may or may not limit you.**

187

188

189

190 Preparing for exercise

191

192 To prepare for exercise (sports, heavy work), the normal protocol with a pump or hybrid closed loop
193 is to take actions that reduce insulin on board prior to exercise

194 With your full closed loop, the algorithm is tuned to detect meals and to give you insulin to counter
195 glucose rises automatically. Setting a high temp. target and lower %profile right away (effective al-
196 ready around meal start) could be a problem.

197 Unusual activity levels therefore likely require **disciplined preparation** (especially **if you want to**
198 **keep the need to snack during sports low**)

199 In autoISF you can reduce this problem to some extent via two or three keystrokes on your
200 AAPS home screen. While certainly being a deviation from the FCL idea(l), this would be
201 one of the exceptional situations where you better “flick a lever” from your “FCL cockpit” to
202 keep iob low (example see [case study 6.2](#)).

203

204 Extra hurdles to establish FCL for kids

205

206 To establish and maintain a FCL for kids brings about some extra challenges if:

- 207
- Lyumjev is not available or well tolerated
 - Hourly basal rate is very low, providing a poor basis for big SMBs
 - Diet is rich in sweet components. With the typical low blood volume of a small body, strong
210 tendency towards very high bg spikes!
 - Going through marked changes of insulin sensitivity or of circadian pattern makes it difficult
212 to keep the FCL appropriately tuned.

213 This problem is about the same in Hybrid Closed Looping. However, now you might
214 expect miracles from the FCL. This is not going to happen. You still should try to set
215 appropriate temp. changed profiles, that serve also as a basis for your autoISF FCL.

- 216
- Discipline is poor regarding keeping Bluetooth connectivity and infusion sites perfectly run-
217 ning
 - Between kid and supervising parent it must be guaranteed, especially in the initial weeks,
218 that an eye is kept on whether the FCL is working about as to be expected.
219

220 More see [section 7](#).

221

222

223

224

1.7 Time required for setting-up

Lastly, before enjoying a functioning full closed loop you need to have a period of a some weeks with some free time and „free head“ for set-up –. Can you get, in the time you are willing to invest, to a result that you consider good-enough is really the question. Depending on your „habits“, and which – if any - compromises (like doing cannula/pod changes more often, never starting meals when bg sits high ...) are you willing to make (and everyday able to stick to), for the ease of not having to deal with assessing meals and bolussing for them?

While setting up your personal FCL using autoISF is a substantial project, there is no need to implement it fully in one step. There is nothing wrong to go in your well running Hybrid Closed Loop mostly, while switching to FCL only for dinners, for instance, or only for weekend lunches, as a start. Once you found feasible settings, you can expand to other meal times, and lastly towards figuring out your best strategies for challenges outside of meal windows, as we shall discuss in [sections 5. and 6.](#)

There are alternatives to using autoISF for FCL, as well. See [section 13.](#) for more info.

Notably [13.1](#) FCL using AAPS Master and Automations could be a much easier and more error-tolerant way of stepping into FCL. In a clinical study with 16 participants about 80% TIR was achieved without much tuning effort (source: see in [section 13.1](#)).

To close the circle to where we had started ([section 1.1](#)): **A very time consuming pre-requisite might actually be to *first sort out your Hybrid Closed Loop*, so your profile parameters are set „right“, and your “old” data really can serve as a blueprint for what, now, you would like *your loop to do* in FCL mode**

Note that if you had used dynamic parameters or special Automations („loops inside the loop“) this might have balanced some principal errors, but leaves you now without a good starting point as you must get rid of these over-patches (see also warnings at start of [section 4](#))..

You will see also success stories of loopers who just jump into using more powerful tools, in kind of a trial and error mode, and frequently add the latest add-on, or self-constructed patch (often in form of an Automation) to counter-balance problems.

So, yes, you can also continue in that spirit. Resulting solutions may be good-enough. But they tend to be unstable and not well-understood. That is a poor basis for managing arising problems (-> fine tuning), and for adjusting to special situations (-> which setting to temporarily change). But it certainly is an alternative avenue for the impatient, less analytically, and more adventurous inclined.

264 In any case, PLEASE always observe the safety settings/instructions coming with the DIY
265 dev- variant of software you select.
266

267 One key safety measure every AAPS user going towards FCL should have in place is to set an **iob**
268 **threshold** (iobTH; size a bit below what you used as a bolus for bigger meals in HCL) above which
269 no more SMBs can be given by your FCL.
270

271 This is an integrated feature of autoISF, from 3.0 version onwards (see [section 2.4](#)).
272

273 Other FCL methods may require to set up an Automation for a temporary iob threshold that
274 blocks SMBs from being delivered, see e.g. here for AAPS FCL w/Automations:
275 <https://androidaps.readthedocs.io/de/latest/Usage/FullClosedLoop.html#iob-threshold>