

## 5. Modulation of autoISF aggressiveness.

V 2.61

**Please note that with autoISF 3.0 you are in an early-dev. environment,** where the user interface is **not optimized for safety** of users who stray away from intended ways to use. Good safety features exist, but these are only as good as the development-oriented user understands and implements them. This is not a medical product, refer to disclaimer in [section 0](#)



### 5.1 Automatic modulation of loop aggressiveness

- 5.1.1 “autoISF off” outside of meal times
- 5.1.2 SMB off @ odd profile target
- 5.1.3 SMB off @ odd temp. target
- 5.1.4 diff. of FCL aggressiveness via Automations
- 5.1.5 diff. of FCL aggressiveness via Activity Monitor
- 5.1.6 Pro/con completely hands-off FCL

### 5.2 Manual modulation of FCL aggressiveness (DIY cockpit)

- 5.2.1 Status recognition
- 5.2.2 Manual interventions from DIY cockpit
  - 5.2.2.1 Temp. %profile or TT settings
  - 5.2.2.2 Temp. settings in /preferences
  - 5.2.2.3 Grey DIY cockpit buttons for FCL responses
- 5.2.3 Temporary exit from FCL

### 5.3 Manual modulation via improved cockpit

- 5.3.1 Violet FCL icon and underlying buttons
- 5.3.2 Bottom buttons “insulin” etc.
- 5.3.3 Top three fields
  - 5.3.3.1 TT dialogue field
  - 5.3.3.2 Exercise button / dialogue field
  - 5.3.3.3 Profile dialogue field

### 5.4 Recognizing loop state from the AAPS home screen

- 5.4.1 Color scheme (grey/yellow) of the top 3 fields
- 5.4.2 Info on the top 3 fields (profile, exercise, TT)
- 5.4.3 FCL related indicator fields
- 5.4.4 Overall AAPS home screen

[Available \(related\) case studies:](#)

Case study 5.2: Sweet snack, juice snacks ...../ or Glühwein (5.2.2.3)

Once the initial tuning according to [section 4](#), is done, you are ready to use autoISF for your fully automated meal management.

You will have three major *other* challenges to manage:

- recognize and manage (partial) occlusions, or other technical (CGM or BT related) obstacles (see [section 1](#) on pre-requisites of FCL, and related case studies)
- deal with times when insulin given by the loop must be restricted (e.g. a snack could be “misinterpreted” as a meal)
- deal with times when the loop should be set “milder” as a precaution (e.g. nights; or in an exercise context).

46 How big the remaining challenge is depends very much on your individual lifestyle.

47 [Sections 5](#) and [6](#) discuss this in more detail.

48 **In order to run the loop fully automatically around the clock, the times *outside* the meal**

49 **blocks must also be precisely analyzed, and solutions to problems (if any) must be sought.**

50

51 It is up to every user to decide where to draw the line.

52 • With a technically well functioning system, moderate meals, moderate or no exercise,  
53 moderate %TIR expectations and a bit of mindfulness it should be possible to go into **Full**  
54 **Closed Loop 24/7**, after working through, and observing, [sections 1-4](#).

55 • Especially if you are a bit shy of using the emulator for really detailed analysis, it is likely  
56 that you will not hit *one* real good system calibration ([section 4](#)) for your *entire range* of  
57 diets.

58 In that case you will occasionally run out of range, and your options to prevent, react, or  
59 improve are

60 ○ accepting a few % higher time outside range for that day (and, if feasible, in the  
61 future avoiding what seemed to have caused it)

62 ○ taking a snack (whenever you tend to go low from the “tails” of insulin activity that  
63 was required to fight a peak)

64 ○ doing a manual override (if you can think of one in time, to manage the problem  
65 manually)

66 ○ temporarily resorting to the well-known hybrid closed loop.

67

68 Instead of accepting such instances, you could launch “improvement projects”

69 • that refine your initial tuning ([section 4](#). and [sections 8-9](#))

70 • that make you and your FCL loop fit to manage an increasing number of disturbances  
71 either automatically, or via an “informed”, maybe pre-programmed, user intervention  
72 (notably, an exercise “announcement”) ([sections 5-6](#)).

73  
74 To tailor the loop's response to disturbances *other-than* your major meals probably will require  
75 specific **modulation of the aggressiveness** that you set according to [section 4](#) for your meal  
76 management.

77  
78 There are many avenues to achieve this. The main ones, that are also easy accessible via  
79 Automations in AAPS, are:

- 80
- 81 • temporary shut-off SMBs (odd-numbered target)
  - 82 • temporary reduce bgAccel\_ISF-weight
  - 83 • temporary reduce iobTH\_percent
  - 84 • temporary reduce set %profile
  - 85 • temporary set higher TT (especially in connection with exercise mode)t

86  
87 After set up of your core FCL for fully automatic meal management according to [section 4](#), and  
88 some performance monitoring ([section 8](#)) and trouble shooting ([section 9](#)) you probably will identify  
89 areas that you like to further improve (notably to cut down the occasional need for snacks in an  
90 exercise context).

91  
92 In setting up your FCL, you then have another difficult and time-consuming job at hand, to define  
93 solutions for any of your „other“ situations (outside of meal management) that tend to drive glucose  
94 outside of the desirable range.

- 95
- 96 • In [section 5.1](#) we explore avenues towards *fully automated* management that in daily life  
97 will require no user intervention at all.
  - 98 • In [section 5.2](#) and [5.3](#) we will look at solutions that involve an easy *user interaction like a*  
99 *data entry or button push*.

100  
101  
102  
103  
104  
105  
106  
107  
108  
109  
110  
111  
112  
113  
114  
115  
116  
117  
118  
119  
120  
121  
122  
123  
124  
125  
126  
127  
128  
129  
130  
131  
132  
133  
134  
135

## 5.1 Fully automatic modulation of FCL aggressiveness

The following subchapters describe set-ups you may want to use for allowing **completely hands-off FCL in as many daily situations as possible.**

### 5.1.1 autoISF ISF adaptations generally switched off outside of meal-time windows

If, aside from having to bolus for meals, your hybrid closed loop was running pretty well *without* other interventions from your side, you could continue to run in that mode, and just focus your new autoISF FCL on management of meals.

In your initial transitioning phase this approach makes a lot of sense, and even by focusing autoISF on just a sub-set of them, like only dinners.

Also in the long run this avenue is taken by many FCL users for the night times, “hanging on” to their well performing hybrid closed loop with standard `oref(1)` SMB+UAM

For this, you define Automations

- that set meal time windows in which “Enable ISF adaptation by glucose behavior” (autoISF) is turned on in AAPS preferences/OpenAPS SMB
- or: that turn *all* autoISF’s ISF modulations (*or* just `bgAccel_ISF` ) off in time windows in which surely no meal occurs. For instance, you can go for all nights back into your Hybrid Closed Loop, as you had before.

Other early DEV AAPS variants (see [section 13.3](#)) all work with meal-time windows. The window is either set by time of day in the settings, or it always must be „set“ by the user via giving a mandatory small pre-bolus before any meal starts. **Outside** of these time windows, these loops then runs with less aggressive SMBs like `oref(1)` SMB+UAM in AAPS Master. This mode is not really FCL, but an advance over traditional HCL that often achieves satisfying degrees of automation and performance.

The term **Meal Announcement** (MA) is often used to label this closed looping mode. Trigger to set a meal time window could also be a pre-bolus given by the user, a carb entry made, an EatingSoonTT set, or a meal announcement button pushed.

136 Note: Outside of the meal time windows you would be in hybrid closed loop. To the extent you  
137 rarely face disturbances (aside from meals), you could be looping in full automatic mode around  
138 the clock,

139

140 Your temp. "autoISF shut-down" (exiting autoISF FCL = shutting off "Enable ISF adaptation by  
141 glucose behaviour") is meant to prevent problems from the loop *over-reacting* to bumps in the  
142 glucose curve in times of day (night) when standard oref(1) performance is sufficient.

143

144 A very good alternative to fully resorting to Hybrid Closed Loop is "taming" the FCL via a night time  
145 SMB shut-off (see next [section 5.1.2](#)).

146

#### 147 5.1.2 Odd-numbered *profile* targets used to block SMBs

148

149 An alternative route of preventing the FCL loop from over-reacting to bumps in the glucose curve  
150 would be to make use of the option to temporarily shut down SMBs

151

152       Ensure the even/odd logic in the settings is toggled on in Preferences> openAPS SMB>  
153       autoISF settings> smb delivery settings>: "Enable alternative activation of SMB depending  
154       on profile target" ON.

155

156 In time blocks with an odd-numbered profile target you can prevent any SMBs being given by your  
157 loop. The (unchanged) aggressive settings then can only translate within the limits set by %TBR  
158 possible.

159

160 This will very much slow down any more insulin being given, and is an excellent solution for night  
161 times, especially if you occasionally experience compression lows.

162

163       Alternatively, you could use the new included options for Automation Conditions and  
164       temporarily tune your bgAccel\_ISF\_weight much lower ([section 5.1.4](#)).

165

166       Yet another alternative was already presented ([section 5.1.1](#)) = to go into hybrid closed  
167       loop for the night.


168       That is possible to do with SMBs available (without them getting boosted via autoISF), and,  
169       for a long time, was the author's favoured solution for the nights.

170

171       But, my current favourite builds on the method of this section (5.1.2, odd profile target  
172       provides SMB shut off), but then allowing some, automatically triggered when needed:

173

174 In case you occasionally do have nights that would benefit from a couple of SMBs (to treat temp.  
 175 highs from a late fatty pizza, raclette and such): Define suitable Automations like the two „night“  
 176 ones in this *example*:

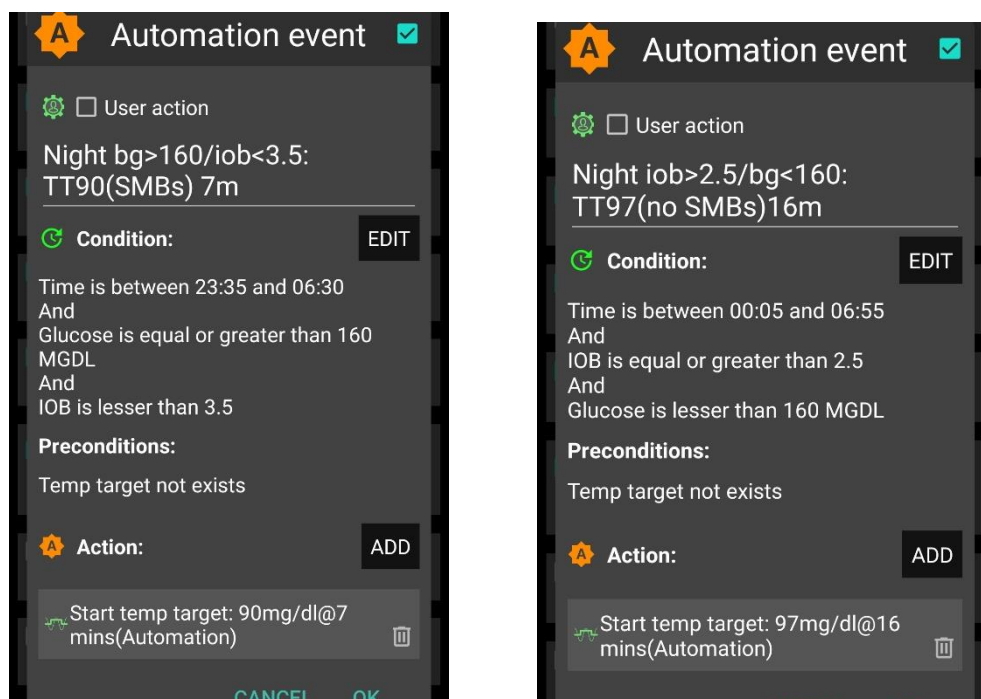

177  
 178  
 179 Never underestimate the „trickyness“ of getting your Automations „right“.  
 180 With your thought-out Automations in place, night data need to be analyzed to see  
 181 • whether the bg and iob limits defined in the given example work sensibly for your data  
 182 pattern  
 183 • whether the TT duration is chosen appropriately  
 184 • how swapping the sequence in which the automations appear in the Automation list would  
 185 lead to different SMB impacts.

### 187 5.1.3 Odd-numbered *temp.* targets (TT) set via Automation used to block SMBs

188  
 189 A widely used Action that strongly modifies how fast your FCL can add more iob is setting an **odd-**  
 190 numbered **temp. glucose target** which makes the loop operate without giving any SMBs (%TBR  
 191 modulation only).

192 Ensure the even/odd logic in the settings is toggled on in Preferences> openAPS SMB>  
 193 autoISF settings> smb delivery settings>: "Enable alternative activation of SMB depending  
 194 on TempTarget" ON.

195  
 196 So, from patterns you find in YOUR data, at times where you want your loop act differently, you  
 197 need to carve out Conditions that describe the respective situations (and either *for how long* it  
 198 typically lasts, or at which *other* Conditions you want your loop get back to default FCL operation).

199 An odd TT is often set for an *anti-hypo* snack or *sports* snack. In both instances, you do not want  
200 SMBs to quickly counter act.

201 In case of *sweet “fun”* snacks, this is entirely different -> [section, 5.2.1](#) or for regular snacks  
202 (*e.g. at school break*) see next [section 5.1.4](#)

203

#### 204 5.1.4 Automatic differentiation of FCL aggressiveness using Automations

205

206 **Personalized Automations** tailor the loop exactly to **YOUR** data so **fully automated** handling of  
207 situations with **different aggressiveness** of the loop can be made.

208

209 From, autoISF 3.0 onwards, also the following parameters are provided as Condition and/or as  
210 Action for defining YOUR Automations:

- 211 • Enable ISF adaptations by glucose behaviour => Allows temp. ON/OFF for the key ISF  
212 modulation parts of autoISF (and, as a result, will usually decrease loop aggressiveness)
- 213 • Trigger/set iobTH percent => Keeps default aggressiveness, but only until a iob threshold  
214 (that your Automation modifies) is surpassed (which is when any further SMBs will be  
215 blocked blocked)
- 216 • Trigger/set bgAccel\_ISF\_weight => Modifies the default aggressiveness of just the  
217 acceleration component

218

219 To set up suitable Automations, you first must **analyze patterns** you find in **YOUR data**, at times  
220 (or geo-locationa, or bg and iob patterns that point to a problem ...) **where you want your loop**  
221 **act differently**, to carve out Conditions that describe the respective situations (and either for how  
222 long it typically lasts, or at which *other* Conditions you want your loop get back to default FCL  
223 operation).

224

225 A variant of this mode is to define several windows in which autoISF aggressiveness  
226 (bgAccel\_ISF\_weight) and/or iobTH are automatically set differently

- 227 • for **different meal time slots** of your day –

228 (*Breakfast at home, school lunches, school intermission snacks, dinners at home* could for  
229 example all deserve special settings regarding ISF\_weights and iobTH).

- 230 • or even for a geo-location etc –


231 (*School lunches, or mother-in-law visits, would be examples*).

232 An example for this was given in section 3 already:



Here is an example set of automations to alternate between two values of iobTH:

I use two different values of *iob\_threshold\_percent* during a normal day. It is 40% for lunch time and 60% for dinner time. I have these two rules to switch by time of day and only if the current value equals the value from the earlier shift. Any other value is treated as a manual override for special occasions until I manually set it to its regular value. The time windows for switching are long enough to catch an opportunity to be processed and do not need to be actioned half a day each.

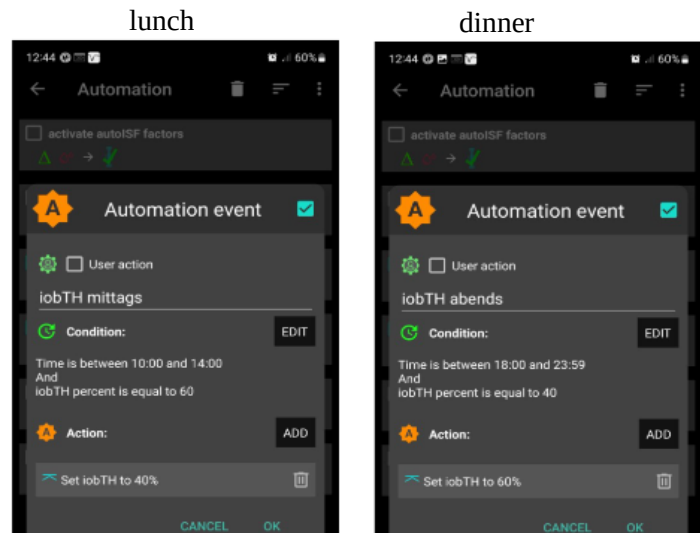


Unless your meals differ vastly in size and in fast carb content all this may not be needed.

Still, personalized Automations might help ease your initial job of setting the various ISF\_weights, and a best-suitable iob\_theshold\_percent that would work “always”.

**Caution:** If (as in autoISF 3.0) setting a different iobTH or bgAccel\_ISF\_weight can not be done temporarily (i.e. with a duration attached) you **must** define a suitable **additional Automation**, that must be active in tandem, that **restores the default set iobTH or bgAccel-ISF\_weight again**. Else, once your Automation set in, it will *forever* shift this important parameter setting!

If for instance you have several Automations that, in combination with a set elevated TT also set a lower iobTH: Don't be fooled, the duration only applies to the TT. You need an extra Automation for all of them. I picked out the highest of the altered iobTH values that these Automations can set (45\_percent), and then I can automatically restore my default desired 60% via this one Automation (see screenshot - - > )

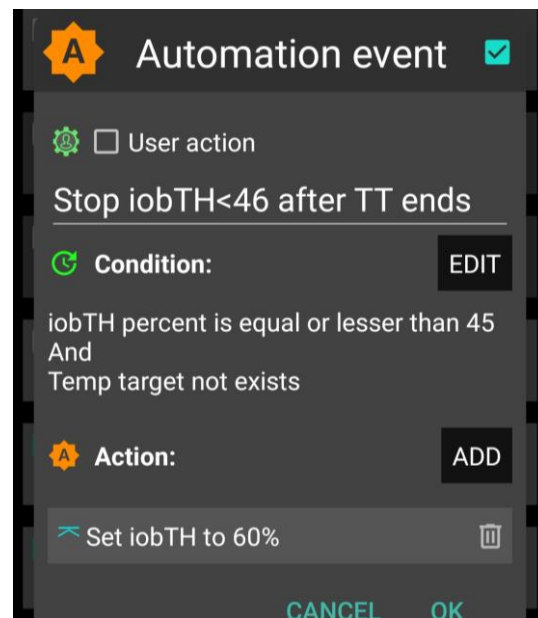

#### 5.1.5 Automatic adjustment of FCL aggressiveness via the Activity Monitor

With the autoISF variant of AAPS you can make use of your smartphone's **stepcounter** and use it to fully automatically adjust insulin sensitivity ratio to **activity level in the past minutes to one hour** time frame.



260 This feature comes with yet another little tuning opportunity, in which you study your body's  
261 response to light exercise (like walking) or to not moving at all (like desk, couch), and select  
262 appropriate settings which, in the future, will automatically adjust insulin delivery to suit activity  
263 state of the past minutes (up to 1 hour).(AAPS Preferences/OpenAPS SMB/Activity modifies  
264 sensitivity/ -> set **two scaling factors**.)

265  
266 This autoISF feature (new since V.3.0) is much quicker responding than Autosens or dynamicISF  
267 to adjust insulin sensitivity to your current „lifestyle state“.

268  
269 For loopers who do not have huge variations in exercise levels in their everyday lives, this feature  
270 might fairly much close the gap towards being able to do a 24/7 hands-off FCL.

271  
272 [Sections 3.5](#) and [6.5](#) describe the Activity monitor in more detail.

273  
274 [5.1.6 Pro/con completely hands-off Full Closed Loop](#)

275  
276 To stay 24/7 in a completely „hands-off“ FCL can be a realistic goal with autoISF 3.0 if besides  
277 meals also some special challenges, as discussed in this [section 5.1](#), were analyzed and could be  
278 addressed.

279  
280 Clearly it depends very much on your lifestyle, and how interested, willing, and capable you are to  
281 recognize, deal with, (and in the future avoid?) situations that get you outside of your desired %TIR  
282 on occasion.

283         So, this is also about what %TIR you are aiming at, and can accept, as it averages out for  
284         the week, for instance.

285  
286 Everybody must weigh for her/himself

- 287         • how much **upfront effort** to put into the setting up process for getting it all 100% automatic
- 288         • **or whether to take an easier start, with a couple of situations left to take care of when**
- 289         **and as they arise in daily life**

290  
291 Also, even if a principal capability for a fully automatic running FCL is given, this still  
292 means that

- 293         • the user should be knowledgeable about what exactly is going on, and
- 294         • have a principal capability to „nudge“, or to take over in a manual mode.

295

296 In the sections that immediately follow, we present the options to nudge or temporarily take over  
297 from the AAPS home screen which will be serving as your **FCL cockpit**:

- 298 • [Section 5.2](#) describes how you can use available “buttons” from your AAPS home screen,  
299 and how to complete it towards a suitable DIY FCL cockpit, for an even better FCL  
300 experience.
- 301 • [Section 5.3](#) describes how you might be able to manage “disturbances” even better (with  
302 more convenience) with an **improved FCL cockpit in the future**. **(No need to read any of**  
303 **the green lines, unless you are interested in contributing to define/design/program further**  
304 **improvements)**

305

## 306 5.2 Modulating aggressiveness manually, from the DIY-FCL-Cockpit\*

307

308 \* Like in the airplane cockpit: Cruising in full auto mode should involve having an eye on the  
309 instruments, and on potential disturbances ahead in the environment.

310

311 In [section 4](#), we dealt with major meals. In [section 5.1](#) we looked into fully automatable manage-  
312 ment of other situations. Life in Full Closed Loop can become extremely easy then...

313

314 However: Other **disturbances** might come up, that:

- 315 • are not noticeable in-time, or foreseeable, by the loop (*e.g. your plan to start exercise in an*  
316 *hour or two*), but **that influence sensitivity dramatically** and therefore require temporary  
317 non-default settings in order to remain in-range, and/or
- 318 • **require** a different “starting point” regarding iob and bg, which translates into **a different**  
319 **iobTH** that should **temporarily** be set much lower (*in case of exercise*) or noticeably higher  
320 (*e.g. with very fast absorbing carbs in a sweet snack “sin”*) .

321

322 In [section 5.1](#) we looked into ways to automate also a modified loop response to *foreseeable* situa-  
323 tions (tied to a time of day, geo-location etc), or to those *the loop could recognize* (with enough  
324 time to react).

325

326 Other „disturbances“ might come up, and you must find an easy way to

- 327 • call up a pre-programmed routine for automatic management, with adjusted  
328 aggressiveness, or:
- 329 • manually tweak a setting or two, to temporarily adjust the aggressiveness

330 • There may also arise a desire to just exit the FCL mode, and be your own captain for  
331 mastering a special situation.

332 For peace of mind, to learn, and to stay informed (especially so in your initial tuning phase, or  
333 when your glucose curve goes in unexpected ways) we also must be able to

334 • find the key parameters that frame and drive the recent and upcoming loop decisions.

335  
336 All this is facilitated within seconds right from the AAPS home screen, serving as a **FCL cockpit**  
337 after you built a couple of DIY cockpit features via Automations (as described below and in [case](#)  
338 [studies 5.2](#) and [6.2](#)):

339  
340 Thoughts went also into [how to improve the cockpit in future releases](#), see [section 5.3](#)

341

342

### 343 5.2.1 Status recognition

344

345 Before considering any manual interventions into the ongoing FCL, you should be aware what the  
346 current mode of action is, and hence which button eventually to fine-tune or lever to switch, in or-  
347 der to adjust to the disturbance you see coming up.

348

349 See [section 5.4](#)

350

### 351 5.2.2 Manual interventions from the (DIY-) FCL cockpit

352

353 Trouble with all these is, not to forget to set back manually, too (=> better solutions in 5.3)

354

#### 355 5.2.2.1 Temporary tuning of FCL aggressiveness via temp. %profile or TT settings

356

357 The set **% profile** multiplies with both, the ISF resulting from autoISF, and also with the default  
358 iobTH you have set, so both are nicely modulated in a linear way with the % temporarily chosen

359

360 Just taking profile e.g. to 110% for an afternoon might be an easy way to explore whether  
361 you might benefit from 10% more “aggressiveness” in your core settings for lunches (like  
362 bgAccel\_ISF\_weight). Make sure, though, that the extra 10% are not cut away by set safety  
363 limits.

364

365 A lowered (relative to profile glucose target) temporary **bg target (TT)** signals lowered sensitivity  
366 (more insulin need), and an elevated TT (as often used with exercise) increases sensitivity and  
367 hence works in the direction of a lowered % profile to also reduce insulin given by the loop.

368

369 Moreover, the **exercise button** (top center on your AAPS home screen) can be activated (turns  
370 yellow, then). This will **further boost** how your set TT elevates the resulting ISF, and sharply  
371 lowers iobTH, as often desired for sports. See [section 6.1](#)).

372

#### 373 5.2.2.2 Making temporary changes in settings made in AAPS/preferences/Open APS SMB

374

375 Going into AAPS/preferences/Open APS SMB allows to:

- 376 - set milder or stronger ...\_ISF\_weights
- 377 - set different iob\_threshold\_percent (or iobMAX)
- 378 - elevate or lower the SMB\_delivers\_ratio
- 379 - limit or expand max. allowed SMB size
- 380 - change the the even <-> odd logic for SMB on/off

381

382 Doing temporary changes in AAPS/preferences should be the exception because

- 383 - they require multiple steps, including entering a password
- 384 - you will often forget to set everything back to default a couple of hours, or minutes, later

385

#### 386 5.2.2.3 Triggered Automations: Grey extra DIY cockpit buttons for pre-programmed “responses”

387

388 Recognizing conditions for fully automatic handling by the loop may not be not possible, or come  
389 too late for the loop to act on. Examples would be

- 390 • *exercise*: Minimum an hour before starting “the loop should know” to be able to lower iob  
391 and elevate bg by the time exercise starts.
- 392 • *snacks*: High carb snacks, sweets, consuming ice cream or having a sweet drink comes  
393 with the problem of even steeper glucose rises, but overall a lesser insulin need, compared  
394 to major meals (for which we tuned our FCL according to [section 4](#)).

395

396 This not necessarily implies that snacks need different settings than a meal. After all, autoISF  
397 was designed to act to all available data, especially to where the developing glucose curve is  
398 headed. So, depending on your effort to set parameters for a broad variety of meals (notably:

399 how well you avoid to invariably bounce fast against your iobTH), you might be able to accom-  
400 modate low carb, snack, and major meals with *one* set of settings.

401

402 In FCL autoISF, this is a bit more difficult than in HCL autoISF applications, because FCL  
403 involves revving up iob supply (largely via big bgAccel\_ISF-weights) often too much to be  
404 balanced by just a snack getting absorbed.

405

406 For that reason, or just for increased comfort and safety, you might want to differentiate, and make  
407 use of what follows for the *sweet snack* example, [case study 5.2](#).

408

### 409 **Tuning aggressiveness**

410 Key is that a sweet snack likely benefits from even more aggressive initial FCL  
411 performance than the meals in your normal spectrum of diets require.

412 Therefore, you could set

- 413 • a higher **temp. profile%** and/or
- 414 • a temp.elevated **bgAccel\_ISF-weight** (see screenshot of my Automation).
- 415 • a **low temp. target** (76 for instance; this additionally helps maximize the first SMBs  
416 that will automatically be triggered at detection of acceleration)..

417

418 When first defining and testing this Automation, also check:

- 419 • that the safety limits as discussed in [section 2](#) will not block the intended elevated  
420 aggressiveness
- 421 • SMBs will not get outrageously big and iobTH sometimes exceeded by too much

422 Note that “the last SMB” is allowed to overshoot the effective iobTH by 30%


423

### 424 **Limiting iob**

425 For “just a snack”, total insulin need will be lower than for a meal.

426 If you would just have your sweet drink, and your meal-oriented FCL would “attack”,  
427 iob likely would become too high, and a glucose rollercoaster would start, with you  
428 needing to consume more =>

429 If you just have a snack, or drink a glass of juice, you can lower the **iobTH\_percent**  
430 accordingly.

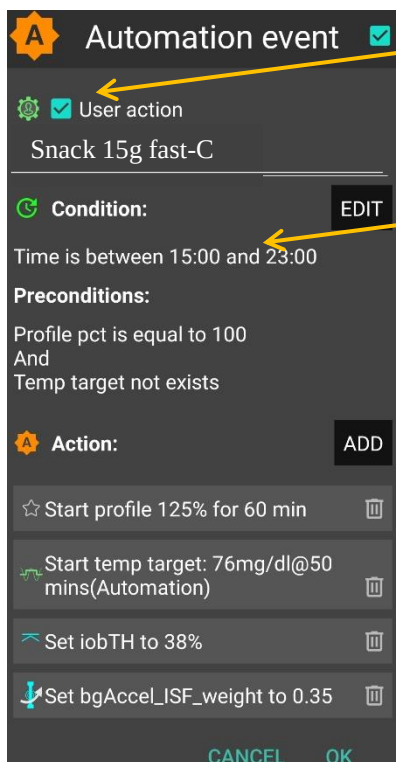


„User action“ is always ticked-on

This will, in the defined time space ..

..offer the “DIY cockpit” button..

..which I must press any time (~90...30 minutes) before my snack.

Note: Make sure that iobTH and bgAccel\_ISF\_weight revert to default afterwards

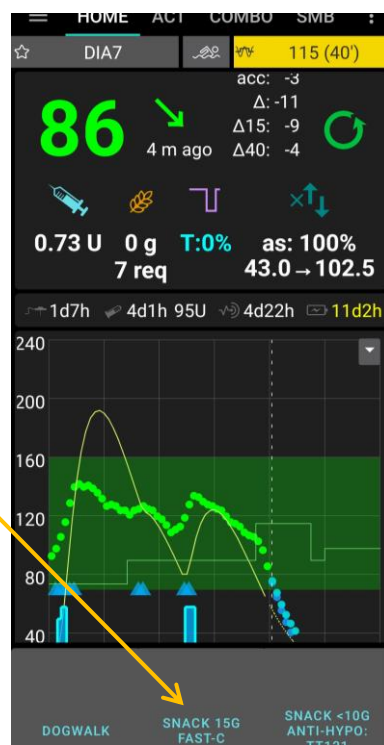

431

432

433 So, this can be a little extra “project” when setting up your FCL.

434 You need to research your snack habits (if any), and over time find out which settings in  
435 the snack-related Automations work well.

436

437 In everyday life you then just must press the related button in your cockpit (which is  
438 not time critical at all, except it should be clicked *latest* a couple of minutes after you took  
439 the drink or snack).

440

441 If you consume more, and also eat something with your sweet drink, this will more  
442 resemble a full meal... however, with unusual amounts of fast carbs.

443

444 Note: Pressing your snack button a *second time* would **not** help because the  
445 lowered iobTH does not allow iob going high enough. So you are better off just  
446 letting your normal FCL meal routine run, after your snack mode expired.


447

448 Other options when snacks keep extending would require a manual modulation  
449 regarding %profile and/or bgAccel\_ISF, but keeping the full default set  
450 iobTH\_percent, or even elevating it (refer to [section 5.2.3](#)). If that happens often,  
451 define for yourself an extra User action Automation for a bigger snack (= another  
452 grey DIY cockpit button).

453

454 **Caution:** If (as in autoISF 3.0) setting a different iobTH or bgAccel\_ISF\_weight can not be done  
455 temporarily (i.e. with a duration attached) you **must** define a suitable **additional Automation**, that  
456 must be active in tandem, that **restores the default**  
457 **set iobTH or bgAccel-ISF\_weight again**. Else, once  
458 your Automation set in, it will *forever* shift this  
459 important parameter setting!

460 If for instance you have several Automations that, in  
461 combination with a set elevated TT also set a lower  
462 iobTH: Don't be fooled, the duration only applies to  
463 the TT. You need an extra Automation for all of them.  
464 I picked out the highest of the altered iobTH values  
465 that these Automations can set (45\_percent), and  
466 then I can automatically restore my default desired  
467 60% via this one Automation (see screenshot - - > )

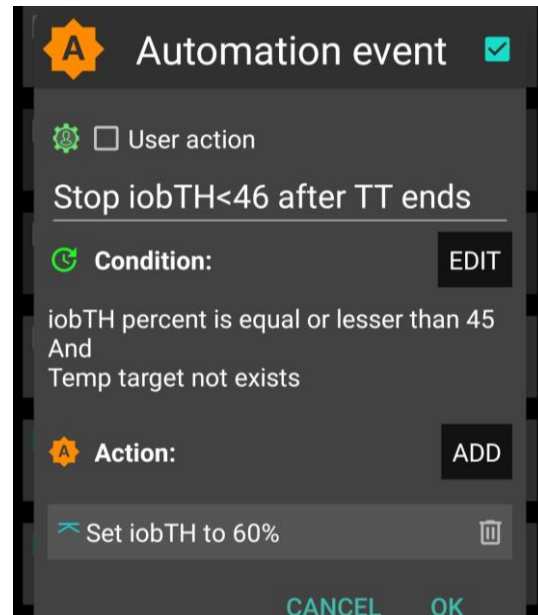

468

### 469 Installing the DIY cockpit button

470 In the related Automation, just keep the "User action" box clicked at all times, and define in the  
471 Conditions when you want to see that button available for cockpit use (see screenshot above) =>  
472 you will see that button offered.

473

474 Besides snacks, also any **other recurring special situations can be addressed via a DIY**  
475 **cockpit button, and receive different aggressiveness up to a suitable iobTH level.**

476

477 Over time you can have a big number of User action Automations, and keep them "shelved" rather  
478 invisibly (clicked in-active, top left box) in your long list of potential Automations. Even when active,  
479 they only show in your cockpit (bottom grey field of your AAPS home screen) in the time slot you  
480 assigned as potentially relevant.

481

482 In the future you might be able to set the stage for a snack and other "disturbances" also via  
483 an extended menu behind the TT button on the AAPS home screen, see [section 5.3.3.1](#)

484

### 485 Discussion

486 In case you do have a snack habit and

487 • can not find settings, as in [section 4](#). defined for your meals, also suit your snacks

488 • can not pin a as in [section 5.1.4](#)



489 you minimum need a “snack announcement” for which the extra button in your DIY cockpit  
490 provides a time-uncritical 1-button-push solution.

491 This could be a good solution for kids in kindergarten, too. Make sure caregivers  
492 understand to use it only once for one snack. Continued snacking would require iob as for a  
493 meals. This is what the FCL loop takes care of automatically, - while using the snack button  
494 several times in a row would limit iobTH at a too-low level!

495 In a software update, we might try to automatically block usage of that type of  
496 Automation for 2 hours, after it was once used.

497

### 498 5.2.3 Temporarily exiting the FCL

499

500 The “last resort” alternative always is to **temporarily** leave the FCL mode, and handle any  
501 disturbance “the traditional way” in **hybrid closed loop**. For this, we switch the automatic  
502 aggressive adaptations of ISF to the bg curve off that are only needed in FCL ....

503 (if in hybrid closed loop you like e.g. the dura\_ISF adaption still, you alternatively could elect  
504 to just set bgAccel\_ISF\_weight temp. to zero, instead)  
505 ... and now around meal starts giving a bolus will be necessary.

506

507 The suggested improved FCL cockpit user interface with an extra version of violet loop on the  
508 AAPS home screen ([section 5.3.1](#)) would facilitate this transition FCL < - > HCL, including  
509 automatic removal and re-appearance of the insulin button at the bottom of the APS home screen.

510

511 In case [this feature](#) is not yet available, you must:

512 Exit the FCL mode by going to AAPS/preferences/put in your password/OpenAPS SMB/scroll down  
513 to autoISF settings and switch “Enable ISF adaptation..” OFF

514 (or, alternatively, set bgAccel\_ISF\_weight to zero).

515

516 **Caution:** Unfortunately, there is no way yet for your full closed loop with ISF adaptations to come  
517 automatically back on, after a selected time for instance. So do not forget to switch your autoISF  
518 fully back on, later.

519

520 As this will often be forgotten, it may be worth doing a “User action” Automation, for a “temp.  
521 FCL OFF” grey button (see [section 5.2.2.3](#)).

522 Caution though, there is very limited experience with this brand new feature. Make sure your  
523 Automation definition really applies a duration (or other condition) that will automatically  
524 terminate all non-default settings it made. As we have seen e.g. in [section 5.1.4](#), this is not  
525 always the case.

526

527 To recognize whether autoISF currently runs with ISF adaptation or not, you must consult the  
528 profile\_sens -> actual\_sens indicator below the Autosens%. However, this gets also modified  
529 by %profile switches or TT +/- exercise mode. So it is not as easy as it would be with the “violet  
530 loop” proposal mentioned already above.

531 Ultimately, you can of course study the SMB tab to find out what is going on.

### 532 533 534 5.3 Modulating aggressiveness manually from the improved FCL-cockpit 535 536

537 autoISF 3.0 is an early dev variant of AAPS, and as user you are participating in an on-going  
538 development. Of note, autoISF 3.0 is first launched without many of the cockpit features that are  
539 presented below in this font color.

540  
541 Only what is written in black is at this point of some relevance for using autoISF 3.0.

542 **No need to read any of the green lines, unless you are interested in contributing to**  
543 **define/design/program further improvements.**

544 This is also an open invitation for you to contact us in case you could help program a  
545 module for one of the suggested user interface extras.

546 For future integration into AAPS Master, an eye should be kept also on the question which  
547 other modes (like FCL using Automations and others mentioned in section 13; and maybe  
548 also HCL) might benefit from some of the extra features.

549 For the time being, multi-step work-arounds may become necessary

- 550 • In many cases, going into AAPS Preferences and changing settings would be needed  
551 (...plus not forgetting to change these settings back, afterwards).
- 552 • Automations allow a DIY FCL cockpit, see section 5.2 and case studies 5.2 and 6.2

553  
554 Keep in mind, though, that the **goal should be to interfere with the loop as little as possible.**

555 Under the described conditions it can run **fully automatically** without any user interaction ( = after  
556 the initial tuning phase, and related settings made in AAPS /preferences/SMB/autoISF. See section  
557 4. and 5.1).

558  
559 However, just like in the airplane cockpit: Cruising in full auto mode should involve having an eye  
560 on the instruments, and on potential disturbances ahead in the environment.

561 *E.g.: storm ahead => instruct your plane to climb to another flight height.*

562 *Analogy: exercise ahead => setting an exercise TT, or => pressing a button that activates a*  
563 *sequence of instructions (some of them probably hinging on conditions, like actual iob) how*  
564 *to manage through that exercise situation).*

565

566 So, for the occasional „disturbance“ coming up, you should find an easy way to

- 567 • call up a pre-programmed routine for automatic management, with auto-adjusted
- 568 aggressiveness, or:
- 569 • tweak a setting or two, to temporarily adjust the aggressiveness
- 570 • There may also arise a desire to just exit the FCL mode, and be your own captain for
- 571 mastering a special situation.

572 All this is facilitated within seconds right from the AAPS home screen's **cockpit features** to the

573 **extent they are already incorporated**, or to the extent you can build alike DIY cockpit features via

574 Automations, as described in [section 4.1.3](#) and [case studies 5.2](#) and [6.2](#)):

575

- 576 • The button that is integrated into the **violet FCL icon** serves as emergency off button, to
- 577 quickly stop FCL, or to at least to immediately stop any more SMBs (...just for a couple of
- 578 minutes, or for the remaining meal time: pick from the options offered with just one
- 579 keystroke).

580 Via the violet FCL icon on your AAPS home screen, you also can access a temp. switch-off

581 button for SMBs (see section that next follows below).

582

- 583 • The **three top fields** (%profile, exercise, TT) provide access to temp. tuning of core
- 584 parameters, and/or to some pre-programmed routines.

585 Taken together with some **new indicator fields** about your loop state ([section 5.4.3 – 5.4.4](#)), and

586 the **grey DIY cockpit buttons** ([section 5.2.2.3](#)) this makes the AAPS home screen your **cockpit**

587 for Full Closed Looping.

588

589 Let us look on each of these cockpit elements in some detail:

590

### 591 5.3.1 Violet FCL icon and underlying buttons

592

593 Novices to FCL, or really anyone running into a very special situation, may appreciate that the new

594 closed loop icon on the AAPS home screen in pink (for FCL) has buttons to quickly shut off getting

595 more SMBs (1st row), or to enter other loop modes (second row).

596

597 It functions very much as the other ones that you know from HCL already, and in fact you

598 get offered some of the same options (for instance, to switch the (full) closed loop off for 15

599 minutes for going to take a shower)


600

601 Note that in FCL you leave all BG regulation, notably against meal spikes, to the loop. So, try not to  
602 disconnect in phases when your FCL must ramp up your job.

603 The required insulin would still be supplied *after* you reconnect. However, without the user  
604 pre-bolussing, the delay would be more of an issue in FCL than it had been in HCL.

605

606 Just pressing on the FCL icon, a dialogue box comes up:

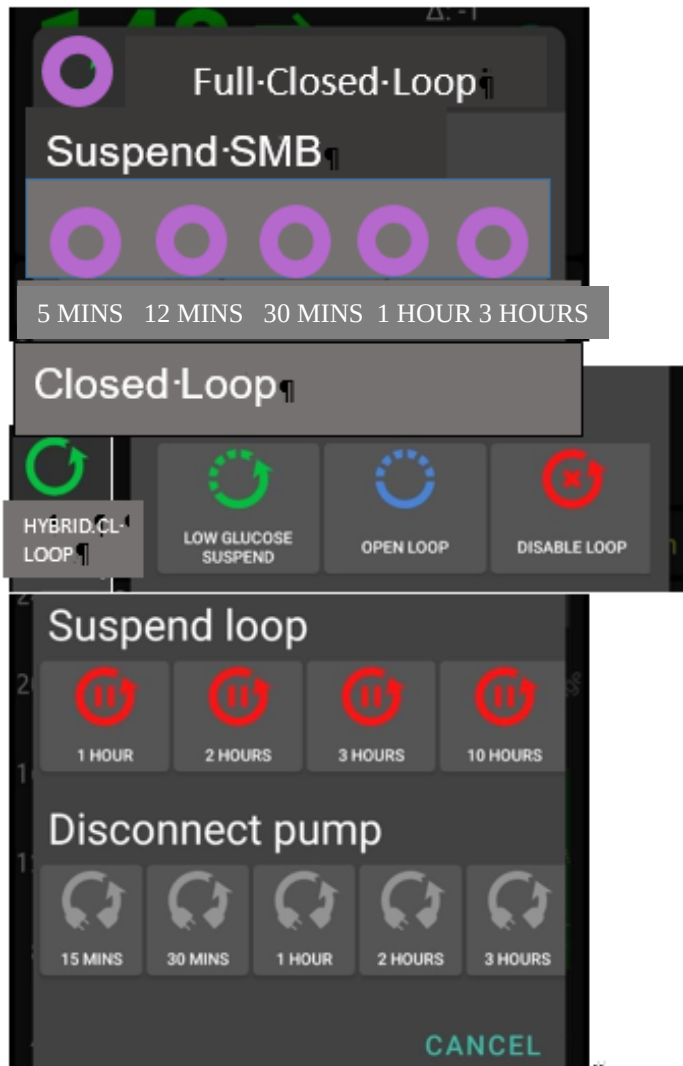

```
<-add-arrow, form-like-like-CL-Loop...¶
```

```
<-dotted+arrow, form-like-like-LGS...¶
```



**DEV:** when in Hybrid-Closed-Loop = top headline + green circle, there is no Suspend SMB part but it starts right with „Closed-Loop“. In that screen, there the 1st element must be option to pick ¶

.....pink-circle, FULL-CL-LOOP ¶

where here the green HYBRID-CL-LOOP stands ¶

For loopers who did not set up FCL, a feedback must come up ~ „FCL not installed“ if they press on that button ¶

607

608 Pressing „**Suspend SMB**“ provides fast and easy „emergency braking“ regarding delivery of more  
609 SMBs:

610 Select the one with the desired number of minutes: 5 or 12 for just blocking the potential next  
611 SMB(s), and up to 3 hours to manage the entire rest of this meal with %TBR from then on.

612

613 Whenever, and why-ever, your FCL is in „no SMBs allowed“ mode (e.g. automatically after  
614 surpassing an iobTH, or triggered by a set odd TT), the FCL icon will turn into a dotted one.

615 Instead of remaining **duration to end time** it indicates in the middle „the condition“, „**iob**“ or „**TT**“

616 Add an indication **if** suspend SMB comes from an Automation, e.g. add an „**(A)**“, **underneath** the

617 #minutes, iob, or TT in the middle of the dotted violet field.

618 So, as in other (already in HCL existing) cases, those icons show in the middle the minutes left that  
619 they will be running, or the condition which would have to go away for this temp. setting to stop.  
620 It always auto-reverts into the FCL state and FCL icon, when time (or other condition) has elapsed.

621  
622 Pressing „**HYBRID CL. LOOP**“ or other buttons from the 2nd row provides fast and easy  
623 „emergency exit“ into other modes.

624 This enables beginners an easy „temp. escape“ into their well-known HCL (green) at any  
625 point of time. bgAccel\_ISF\_weight is set to zero when going FCL->HCL. HCL can run with  
626 autoISF (for instance dura\_ISF) uninhibited otherwise. (check implications for HCL users of  
627 autoISF ?? ).

628 Note: These options from row 2 have no time limit. Loop will not by itself go back to FCL. You see  
629 the different loop icon as a reminder to manually revert, when ready.

630

631

632 5.3.2 Buttons „Insulin“, „Calculator“ etc at bottom of AAPS home screen

633

634 These buttons are **not useful any longer in FCL**, and automatically disappear whenever in FCL  
635 mode (also in Suspend SMB state), and re-appear when leaving FCL. This applies also when an  
636 Automation or technical system failure shut off FCL.

637 Users who, maybe in the beginning phase, feel better having those buttons, can override  
638 the removal (of the insulin button, or any other) by going into /preferences/overview/buttons  
639 and forcing them on. They only remain on until the next re-entry into FCL mode, when auto-  
640 off happens again.

641 The reason why we do this: It really is important to let the loop loop, and not interfere more  
642 than absolutely needed. Any bolus the user gives will sure distort the bg curve, on which  
643 autoISF, especially when aggressively tuned for FCL, builds a lot of its decisions!

644

645 5.3.3. Three top fields (%profile, exercise, TT)

646

647 Depending on the variedness of lifestyle, the desired %TIR, and the initial tuning effort put in, the  
648 user may want occasionally to „tweek“ the aggressiveness of her/his FCL.

649

650 The top 3 fields (grey in default mode, **yellow when temp. in mode with changed**  
651 **aggressiveness**) serve as quick and easy entry points to make temp. switches (as users will be  
652 used to for %profile switches, or for setting an EatingSoonTT in HCL, .. which they still can do in  
653 FCL ... but more:)

654

655 Expert FCL users might need this feature rarely, but probably at least to manage activity after  
656 meals: Each require opposite aggressiveness, and the switch has to come in a certain point in  
657 time that would be difficult to capture. (More see [section 6.4](#) )

658  
 659 5.3.3.1 TT dialogue field ( Currently not available in the pictured form and function ! )  
 660  
 661 The TT field (top right of AAPS home screen) is a primary daily interface, and a dialogue field  
 662 opens when pressing on it

Initially, the form can be entirely empty re. TT inputs and just show the default iobTH and bGAccel parameters.

Two lines appear for target and get labeled ES, AC or HY depending on what was pushed at the bottom buttons. (HY => red frame; evtl. 1. empty => grey)

TT entries can be made or overridden.

iobTH calculates from Target AC and other settings shifts away from default set in /preferences (here 0.7) but could be overridden here 0!..200% is allowed

bGAccel ISF\_weight can be modulated here, too. Note: it can change again if % profile is also changed. -- 0...200% is allowed

CANCEL allows to start fresh (select one or two of the square buttons, ES, AC or HY)

OK needed to use the settings

Duration input is made in minutes. In the exceptional case that both, ES and AC targets are defined, the duration input is pr-AC and framed blue. (This is because the preceding AC mode is automatically determined in length by the loop observing when iobTH is exceeded)

663  
 664  
 665 This looks complicated but only because it allows 4 different modes of use. Each user will primarily  
 666 use her/his preferred one.  
 667 (1) Who is happy with the initially well tuned FCL and does not have huge variations in daily eating  
 668 and moving around, will **not use** the TT **at all**. FCL is possible without an intervention via the



669 TT button in your cockpit. Actually 4 of 8 modes (GGG ...YYY permutations, list see [section](#)  
670 [5.4.1](#)) are not making use of TT.

671 (2) Super easy is also, to just input **any odd-TT** (odd-numbered temporary target) that will shut out  
672 any SMBs for the set duration. *That can be a good idea when having a small snack, for*  
673 *instance.*

674 Super quick access to stop SMBs is possible also via the loop icon ([section 5.3.1](#)).

675 Specifically, an **EatingSoon TT** can be activated here (limited relevance see [section 2.5](#)). It is  
676 time-un-critical, can be manually set, or come up via an Automation.

677 The cockpit enables you to set the iobTH differently (override) for the current meal.

678 Alternatively, iobTH can be temporarily changed in /preferences or using an Automation.  
679

680 Temp. iobTH will always revert to default when the TT expires. If another TT immediately  
681 follows, like in the example of the screen above, it will calculate, (then) show and use a new  
682 temp. iobTH.

683 (3) The third way is to **use the input mask** (if already ncluded in your software version see picture  
684 above) **to freely modulate the loop aggressiveness** for a declared number of minutes. Click  
685 the bottom big square(s): Either HYPO, or ACTIVITY, or EATING SOON, or ACTIVITY and  
686 EATING SOON (example in the pictured screen above). Make or override entries in the offered  
687 fields. Press OK.

688 (4) The fourth way is to exclusively use one of the 4+4+2 little buttons seen in the bottom part of  
689 the TT dialogue box (if already included in your software version). They provide a set of  
690 settings (as will immediately show in all input fields above) that the user has set up in  
691 Preferences/SMB/autoISF/FullLoop (refer to [section 6.3](#)), and can freely label there. *For*  
692 *instance „hiC“ at high carb EatingSoon, „piz“ for Pizza/fatty meals, „grd“ for garden work,*  
693 *„mtb“ for mountain biking ...*

694 **Capturing good settings for not-everyday situations in /preferences** (if already included)  
695 **allows calling them up within 1 second**, from your cockpit on the AAPS home screen (...and  
696 won't ruin the FCL experience at all , especially because in most cases it is not time-critical,  
697 how long before the intended exercise the buttons are pressed).

698 [Case study 6.2](#) demonstrates that nearly the same performance and comfort can be reached  
699 via the **DIY FCL cockpit** with the grey extra buttons appearing at the bottom of the AAPS home  
700 screen, based on Automations with User action (see also [section 5.2.2.3](#)).

701 The example picture given above, and also [case study 6.2](#), is the most complicated (but also most  
 702 useful) case, **when exercise follows after a sizeable meal**. It is then that you need (a) aggressive  
 703 FCL initial performance at the meal, but, *exactly when (!)* a (for the intended sport already  
 704 temp.lowered) *iobTH* is exceeded, you need (b) to have SMBs automatically switched off and go  
 705 into the „milder“ mode, as defined for the exercise (with *high* instead of the immediately prior  
 706 *low*TT, that automatically significantly reduces iobTH again, and insulin sensitivity(resistance)  
 707 settings too).

708  
 709 Pressing exercise related buttons will automatically also light the **exercise button** on the main  
 710 screen yellow.

711  
 712 To summarize, the TT dialogue field offers easy but powerful ad-hoc [modulation of loop](#)  
 713 [aggressiveness](#) for FCL (if already included).

714

715

716 5.3.3.2 Exercise button (see more in [section 6](#).)

717

718 The exercise button automatically lights yellow when exercise related TTs are activated [in the TT](#)  
 719 [dialogue box](#).


720 4 of 8 principal FCL modes ([section 5.4.1](#)) are making use of the exercise button.

721

722 If pressing on the exercise button, a dialogue box appears (*if extended design for FCL cockpit is*  
 723 *already launched*) with info on exercise setting first (and opportunity to override), plus below the  
 724 activity monitor (experimental for auto-tracking of lighter movement during the day, and effects on  
 725 sensitivity that may have. See [section 4.5](#)).

726

727 So, first the exercise settings (as set under TT) are there to read. Example :

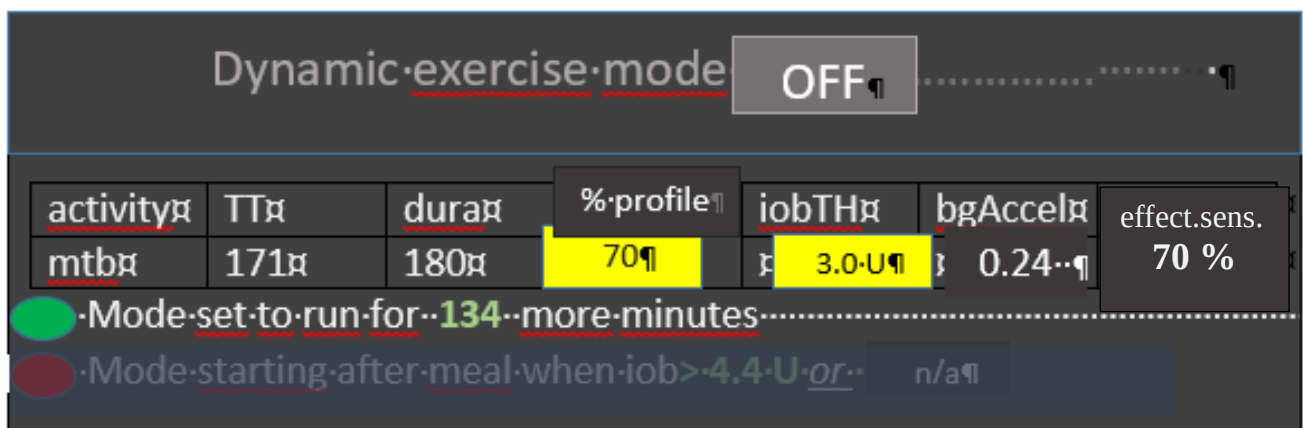

728

729

730 The exercise (here mtb) is selected in the dialogue box of the neighboring TT field, and there auto-  
731 filled with settings made in the set-up and tuning stage by the user under preferences. They are  
732 reported also under the exercise button here, and TT, duration, and % sens (the temp. profile  
733 sensitivity that also shows on the %profile field on the left side of the exercise button) can be temp.  
734 changed there.

735 iobTH, bgAccel\_ISF and overall resulting effective sensitivity ratio (effect.sens. %) is given in the  
736 other fields.

737 The **middle field** of the table, „% **profile**“ either picks up the % set under the %profile button, or  
738 an input can be made here, in the exercise button domain, which will:

- 739 • turn the neighboring %profile button on yellow and show that inputted % on it, too
- 740 • be multiplied with the result from the exercise mode settings per se, and change the %  
741 overall, accordingly.

742 So, if this middle field of above table (dialogue box of sports button) contains a figure other than  
743 100, input field becomes yellow, and you are operating with a combination of traditional PLUS new  
744 exercise mode (with all three top buttons of your FCL cockpit yellow). This maximally will soften  
745 aggressiveness, for which you get an idea by the last calculated figure.

746

747 The mode is either running already (for another number of minutes, as probably also shown in the  
748 yellow TT field anyways). Or it is scheduled to run, after insulination for a started meal reaches  
749 iobTH (as in table). Or, no exercise is scheduled (both points red, no entries).

750

751 The lower part of the exercise dialogue box (not pictured above, but see in [section 6.5](#)) is  
752 dedicated to the Activity Monitor

753

#### 754 5.3.3.3 Profile button

755 The profile button can still be used to set a different profile, or profile%, for instance to adjust for  
756 days with sickness (as you are used to from hybrid closed looping). 4 of 8 modes are not making  
757 use of the profile button.

758

759 Any inputs made here will be used to modify profile\_ISF on which all further changes are made on  
760 (multiplied with).

761

762 The profile field remains grey if standard profile is applied.

763 It turns yellow, displaying a %figure relating to any altered loop overall aggressiveness:

- 764 • When no inputs (changes from 100% profile) are made here, but inputs in the TT field,  
765 e.g. for exercise, automatically lead to different effective sensitivity ratio

766 • when% is changed by input in the profile button itself, it will be multiplied with with  
767 profile\_ISF and be used in place of profile\_ISF *by the algorithm*.

768 However, for exercise (sports) you no longer must make an entry here, because  
769 reasonable %reductions should be automatically provided, driven by your set TT (and half-basal  
770 exercise target), see [section 6](#).

771

## 772 5.4 Recognizing your loop state in the AAPS home screen

773

### 774 5.4.1 Color scheme of top cockpit buttons tells kind of closed loop that is running

775

776 3 Buttons (%profile; exercise; TT) each in 2 states (yellow Y, or grey G) make  $2 \times 3 =$   
777 **eight principal FCL modes** possible:

778

779 GYY = dynamic exercise mode

780 YGY = not-dynamic „traditional“ exercise mode (if <100%) or hypo mode (if >100%)

781 GYG = **basic closed loop with Activity Monitor running**

782 (Note: activity monitor on/off set in preferences will **not** affect the button color in **autoISF V.3.0**; you may  
783 recognize Activity monitor is running by the indicated adaptation of sensitivity despite no TT or temp% are set)

784 GGG = basic closed loop (FCL or HCL) without any altered sensitivities etc

785 YGG = basic closed loop but with a „long wave“ sensitivity shift (e.g. sickness)

786 GGY = temp. target like e.g. EatingSDoonTT is set; or Hypo mode

787 YYG = closed loop with „long wave“ sensitivity adjustment and Activity Monitor running

788 YYY = dynamic exercise mode in time with additional „long-waved“ sensitivity shift

789

### 790 5.4.2 Information printed on the top buttons

791

792 The yellow TT field shows **the currently valid TT** (and further duration):

793 (profile) stands for the abbreviation you labeled your selected running profile



795 In the special case of settings for meal preceding sports, the field will look slightly differently:



797 ...and ...

797 ... when iobTH is first time exceeded, this automatically switches to:



798  
799

800 Likewise, if on the AAPS main screen just an **EatingSoonTT** is set (e.g.72), this is entered with the  
801 desired duration. Afterwards, it automatically reverts to profile target and the display turns grey  
802 again there with e.g. 90 on it (and no time limit).  
803 Without sports context, the middle field remains grey.



804  
805

806 Independently from setting a TT, the user can choose to set a **%profile in the left top field**, for an  
807 independent number of minutes, e.g. 70% in this screen example: Also, or additionally, this will  
808 influence the resulting ISF and sensitivity%

809



810  
811

812 The % might change and turn yellow also in context of making TT inputs in the related dialogue  
813 box (see chapter TT dialogue field, above). Still, the % (or the length of time the profile switch shall  
814 be active) can be independently overridden in the top left field, if so desired.

815

816 If an **Automation** sets a %profile, and/or a TT (e.g. *automatic detection of meal start at condition*  
817 *e.g. when delta >10*) , this would automatically show in respective field(s) turning yellow and  
818 showing the temp. setting. To show the set parameter comes from an Automation, „ **(A)** „ is added  
819 in the end of button text.

820 Note that an Automation is **usually/ always (?)** only permitted to temp. change default profile  
821 settings, not other pre-existing temp. settings. This is for a good reason : Why should a  
822 sometimes in the past thought-out Automation supersede your - just for the occasion  
823 specified – temp.settings that you consciously activated for the day?

824 Advice: Try to stay away from Automations that also aim at temp. modifying  
825 aggressiveness. For the reason just given in above note, they often will not kick in anyways.  
826 Generally, it also is no good idea to double up sub-algorithms for tweaking loop behaviours.

827 Also, as pointed to above twice already:

828

829 Try to **keep things as simple and clear as possible**.

830

831 That said, a limited number of Automations can be of help in distinct scenarios (that differ in  
832 purpose and in applicable time of day).

833

834 A good one could be for night time, when your odd profile TT has SMBs shut off, but your  
835 experience after pizza nights tells you that, under certain condition patterns (bg, iob), an  
836 SMB or two should be „allowed in“ (see example given in [section 5.1.2](#); used also in [case  
837 study 4.3](#) ).

838 Another good example, if you go usually FCL without any use of the TT button (which would  
839 be a meal announcement of sorts), is to define an Automation that, after detecting a meal  
840 start, automatically sets a low TT to get maximally aggressive first SMBs (as is the author's  
841 preferred way, mentioned already in [section 2.5](#), used also in [case study 4.3](#) ).

842

#### 843 5.4.3 FCL related indicator fields in the AAPS home screen

844

845 In extra data fields of the AAPS main screen you can always see (not change) the key  
846 „aggressiveness“ parameters your loop currently operates with (see also home screen  
847 example below):

848 • how profile sensitivity (**ISF**) adjusts by the %profile input, by autoISF, and/or a set  
849 exerciseTT, resulting in an effective sensitivity (ISF that is used to determine  
850 insulinRequired. Details for every loop decision see result/debug section of the SMB tab).

851 • next to current available iob number is an indication of your **valid iobTH** (the iob above  
852 which no more SMBs will be given)

853 • The AAPS home screen additionally shows, above the deltas, the current **acceleration**

854 Having a look at that can be valueable. For instance, when glucose is relatively low and still  
855 falling, a positive (and getting more positive) acceleration indicates that bg will swing back  
856 up, rather than crash low. This will give info about necessary snack size, and hence help  
857 avoid both, unnecessary calories, and going on a bg roller coaster.

858  
859  
860 5.4.4 Overall home screen:  
861

Overall home screen:



862  
863  
864