

## 5. Modulation of autoISF aggressiveness.

V 4.3.2

Please note that with autoISF you are in an early-dev. environment, where the user interface is **not optimized for safety** of users who stray away from intended ways to use. Good safety features exist, but these are only as good as the development-oriented user understands and implements them. This is not a medical product. Refer to disclaimer in [section 0/readme.md](#)



### 5.1 Automatic modulation of loop aggressiveness

- 5.1.1 “autoISF off” outside of meal times
- 5.1.2 SMB off @ odd profile target
- 5.1.3 SMB off @ odd temp. target
- 5.1.4 diff. of FCL aggressiveness via Automations
- 5.1.5 diff. of FCL aggressiveness via Activity Monitor
- 5.1.6 Pro/con completely hands-off FCL

### 5.2 Manual modulation of FCL aggressiveness (DIY cockpit)

- 5.2.1 Status recognition
- 5.2.2 Manual interventions from DIY cockpit
  - 5.2.2.1 Temp. %profile or TT settings
  - 5.2.2.2 Temp. settings in /preferences
  - 5.2.2.3 Grey DIY cockpit buttons for FCL responses
- 5.2.3 Temporary exit from FCL

### 5.3 Recognizing loop state from the AAPS home screen

- 5.3.1 Modulated loop aggressiveness (3 top buttons)
- 5.3.2 Color scheme of the top 3 buttons
- 5.3.3 Active Automations
- 5.3.4 FCL related indicator fields
- 5.3.5 Overall AAPS home screen
- 5.3.6 Info given every 5 minutes in the SMB tab
- 5.3.7 Info about last 15 autoISF decisions
- 5.3.8 SMB tab info when operating 1-minute Libre
- 5.3.9 Summary: Your personal FCL cockpit

### 5.4 Manual modulation via improved cockpit

- 5.4.1 Violet FCL icon and underlying buttons
- 5.4.2 Bottom buttons “insulin” etc.
- 5.4.3 Top three fields
  - 5.4.3.1 TT dialogue field
  - 5.4.3.2 Exercise button / dialogue field
  - 5.4.3.3 Profile dialogue field

[Available \(related\) case studies:](#)

Case study 5.2: Sweet snack.

Case study 5.3: Compression low

*Skip what is in green writing:*

= Drafted fragments or not implemented ideas. Read only if you like to participate in discussion – Suggestions will probably be radically reduced because FCL *can and should* be done *without* most discussed extra features

Once the initial tuning according to [section 4.](#) is done, you are ready to use autoISF for your fully automated meal management.

## Challenges beyond managing main meals

There are up to four other major challenges you might have to manage:

1. Secure the technical pre-requisites for your FCL remain a given (see [section 1](#)):
  - Prevent, recognize and manage (partial) occlusions (=>frequent cannula changes)
  - Secure “meaningfulness” of 100% all your CGM values

- 49           ○ Secure 24/7 BT connectivity of all your components
- 50    2. Deal with times when insulin given, after meal detection, must be restricted, e.g.:
- 51           ○ a snack should not be “treated” as a meal
- 52           ○ meal before exercise should not be treated “all the way” towards high iob.
- 53    3. Deal with times when the loop should be set "less aggressive" as a precaution, e.g.
- 54           ○ during nights
- 55           ○ in an exercise context
- 56    4. Deal with times when the loop should be set "more aggressive" e.g.
- 57           ○ at sickness
- 58           ○ at tempting hotel breakfast buffet

59 How big the remaining challenge really is, depends very much on your individual lifestyle.

60

#### 61 Hands-off FCL – an elusive goal?

62

63 To run the FCL “hands-off” around the clock, you now must analyze also **your** data in the times  
64 *outside* the meal blocks, and seek solutions to problems (if any).

65 With a technically well-functioning system, moderate meals, moderate or no exercise, moderate %TIR  
66 expectations and a bit of mindfulness it should be possible to go into **Full Closed Loop 24/7**, after  
67 working through, and observing, [sections 1-4](#). See [case study 4.3](#).

68

- 69 • In [section 5.1](#) we explore avenues towards **fully automated** management that in daily life will  
70 require no user intervention at all.

71 Often it will be **your** choice whether you want to bother with researching in **your** data for, and  
72 defining, dedicated (“personalized”) Automations (see also [section 5.1.6](#)). Or whether you accept  
73 instances, where you do manual steps:

74

- 75 • In [section 5.2](#) (and later, related to exercise, in [section 6](#).) we will look at solutions that involve  
76 an easy **user interaction** like a data entry, or a button push.

77

#### 78 Avenues for temporary modulation of the FCL aggressiveness

79

- 80    1. temporary shut-off SMBs (odd-numbered target)
- 81    2. temporary change bgAccel\_ISF-weight

- 82 3. temporary change iobTH\_percent
- 83 4. temporary change the set %profile
- 84 5. temporary set different bg target (especially in connection with exercise mode)t

85 All of these are also easy accessible via Automations in AAPS

86 This means: If, in **your** data, you find sets of conditions where any of the actions 1.-5. could  
87 help, you can make this aspect integrated part of **your** hands-off FCL!

88

89 Note that in Automations you can go beyond above mentioned 5 avenues, and come up with a  
90 nearly surgical definition as to from which iob on, and many other criteria, some temporary  
91 modulation of aggressiveness shall happen.

## 92

### 93 5.1 Fully automatic modulation of FCL aggressiveness

#### 94

95 The following subchapters describe set-ups you may want to use for allowing **completely hands-**  
96 **off FCL** in as many daily situations as possible (and potentially all the time, as in [case study 4.3](#))

#### 97

##### 98 5.1.1 All autoISF ISF adaptations switched off outside of meal-time windows

#### 99

100 If, aside from having to bolus for meals, your hybrid closed loop was running pretty well *without*  
101 other interventions from your side, you could continue to run in that mode, and just focus your new  
102 autoISF FCL on management of meals.

103

104 In your initial transitioning phase, this approach makes a lot of sense, and even by focusing  
105 autoISF on just a sub-set of meals, like only dinners.

106

107 Also in the long run this avenue is taken by many FCL users for the night times, “hanging  
108 on” to their well performing hybrid closed loop with standardoref(1) SMB+UAM

109

110 For this, you define Automations

111

- 112 • that set meal time windows in which “Enable ISF adaptation by glucose behavior” (autoISF) is  
113 turned on in AAPS preferences/OpenAPS SMB
- 114 • or: that turn *all* autoISF's ISF modulations (*or just bgAccel\_ISF* ) off in time windows in which  
115 surely no meal occurs.

116 For instance, you can go *for all nights* back into your Hybrid Closed Loop, as you had before.

117

118 Your temp. "autoISF shut-down" (exiting autoISF FCL = shutting off "Enable ISF adaptation by  
119 glucose behavior" in AAPS Preferences / SMB) is meant to prevent problems from the loop over-  
120 reacting to bumps in the glucose curve in times of day (night), when standard oref(1) performance  
121 is sufficient.  
122 A very good alternative to defining "meal-time windows", or to resorting to night-time Hybrid Closed  
123 Loop, is letting the autoISF FCL run 24/7, and "taming" the FCL via a night time SMB shut-off (see  
124 next [section 5.1.2](#)).

125

126 *Off topic: **Other "early dev" AAPS variants** (see [section 13.3](#)), do require working with meal-time*  
127 *windows. The window is either set by time of day in the settings, or it always must be „set“ by the user.*  
128 *Trigger to set a meal time window could be a pre-bolus given by the user, a carb entry made, an*  
129 *EatingSoonTT set, or a meal-announcement-button pushed (none of these things are required in*  
130 *autoISF FCL).*  
131 *Outside of these time windows, these loops then run with less aggressive SMBs, just like oref(1)*  
132 *SMB+UAM in AAPS Master, with or without (modified versions of) dynamicISF. So, outside of the meal*  
133 *time windows, you would be in "your normal" hybrid closed loop.*  
134 *The term **Meal Announcement** (MA) is often used to label this closed looping mode.*  
135 *It is not really FCL, but an advance over traditional HCL.*

136

### 137 [5.1.2 Odd-numbered profile targets, to block SMBs](#)

138

139 An alternative route of preventing the FCL loop from over-reacting to bumps in the glucose curve  
140 would be to make use of the option to temporarily shut down SMBs

141

142 Ensure the even/odd logic in the settings is toggled on in Preferences> openAPS SMB>  
143 autoISF settings> smb delivery settings>: "Enable alternative activation of SMB depending  
144 on bg target": ON.

145

146 In time blocks with an odd-numbered profile target you can prevent any SMBs being given by your  
147 loop. The (unchanged) aggressive settings then can only translate within the limits set by %TBR  
148 possible.

149

150 This will very much slow down any more insulin being given, and is an excellent solution for night  
151 times, especially if you occasionally experience compression lows.

152

153 Alternatively, you could use the new included options for Automation Conditions and temporarily  
154 tune your bgAccel\_ISF\_weight much lower (see [section 5.1.4](#)).

155

156 The same situation can be achieved if you generally operate with a mild bgAccel\_ISF, and make  
157 your autoISF only really aggressive for meal-time slots (if you have similar enough times every day,  
158 or also can “employ” geo-fencing in your Automation (or middleware, in iAPS) conditions).

159 In these cases you would not need to have night profiles that disable SMBs: - Which is the better  
160 way would depend on a lot of personal factors relating to how high-carb the diet is, regularity of  
161 meals, snacking habit, CGM quality and incidence of compression lows, and probably more. - I  
162 would try both routes, or, as this is fairly complex to tune, just one, and stick with what is working  
163 good enough.

164 Yet another alternative was already presented ([section 5.1.1](#)) = to go into hybrid closed loop for  
165 the night.



166 That is possible, with SMBs available (without them getting boosted via autoISF), and, for a  
167 long time, was the author’s favored solution for the nights.

168 This solution is evidently similar to the prior discussed one, of having a mildly tuned autoISF  
169 24/7, boosted to high aggressiveness only in meal-time slots.

170  
171 **Adjunct Automations to allow a few SMBs, in nights with odd profile target SMB shut-out**

172  
173 My current **favorite** builds on the **method** of this section (5.1.2, odd profile target provides SMB  
174 shut out), but then allowing some, automatically triggered SMBs, when needed:

175  
176 In case you occasionally do have nights that would benefit from a couple of SMBs (to treat temp.  
177 highs from a late fatty pizza, raclette and such): Define suitable Automations like the two „night“  
178 ones in this *example*:

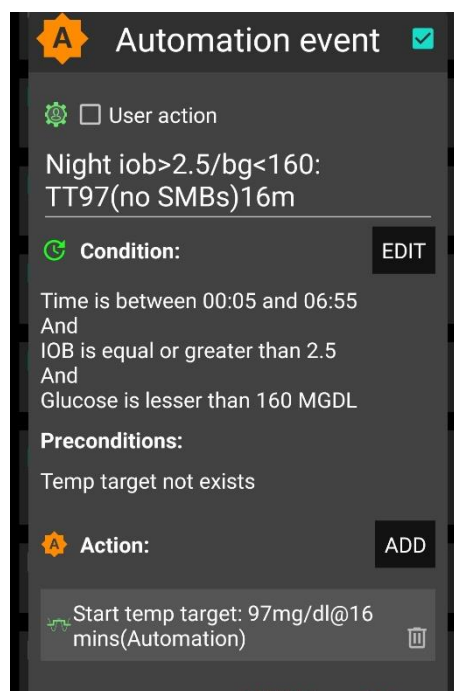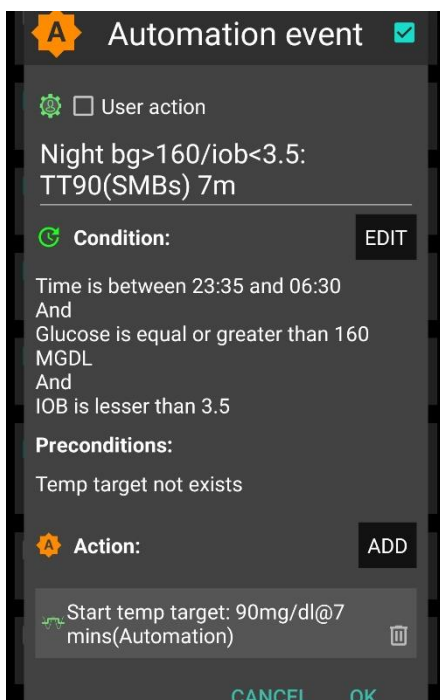

179  
180  
181 Caution: Never underestimate the „trickiness“ of getting your Automations „right“.

182

183 With your thought-out Automations in place, night data need to be analyzed to see

- 184 • whether the bg and iob limits, as defined in the given example, work sensibly four your data
- 185 pattern
- 186 • whether the TT duration is chosen appropriately
- 187 • how swapping the sequence in which the automations appear in the Automation list would
- 188 lead to different SMB impacts.

189

### 190 5.1.3 Odd-numbered temp. targets (TT) set via Automation, to block SMBs

191

192 A widely used Action that strongly modifies how fast your FCL can add more iob is setting an **odd-**  
193 numbered **temp. glucose target** which makes the loop operate without giving any SMBs (%TBR  
194 modulation only).

195 Ensure the even/odd logic in the settings is toggled on in Preferences> openAPS SMB>  
196 autoISF settings> smb delivery settings>: "Enable alternative activation of SMB depending  
197 on bg target": ON.

198

199 Then, from patterns you find in *your* data, at times where you want your loop act differently, you  
200 need to carve out Conditions that describe the respective situations (and either *for how long* it  
201 typically lasts, or at which *other* Conditions you want your loop get back to default FCL operation).

202

203 An odd TT is often set for an *anti-hypo* snack or *sports* snack. In both instances, you do not want  
204 SMBs to quickly counter-act.

205

206 In case of *sweet "fun"* snacks, this is entirely different -> [section, 5.2.1](#) or for regular snacks  
207 (*e.g. at school break*) see next [section 5.1.4](#)

208

### 209 5.1.4 Automatic differentiation of FCL aggressiveness using Automations (or middleware)

210

211 **Personalized Automations** tailor the loop exactly to **your** data so **fully automated** handling of  
212 situations with **different aggressiveness** of the loop can be made.

213

214 Automations are an integrated and very easy-to-use feature in AAPS.

215 (The i-Phone platforms Trio or iAPS lack this feature. However, so-called **middleware** has been  
216 developed as add-in to your code, see: <https://github.com/macconnellk/RoboSurfer/tree/main> )

217 .

218 From, autoISF 3.0 onwards, also the following parameters are provided as Condition and/or as  
219 Action for defining YOUR Automations:

- Enable ISF adaptations by glucose behavior => Allows temp. ON/OFF for the key ISF modulation parts of autoISF (and, as a result, will usually decrease loop aggressiveness)
- Trigger/set iobTH percent => Keeps default aggressiveness, but only until a iob threshold (that your Automation modifies) is surpassed (which is when any further SMBs will be blocked blocked)
- Trigger/set bgAccel\_ISF\_weight => Modifies the aggressiveness of just the acceleration component

To set up suitable Automations, you first must **analyze patterns** you find **in your data**, at times (or geo-location, or bg and iob patterns that point to a problem ...) **where you want your loop act differently**, to carve out Conditions that describe the respective situations (and either for how long it typically lasts, or at which *other* Conditions you want your loop get back to default FCL operation).

A variant of this mode is to define several windows in which autoISF aggressiveness (bgAccel\_ISF\_weight) and/or iobTH are automatically set differently

- for **different meal time slots** of your day –

(*Breakfast at home, school lunches, school intermission snacks, dinners at home* could for example all deserve special settings regarding ISF\_weights and iobTH).

Note: Circadian differences in insulin sensitivity between meal times are included via your ISF profile and should not be a reason for different \_weights needed between meals!

- or even for a geo-location etc –


(*School lunches, or mother-in-law visits, would be examples*).

An example for this was given in section 3 already:


Here is an example set of automations to alternate between two values of iobTH:

I use two different values of *iob\_threshold\_percent* during a normal day. It is 40% for lunch time and 60% for dinner time. I have these two rules to switch by time of day and only if the current value equals the value from the earlier shift. Any other value is treated as a manual override for special occasions until I manually set it to its regular value. The time windows for switching are long enough to catch an opportunity to be processed and do not need to be actioned half a day each.

lunch

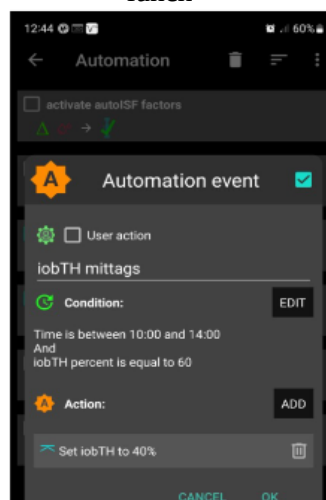

dinner

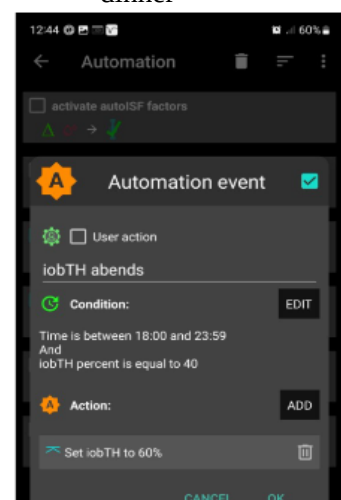



244 Unless your meals differ vastly in size and in fast carb content **all this may not be needed.**

245

246 Still, personalized Automations might help ease your initial job of setting the various ISF\_weights,  
247 and a best-suitable iob\_threshold\_percent that would work “always”.

248 Note: Only the main two parameters (bgAccel\_ISF\_weight for “initial aggressiveness”, and  
249 iobTH\_percent for “where SMBs stop”) are available in Automations. So, finding your  
250 parameter sets *for each of* the time slots, will not be trivial. => **Spending more effort to set**  
251 **the ..\_weights so they accommodate just one, broader spectrum (section 4.) should be**  
252 **the first, and standard, approach.**

253 An intermediate (maybe only temporary) approach could be to use a profile switch (for low  
254 carb meal, or eating half, setting %profile to 60% for instance, and only for the brief, less  
255 than an hour, initial meal period). See [section 5.2.2](#) and [5.4.5](#)

256


257 **Caution:** Setting a different iobTH% or bgAccel\_ISF\_weight can probably not be done with a  
258 duration attached. Then you **must** define a suitable **additional Automation** that must be

259 active in tandem, to **restore the values you had set**  
260 **in /Preferences for your iobTH% or bgAccel-**  
261 **ISF\_weight.** Else, once your Automation set in, it will  
262 *forever* shift these important parameter settings!

263

264 If for instance you have several Automations that, in  
265 combination with a set elevated TT also set a lower iobTH:  
266 Don't be fooled, the duration only applies to the TT. You  
267 need an extra Automation for all of them.

268 I picked out the highest of the altered iobTH values that  
269 these Automations can set (45\_percent), and then I can  
270 automatically restore my default desired 60% via this one  
271 Automation (see screenshot - - > )

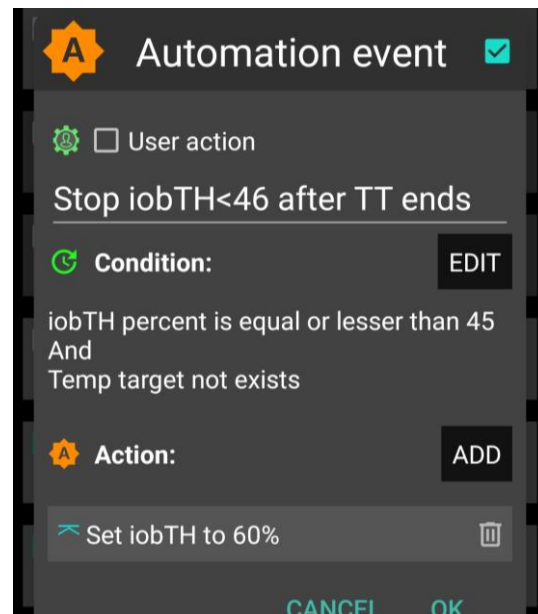

271

#### 272 5.1.5 Different FCL aggressiveness set by the Activity Monitor

273

274 With the autoISF variant of AAPS you can make use of your smartphone's **stepcounter** and use it  
275 to fully automatically adjust insulin sensitivity ratio to **activity level in the past minutes to one**  
276 **hour** time frame.

277

278 This feature comes with yet another little tuning opportunity, in which you study your body's  
279 response to light exercise (like walking) or to not moving at all (like desk, couch), and select  
280 appropriate settings (in AAPS Preferences/OpenAPS SMB/Activity modifies sensitivity/ -> set **two**  
281 **scaling factors**).



282 This **will automatically adjust insulin delivery** (basal, ISF, and iobTH; see 1<sup>st</sup> screen of AAPS  
283 SMB tab (example in [section 5.4.5](#))), to suit activity state of the past minutes (up to 1 hour).

284

285 This autoISF feature (new since V.3.0) is much quicker responding than Autosens or dynamicISF  
286 to adjust insulin sensitivity to your current „lifestyle state“ (if it is largely (in-)activity related, which  
287 often is the case).

288         For loopers who do not have huge variations in exercise levels (for which they would need the  
289         exercise mode) in their everyday lives, the Activity Monitor might fairly much close the gap towards  
290         being able to do a 24/7 hands-off FCL.

291

292 The Activity monitor is described in more detail on page 9 of the [autoISF Quick guide](#)  
293 (<https://github.com/ga-zelle/autolSF>) and in [section 6.6](#)

294

295 While the Activity Monitor takes automatic care of *light* deviations from an average activity level,  
296 exercise enthusiasts, or heavy workers, should make use of the “*heavier tools*” (which  
297 automatically supersede (shut off) Activity Monitor: *Exercise mode*, discussed in [section 6.](#))

298

#### 299 [5.1.6 Pro/con completely hands-off Full Closed Loop](#)

300

301 To stay **24/7** in a **completely „hands-off“ FCL** can be a realistic goal with autoISF, if besides  
302 meals also some special challenges, as discussed in this [section 5.1](#), were analyzed and could be  
303 addressed. Example see in [Case study 5.3](#).

304

305 Clearly it depends very much on your lifestyle, and how interested, willing, and capable you are to  
306 recognize, deal with, (and in the future avoid?) situations that get you outside of your desired %TIR  
307 on occasion.

308         So, this is also about what %TIR you are aiming at, and can accept, as it averages out for  
309         the week, for instance.

310

311 Everybody must weigh for her/himself

- 312         • how much **upfront effort** to put into the setting up process, for getting it all 100% automatic
- 313         • **or whether to take an easier start, with a couple of situations left to take care of when**
- 314         **and as they arise in daily life**

315

316 Even if a principal capability for a fully automatic “hands-off” running FCL is given, this still  
317 means that

- 318         • the user should be knowledgeable about what exactly is going on, and

319 • have a principal capability to „nudge“, or even to completely take over, in a manual mode.

320 In the sections that immediately follow, we present the options to nudge or temporarily take over  
321 from the AAPS home screen which will be serving as your **FCL cockpit**:

322  
323 [Section 5.2](#) describes how you can use available “buttons” from your AAPS home screen to  
324 manually “nudge” aggressiveness of your FCL.

325  
326 [Section 5.3](#) will show how you recognize FCL state (aggressiveness) and understand recent loop  
327 decisions

328  
329 **No need to read section 5.4, unless you are interested in contributing to define/design/program**  
330 **further improvements). Ideas for an improved FCL cockpit in the future (probably an over-design; I**  
331 **intend to reduce to a few “really nice to-have” features)**

332 •

## 333 5.2 Modulating aggressiveness **manually**, from the DIY-FCL-Cockpit\*

334

335 \* Like in the airplane cockpit: Cruising in full auto mode should involve having an eye on the  
336 instruments, and on potential disturbances ahead in the environment.

337

338 In [section 4](#), we dealt with major meals.

339

340 In [section 5.1](#) we looked into ways to **automate** also a modified loop response to *foreseeable* situ-  
341 ations (tied to a time of day, geo-location etc), or to those *the loop could recognize* (with enough  
342 time to react).

343

344 However: *Other disturbances* might come up, that:

- 345 • are not noticeable in-time, or foreseeable, by the loop (*e.g. your plan to start exercise in an*  
346 *hour or two*), but **that influence sensitivity dramatically**, and therefore require temporary  
347 modified settings in order to remain in-range, and/or
- 348 • **require** a different “starting point” regarding iob and bg, which translates into **a different**  
349 **iobTH** that should **temporarily** be set much lower (*in case of exercise*) or noticeably higher  
350 (*e.g. with very fast absorbing carbs in a sweet snack “sin”*) .

351

352 *In these scenarios*, you must find an easy way to

- 353 • call up a pre-programmed routine for automatic management, with adjusted  
354 aggressiveness, or:

- manually tweak a setting or two, to temporarily adjust the aggressiveness
- There may also arise a desire to just exit the FCL mode, and “be your own captain” for mastering a special situation.

Lastly, for peace of mind, to learn, and to stay informed (especially so in your initial tuning phase, or when your glucose curve goes in unexpected ways) we also must be able to...

- ...find the key parameters that frame and drive the recent and upcoming loop decisions.

All this is facilitated within seconds right from the AAPS home screen, serving as a **FCL cockpit** after you built a couple of DIY cockpit features via Automations (as described below and in [case studies 5.2](#) and [6.2](#)):

### 5.2.1 Status recognition

Before considering any manual interventions into the ongoing FCL, you should be aware what the current mode of action is (refer to [section 5.3](#)), and hence how you might be able to “nudge” your loop in order to adjust to the disturbance that you see coming up.

### 5.2.2 Manual interventions from the (DIY-) FCL cockpit

Caution: Trouble with some of these is, not to forget to stop and set them back, manually, too.

#### 5.2.2.1 Temporary tuning of FCL aggressiveness via temp. %profile or TT settings

##### %profile set

The set **% profile** multiplies with both, the ISF resulting from autoISF, and also with the default iobTH you have set, so both are nicely modulated in a linear way with the % temporarily chosen

Just taking profile e.g. to 110% for an afternoon might be an easy way to explore whether you might benefit from 10% more “aggressiveness” in your core settings for lunches  
Make sure, though, that the extra 10% are not cut away by set safety limits.

##### Elevated or lowered temp. bg target (TT)

A lowered (relative to profile glucose target) temporary **bg target (TT)** signals lowered sensitivity (more insulin need), and

391 An elevated TT (as often used with exercise) increases sensitivity and hence works in the direction  
392 of a lowered % profile to also reduce insulin given by the loop.

393 Caution: In preferences/SMB..., make sure you set “High TT raises...” and “Low TT lowers  
394 sensitivity”.

395

396 Even or odd numbered temp. bg target (TT) (of similar level)


397

398 As we learned in section 5.1.3 already, switching between an even or odd bg target can be used to  
399 strongly modify loop aggressiveness by switching SMBs on/off temporarily.

400

401 Note this is only an on/off choice, so not as finely tune-able as the other presented routes. (Only by  
402 choosing very short vs. a bit longer run time, there is “a bit of tuning” possible).

403

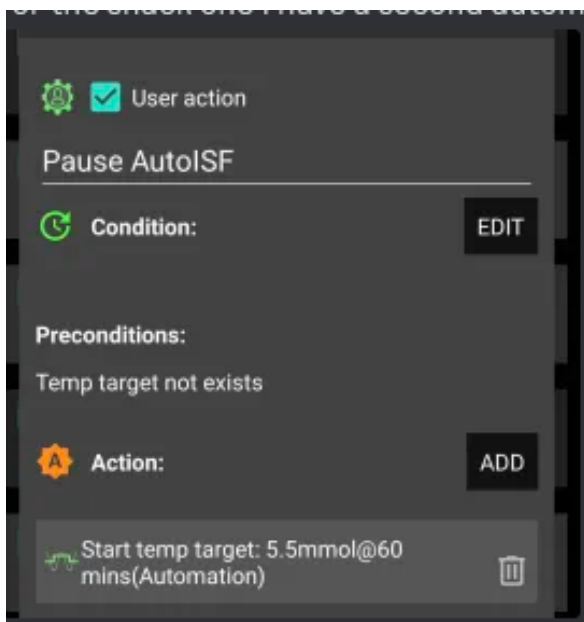

404

405

406 Exercise button

407

408 Moreover, the **exercise button** (top center on your AAPS home screen) can be activated (turns  
409 yellow, then). This can **\*) further boost** how your set TT elevates the resulting ISF, and sharply  
410 lowers iobTH, as often desired for sports. **\*) see below in [section 5.3.2](#)**.

411

412 5.2.2.2 Making temporary changes in settings made in AAPS/preferences/Open APS SMB

413

414 Going into AAPS/preferences/Open APS SMB allows to:

415 - set milder or stronger ...\_ISF\_weights

416 - set different iob\_threshold\_percent (or iobMAX)

It is easy to use the **top right “TT” button** in the AAPS screen to switch from an even target to an odd numbered one, in order to shut off SMBs for a time you can set.

<- An even easier to operate alternative is setting up this **User Action Automation** called “Pause autoISF”

(suggested by Slip, Discord Dec.2024).

*(Actually, autoISF keeps running, but it seems paused because it can't translate big calculated insulin requirements into SMBs)*

417 - elevate or lower the SMB\_delivery\_ratio

418 - limit or expand max. allowed SMB size

419 - change the the even <-> odd logic for SMB on/off

420 Doing temporary changes *in AAPS/preferences* should be the exception because

421 - they require multiple steps, including entering a password

422 - you will often forget to set everything *back to original* settings, a couple of hours (or already  
423 minutes) later.

424

425 5.2.2.3 Triggered Automations: Grey extra DIY cockpit buttons for pre-programmed “responses”

426

427 Recognizing conditions for fully automatic handling by the loop may not be not possible, or come  
428 too late for the loop to act on. Examples would be

429

430 • *exercise*: Minimum an hour *before* starting exercise, “the loop should know” to be able to  
431 lower iob and elevate bg by the time exercise starts.

432

433 • *snacks*: High carb snacks, sweets, consuming ice cream or having a sweet drink comes  
434 with the problem of even steeper glucose rises, but overall a lesser insulin need, compared  
435 to major meals (for which we tuned our FCL according to [section 4](#)).

436

437 This not necessarily implies that snacks *need* different settings than a meal. After all, autoISF  
438 was designed to react to all available data, especially to where the developing glucose curve is  
439 headed. So, depending on your effort to set parameters for a broad variety of meals (notably:  
440 how well you avoid to invariably bounce fast against your iobTH), you might be able to accom-  
441 modate low carb, snack, and major meals with *one* set of settings.

442

443 In FCL autoISF, this is a bit more difficult than in HCL autoISF applications, because FCL  
444 involves revving up iob supply (largely via big bgAccel\_ISF-weights) sometimes too much,  
445 to be balanced by just a snack getting absorbed.

446 In case a snack did trigger a “full meal response”: (1) You probably must continue snacking  
447 to prevent a hypo from your initial FCL over-reaction. (2) For future days, analyze your data  
448 (and snacking habit) to define how to prevent this from happening often.

449

450 For increased comfort and safety, you might have to differentiate, and make use of what follows for  
451 the *sweet snack* example, [case study 5.2](#).

452

453 Note that in the iPhone versions of autoISF (**Trio** and **iAPS**) there are no Automations . Instead you  
454 need so-called **Middleware**, like for instance suggested for %sensitivity (profile ISF) adaptation by  
455 one user here: <https://discord.com/channels/953929437894803478/1025731124615458848/1238099464531611668>  
456

457

## 458 **Tuning aggressiveness**

459

460 A sweet snack likely benefits from even more aggressive initial FCL performance than set for the  
461 meals in your normal spectrum of diets.

462 Therefore, you could set

- 463 • a higher **temp. profile%** and/or
- 464 • a temp.elevated **bgAccel\_ISF-weight** (see screenshot of my Automation).
- 465 • a **low temp. target** (76 mg/dl for instance; this additionally helps maximize the first SMBs  
466 that will automatically be triggered at detection of acceleration)..

467

468 When first defining and testing this Automation, also check:

- 469 • that the safety limits as discussed in [section 2](#) will not block the intended elevated aggressiveness  
470
- 471 • SMBs will not get outrageously big, and iobTH sometimes exceeded by too much.

472 Note that “the last SMB” is allowed to overshoot the effective iobTH by up to 30%, where it  
473 will be cut ~~(or by up to 20% at even target > 100 mg/dl).~~

474


## 475 **Limiting iob**

476

477 For “just a snack”, total insulin need will be lower than for a meal.

478 If you would just have your sweet drink, and your meal-oriented FCL would “attack”, iob  
479 likely would become too high, and a glucose rollercoaster would start, with you needing to  
480 consume more =>

481 If you just have a snack, or drink a small glass of juice, you can lower the **iobTH\_percent**  
482 accordingly.

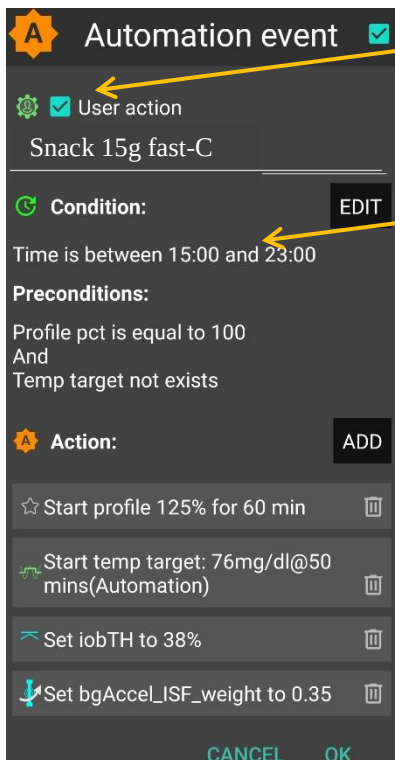


„User action“ is always ticked-on

This will, in the defined time space ..

..offer the “DIY cockpit” button..


..which I must press any time (~90...30 minutes) before my snack.

Note: Make sure that iobTH and bgAccel\_ISF\_weight revert to default afterwards

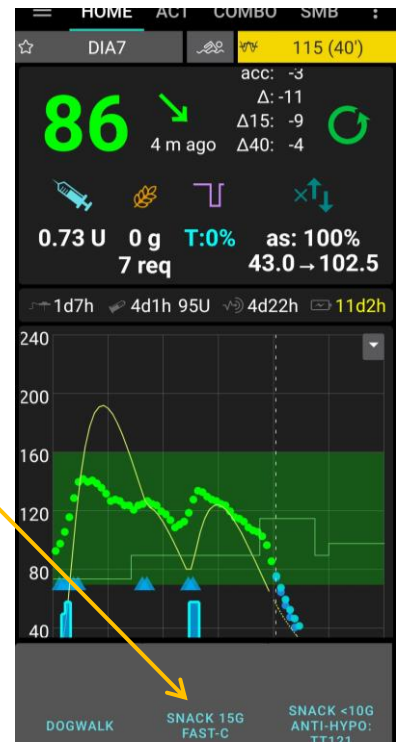

483

484

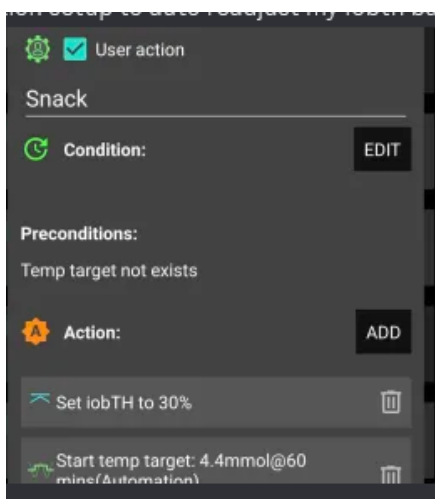

<.- A simplified way was proposed in Discord (Slip, Dec.2024 )

Again, note: Make sure that iobTH% reverts to default afterwards (ideally, automatically –see below ~ line 524)

485

486

487 So, this can be a little **extra “project” when setting up your FCL.**

488 You need to research your snack habits (if any), and over time find out which settings in  
489 the snack-related Automations work well.

490

491 **In everyday life you then just must press the related button in your cockpit** (which is  
492 not time critical at all, except it should be clicked *latest* a couple of minutes after you took  
493 the drink or snack).

494

495 If you consume more, and also eat something with your sweet drink, this will more  
496 resemble a full meal... however, with unusual amounts of fast carbs.



497 **Caution:** Pressing your snack button a *second time* would **not** help because the lowered  
498 iobTH does not allow iob going high enough. So you are better off just letting your *normal*  
499 FCL meal routine run, after your snack mode expired.

500

501 Developers feedback re. the **option to automatically block out** using the same User

502 Action Automation a second time within, say, 2 hours:

503 1) In the tandem Automation, line 524, add an ACTION for 120 minutes setting ("value")  
504 for parameter ("whatever") (= whatever does not hurt your looping, or block any of your  
505 other potentially important other Automations)

506 **Caution:** Better **not** use a TT for your ("whatever") 1) because of the even / odd SMB  
507 shut-out 2) in case you have Automations that demand a clean slate regarding no TT  
508 set, better use any other parameter from the long list that is given to choose from, when  
509 you set up Automations.

510 Setting a 105% profile, for 2 hours, might work, for instance. **Note:** Whether this really  
511 is a good idea also for you to pick for ("whatever"), depends very much on your other  
512 Automations, and loop settings.

513 2) To the User Action Automation, line 465, add as an additional CONDITION: ...  
514 ("whatever"... ) is NOT ("value" what you set it to block out 2 hours

515 In case you chose 105% profile for your ("whatever"): additional CONDITIONS would  
516 be: ... AND Profile pct is lesser than 105 OR Profile pct is greater than 105

517 ⇒ At 105% profile, which gets set (for the 2 hours block-out period) by your other  
518 Automation, you cannot activate a (would-be 2nd) User Action "snack" Automation.

519 Other options (when you just can't stop snacking) would require a manual modulation  
520 regarding %profile and/or bgAccel\_ISF, but keeping the full default set iobTH\_percent, or even  
521 elevating it (refer to [section 5.2.3](#)). If that happens often, define for yourself an extra User action  
522 Automation for a bigger snack (= another grey DIY cockpit button).


523

524 **Caution:** Setting a different iobTH or bgAccel\_ISF  
525 \_weight can not be done with a duration attached.

526

527 Hence you **must** define a suitable **additional**  
528 **Automation, that** must be active in tandem, and  
529 **restores the iobTH or bgAccel-ISF\_weight** in  
530 AAPS/Preferences. Else, once your Automation set in,  
531 it will *forever* shift these important parameter settings!

Restoring setting in /preferences :TandemAutomation:

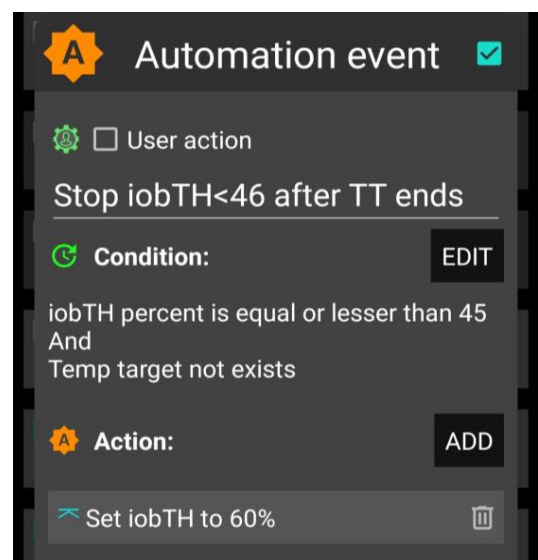

532 In my TandemAutomation example given above, I picked out the *highest* of the altered iobTH values that my  
533 iobTH-reducing Automations can set (45\_percent), and then I can automatically restore my default desired  
534 60% via this one Automation (see screenshot given above )

535

536 In the related Automation, just keep the “User action” box clicked at all times, and define in the  
537 Conditions when you want to see that button available for cockpit use (see screenshot above) =>  
538 you will see that button offered.

539

540 Besides snacks, also any **other recurring special situations can be addressed via a DIY**  
541 **cockpit button, and receive different aggressiveness up to a suitable iobTH level.**

542

543 Over time you can have a big number of User action Automations, and keep them “shelved” rather  
544 invisibly (clicked in-active, top left box) in your long list of potential Automations. Even when active,  
545 they only show in your cockpit (bottom grey field of your AAPS home screen) *in the time slot you*  
546 *assigned* as potentially relevant (and if all other conditions are met).

547

548

549 Discussion

550

551 In case you do have a snack habit and ...

- 552 • can not find settings, as in [section 4](#). defined for your meals, also suit your snacks
- 553 • can not pin a time slot or other Condition to it for programming an Automation response
- 554 as in [section 5.1.4](#)

555 .. then you minimum need a “snack announcement” for which the extra button in your DIY cockpit  
556 provides a time-uncritical 1-button-push solution.

557

558 This could be a good solution for kids in kindergarten, too. Make sure caregivers  
559 understand to use it *only once* for *one* snack. Continued snacking would require iob as for a  
560 meals. This is what the FCL loop takes care of automatically; **using the snack button**  
561 **several times in a row would limit iobTH at a too-low level!**

562 In a software update, we might try to automatically block usage of that type of  
563 Automation for 2 hours, after it was once used.

564

565 5.2.3 Temporarily exiting the FCL

566

567 The “last resort” alternative always is to **temporarily** leave the FCL mode, and handle any  
568 disturbance “the traditional way” in **hybrid closed loop**.

569

570 For this, we switch the automatic aggressive adaptations of ISF to the bg curve OFF that are only  
571 needed in FCL ....

572 (if in hybrid closed loop you like e.g. the dura\_ISF adaptation still, you alternatively could elect  
573 to just set bgAccel\_ISF\_weight temp. to zero, instead)

574

575 Do not forget that, before meal starts, giving a bolus will then be necessary again.

576

577 The suggested improved FCL cockpit user interface with an extra version of violet loop on the  
578 AAPS home screen ([section 5.4.1](#)) would facilitate this transition FCL < - > HCL, including  
579 automatic removal and re-appearance of the insulin button at the bottom of the APS home screen.

580

581 In case [this feature](#) is not yet available, you must:

582 Exit the FCL mode by going to AAPS/preferences/[put in your password](#) \*)/OpenAPS

583 SMB/scroll down to autoISF settings and switch "Enable ISF adaptation.." OFF

584 (or, alternatively, set bgAccel\_ISF\_weight to zero).

585 *\*) if you set a short password (recommended!) to avoid accidental clicks*

586 **Caution:** Unfortunately, there is no way yet for your full closed loop with ISF adaptations to come  
587 automatically back on, after a selected time for instance. So **do not forget to switch your**  
588 **autoISF fully back on, later.**

589

590 As this will often be forgotten, it may be worth doing a "User action" Automation, for a "temp.  
591 FCL OFF" grey button (see [section 5.2.2.3](#)).

592 Caution though, there is very limited experience with this brand new feature. Make sure your  
593 Automation definition really applies a duration (or other condition) that will automatically  
594 terminate all settings changes it made. As we have seen e.g. in [section 5.1.4](#), this is not always  
595 the case.

596

597 To recognize whether autoISF currently runs with **ISF adaptation** or not, please consult the  
598 **"ai: %" indicator** below the Autosens% **on the AAPS home screen.**

599

600 From autoISF 3.0.1 onwards, there is also a very easy way to see effective ISF and effective iobTH  
601 in the 1<sup>st</sup> screen of the **SMB tab**. At the same time, there you see the adaptation of sensitivity to:

602 • a set %profile change (or effect from Autosens, in case you have that activated)

603 • a set temporary target

604 • the Activity Monitor

605 • +/- exercise mode

606 So, occasionally (especially in your early set-up phase, after starting of a meal) it is a great idea to study the  
607 SMB tab to find out what is going on. See example given in [section 5.3.6](#)

608

## 609 5.3 Recognizing your loop state in the AAPS home screen

610

### 611 5.3.1 Modulated loop aggressiveness via the 3 top buttons in the AAPS home screen

612

613 3 Buttons (%profile; exercise; TT) allow temporary modulation of aggressiveness (sensitivity) as  
614 may be needed e.g.

- 615 • going into meals (EatingSoonTT => providing lower bg “starting point” and more pos. iob),
- 616 • or doing exercise (Exercise mode, limiting basal, iob etc for hypo prevention).

617 ,

618 All buttons grey (GGG) indicates: The loop is running using un-modulated profile parameters

- 619 • Note that features like autoISF might still build on profile (regardless whether original, or a  
620 temp. modulated profile) for further modifications (See calculation cascade for used  
621 sensitivity (sens) in the SMB tab. Example see below in [section 5.3.6](#) ).

622 Any **yellow** (y or Y) button indicates you are running with a modified (*e*levated, or **l**owered) **profile**  
623 **sensitivity**.

624 *Actually, the middle button could be yellow without being active, see table below*

625

	<b>Less aggressive</b> loop = more sensitive user needs temp. less insulin	Standard loop using the set profile	<b>More aggressive</b> loop = less sensitive user needs temp. more insulin
% profile	under 100 % (y)	= 100% (G)	above 100% (Y)
Exercise button	ON (y)  Note that it will only work in combination with elevated TT	OFF (G)  <i>If ON (y), it is <u>not</u> active but ready to automatically become active whenever any TT above target is set (e.g. by any Automation)</i>	OFF (G)
TT	above profile target (y)	= profile bg target (G)	below profile target (Y)
Example	Exercise mode with weakened %profile, yyy		EatingSoonTT, GGY
..used e.g. in :	Case study 6.2 (yyy)		Case study 6.2 (G_Y via Automation) Case study 5.2 (YGY)

626 [Case study 6.2](#) uses extremes from both ends of this sensitivity modulation spectrum:

- 627 • yyy for reducing basal and iob, and getting weaker ISF, during exercise
- 628 • GGY for temporary “EatingSoonTT” boost when, in FCL, a high carb meal had

- The even stronger YGY case is contained in [Case study 5.2](#), where for a high carb snack an extra boost, using >100% profile sensitivity, is temp. added

### 5.3.2 Color scheme of top cockpit buttons tells kind of closed loop that is running

3 Buttons (%profile; exercise; TT), each in 2 states, yellow (for modified), else: grey (G)

Note that yellow could be less (y) or more (Y) aggressive than standard profile

Any yellow setting modulates sensitivity *or attempts to modulate, whenever another required condition may become true (e.g. also via an Automation)*

**Note, though, there are EXCEPTIONS (see below, under GGG, and under Gyy) where sensitivity is auto-adjusted without the user, or an Automation, making a different setting on any of the 3 top buttons. Autosens, or the Activity Monitor, can adjust the %profile sens without the top left button turning yellow or showing that %value. (See sens field in AAPS main screen, or SMB tab!)**

- GGG is the loop running with set profile.

**EXCEPTION:** If the Activity Monitor or Autosens are running, profile sensitivity could be adjusted to (in)activity despite the %profile button remaining grey (see top of SMB tab text from AAPS main screen).

Liikewise, Autosens can temp. tweak the set profile, with the button still remaining grey,

Overall, 2 exp 3 = **eight principal FCL modes** are possible. They need to be **differentiated further**, based on whether the modulation...:

- ...goes into the **more aggressive (Y, higher % and/or lower TT)**, *or...*
- ...the *less aggressive (y, lower % and/or higher TT and evtl exercise)* direction.
- Note the **Exercise button** works only one way, to make the loop less aggressive (“y”). However, whether the dynamic Exercise Mode (which strongly can adjust sensitivity, see [section 6.1.3.1](#)) is at work, **can~~not~~** be recognized by the color of the Exercise button:

Dynamic Exercise Mode (dyn.EM) ?	bg target as in profile: grey	TT > profile target set (Yellow)
Exercise : grey	dyn.EM off	dyn.EM off
Exercise: yellow	dyn.EM <b>off (!) but “on in waiting” = after</b> an Automation sets a TT>profile target	dyn.EM on

657

658

DEV idea: Should we make the **Y** to violet or orange (more aggressive) and keep the **y** yellow (less aggressive

659

than profile)? Then beginners would easier recognize whether they go more or less aggressive ...

Color combinations	Lower aggressiveness	.....	Higher aggressiveness
GGG		profile	
Gyy .... GyY	dynamic exercise mode (if TT > profile target)		EatingSoonTT (as TT < profile target: middle y inactive, GyY = GGY 2 )
	<b>EXCEPTION:</b> Profile sensitivity is auto-adjusted in dynamic exercise mode, but this leaves the 1 <sup>st</sup> button grey (G)		
yGy-yGY- YGy-YGY	„traditional“ exercise mode (if <100% sens AND TT > profile target)	<b>mixed cases 1)</b>	boosted EatingSoonTT (if >100% sens AND TT < profile target)
GyG		GyG = GGG 2)	
yGG .... YGG	Loop adjusted for period of lower insulin need (e.g. 80% for night after exercise)		Loop adjusted for period of increased insulin need (e.g. 120% during sickness)
	<b>The set %profile might be further adjusted by Activity Monitor or by Autosens</b>		
GGy ..... GGY	Loop running with temp elevated target e.g. for enhanced safety against going low (HypoTT)		EatingSoonTT (if TT < profile target)
yyG ..... YyG	same as yGG	y inactive 2)	same as YGG
yyy - Yyy- yyY- YyY	dynamic exercise mode further modulated by <100% profile	<b>mixed cases 3)</b> yyY = yGY, 2)	Boosted EatingSoonTT (if >100% sens AND TT < profile target)( YyY = YGY, 2)

660

661 1) **yGY** would be the combination of setting a milder % profile, yet orienting the loop to a lower  
 662 target. Usually we do not use this combination, but it would be OK for an instance where we  
 663 want to start the next meal at a low target, but we want our loop to be careful and not rush  
 664 getting there.

665 **YGy** would be the combination of setting a more aggressive % profile, yet orienting the loop to  
 666 a higher target. Usually we do not use this combination, but it would be OK if we believe we are  
 667 in a state of enhanced insulin need, and want to get to a set elevated TT (which could be what  
 668 we want to go max down to to have room for exercise, or just to be sure not to go into a  
 669 hypoglycemia if we misjudged the symptoms about enhanced insulin need.

670 2) GyG is the same as GGG, differing only *if* (e.g. an Automation) sets a *TT>profile target*. The  
671 GyG setting *only then* provides softer loop response, due to the exercise mode kicking in.

672 Generally, a yellow middle button y works really only in the combination ..yy. Any combination ending yG  
673 or yY will work the same as if it were ending GG or GY. Although the yellow middle button does not  
674 make a difference in these cases, it might be “justified” by the loop showing its readiness to employ the  
675 exercise mode, whenever a TT > profile target comes around, e.g. via an Automation.

676  
677 3) yyY would essentially working like yGY (discussed under 1), because the middle button,  
678 whether y or G, makes no difference with an elevated TT.

679 Again, this the combination of setting a milder % profile, yet orienting the loop to a lower target.  
680 Usually we do not use this combination, but it would be OK for an instance where we want to  
681 start the next meal at a low target, but we want our loop to be careful and not rush getting  
682 there.

683 Yyy would be the combination of setting a more aggressive % profile, yet orienting the loop to a  
684 higher target. As opposed to the YGY case (discussed under 1), in Yyy the Exercise mode is  
685 active and strongly reduces aggressiveness. It is hard to imagine, why this combination this  
686 combination of settings could make sense; maybe if a “long waved” state of reduced insulin  
687 sensitivity meets a case where normally Gyy would be sought (see there).

688

689

690 5.3.3 Active Automations (or middleware in Trio, iAPS)

691

692 **Active** Automations can temporarily set different aggressiveness will automatically adjust to the  
693 color scheme discussed in the previous chapter.

694

695 It is important to be aware also, which of your Automations (see listed in your AUTO tab, or via the  
696 top left burger menu, in AAPS) is **active to eventually “kick in”**.

697 The ones where you clicked “User Action” should also show as an **extra grey button** on the  
698 bottom of your AAPS screen (only in the suitable time-of-day bracket, if you assigned one).

699

700 Note that an **Automation might not be permitted** to change settings if still **another Automation**  
701 **is running**. (Always consider that, **try to use short durations** in your Automations!).

702 For instance, you cannot switch from 130% profile to 110%. Either the 130% times out, or you

703 **need an extra “in-between” Automation that terminates** the 130% under described conditions

704 (example see in [Case study 6.2](#)).

705 This “design” is for a good reason: The assumption here is, that your 1<sup>st</sup> Automation (the

706 130% in the example) is designed well and runs for a reason. It should either “get finished”



when the job might be done (and kick in again, if not), or, in exceptional cases, it should be consciously terminated by another well thought through 2<sup>nd</sup> Automation (describing the conditions in which you would find that other Automation more important than “finishing up” the one that was already running). That “in-between” Automation makes the loop return to base profile, which is a signal *to all Automations*, to now check whether any conditions exist, to activate a 3<sup>rd</sup> Automation (as in example of [Case study 6.2](#)).

Advice: 1) Do some “house cleaning”: Occasionally check which of your Automations might work with shorter durations assigned. Reduce your long list of Automations, or at least de-activate those that will not be needed.

In case you have many Automations, you could make it an evening routine to activate only those Automations that you might need next day.

2) Try to stay away from Automations that also aim at temp. modifying aggressiveness (especially if triggered by bg level). Often, they will not kick in anyways. And generally, it is not a good idea, to “double up” sub-algorithms for tweaking loop behaviors (“loop inside a loop”).

#### 5.3.4 FCL related indicator fields in the AAPS home screen

In extra data fields of the AAPS main screen you can always see (not change) the key „aggressiveness“ parameters your loop currently operates with (see also home screen example below):


- To recognize whether autoISF currently runs with **ISF adaptation** or not, please consult the “**ai: %**” **indicator** below the Autosens% **on the AAPS home screen**.
- Details for every loop decision see result/debug section of the **SMB tab**.
- The AAPS home screen additionally shows, above the deltas, the current **acceleration**

Having a look at that can be valuable. For instance, when glucose is relatively low and still falling, a positive (and getting more positive) acceleration indicates that bg will swing back up, rather than crash low. This will give info about necessary snack size, and hence help avoid both, unnecessary calories, and going on a bg roller coaster.

740 5.3.5 Overall home screen: (do a new picture

741

Overall home screen:

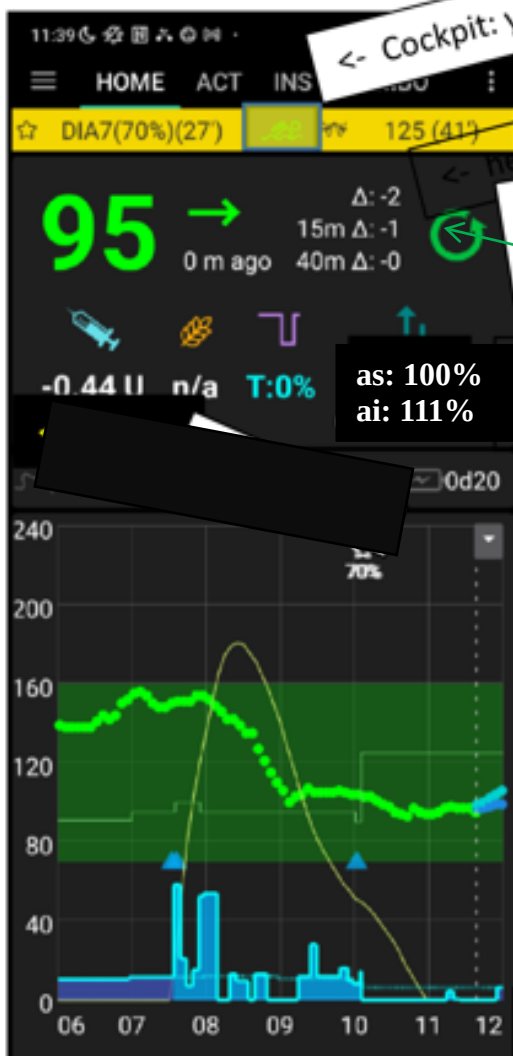

<- Cockpit: yellow fields=>temp. modulated sens.

<- here additionally: acceleration-factor

as: 100%  
ai: 111%

<- Autosense status  
<- Factor resulting from autoISF

<- buttons „bolus“ „carbs“ etc. eliminated

742

743

744

### 5.3.6 Info given every 5 minutes in the SMB tab

(If **tab** not shown -> Config.builder/APS: tick box on the right)

When clicking on the SMB tab, you see how your standard and temporary settings, as well as the latest bg and iob status, influenced the last decision of your FCL.

Example 1: A 80% temp. profile modulates 60% iobTH to 80% of 60% = 48%:



The profile ISF of 41 mg/dl/U got modified by the set 80% temp. profile to  $41/0.8 = 51.3$  mg/dl/U, called “ISF unchanged” (before “start autoISF”).

autoISF applies the 4 sub-categories (acce, bg, pp and dura\_ISF), and depending on the bg curve form suggests various ISF changes.

The final resulting factor “sens” (see flowcharts in [section 3](#)) is 1.11 (in our case, driven by bgAccel\_ISF). This changes the 51.3 “unchanged” ISF to  $51.3/1.11 = 42.4$  mg/dl/U

Further down in the SMB tab, you can see how this ISF is applied to define the SMB size to be given, and whether any limitations – notably by autoISFmax, max possible SMB size, or maxIOB – cut the amount.

```

COMBOV2 SMB AUTO NSI
Last run : 3/30/24 22:18

Result

Script debug : d:
  Activity monitor disabled
  inactivity detection: sleeping
  hours
  Autosens ratio: 1;
  Basal unchanged: 0.48;
  ISF unchanged: 51.3
  CR: 10

  -----
  start autoISF 3.0.1
  -----
  User setting iobTH=60%
  modulated to 48% or 4.8U
  due to profile % and/or
  exercise/activity mode
  SMB enabled; current target
  90 is even number
  Loop allows maximum power
  acce_ISF adaptation is 1.11
  bg_ISF adaptation is 1
  pp_ISF adaptation is 1.09
  dura_ISF adaptation is 1.03
  because ISF 51.3 did not do it
  for 15 m
  final ISF factor is 1.11
  -----
  end autoISF
  -----
  currenttemp: 0.34
  lastTempAge: 0 m
  tempModulus: 28 m
  
```

Message	Condition	What does it affect?
Loop allows maximum power	even target < 100	increase in bg limited to 30%, otherwise no SMB; actual SMB delivery ratio is max of fixed smb_delivery_ratio and linearly growing ratio
Loop allows medium power	even target >= 100	increase in bg limited to 20%, the AAPS default, otherwise no SMB; actual SMB delivery ratio is either fixed smb_delivery_ratio or linearly growing ratio
Loop allows minimal power	odd target	no SMB, only TBR available for action
Loop power level temporarily capped	IOB > effective iobTH	Last SMB capped to stay below iob threshold + 30% overrun; IOB getting above user defined iobTH, potentially modulated by exercise mode, activity monitor and profile percent
Loop allows AAPS power level	no even/odd target option active	SMB enabled/disabled according to standard AAPS rules and settings; no iobTH threshold is active

Note that in the SMB tab you can– in “real time” - capture and analyze *only one* decision.

### 771 5.3.7 Info about last 15 autoISF decisions

772  
773 Refer to [section 11.2.2](#) for an option that enables extended analysis of the on-going ISF modul-  
774 ations from autoISF. (To do this on your **Android** loop phone requires QPython and a logfile  
775 emulator).

776

777 You get tables like this ...

17:05

...

<

778  
779 ... which gives you a quick orientation about recent loop decisions, and relative contributions of the  
780 various autoISF contributors (ace, pp, bg, dura).

781

782 A similar table is available also on the **iPhone** for Trio and iAPS users of their autoISF variants.

783

### 784 5.3.8 SMB tab info when operating in 1-minute mode with Libre3

785

786 Users: Please supply text and screenshot

787

788

789

### 790 5.3.9 Summary: Your personal FCL cockpit (for maneuvering through disturbances)

791

792 A lot of avenues were shown that could help you or your loop maneuver through a variety of  
793 “disturbances”. You should not have to try out many of them, and (like the author), you should find  
794 a way to narrow it down to what really helps in **your** everyday T1D management.

795

796 Try to **keep things as simple and clear as possible**.

797

798 Especially, do not pre-maturely rush into setting up Automations as quick over-patches for what  
799 you may not like to see. Limit the number of Automations, and further limit which ones from your  
800 list are every-day active (vs. switched off, and only ticked active for special days).

801

802

## 5.4 Modulating aggressiveness manually from the improved FCL-cockpit

***Skip this section 5.4 (next ~8 pages. Continue with section 6)** unless you are deeper interested in discussing further user interface upgrades. Actually, some suggestions made are probably an "over-design". After trying a lot of options for refinements out, the author returned pretty much to a "keeping-it-simple" route.*

*My main suggestion is to get that violet loop button (sections 5.4.1-5.4.2), something I think many would use - very handy certainly in the setting-up stage, too, for easy switching between the "old" HCL, and new territory in FCL.*

autoISF is an early dev variant of AAPS, and as user you are participating in an on-going development. Of note, autoISF 3.0.x was launched without many of the cockpit features that are suggested below in green font color.

Only what is written in black is currently of some relevance for using autoISF.

**No need to read** any of the green lines, unless you are interested in contributing to define/design/program further improvements.

This is also an open invitation for you to contact us in case you could help program a module for one of the suggested user interface extras.

For future integration into AAPS Master, an eye should be kept also on the question which other modes (like FCL using Automations and others mentioned in [section 13](#); and maybe also HCL) might benefit from some of the extra features.

Keep in mind, that the **goal should be to interfere with the loop as little as possible**. Under certain conditions, it can run **fully automatically** without any user interaction, as described in the preceding [section 4](#). + [section 5.1](#).

Just like in the airplane **cockpit**, also cruising in full auto mode should involve having an eye on the instruments, and on potential disturbances ahead in the environment.

*E.g.: storm ahead => instruct your plane to climb to another flight height.*

*Analogy: exercise ahead => setting an exercise TT, or => pressing a button that activates a sequence of instructions (some of them probably hinging on conditions, like actual iob), how to manage through that exercise situation).*

So, **for the occasional „disturbance“ coming up**, you should find an easy way to

- call up a pre-programmed routine for automatic management, with auto-adjusted aggressiveness, or:
- tweak a setting or two, to temporarily adjust the aggressiveness

844 • There may also arise a desire to just exit the FCL mode, and “be your own captain” for  
845 mastering a special situation.

846 All this is facilitated within seconds right from the AAPS home screen's **cockpit features to the**  
847 **extent they are already incorporated, or** to the extent you can build **alike** DIY cockpit features via  
848 Automations, as described in [section 4.1.3](#) and [case studies 5.2](#) and [6.2](#)):

849

850 • The button that is integrated into the **violet FCL icon** serves as emergency off button, to  
851 quickly stop FCL, or to at least to immediately stop any more SMBs (...just for a couple of  
852 minutes, or for the remaining meal time: pick from the options offered with just one  
853 keystroke).

854 Via the violet FCL icon on your AAPS home screen, you also can access a temp. switch-off  
855 button for SMBs (see section that next follows below).

856

857 • The **three top fields** (%profile, exercise, TT) provide access to temp. tuning of core  
858 parameters, and/or to some pre-programmed routines.

859 Taken together with some **new indicator fields** about your loop state ([section 5.3.3](#) and [5.3.4](#)),  
860 and the **grey DIY cockpit buttons** ([section 5.2.2.3](#)) this makes the AAPS home screen your  
861 **cockpit** for Full Closed Looping.

862

863 Let us look on each of these **suggested** cockpit elements in some detail:

864

#### 865 5.4.1 Violet FCL icon and underlying buttons

866

867 Novices to FCL, or really anyone running into a very special situation, may appreciate that the new  
868 closed loop icon on the AAPS home screen in pink (for FCL) has buttons to quickly shut off getting  
869 more SMBs (1st row), or to enter other loop modes (second row).

870

871 It functions very much as the other ones that you know from HCL already, and in fact you  
872 get offered some of the same options (for instance, to switch the (full) closed loop off for 15  
873 minutes for going to take a shower)

874


875 Note that in FCL you leave all BG regulation, notably against meal spikes, to the loop. So, try not to  
876 disconnect in phases when your FCL must ramp up your job.

877 The required insulin would still be supplied *after* you reconnect. However, without the user  
878 pre-bolussing, the delay would be more of an issue in FCL than it had been in HCL.

879



880 Just pressing on the FCL icon, a dialogue box comes up:

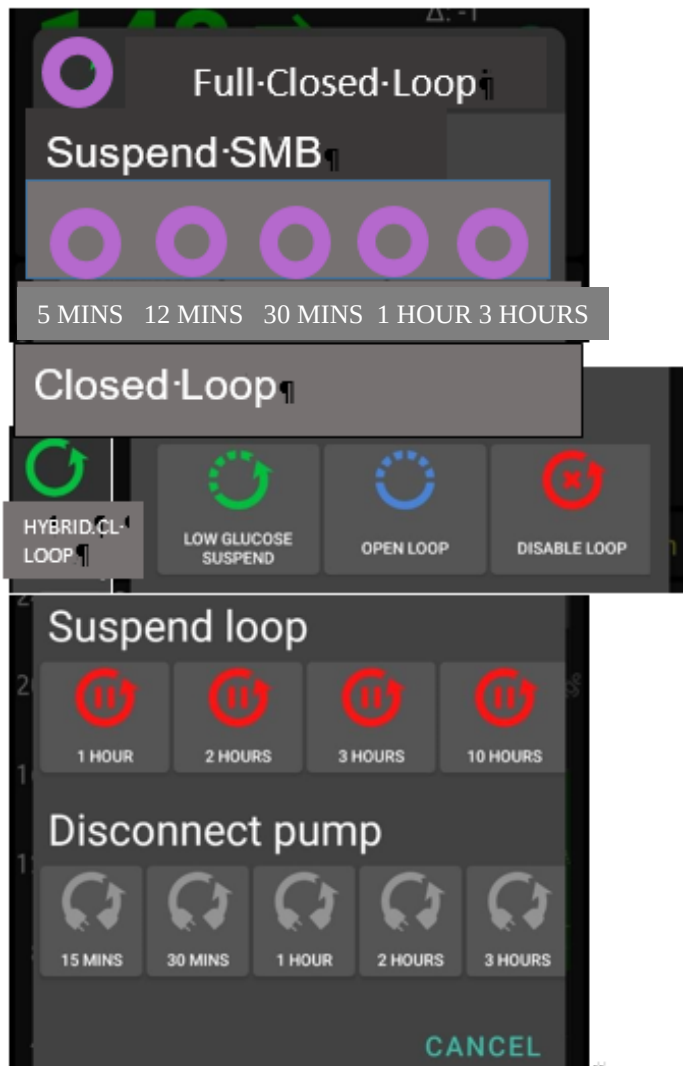

<-add-arrow, form-like-like-CL-Loop...¶

<-dotted+arrow, form-like-like-LGS-¶

DEV: when in Hybrid-Closed-Loop = top headline + green circle, there is no Suspend-SMB part but it starts right with „Closed-Loop“. In that screen. There, the 1st element must be option to pick ¶

.....pink circle, FULL-CL-LOOP-¶

where here the green-HYBRID-CL-LOOP stands-¶

For loopers who did not set up FCL, a feedback must come up ~ „FCL not installed“ if they press on that button ¶

Skip:

= This User Interface is a not implemented idea only.  
Please contribute, or wait for update with the missing info

881

882 Pressing „**Suspend SMB**“ provides fast and easy „emergency braking“ regarding delivery of more  
883 SMBs:

884 Select the one with the desired number of minutes: 5 or 12 for just blocking the potential next  
885 SMB(s), and up to 3 hours to manage the entire rest of this meal with %TBR from then on.

886

887 Whenever, and why-ever, your FCL is in „no SMBs allowed“ mode (e.g. automatically after  
888 surpassing an iobTH, or triggered by a set odd TT), the FCL icon will turn into a dotted one.

889 Instead of remaining **duration to end time** it indicates in the middle „the condition“, „**iob**“ or „**TT**“

890 Add an indication if suspend SMB comes from an Automation, e.g. add an „**(A)**“, **underneath** the  
891 #minutes, iob, or TT in the middle of the dotted violet field.

892 So, as in other (already in HCL existing) cases, those icons show in the middle the minutes left that  
893 they will be running, or the condition which would have to go away for this temp. setting to stop.

894 It always auto-reverts into the FCL state and FCL icon, when time (or other condition) has elapsed.

895

896 Pressing „**HYBRID CL. LOOP**“ or other buttons from the 2nd row provides fast and easy

897 „emergency **exit**“ into other modes.



898 This enables beginners an easy „temp. escape“ into their well-known HCL (green) at any  
899 point of time. bgAccel\_ISF\_weight is set to zero when going FCL->HCL. HCL can run with  
900 autoISF (for instance dura\_ISF) uninhibited otherwise. (check implications for HCL users of  
901 autoISF ?? ).

902 Note: These options from row 2 have no time limit. Loop will **not** by itself go back to FCL. You see  
903 the different loop icon as a reminder to manually revert, when ready.

904

905

906 5.4.2 Buttons „Insulin“, „Calculator“ etc at bottom of AAPS home screen

907

908 These buttons are **not useful any longer in FCL**, and automatically disappear whenever in FCL  
909 mode (also in Suspend SMB state), and re-appear when leaving FCL. This applies also when an  
910 Automation or technical system failure shut off FCL.

911 Users who, maybe in the beginning phase, feel better having those buttons, can override  
912 the removal (of the insulin button, or any other) by going into /preferences/overview/buttons  
913 and forcing them on. They only remain on until the next re-entry into FCL mode, when auto-  
914 off happens again.

915 The reason why we do this: It really is important to let the loop loop, and not interfere more  
916 than absolutely needed. Any bolus the user gives will sure distort the bg curve, on which  
917 autoISF, especially when aggressively tuned for FCL, builds a lot of its decisions!

918

919

920 5.4.3. Three top fields (%profile, exercise, TT)

921

922 Depending on the variedness of lifestyle, the desired %TIR, and the initial tuning effort put in, the  
923 user may want occasionally to „tweak“ the aggressiveness of her/his FCL.

924

925 The top 3 fields (grey in default mode, **yellow when temp. in mode with changed**  
926 **aggressiveness**) serve as quick and easy entry points to make temp. switches (as users will be  
927 used to for %profile switches, or for setting an EatingSoonTT in HCL, .. which they still can do in  
928 FCL ... but more:)

929

930 Expert FCL users might need this feature rarely, but probably at least to manage activity after  
931 meals: Each require opposite aggressiveness, and the switch has to come in a certain point in  
932 time that would be difficult to capture. (More see [section 6.4](#) )

933

934 Information printed on the top buttons

935

936 The yellow TT field shows the currently valid TT (and further duration):

937 (profile) stands for the abbreviation you labeled your selected running profile



939 In the special case of settings for meal preceding sports, the field will look slightly differently:



...and ...

941 ... when iobTH is first time exceeded, this automatically switches to:



943

944 Likewise, if on the AAPS main screen just an **EatingSoonTT** is set (e.g.72), this is entered with the  
945 desired duration. Afterwards, it automatically reverts to profile target and the display turns grey  
946 again (can also turn green first, for a short transition period), with the profile target (e.g. 90) on it.  
947 (No time limit, then, for the profile value as set in preferences).

948 Without sports context, the middle field remains grey.



950

951 Independently from setting a TT, the user can choose to set a **%profile in the left top field**, for an  
952 independent number of minutes, e.g. 70% in this screen example: Also, or additionally, this will  
953 influence the resulting ISF and sensitivity%

954



956

957 The % might change and turn yellow also in context of making TT inputs in the related dialogue  
958 box (see chapter TT dialogue field, above). Still, the % (or the length of time the profile switch shall  
959 be active) can be independently overridden in the top left field, if so desired.

960

961 If an **Automation** sets a %profile, and/or a TT (e.g. *automatic detection of meal start at condition*  
962 *e.g. when delta >10*) , this would automatically show in respective field(s) turning yellow and  
963 showing the temp. setting. To show the set parameter comes from an Automation, „ **(A)** „ is added  
964 in the end of button text.

965

966  
 967 5.4.3.1 TT dialogue field ( Currently not available in the pictured form and function ! )  
 968  
 969 The TT field (top right of AAPS home screen) is a primary daily interface, and a dialogue field  
 970 opens when pressing on it

Initially, the form can be entirely empty re. TT inputs and just show the default iobTH and bGAccel parameters.

Two lines appear for target and get labeled ES, AC or HY depending on what was pushed at the bottom buttons. (HY => red frame; evtl. 1 empty => grey)

TT entries can be made or overridden.

iobTH calculates from Target AC and other settings shifts away from default set in /preferences (here 0.7) but could be overridden here 0!..200% is allowed

bGAccel ISF\_weight can be modulated here, too. Note: it can change again if % profile is also changed. -- 0...200% is allowed

CANCEL allows to start fresh (select one or two of the square buttons, ES, AC or HY)

OK needed to use the settings

Skip: = This User Interface is a not implemented idea only. Please contribute, or wait for update describing the eventual implementation

Duration input is made in minutes. In the exceptional case that both, ES and AC targets are defined, the duration input is for AC and framed blue. (This is because the preceding AC mode is automatically determined in length by the loop observing when iobTH is exceeded)

971  
 972  
 973 This looks complicated but only because it allows 4 different modes of use. Each user will primarily  
 974 use her/his preferred one.  
 975 (1) Who is happy with the initially well tuned FCL and does not have huge variations in daily eating  
 976 and moving around, will **not use** the TT **at all**. FCL is possible without an intervention via the

977 TT button in your cockpit. Actually 4 of 8 modes (GGG ...YYY permutations, list see [section](#)  
978 [5.4.1](#)) are not making use of TT.

979 (2) Super easy is also, to just input **any odd-TT** (odd-numbered temporary target) that will shut out  
980 any SMBs for the set duration. *That can be a good idea when having a small snack, for*  
981 *instance.*

982 Super quick access to stop SMBs is possible also via the loop icon ([section 5.3.1](#)).

983 Specifically, an **EatingSoon TT** can be activated here (limited relevance see [section 2.5](#)). It is  
984 time-un-critical, can be manually set, or come up via an Automation.

985 The cockpit enables you to set the iobTH differently (override) for the current meal.

986 Alternatively, iobTH can be temporarily changed in /preferences or using an Automation.  
987

988 Temp. iobTH will always revert to default when the TT expires. If another TT immediately  
989 follows, like in the example of the screen above, it will calculate, (then) show and use a new  
990 temp. iobTH.

991 (3) The third way is to **use the input mask** (if already ncluded in your software version see picture  
992 above) **to freely modulate the loop aggressiveness** for a declared number of minutes. Click  
993 the bottom big square(s): Either HYPO, or ACTIVITY, or EATING SOON, or ACTIVITY and  
994 EATING SOON (example in the pictured screen above). Make or override entries in the offered  
995 fields. Press OK.

996 (4) The fourth way is to exclusively use one of the 4+4+2 little buttons seen in the bottom part of  
997 the TT dialogue box (if already included in your software version). They provide a set of  
998 settings (as will immediately show in all input fields above) that the user has set up in  
999 Preferences/SMB/autoISF/FullLoop (refer to [section 6.3](#)), and can freely label there. *For*  
1000 *instance „hiC“ at high carb EatingSoon, „piz“ for Pizza/fatty meals, „grd“ for garden work,*  
1001 *„mtb“ for mountain biking ...*

1002 **Capturing good settings for not-everyday situations in /preferences** (if already included)  
1003 **allows calling them up within 1 second**, from your cockpit on the AAPS home screen (...and  
1004 won't ruin the FCL experience at all , especially because in most cases it is not time-critical,  
1005 how long before the intended exercise the buttons are pressed).

1006

1007 [Case study 6.2](#) demonstrates that nearly the same performance and comfort can be reached  
1008 via the **DIY FCL cockpit** with the grey extra buttons appearing at the bottom of the AAPS home  
1009 screen, based on Automations with User action (see also [section 5.2.2.3](#)).

1010 The example picture given above, and also [case study 6.2](#), is the most complicated (but also most  
 1011 useful) case, **when exercise follows after a sizeable meal**. It is then that you need (a) aggressive  
 1012 FCL initial performance at the meal, but, *exactly when* (!) a (for the intended sport already  
 1013 temp.lowered) *iobTH* is exceeded, you need (b) to have SMBs automatically switched off and go  
 1014 into the „milder“ mode, as defined for the exercise (with *high* instead of the immediately prior  
 1015 *low*TT, that automatically significantly reduces iobTH again, and insulin sensitivity(resistance)  
 1016 settings too).

1017  
 1018 Pressing exercise related buttons will automatically also light the **exercise button** on the main  
 1019 screen yellow.

1020  
 1021 To summarize, the TT dialogue field offers easy but powerful ad-hoc [modulation of loop](#)  
 1022 [aggressiveness](#) for FCL (if already included).


1023  
 1024  
 1025 [5.4.3.2 Exercise button](#) (see more in [section 6](#).)

1026  
 1027 The exercise button automatically lights yellow when exercise related TTs are activated [in the TT](#)  
 1028 [dialogue box](#).

1029 4 of 8 principal FCL modes ([section 5.4.1](#)) are making use of the exercise button.

1030  
 1031 If pressing on the exercise button, a dialogue box appears (*if extended design for FCL cockpit is*  
 1032 *already launched*) with info on exercise setting first (and opportunity to override), plus below the  
 1033 activity monitor (experimental for auto-tracking of lighter movement during the day, and effects on  
 1034 sensitivity that may have. See [section 4.5](#)).

1035  
 1036 So, first the exercise settings (as set under TT) are there to read. Example :

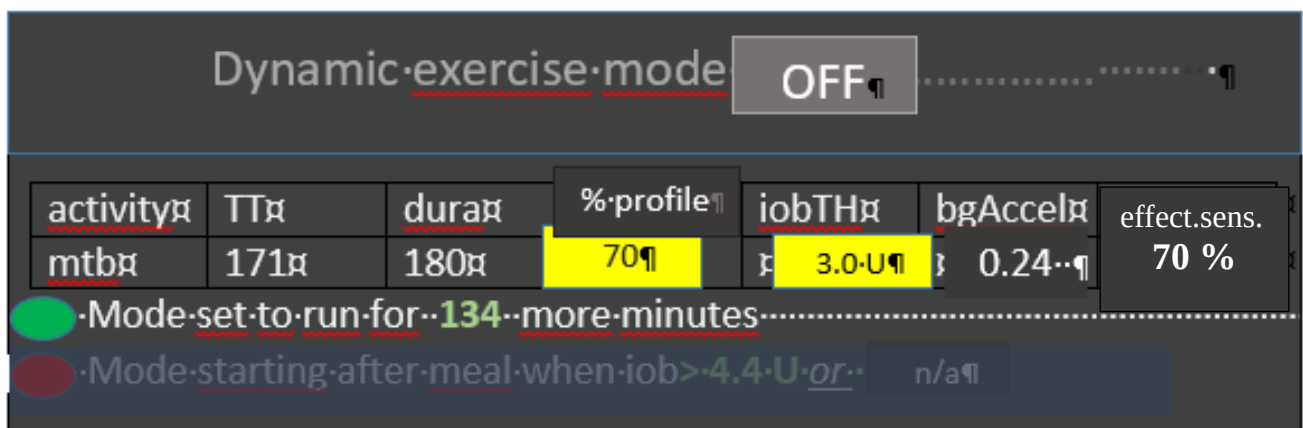

1037  
 1038

1039 The exercise (here mtb) is selected in the dialogue box of the neighboring TT field, and there auto-  
1040 filled with settings made in the set-up and tuning stage by the user under preferences. They are  
1041 reported also under the exercise button here, and TT, duration, and % sens (the temp. profile  
1042 sensitivity that also shows on the %profile field on the left side of the exercise button) can be temp.  
1043 changed there.

1044 iobTH, bgAccel\_ISF and overall resulting effective sensitivity ratio (effect.sens. %) is given in the  
1045 other fields.

1046 The **middle field** of the table, „% **profile**“ either picks up the % set under the %profile button, or  
1047 an input can be made here, in the exercise button domain, which will:

- 1048 • turn the neighboring %profile button on yellow and show that inputted % on it, too
- 1049 • be multiplied with the result from the exercise mode settings per se, and change the %  
1050 overall, accordingly.

1051 So, if this middle field of above table (dialogue box of sports button) contains a figure other than  
1052 100, input field becomes yellow, and you are operating with a combination of traditional PLUS new  
1053 exercise mode (with all three top buttons of your FCL cockpit yellow). This maximally will soften  
1054 aggressiveness, for which you get an idea by the last calculated figure.

1055

1056 The mode is either running already (for another number of minutes, as probably also shown in the  
1057 yellow TT field anyways). Or it is scheduled to run, after insulination for a started meal reaches  
1058 iobTH (as in table). Or, no exercise is scheduled (both points red, no entries).

1059

1060 The lower part of the exercise dialogue box (not pictured above, but see in [section 6.5](#)) is  
1061 dedicated to the Activity Monitor

1062

#### 1063 5.4.3.3 Profile button

1064 The profile button can still be used to set a different profile, or profile%, for instance to adjust for  
1065 days with sickness (as you are used to from hybrid closed looping). 4 of 8 modes are not making  
1066 use of the profile button.

1067

1068 Any inputs made here will be used to modify profile\_ISF on which all further changes are made on  
1069 (multiplied with).

1070

1071 The profile field remains grey if standard profile is applied.

1072 It turns yellow, displaying a %figure relating to any altered loop overall aggressiveness:

- 1073 • When no inputs (changes from 100% profile) are made here, but inputs in the TT field,  
1074 e.g. for exercise, automatically lead to different effective sensitivity ratio



1075       • when% is changed by input in the profile button itself, it will be multiplied with with  
1076       profile\_ISF and be used in place of profile\_ISF *by the algorithm*.

1077 However, for exercise (sports) you no longer must make an entry here, because  
1078 reasonable %reductions should be automatically provided, driven by your set TT (and half-basal  
1079 exercise target), see [section 6](#).

1080

#### 1081 5.4.4 FCL related indicator fields in the AAPS home screen

1082

1083 In extra data fields of the AAPS main screen you can always see (not change) the key  
1084 „aggressiveness“ parameters your loop currently operates with (see also home screen  
1085 example below):

1086       • how profile sensitivity (**ISF**) adjusts by the %profile input, by autoISF, and/or a set  
1087       exerciseTT, resulting in an effective sensitivity (ISF that is used to determine  
1088       insulinRequired. Details for every loop decision see result/debug section of the SMB tab).


1089       • next to current available iob number is an indication of your **valid iobTH** (the iob above  
1090       which no more SMBs will be given)

1091       • The AAPS home screen additionally shows, above the deltas, the current **acceleration**

1092       Having a look at that can be valueable. For instance, when glucose is relatively low and still  
1093       falling, a positive (and getting more positive) acceleration indicates that bg will swing back  
1094       up, rather than crash low. This will give info about necessary snack size, and hence help  
1095       avoid both, unnecessary calories, and going on a bg roller coaster.

1096  
1097 5.4.5 Overall home screen with discussed further features:  
1098  
1099

Overall home screen:

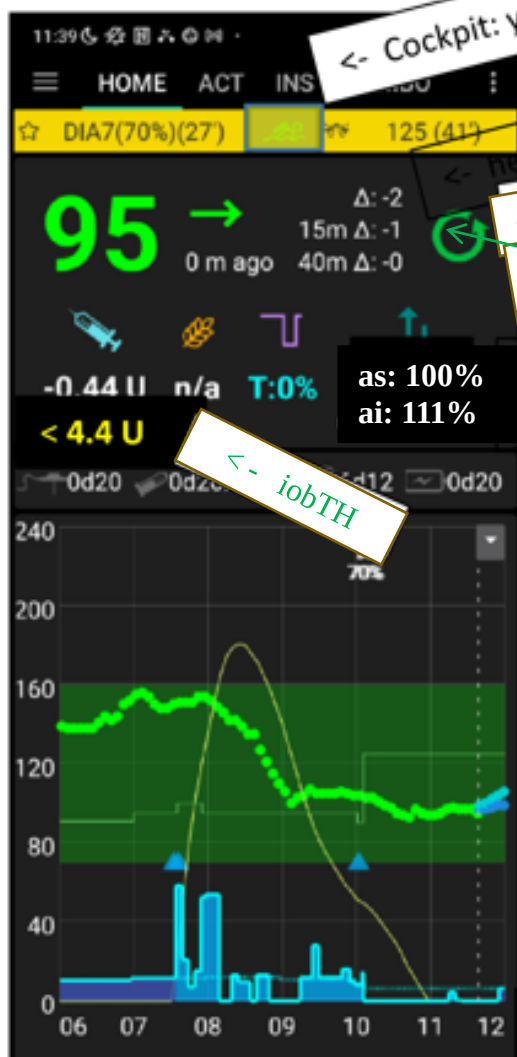

<- Cockpit: yellow fields=>temp. modulated sens.

here additionally: acceleration-factor

<- Violet <-> green circle for FCL <-> HCL dotted if SMB off  
in the middle: minutes counting down if temp.set;  
below: (A) if coming from an Automation that is running

as: 100%  
ai: 111%  
<- Autosense status  
<- Factor resulting from autoISF

<- iobTH

Note: iobTH is currently not shown here, but can  
easy be found in SMB tab, see next page

<- buttons „bolus“ „carbs“ etc. eliminated  
(auto- re-appearing when violet -> green loop)

1100  
1101