

Case Study 1.4: Importance of having proper BT connections, notably after a meal started

V.2.0




Here an incidence report that emphasizes the high importance of having proper Bluetooth connectivity with the pump (and with the CGM), notably after starting a

Glucose going higher than usual after starting a meal

After a ~12:30 h lunch, I saw glucose rise, by 13:40 h, to about 230 mg/dl, which is an unusually high level for me.

My AAPS home screen ([graph 1](#)) reveals a super low iob (0.01 U) and the yellow insulin activity curve seems missing. Also, there are no blue triangles from SMBs. Nearly no insulin had been given despite glucose rising sharply to over 200 mg/dl:

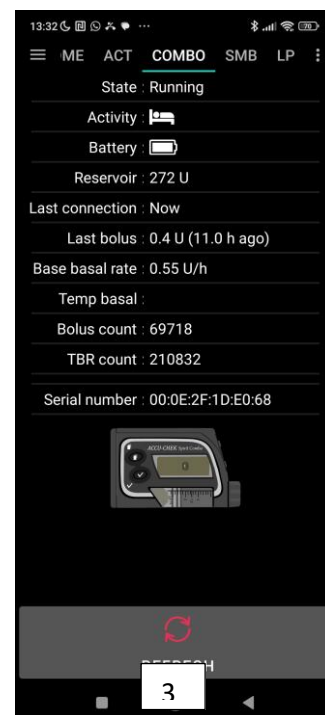


14

The dark blue block (instead of jittery TBR curve) for basal hinted at a pump problem, which I could confirm looking into the pump tab ([graph 2](#)).


By pressing Bluetooth on my pump, and “refresh” at the bottom of the pump tab in AAPS, the pump connection could immediately be re-established ([graph 3](#)).

22 The SMB size was surprisingly low, though. Graphs 4 – 7 show why that is:

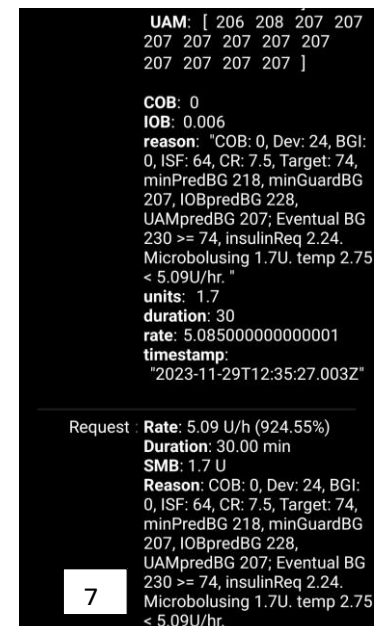


33 The meal had been such that 2 solid hours (~ 13 to 15 h), worth of 30 g max carb absorption
34 each, had to be expected. Therefore, about 45 minutes after my given 4.2 U bolus, the

glucose curve turned a bit upwards again. Based on the detected significant acceleration, autoISF finally was in a position to issue a significant SMB of 3.2 U at 14:20 (02:20 PM). Graphs 8 and 9 show that this SMB at 02:20PM was important to get glucose back into range.

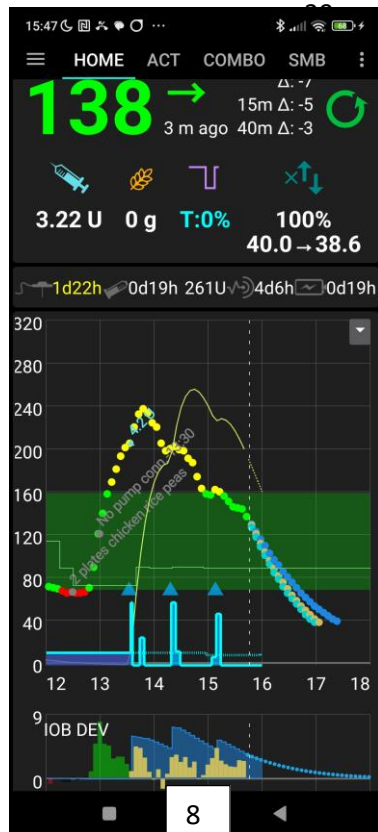

8

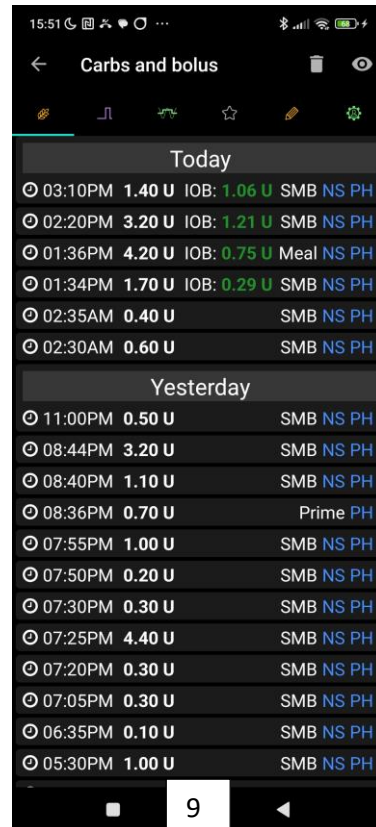

9

52

Conclusions

A no-bolus FCL absolutely requires intact pump connection, notably in the hour after any meal start.

Same of course applies to needing un-interrupted glucose values. But these are less likely to go unnoticed for a while, because lacking glucose values are easy recognized from just glancing on screen or watch occasionally.

autoISF „designs“ SMB response to be aggressive only when acceleration and big glucose deltas are seen. However, this phase might be missed while the pump had lost connection to AAPS (and gave only basal).

Upon re-connecting, the loop response might be significantly weakened by falling into a decelerating phase of the bg curve.

In this situation, to issue a user bolus (and immediately after to go into FCL again) can help get back on track (into range) faster.

Avoid giving a bolus just under the valid iobTH, though: Your relative big „user bolus“ will, even with Lyumjev, take a couple of minutes to have any effect. So, in the meantime, you probably will see further bg rise, with a chance that a fairly big additional SMB might be triggered by the loop and significantly surpass iobTH.

Avenues to avoid the problem

Avenues to avoid problems from temporary loss of BT-connection to pump could be:

- It is probably worth defining an alarm for lost connectivity to the pump during daytime (at night, an on-going basal should be OK, too).
This could be done via an Automation: Conditions: Time between 07 and 22 h AND Last connection to pump greater than 12 min ago ACTION: Notification and/or Alarm
- The easiest alternative is to **make use of the „safety“alarm setting in AAPS.../Preferences /Local alerts** (one of the last headlines)

Besides a pump unreachable after (your input) minutes without connection, you can define an alarm for lacking CGM values there, as well.
- An Automation could build also on the observation $bg > 140$ and iob remained near zero ...