

1. Pre-requisites for Full Closed Loop

V 3.7

Please note that with autoISF you are in an early-dev. environment, where the user interface is **not optimized for safety** of users who stray away from intended ways to use. Good safety features exist, but these are only as good as the development-oriented user understands and implements them. This is not a medical product, refer to disclaimer in [section 0](#)



- 1.1 Well tuned hybrid closed loop
- 1.2 Fast insulin
- 1.3 Reliable insulin delivery from pump and cannula
- 1.4 Excellent CGM
- 1.5 Meal-related limitations?
- 1.6 Lifestyle-related limitations?
- 1.7 Time required for setting-up

Available related case studies:

- Case study 1.1: Occlusion
- Case study 1.2: Comparing insulins for FCL
- Case study 1.3: Jumpy CGM
- Case study 1.4: Lost pump connection
- Case study 1.5: Overlapping 2 x G6
- Case study 1.6: When catastrophe strikes
- Case study 1.7: MDI when the loop dies
- Case study 1.8: Libre3 (1 minute) placeholder

1.1 Well-tuned hybrid closed loop (HCL)

It is advisable to first establish a well-tuned **hybrid closed loop** before considering the transition to FCL. Best if you achieved good HCL performance **without** using Autotune or dynamicISF (which can introduce, or cover up, problems that would get exposed in your transition to Full Closed Loop (FCL); more see at beginning of [section 4](#)).

There are three important reasons for starting on a solid basis (profile):

- The UAM full closed loop requires a highly personalized (individual) tuning of settings, so the loop will give insulin mimicking YOUR successful hybrid closed loop mode.
- Key autoISF tuning parameters multiply with your profile_ISF. You could be “big time starting on the wrong foot”, and potentially waste a lot of tuning time with your autoISF! See for example: red box in [section 4.3.2](#)
- The UAM full closed loop comes with new parameters to be set and tuned. It would be problematic to set and tune several new parameters before the basics were tuned “right”. Errors could easily be balanced with counter-errors. This can work in single scenarios, but would create a highly unstable system. It would be very hard to re-calibrate better later, too.

1.2 Fast insulin (Lyumjev, Fiasp, Apidra?)

It should be clear without saying, but it is absolutely necessary to feed your loop with correct time-to-peak and DIA for your insulin, so it has any chance to know how the job will get active in your body. See “Insulin_DIA...pdf” in: <https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings>

39 If the user does not bolus for meals, clearly a very fast insulin is needed so, upon realization of a
40 starting meal-related glucose rise, the loop has any chance to eventually keep glucose in range (by
41 common definition, under 180 mg/dl (10 mmol/l))

42

43 A modelling study..

44 key findings are summarized in initial section of [case study 1.2](#); for more see: “The artificial
45 pancreas and meal control...pdf” in: [https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-](https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings/blob/HCL-.settings-main-repo-(pdf)/)
46 [settings/blob/HCL-.settings-main-repo-\(pdf\)/](https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings/blob/HCL-.settings-main-repo-(pdf)/)

47 ...can show in quantitative terms, that **faster insulins**:

- 48 • will result in significantly **lower** glucose **peaks** than slower insulins
- 49 • **tolerate** a couple of minutes **delayed** first meal bolus while not incurring unacceptable
50 height of peaks
- 51 • minimize the effect on glucose peak from **different** carb load (**meal sizes**).

52 In conclusion, do not attempt FCL with other insulin than Lyumjev® or Fiasp®.

53 Potential exceptions:

- 54 ○ Being consistently on low carb, on a GLP-1 drug, or with gastroparesis -all of which ease the
55 job of a FCL significantly.
- 56 ○ According to [case study 1.2](#), Apidra® might work, too, but Humalog® would not work well
57 (with “normal diet”).

58

59 1.3 Reliable insulin delivery from the used pump/cannula/insulin system

60

61 **Good tolerance of Lyumjev (or Fiasp): Occlusions threaten the function of the full closed loop.**

62

63 It is very important to have an eye on the time a **cannula (or pod)** is in use (many find **48 hrs** to be
64 the **limit**), and whether hard-to-explain glucose rises happen at ever increasing „fake“ iob (even
65 before a 48 hr routine replacement). (You easily lose 25% TIR that day, see [case study 1.1](#))

66 It is absolutely contra-indicated to attempt FCL coming from leaking pods and associated erratic
67 sensitivity swings that may or may not have been somewhat controlled and tolerable by
68 dynamicISF or other measures when you were Hybrid Closed Looping,


69

70 **Stable pump connection** (see [case studies 1.4](#) and [1.6](#))

71 In FCL you absolutely rely on your pump delivering,
72 without any further delay, the much needed insulin,
73 after any meal start.

74 Hence it is absolutely essential to avoid any
75 problems from a lost Bluetooth connection. In AAPS
preferences / Local alerts, switch alert on!

An Automation similar to the one pictured here →
could also help recognizing eventual problems.

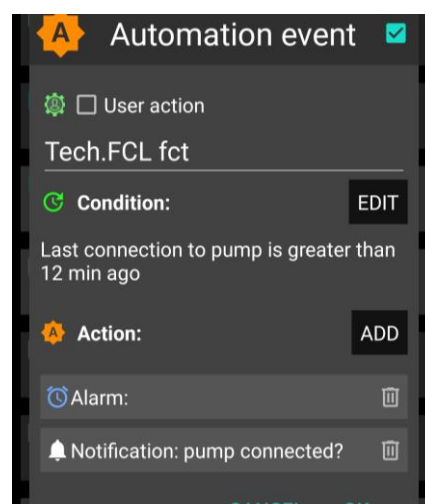

1.4 Excellent CGM

You do not give a meal size-related bolus any longer. That leaves all insulination jobs to the algorithm! Around meals, a **stable Bluetooth** connectivity is absolutely essential, too, so CGM, loop, and pump can do their job without losing more valuable time (see [case study 1.4](#)).

As glucose values are the very basis for your autoISF loop, please **inform yourself well about your CGM:**

- How it principally performs (e.g. you absolutely must be “SMB always-ready” at cob=0)
- Whether you are using the best suited intermediate app that reports the “raw” value from the CGM transmitter into AAPS
- Specifically, how and where any smoothing is done, and what this might imply for the ISF boosting method you will be using See for instance here:
<https://androidaps.readthedocs.io/en/latest/Usage/Smoothing-Blood-Glucose-Data.html>
- Go through your data (in *all* day and also night *times*) to see whether your CGM produces any artefacts (jumpy values; see [case study 1.3](#)) that the loop could **misinterpret** as sign of a starting meal.
For some of these problems, e.g. “jumps” associated with nighttime compression lows, there are options to mitigate (see [section 5.1.2./3.\)](#) . See also the User Action Automation discussed about 2 pages below ([line 149ff](#)).
- In case your CGM requires calibrations: Note that calibrations often produce jumps. In that case, be prepared to do, in the future, an extra handling step to protect from your FCL reacting harshly.

autoISF has also a couple of in-built checks on the quality of the recent CGM values.

They are particularly important in the acceleration detection.

You will learn more about how a **“lousy” CGM performance can badly interfere with your tuning**, and resulting autoISF loop performance, in [section 4.2.4](#) (see red warning box there).

To summarize, a CGM with more scatter will make the loop lose more time, and lead to higher peaks and lower %TIR.

111 So, if you are unhappy with a slow reaction of your loop it could be because the loop is
112 unhappy with your CGM.

113 Consult the detail info given (at the time) in your AUTO ISF tab, or look it up later in the logfiles
114 (using the Emulator, [section 10](#), eventually).

115

116

117 1.4.1 Dexcom G6 and other 5-minute CGMs

118

119 The best proven way to stay out of trouble currently is to use Dexcom **G6**, and to ensure via
120 **overlapping** right and left arm sensor and transmitter utilization always good quality values, that
121 can be used by the Full Closed Loop ([case study 1.5](#)).

122

123 Other ways (using values from just *one* G6, or Dexcom ONE, G7, Libre2, or other new AAPS
124 integrated methods) are possible, but come with a lot of monitoring effort (best via watch), and
125 occasional time-outs for your FCL.

126

127 **Discontinuation** of the **G6** CGM is announced for mid 2026 by the manufacturer.

128 Fortunately, recently several new companies started to offer alternatives, and we do have an
129 analysis of 10 available CGMs by an autoISF FCL user:

130 (if you are in Facebook:) <https://www.facebook.com/groups/AndroidAPSUsers/permalink/4351134601774590/?app=fbl>

131 (else try:)

132 https://cdn.discordapp.com/attachments/968203989546070067/1455662683771895859/Lukas_Trcka_10_CGM_test_RESULTS_2025_12_24.pdf?ex=69578544&is=695633c4&hm=6ef57b3a7e0f51dbd659d938bd709fdb4355d9e108eaa8dbad1bce3fba00621d&

133 [12_24.pdf?ex=69578544&is=695633c4&hm=6ef57b3a7e0f51dbd659d938bd709fdb4355d9e108eaa8dbad1bce3fba00621d&](https://cdn.discordapp.com/attachments/968203989546070067/1455662683771895859/Lukas_Trcka_10_CGM_test_RESULTS_2025_12_24.pdf?ex=69578544&is=695633c4&hm=6ef57b3a7e0f51dbd659d938bd709fdb4355d9e108eaa8dbad1bce3fba00621d&)

134

135 According to his elaborate data, **Syai Ultra** (Syai Health, Singapore) or **Caresens Air** (i-Sens, S-
136 Korea) might be the best alternatives to switch to.

137 Watch for more reports in the Loopers' communities, and have a look also into Tim Street's
138 frequent tests on newly available options.

139

140 Note that, as a FCL user, you must put an elevated emphasis on avoiding jumpy values that could
141 be mistaken by your FCL as start of a meal.

142

143 To secure top performance, you might need overlapping "both arms" sensor use also with the new
144 system. It could be therefore wise to prepare for the transition from G6 to your next system by

145 • Building a bit of G6 sensor stock (be it via saving some, using the Anubis extention; or by
146 collecting some from others who discontinue usage)

147 • Keeping (or acquiring, from those who discontinue G6) an Anubis

- Start using the new sensor system in a fashion that every 2nd sensor is your G6/Anubis. This will allow the usual overlap, while getting to know the relative performance of your old vs your new system. Also, it will help you build a stock of extra sensors from your new CGM system, giving you the freedom to overlap those in the future.
- *If anyone does the mixed use with G6 and a new system before me, please provide material for an earlier update of [Case Study 1.5](#).*

One safety feature in autoISF is a **blockage of SMB delivery whenever delta bg** (within the last two 5 minute values) is **higher than 30% of that bg** (or higher +20%, at bg targets above 100 mg/dl).

Example: From 74 mg/dl, a jump to 97 (+23 mg&dl = + 31% of 74) or more would not receive SMB “response”. From 100 mg/dl to 131 mg/dl (+31) would neither.

Check in your (HCL or FCL) data whether at meals or sweet drinks with rapid absorbing carbs you could run into the problem that jumps are “too high” and much needed insulin will be blocked (only come via very much smaller portions.

For example, 400%TBR @ 0.6 U/h => 0.2 U in 5 minutes, instead of one ~3 U SMB. The difference (of 2.8 U missed) translates @ ISF~ 40 mg/dl/U into up to + 112 mg/dl higher bg peak! It will not become quite that bad, because the loop will catch up to the insulinRequired with its next couple of decisions.

Instead searching in old data, you can also just have an eye on instances where you think a first SMB was due, *but blocked*. Confirm that (by looking in the AUTO ISF tab) and think about a solution that would not require changing the 30% safety limit in the code.

For instance, not drinking so much juice rapidly around meal start could be a “behavioral” correction to get rid of the problem.

This *blockage* (no SMBs) would likely last only 5 minutes (and go probably unnoticed). However, not only would you lose 5 valueable minutes to get your iob substantially elevated, but also, all following deltas are likely much smaller. As a consequence, if the >30% delta was in fact (largely) due to carb absorption, you would, just when needed most, miss some of the boost sought from bgAccel_ISF.

This example also underscores that the CGM in use cannot be allowed random scatter that leaves no reasonable room for safe detection of (smaller and) bigger “truly carb related” deltas

If or when (like: first half day of a new sensor) you are not sure about sufficient CGM performance you might develop for yourself an Automation with User action ticked (along the lines as used for other purposes in [section 5.2.2.3](#)). It would “ask you” before giving a SMB whether you really want it delivered. That way you could a) have a look on your glucose curve b) .. and on the ai % (underneath Autosens%), which indicates the relative aggressiveness of ISF modulation from “what autoISF sees in your bg curve” c) think about what sense a SMB now makes with respect to your past meal, and the carbs to be still absorbed. Ultimately, you could also d) consult some of the detailed info given (every 5 minutes) in your AUTO ISF tab.

Such User action Automations need not be active at all times, but if you have it for your first half day of a new G6 sensor for instance, you could activate that Automation from your list of Automations; after the values have settled in, you can de-activate (“shelve”) it again.

For a brief period, and if you are tech savvy, another way to deal with uncertainty about CGM would be to employ the emulator method as presented in [section 11](#): Run a “too mildly” tuned FCL, and in parallel run a “what-if” with your more aggressive settings that you really would like to use once you are certain about your CGM.

However, I found it **easiest**, to **lay a solid groundwork by using** 1 Anubis, and **two overlapping G6**, to get rid of most problems (...that I keep seeing in my data, on the worse sensor of the two that often run for some days in parallel; see [case study 1.5](#)).

With a sensible **ioBTH** defined (see [section 2.4](#)), and your standard **alarms** for going towards a hypo not silenced, the worst consequence from any automatically “over-treated” glucose jump should be that you **need an unplanned snack** for the balance of “missing” carbs.

Like you should be used to from anti-hypo snacks also in your Hybrid Closed Looping past, you might also in this case here prevent SMBs for a short while by setting a (here: odd numbered) temporary glucose target (TT).

A disturbing and late sign of dealing with a too unreliable CGM (or too aggressive settings) could be a trend towards increasing TDD and body weight during your early FCL experience!

[1.4.2 CGM with 1 minute values](#) (e.g. Libre3)

Libre 3 showed promising performance too, and using it with **1-minute** values via Jugglucoo is integrated, starting with autoISF version 3.0.1.

In summary the recommended data flow is very simple: Juggluco → AAPS → calibration & 1st order exponential smoothing. Set an SMB delivery ratio below 0.5.

autoISF automatically detects whether values come through Juggluco, and adjusts parabola fit calculations to determine acceleration etc. –

When running on Libre3 / 1 minute, you get many more SMBs that each are, on average, much smaller. This has **implications on** some related **settings** (notably, the smb_range_extention, see [section 2.1](#)).

For more info, please refer to the multi-page section “Using 1 minute data in autoISF” in the autoISF 3.1.0 Quick Guide, <https://github.com/ga-zelle/autoISF>

It is currently too early to tell, but the 1 minute values applied to the parabolic bg curve fit *could* present an avenue to earlier/better acceleration detection than the 5-minute CGMs enable. Actually, 1 minute between two bg readings may sometimes be too short for the algo to go through all calculations + your pump to execute the job at hand:

“My favourite theory for **missed loops** is due to slow pumps. With my original Combo the new V2 driver in AAPS used disconnect/connect between loops which took about 25 sec on average and sometimes much longer before dealing with the real task. With the out-of-date Ruffy driver things were better. So I switched to Dana-i with its much faster communication. But there, the SMB speed is so slow that bottom line the improvement is minor.” (ga-zelle, Discord FCL 6.12.2025)

Tests done prior to introducing the 1 minute Libre3 option (from autoISF 3.0.1 onwards) showed overall comparable but smoother insulin delivery.

Call for a [case study 1.8](#) from a Libre3 user!.

Unfortunately, Libre 3 (1-minute) was not included in the 10 CGM systems comparison we quoted 3 pages above

250 1.5 Meal-related limitations?

251
252 Setting up a full closed loop is relatively easy for people whose diet does not consist **mainly** of
253 components with rapid high effect on blood glucose (more see
254 <https://androidaps.readthedocs.io/en/latest/Usage/FullClosedLoop.html#meal-related-limitations>)
255

256 Meals do not have to be low on carb (provided you use a fast insulin for your FCL)
257 Fat or protein rich diets, or slow digestion/gastroparesis, make things easier rather than harder for
258 the full closed loop because late carbs nicely cover for inevitable “tails” of late action from SMBs
259 needed before or around peak time.

260

261 **Erratic consumption of smaller snacks with fast resorbing carbs can be a problem.**

262 In autoISF you can reduce this problem *to some extent* via one or two keystrokes from your
263 AAPS home screen. While certainly being a deviation from the FCL idea(l), this would be
264 one of the exceptional situations where you better do a quick “nudging” step from your “FCL
265 cockpit”. Details see in [section 5.2.1](#) and [case study 5.2](#)

266

267 **Really, there are no meal limitations.**

268 **The sketched more problematic options just force you to decide which extra “nudging”**
269 **efforts and/or worsened %TIR (and/or “behavioral adjustments”, like less snacking) you are**
270 **willing to accept.**

271

272 1.6 Lifestyle-related limitations?

273

274 Providing a technically stable system

275

276 Full closed looping requires a 24/7 technically stable system, especially regarding

- 277 • reliable CGM signals
- 278 • Bluetooth stability with the pump (see [case study 1.4](#) and [1.6](#))
- 279 • keeping your phone in sufficient proximity at all times
- 280 • avoiding (or at least early recognition of) occlusion (see [case study 1.1](#) and [1.6](#)).

281 This requires a habit (or, unlikely, permanent attention to details) like keeping all components well
282 charged and in close proximity; making cannula (or pod) changes always early enough to lower the
283 risk of occlusion; and having always potentially needed parts with you.

284 Depending on your system, your experience with it, but also on your acceptance and general
285 lifestyle, these aspects may or may not limit you.

286

287 Preparing for exercise

288

289 To prepare for exercise (sports, heavy work), the normal protocol with a pump or hybrid closed loop
290 is to take actions that reduce insulin on board prior to exercise

291 With your full closed loop, the algorithm is tuned to detect meals and to give you insulin to counter
292 glucose rises automatically. Setting a high temp. target and lower %profile, effective already
293 around meal start, could be a problem.

294 Unusual activity levels therefore likely require **disciplined preparation** (especially **if you want to**
295 **keep the need to snack during sports low**)

296

297 In autoISF you can reduce this problem to some extent via two or three keystrokes on your AAPS
298 home screen. While certainly being a deviation from the FCL idea(I), this would be one of the
299 exceptional situations where you better “flick a lever” from your “FCL cockpit” to keep iob low
300 (example see [case study 6.2](#)).

301

302 Extra hurdles to establish FCL for kids

303

304 To establish and maintain a FCL for kids brings about some extra challenges if:

- 305 • Lyumjev is not available or well tolerated
- 306 • Hourly basal rates are very low, providing a poor basis for big SMBs
- 307 • Diet is rich in sweet components. With the typical low blood volume of a small body, strong
308 tendency towards very high bg spikes!
- 309 • Going through marked changes of insulin sensitivity or of circadian pattern makes it difficult
310 to keep the FCL appropriately tuned.

311 This problem is about the same in FCL as it is also in Hybrid Closed Looping. How-
312 ever, now you might expect miracles from the FCL. This is not going to happen. You
313 still must be pro-active by setting suitable temporary % profiles. (These serve as a
314 basis, also for your autoISF FCL, which generally modulates the *at-the-time-valid*
315 profile_ISF, for further enhanced (or reduced) loop aggressiveness).

- 316 • Discipline may be poor regarding keeping 24/7 Bluetooth connectivity
- 317 • Discipline may be poor regarding keeping infusion sites perfectly running.
- 318 • Between kid and supervising parent it must be guaranteed, especially in the initial weeks,
319 that an eye is kept on whether the FCL is working about as expected.

320 More see [section 7](#).

321

322

323 1.7 Time required for setting-up

324

325 Lastly, before enjoying a functioning full closed loop, you need to have a period of a **some weeks**
326 with some free time and „free head“ for set-up –. Can you get, in the time you are willing to invest,
327 to a result that you consider good-enough is really the question. Depending on your „habits“, and
328 which – if any - compromises (like doing cannula/pod changes more often, never starting meals
329 when bg sits high and is not predicted to rush down soon ...) are you willing to make (and
330 everyday able to stick to), for the ease of not having to deal with assessing meals and bolusing for
331 them?

332

333 1.7.1 Recommended structured, step-by-step approach (see also beginning of [section 4](#))

334

335 Setting up your personal FCL using autoISF is a substantial project, for which you should **follow**
336 **the sequence as described** in this e-book.

337

338 But there is **no need to implement it fully in one step**.

- 339 • There is nothing wrong to go in your well running Hybrid Closed Loop mostly, while
340 switching to FCL only for dinners, for instance, or only for weekend lunches, as a start.
- 341 • Once you found feasible settings, you can expand to other meal times...
- 342 • and lastly towards figuring out your best strategies for challenges aside from meal
343 management, as we shall discuss in [sections 5.](#) and [6.](#)

344

345 There are **alternatives** to using autoISF **for FCL**, as well. See [sections 7](#) and [13](#). for more info.

346

347 Notably FCL using AAPS Master and Automations (see in [section 13.1](#)).could be a much
348 easier and more error-tolerant way of stepping into FCL. In a clinical study with 16
349 participants about 80% TIR was achieved without much tuning effort

350

351 Also, watch out for what a 2025 on-going FCL study with AAPS in Australia and NZ will
352 show. See CLOSE IT trial protocol in: *Wilkinson T. et al, BMJ Open 2024;14:e078171. doi:10.1136/bmjopen-2023-078171*

353

354 To close the circle to where we had started; A very time consuming pre-requisite might actually be
355 to **first sort out your Hybrid Closed Loop** ([section 1.1](#)), so your profile parameters are set „right“,
356 and your “old” data really can serve as a blueprint for what, now, you would like *your loop* to do in
357 FCL mode

358 If you feel you better do your homework first there, I highly recommend some of the material in the
359 neighboring HCL repo: <https://github.com/bernie4375/HCL-Meal-Mgt.-ISF-and-IC-settings>

360 1.7.2 “Trial and error” fast track alternative

361

362 Note that if you had used dynamic parameters or special Automations („loops inside the loop“) this
363 might have balanced some principal errors, but leaves you now *without a good starting point*, as
364 now you must get rid of these over-patches (see also warnings at start of [section 4](#))..

365

366 Nevertheless, you will find FCL success stories also from loopers who continue(d) in that spirit, and
367 just jump(ed) into using more powerful tools, in kind of a trial and error mode, frequently adding the
368 latest add-on, or self-constructed patch (often in form of an Automation), to counter-balance
369 encountered problems.

370

371 Resulting solutions may be good-enough.

372

373 However, they tend to be unstable and not well-understood. That is a poor basis for managing
374 arising problems, and for temporarily adjusting to special situations.

375 Nevertheless, this is an alternative avenue for the impatient, less analytically, and more
376 adventurous inclined.

377 Note though, that it is hard to consult (help) someone who, over time, constructed his/her
378 own complicated maze of constructs.

379 Still, I welcome everyone’s creativity. “Not macht erfinderisch” (necessity is the mother of invention)
380 we say in German: Something amazing of broader interest could emerge, and help push
381 innovation for us all.

382

383 Be prepared to needing **a fresh re-start, in case your trial-and-error got you completely lost**.
384 Depending on knowledge level and experience, this could happen on a “fast track” route.

385

386 1.7.3 Safety first

387

388 Regardless which route you choose, PLEASE always observe the **safety** settings/instructions
389 coming with the DIY dev- variant of FCL software you select.

390

391 All FCL methods come with boosted SMBs. So a key safety measure every user going towards
392 any FCL should have in place is to set an **iob threshold** (iobTH; size a bit below what you used as
393 a bolus for bigger meals in HCL) above which no more SMBs can be given by your FCL.

394

- 395 • iobTH is an integrated feature of *autoISF* (see [section 2.4](#)).
- 396 • *Other AAPS-based FCL methods* may require to set up an Automation for a temporary iob
397 threshold that blocks SMBs from being delivered, see e.g. here for AAPS FCL

398 w/Automations; :<https://androidaps.readthedocs.io/de/latest/Usage/FullClosedLoop.html#iob-threshold> .
399

- 400 • *In case there are other methods than autoISF for FCL also on the iAPS or Open iAPS platforms, you*
401 *may have to rely on an adjusted iob_max border, or watch the iob development, and intervene with a*
402 *SMB shut-off, or by opening the loop, when deemed necessary.*

403

404 Also, make use of the easy-to-use feature of **SMB shut-off at odd** profile or temporary **target**. This
405 can **at any time easily** be done manually, via the top right “TT” field in your AAPS screen (set, and
406 time, an odd-numbered target; [section 5.1.3](#)), and can be of enormous help **to temporarily**
407 **safeguard you from aggressive loop actions** (i.e. further growth of iob, no matter how close you
408 already may be to the iobTH)-

409 This is the same concept that you already know from your HCL times, when you wanted to “tame your loop”
410 so it does not “fight” your anti-Hypo Snack with a SMB (An elevated TT> 100 mg/dl was then used, to shut
411 off SMBs in HCL for a while).

412 Lastly, make sure you “train”, for your set-up weeks, how/when to **transition between FCL and**
413 **your prior HCL** mode (refer to [section 5.2.3](#) and [5.1.1](#)).

414